

EUROPEAN SOUTHERN OBSERVATORY

Organisation Européenne pour des Recherches Astronomiques dans l'Hémisphère Austral

Europäische Organisation für astronomische Forschung in der südlichen Hemisphäre

LA SILLA OBSERVATORY

SOFI

User's Manual

Doc. No. LSO-MAN-ESO-40100-0004

Issue 2.5

21/08/2019

Prepared C. Lidman, J-G. Cuby 16/08/2000
Name Date Signature

Revised L. Vanzi; M. Billères 2002
Name Date Signature

Revised V.D. Ivanov 26/04/2006,01/03/2007,05/04/2007,23/05/2011,10/05/2012
Name Date Signature

Revised L.Schmidtobreick; E. Pompei 11/01/2017; 25/08/2019
Name Date Signature

Reviewed I. Saviane 31/08/2019
Name Date Signature

Released I. Saviane 31/08/2019
Name Date Signature

Change Record

Issue/Rev.	Date	Section/Parag. affected	Reason/Initiation/Documents/Remarks
0.9	02/05/98	All	Creation
1.0	14/08/98	All	
1.1	12/11/98	All	New templates
1.2	25/02/99	Some	New Grism
1.3	16/08/00	Some	New IRACE and New Templates
1.4	05/11/02		Some addition and corrections
2.0	26/04/06	All	Merged with the template manual, major addition to the data reduction section
2.1	01/03/07	Some	New Fast Photometry mode, new rotating templates
2.2	23/05/11	All	General updates; general visitor info updated, offset plots redone; all URLs updated.
2.3	10/05/12	Some	updated guiding options, Appendix E, web links, detector non-linearity
2.4	11/01/17	Some	some typos, data backup removed
2.5.	21/08/19	Some	Removed info about decommissioned modes Changes from P2PP to P2

This page was intentionally left blank

Contents

1	Introduction	1
1.1	A First and Final Word	1
1.2	Applicable documents	2
1.3	Abbreviations and Acronyms	2
2	SOFI - Son of ISAAC	4
2.1	Optical Layout	4
2.2	Imaging	4
2.3	Long Slit Spectroscopy	7
2.4	Polarimetry	9
2.5	The DCS - Detector Control System	10
2.5.1	Readout Modes	12
2.5.2	Features on the Detector	13
2.5.3	Windowed Reading	13
2.6	Calibration Unit	14
2.7	Instrument Performance and the Exposure Time Calculator	14
2.8	Instrumental Overheads	15
3	Observing in the IR	16
3.1	The IR Sky	16
3.2	Imaging	16
3.2.1	Selecting the best DIT and NDIT	17
3.2.2	Small Objects or Uncrowded Fields	18
3.2.3	Large Objects or Crowded Fields.	18
3.2.4	Maps of Large Fields	20
3.2.5	Imaging of Moderately Large Object	21
3.2.6	Faint Objects Around Bright Objects	25
3.2.7	Objects with Fast Variability	26
3.2.8	Standard Stars	27
3.3	Polarimetry	29
3.4	Spectroscopy	30
3.4.1	Small Objects and Uncrowded Fields	30
3.4.2	Extended Objects and Crowded Fields	31
3.4.3	Telluric Standards and Flux Calibration	32
3.5	Calibration Frames	33
3.5.1	Darks (Biases)	33
3.5.2	Flat Fields	34
3.5.3	Illumination Corrections	36
3.5.4	Arcs	37
3.6	Finer Points	38
3.6.1	Choosing DIT, NDIT and NINT	38
3.6.2	Autoguiding	39

3.6.3	SOFI Observing Modes	39
3.6.4	Template Parameters - Signatures and Keywords	39
3.6.5	File Naming and Exposure Number	40
3.6.6	Detector Window	41
4	Phase 2 Preparation and Observing with SOFI	42
4.1	General Issues	42
4.2	The VLT environment: P2, BOB, OS, TCS, DCS, ICS	43
4.3	Arriving at the Telescope.	44
4.3.1	Image Analysis	45
4.3.2	Focusing	46
4.4	The SOFI OS GUI Panel	46
4.5	RTD	46
4.6	The Data Flow Path	49
4.6.1	Sofi Pipeline	49
4.6.2	The Archive	50
4.6.3	The Calibration Plan	50
4.7	At the End of the Night and at the End of Your Run	50
5	Data Reduction	51
5.1	Basic Concepts	51
5.2	Imaging	51
5.2.1	Inter-quadrant Row Cross Talk	51
5.2.2	Masking the Bad Pixels	52
5.2.3	Subtracting the Dark/Bias Frame	53
5.2.4	Sky Subtraction	53
5.2.5	Flat Fields and Illumination Corrections	55
5.2.6	Image Alignment and Combination	56
5.3	Long Slit Spectroscopy	57
5.3.1	Inter-quadrant Row Cross Talk	57
5.3.2	Sky Subtraction	57
5.3.3	Flat Fields	57
5.3.4	Removing Slit Curvature	57
5.3.5	Arcs	58
5.3.6	Removing of the Atmospheric Absorption Features and Flux Calibration	58
5.3.7	Alignment and Combination	59
5.4	Polarimetry	59
A	Calibration Arcs	61
B	Atmospheric Absorption	68
C	SOFI Templates: A Reference Guide	70
C.1	General Points	70
C.2	SOFI Imaging Templates	72
C.2.1	SOFI Imaging Acquisition Templates	72
C.2.2	SOFI Imaging Science Templates	74
C.2.3	SOFI Imaging Calibration Templates	88
C.3	SOFI Polarimetric Template:	91
C.3.1	SOFI Polarimetric Acquisition Template	91
C.3.2	SOFI Polarimetric Science Template	91
C.3.3	SOFI Polarimetric Calibration Template	91

C.3.4	SOFI Spectroscopic Templates	94
C.3.5	SOFI Spectroscopic Acquisition Templates	94
C.3.6	SOFI Spectroscopic Science Templates	95
C.3.7	SOFI Spectroscopic Calibration Templates	100
D	Frame Types	105
E	Photometric Standards	106

List of Figures

2.1	Optical layout of SOFI	5
2.2	Orientation of SOFI for rotator angle 0 deg: In imaging the North is to the left and the East to the bottom of the image. The slit for the spectroscopy is oriented North-South. Nota Bene: To align the North-South axis of the field of view (or the slit in case of spectroscopy) long a certain Position Angle, one has to apply in the acquisition template a rotation angle equal to this Position Angle.	6
2.3	SOFI filters. The solid (red) lines indicate the currently available broad-band J , J_s , H and K_s filters. The short-dashed (magenta) lines are the available narrow-band filters. The long-dashed (red) lines show the soon to be commissioned J_s filter (top panel) and the typical broad K filter (Persson et al. 1998, AJ, 116, 2475; see their Table 10). The dotted (blue) line is the atmospheric transmission model for Mauna Kea, for Am=1.0 and water vapour column of 1mm (Lord, S.D. 1992, NASA Technical Memorandum 103957; courtesy of Gemini Observatory). The data for the plot and a SuperMongo script are available from the SOFI web page.	8
2.4	Quantum Efficiency of the SOFI detector at T=78 K. The peak Q.E. is at $1.970\ \mu\text{m}$, and the long wavelength cut-off is at $2.579\ \mu\text{m}$	11
2.5	Deviation from the linearity of the SOFI detector: measured fraction of the expected level for linear response versus the measured illumination level (in log of ADUs) for May 10, 2012.	12
3.1	Example of 4-point observation scheme of a “semi-extended” objects without extra overhead for observation of clear sky. Remember, than although the target is moving on the RTD, it is really the telescope that is moving!	22
3.2	Example of a 6-point observation of “semi-extended” objects alternating between two target fields and one of two different skies.	24
3.3	Examples of Special Dome Flat images. From left to right: lamp off, lamp off with mask, lamp on with mask, lamp on.	36
3.4	Stability of the Special Dome Flat images: accuracy of the photometry as a function of the flat field “age”.	36
4.1	OS of SOFI	47
4.2	Real Time Display, with the “Store Fixed Pattern” option enabled. Note the offset between the negative and the positive star pattern. This offset is necessary to avoid self-cancellation of the objects in the field. The green dot on the top of the image flashes every time when the RTD updates the image. The yellow arrow is a remnant from a previous acquisition – it indicates the direction and the size of the offset that was carried out.	48
A.1	A Xenon arc spectrum taken with the blue grism. The main lines are marked.	62
A.2	A Xenon arc spectrum taken with the red grism. The main lines are marked.	63
A.3	A Xenon and Neon arc spectrum taken with the medium resolution grism at the Z atmospheric window. The main lines are marked.	64

A.4	A Xenon and Neon arc spectrum taken with the medium resolution grism at the J atmospheric window. The main lines are marked.	65
A.5	A Xenon and Neon arc spectrum taken with the medium resolution grism at the H atmospheric window. The main lines are marked.	66
A.6	A Xenon and Neon arc spectrum taken with the medium resolution grism at the K atmospheric window. The main lines are marked.	67
B.1	The atmospheric transmission at a resolution of 8Å.	69
C.1	Positions of the offsets listed in Tables C.42 and C.41 on the detector.	102
E.1	Finding charts for the photometric standards of Persson et al. (1998). I.	109
E.2	Finding charts for the photometric standards of Persson et al. (1998). II.	110
E.3	Finding charts for the photometric standards of Persson et al. (1998). III.	111
E.4	Finding charts for the photometric standards of Persson et al. (1998). IV.	112
E.5	Finding charts for the photometric standards of Persson et al. (1998). V.	113

Chapter 1

Introduction

SOFI or **Son OF ISAAC** is the infrared spectrograph and imaging camera on the NTT. In many ways, it resembles its “parent” ISAAC. They both are focal reducing instruments capable of imaging, spectroscopy and polarimetry.

SOFI offers the following observing modes.

- Imaging with plate scales of 0.273 and 0.288 arc second per pixel in the following modes: Spectroscopic Field and Large Field, respectively; broad and narrow band filters in the wavelength range from 0.9 to 2.5 microns are available.
- Fast imaging and burst mode with detector integration times of order of a few tens of milliseconds, with significantly reduced detector array size (i.e. of order of few tens pixels on the side) via hardware windowing.
- Low resolution $R \sim 600$ (varies across the wavelength range), 0.95-2.52 micron spectroscopy with fixed width slits of 0.6, 1 and 2 arc seconds, and slit length of 4.92 arcmin.
- medium resolution $R \sim 1500$ (varies across the wavelength range) H and K-band spectroscopy, with fixed width slits of 0.6, 1 and 2 arc seconds, and slit length of 4.92 arcmin.
- 0.9-2.5 micron imaging polarimetry with the large field objective (0.288 arcsec per pixel) and the set of filters available in imaging mode.

This manual is divided into several chapters. For proposal writers, it is sufficient to read chapter 2, giving a general description of the instrument. For those who will be observing at the telescope, it is sufficient to read up to the end of chapter 4. Chapter 3 discusses how to use observation templates to set up observations and chapter 4 describes how to use the instrument at the telescope. For those who will reduce data taken with SOFI, it is sufficient to read chapters 2, 3 and 5. Chapter 5 discusses how to reduce data taken with SOFI. At the end, there are appendixes which contain various information useful for observers (detailed template descriptions and examples, fits headers information, calibration data, etc.).

1.1 A First and Final Word

The authors of the manual hope that you find this document useful in writing SOFI proposals or for preparing for your SOFI run. The manual is continually evolving. Some sections still need to be written and some figures still need to be included. If you have any suggestions on how to improve

the manual please contact the La Silla Science Operations department (lasilla@eso.org) or the La Silla coordinator (ls-coord@eso.org).

Nota Bene:

- Some of the SOFI modes have been decommissioned during the life time of the instrument. See the SOFI web page for the latest instrument status. Their description remains in the older versions of the manual for the benefit of users who may download and use old SOFI data from the ESO archive.
- Starting in P105 the seeing constraint used in the proposal preparation will be replaced by a turbulence parameter. For the equivalence, please check the call for proposal and the following link:
<https://www.eso.org/sci/observing/phase2/ObsConditions.html>
- During October 2019, a new tool will replace P2PP for the OB preparation for visitor and service mode runs. This can be found at the following link:
<http://www.eso.org/p2ls>
- The SOFI Template Signature Files (hereafter TSF) Parameters Reference Guide is suppressed by this manual (see Appendix C)!. Those of you who are successful in being awarded SOFI time must read this document and the related P2 Users' Manual before coming to observe with SOFI.
- This document contains a data reduction cook book (Chapter 5) with detailed description of the individual data reduction steps. The users must read it carefully before preparing their observations because the good understanding of the data reduction will help you to plan the observations better and more efficiently.
- There is also a WEB page dedicated to SOFI:
<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/index.html>
It is accessible (go to Instruments) from the main web page of the Science Operations department for La Silla: <http://www.eso.org/lasilla/sciops/>
There you will find the most up-to-date information about the instrument, recent news, efficiency measurements and other useful data that do not easily fit into this manual or a subject of frequent changes. This WEB page is updated regularly.

1.2 Applicable documents

- 1 VLT-MAN-ESO-14100-1510 OS Users' Manual
- 2 VLT-MAN-ESO-14100-1531 DCS Users' Manual
- 3 VLT-MAN-ESO-14100-1094 ICS Users' Manual
- 5 VLT-MAN-ESO-00000-000 SOFI TSF Parameters Reference Guide (obsolete)

1.3 Abbreviations and Acronyms

The abbreviations and acronyms used in this manual are described in Table 1.1.

Acronym	Description
BOB	Broker of Observing Blocks
DCR	Double Correlated Read
DCS	Detector Control System
DEC	Declination
DIT	Detector Integration Time
EFOSC2	ESO Faint Object Spectrograph and Camera 2
EMMI	ESO's Multi-Mode Instrument
ESO	European Southern Observatory
FWHM	Full Width at Half Maximum
ICS	Instrument Control System
IR	Infra-Red
IRACE	InfraRed Array Control Electronics
ISAAC	IR Spectrograph And Array Camera
NDIT	Number of DITs
NDR	Non-Destructive Read
NINT	Number of NDITs
NTT	New Technology Telescope
OB	Observing Blocks
OS	Observing Software
P2PP	Phase 2 Proposal Preparation
P2.	Phase 2: replacement of P2PP
PSF	Point Spread Function
RA	Right Ascension
RON	Read Out Noise
SOFI	Son Of ISAAC
TCS	Telescope Control System
TSF	Template Signature File
USD	User Support Department (formerly User Support Group, USG)
VLT	Very Large Telescope
ZP	Zero Point

Table 1.1: Abbreviations and Acronyms used in this manual.

Chapter 2

SOFI - Son of ISAAC

2.1 Optical Layout

SOFI is mounted on the Nasmyth A focus of the NTT. The light from the tertiary mirror of the telescope enters the front window which has no optical power. Immediately after the front window is a cryogenically cooled mask wheel which coincides with the telescope focus. The mask wheel contains several masks: one for each imaging objective, three for long slit spectroscopy, a pinhole mask and a special mask used with polarimetry.

The optical layout is shown in Fig. 2.1.

The mask wheel is followed by a collimating lens (used to focus the instrument), two filter wheels, a grism wheel, an objective wheel, and then the detector itself. The re-imaged pupil of the telescope, the primary mirror, is located on a slightly undersized stop just before the grism wheel.

The first filter wheel contains the standard broad band near IR filters, several narrow band filters, two order sorting filters for low resolution spectroscopy, an open position and a fully closed position. The second filter wheel contains more narrow band filters, a focus pyramid, an open position and a fully closed position.

The grism wheel contains three grisms for long slit spectroscopy, a Wollaston prism for imaging polarimetry, an open position and a fully closed position.

The objective wheel contains two objectives for imaging at $0''.288$ (large field) per pixel, a spectroscopic objective, an open position and a fully closed position.

2.2 Imaging

SOFI offers imaging at several different pixel scales. The pixel scales and the fields of view are summarised in Table 2.1. The large and small field objectives are used with the corresponding mask in the mask wheel that covers the same field of view. The masks reduce the amount of stray light entering the instrument.

Objective	Pixel Scale	Field of View
Large Field Objective	0.288	4.92' x 4.92'
Spectroscopic Objective	0.273	4.66' x 4.66'

Table 2.1: The fields of view and the pixel scales available with SOFI.

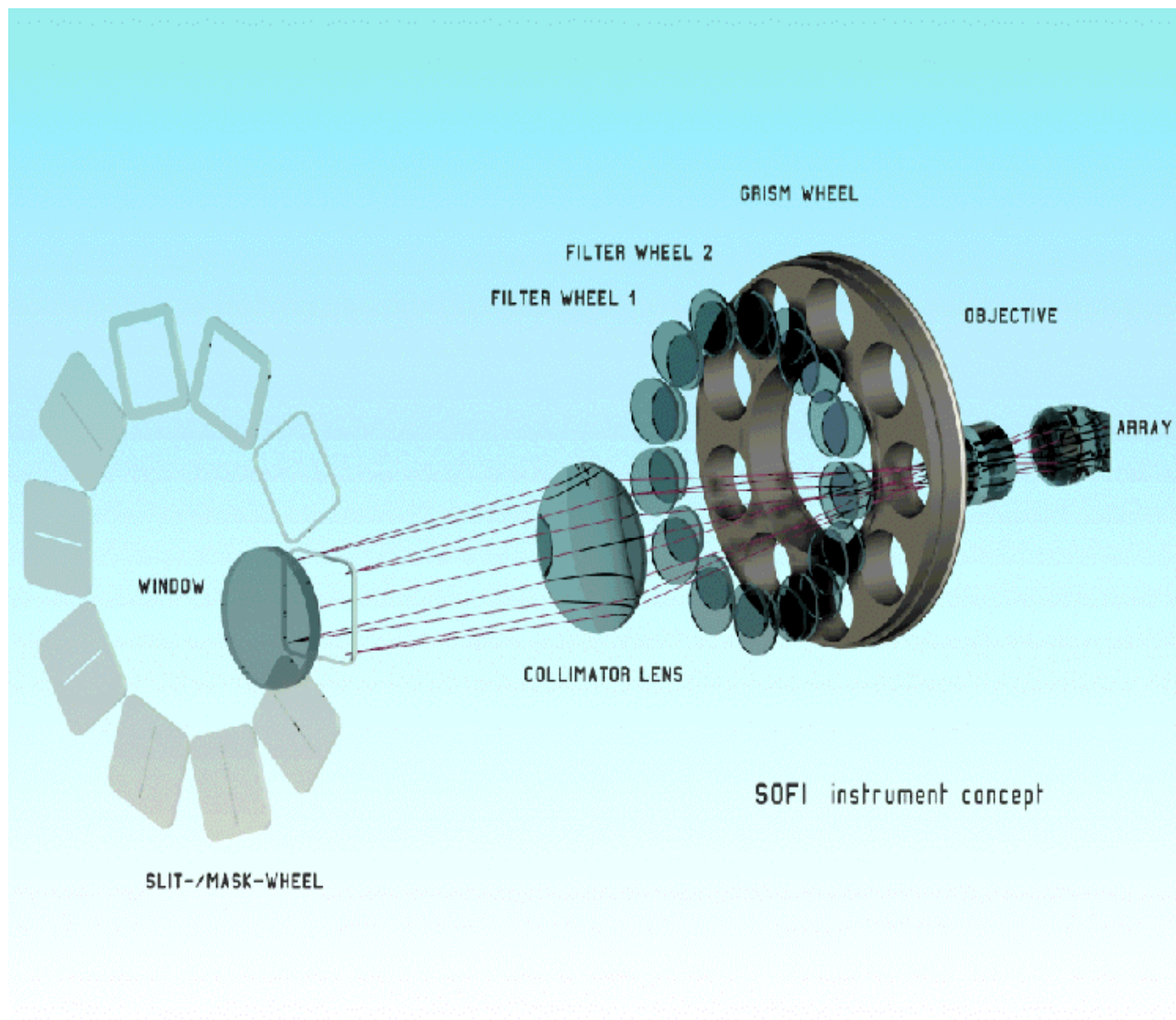


Figure 2.1: Optical layout of SOFI

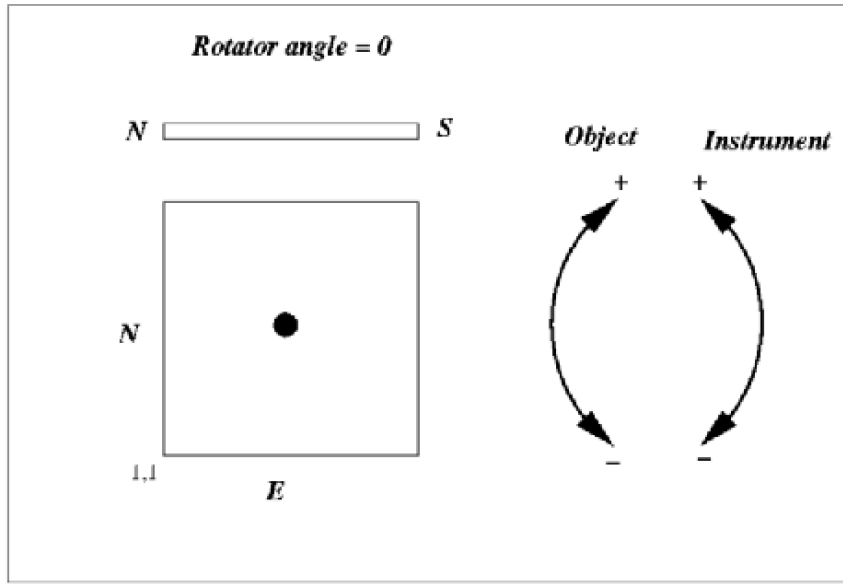


Figure 2.2: Orientation of SOFI for rotator angle 0 deg: In imaging the North is to the left and the East to the bottom of the image. The slit for the spectroscopy is oriented North-South. **Nota Bene:** To align the North-South axis of the field of view (or the slit in case of spectroscopy) long a certain Position Angle, one has to apply in the acquisition template a rotation angle equal to this Position Angle.

For a normal use of SOFI (rotator angle = 0 deg) the orientation is showed in the Fig. 2.2. You could modify the orientation simply using a rotator angle in the template. This option could be useful in spectroscopy if you want to align two objects in the slit (the positive angle is from the north to the east) or in imaging if you want to map some asymmetric target, i.e. an elongated galaxy. A simple rule of thumb says: to align the North-South axis of the field of view (or the slit, in case of spectroscopy) along an axis on the sky with a certain Position Angle, one has to apply in the acquisition template a rotation angle equal to this Position Angle.

Starting from Apr 1, 2007, new templates are available to the users, that allow them to rotate the instrument between the individual pointings in a single template, as opposed to setting one position angle and acquiring all the data in the OB at that orientation. These templates are particularly useful for observations of fields containing bright sources that would cause significant reflections and cross talk artefacts. Due to the rotation, the artefacts will be located on different place of the detector array. Taking data at enough different position angles, derotating/aligning and combining the images will allow the users to reduce the artefacts. The exact number of position angles will depend on the number and the brightness of the bright sources. Note that the movement of the rotator will cause extra overhead.

The spectroscopic objective can be used to take images with a scale that is a bit finer than that of the large field objective; however, this objective is chromatic, as mentioned above. This means that the pixel scale is a function of wavelength. It also has an illumination pattern with a central light concentration.

All the aforementioned objectives suffer various amounts of image degradation over their respective fields of view. We routinely try to improve the image quality during interventions. Please consult the instrument web pages or the instrument scientists for the latest status of the instrument.

Imaging can be done through the standard IR broadband filters: J , J_s , H , K_s , a Z filter which peaks

at 0.9 microns, and many narrow band filters. The filters, together with their central wavelengths and widths at half maximum, are listed in Table 2.2. Listed also are the cut-on wavelengths for the two order sorting filters. These filters are used with the two low resolution grisms. Transmission curves of imaging filters, together with the atmospheric transmission are shown in Fig. 2.3 (see also Appendix B).

Filter Name	Filter Wheel	Central Wavelength (μm)	Width (μm)	Peak Transmission (%)
<i>Z</i>	1	0.9	0.140	
<i>J</i>	1	1.247	0.290	
<i>J_s</i>	1	1.240	0.160	
<i>H</i>	1	1.653	0.297	83
<i>K_s</i>	1	2.162	0.275	88
NB 1.061	1	1.061	0.010	
NB 1.083 HeI J	2	1.083	0.016	61
NB 1.187	1	1.187	0.010	
NB 1.215	2	1.215	0.018	
NB 1.257 [FeII] J	2	1.257	0.019	
NB 1.282 P β	2	1.282	0.019	
NB 1.644 [FeII] H	2	1.644	0.025	
NB 1.710	2	1.710	0.026	
NB 2.059 HeI K	2	2.059	0.028	81
NB 2.090	1	2.090	0.020	
NB 2.124 H ₂ (S1)	2	2.124	0.028	78
NB 2.170 Br γ	2	2.167	0.028	71
NB 2.195	1	2.195	0.030	
NB 2.248	2	2.248	0.030	
NB 2.280	2	2.280	0.030	69
NB 2.336 (CO)	2	2.336	0.031	
GBF	1	0.925 cut-on	—	
GRF	1	1.424 cut-on	—	
Open	1 and 2	—	—	—
Closed	1 and 2	—	—	—

Table 2.2: The broad and narrow band filters available with SOFI.

The *K* short or *K_s* filter is different from both the standard *K* filter and the *K'* filter defined by Wainscoat and Cowie (1992, AJ, 103, 332). The long wavelength edge of the *K_s* filter is similar to that of the *K'* filter, but the short wavelength edge is similar to that of the *K* filter. Thus, the *K_s* filter avoids both the atmospheric absorption feature at 1.9 μm and radiation from the thermal background beyond 2.3 μm . The difference between *K_s* and *K* is given by $K - K_s = -0.005(J - K)$.

Similarly, the *J_s* and *J* filters differ mostly in the long wavelength edge - the transmission of the broader *J* filter is in practice limited by the atmosphere and therefore, the photometry in the *J* filter is less stable than with the narrower *J_s* one. However, the *J_s* filter allows about $\sim 20\%$ less photons to reach the detector than the *J* filter.

2.3 Long Slit Spectroscopy

SOFI offers low and medium resolution, long slit spectroscopy. The large field objective is used in this mode instead of the spectroscopic one because the former is achromatic. There are three grisms

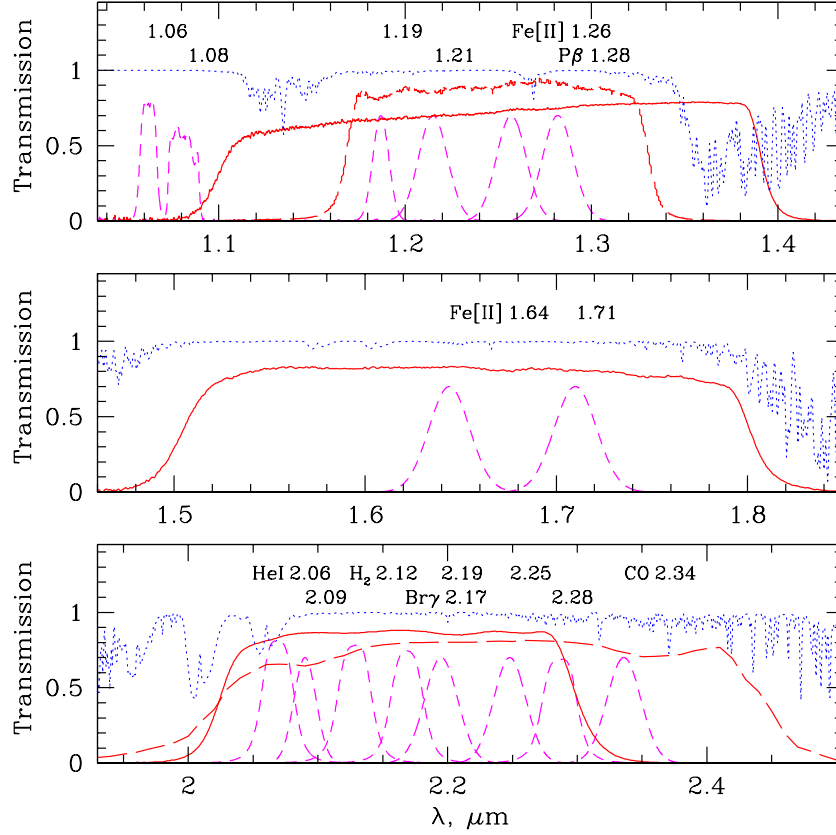


Figure 2.3: SOFI filters. The solid (red) lines indicate the currently available broad-band J , J_s , H and K_s filters. The short-dashed (magenta) lines are the available narrow-band filters. The long-dashed (red) lines show the soon to be commissioned J_s filter (top panel) and the typical broad K filter (Persson et al. 1998, AJ, 116, 2475; see their Table 10). The dotted (blue) line is the atmospheric transmission model for Mauna Kea, for $A_m=1.0$ and water vapour column of 1mm (Lord, S.D. 1992, NASA Technical Memorandum 103957; courtesy of Gemini Observatory). The data for the plot and a SuperMongo script are available from the SOFI web page.

available with SOFI: two low resolution grisms and a medium resolution grism.

Of the two low resolution grisms, one covers the region from 0.95 to 1.64 microns and the second covers the region from 1.53 to 2.52 microns. Both grisms are made of KRS5 (refractive index ≈ 2.44) and the entrance surface of both is inclined to the optical axis. The blaze angle of the grooves is equal to the apex angle of the prism. The wavelength ranges, the resolving powers and the resolutions of the grisms are given in Table 2.3.

The medium resolution grism gives about twice the resolving power of the two low resolution grisms. It is used with the H and K_s filters as order sorting filters in the 3rd and 4th orders to cover respectively the H and K atmospheric transmission windows. The wavelength ranges, the resolving powers and the resolutions of the grisms are given in Table 2.4. The wavelength ranges are defined by the H and K_s filters.

The grism can also be used in the higher orders with the J and Z filters, However, there is significant overlap between these orders, so the useful wavelength range is limited. Furthermore, the line profile degrades in the blue, so that the resolution is not significantly better than that obtained with the low resolution blue grism.

Grism Number	Order Sorting Filter	Wavelength Range (microns)	Resolving Power	Dispersion (\AA)/pixel
Blue	GBF	0.95-1.64	930	6.96
Red	GRF	1.53-2.52	980	10.22

Table 2.3: The wavelength range and the resolution of the low resolution grisms. The resolution is that measured for the 0.6 arc-second slit. The resolution scales inversely with the slit width.

Grism Order Name	Order Sorting Filter	Wavelength Range (microns)	Resolving Power	Dispersion (\AA)/pixel
3	K_s	2.00-2.30	2200	4.62
4	H	1.50-1.80	1500	3.43
5	J	1.20-1.28	1400	2.71
6	J	1.17-1.24	1400	2.22
7	Z	0.89-0.93	1400	1.87
8	Z	0.86-0.95	1400	1.58

Table 2.4: The wavelength range and the resolution of the medium resolution grism. The resolution is that measured for the 0.6 arc-second slit. The resolution scales inversely with the slit width.

Three slits of different fixed widths, 0.6, 1.0 and 2.0 arc seconds, are available. The slit length is 290 arc seconds.

2.4 Polarimetry

SOFI offers imaging polarimetry. A Wollaston prism in the grism wheel splits the incoming parallel beam into two beams that are perpendicularly polarised. The beams are separated by 48 arc seconds. Thus, an image taken with the Wollaston prism will contain two images of every object. To avoid sources overlapping, a special mask, consisting of alternating opaque and transmitting strips, can be inserted at the focal plane. Therefore, in a single exposure, at least half the field will be missing. So three exposures, with telescope offsets in between, are required to cover one field.

The Wollaston prism is not achromatic, so the exact separation between the two beams is a function of wavelength. At J , the separation is 48.3 arc seconds, while at K_s , the separation is 47.4 arc seconds. The beam separation is also a function of position.

To measure the Stokes parameters and hence the degree and position angle of polarisation a second set of images with the Wollaston prism rotated 45 degrees with respect to the first pair are required. This is achieved by rotating the entire instrument. The Stokes parameters are then determined as follows:

$$I = i(90) + i(0) = i(45) + i(135)$$

$$Q = i(0) - i(90)$$

$$U = i(45) - i(135)$$

where $i(\alpha)$ is the intensity of the source which transmits light that is polarised at angle α . We have assumed that the rotator is at a position angle of 0 degrees for the first measurement. This need not

be the case.

The degree of linear polarisation and the polarisation angle are given by;

$$P = \frac{\sqrt{U^2 + Q^2}}{I}$$

$$\theta = 0.5 \times \tan^{-1} \frac{U}{Q}$$

To derive the correct value of θ , attention needs to be paid to the signs of U and Q .

This algorithm neglects instrumental polarisation. Preliminary measurements indicate that the instrument polarisation is 2%. As it is caused by the tertiary mirror, the vector defining the instrument induced polarisation will rotate relative to the sky. A method to eliminate the instrumental polarisation is outlined by Sperello di Serego Alighieri (1989, Proceedings of 1st ESO/ST-ECF Data Analysis Workshop).

Technical reports (Wolf, Vanzi & Ageorges 2002) have been released, which describe in details this mode and its operation. They are available on the SOFI web page:

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/report.ps>

2.5 The DCS - Detector Control System

The DCS is made up of the detector, the front-end electronics, and the controller (InfraRed Array Control Electronics or IRACE).

The IRACE controller controls the detector front-end electronics and manages pre-processing of the data before sending it to the SOFI workstation. It includes an embedded Sparc and two sets of transputers. Further details on the IRACE system can be found at:

<http://www.eso.org/sci/facilities/develop/detectors/controllers/irace.html>

The amount of data pre-processing depends on the readout mode. The readout modes available with SOFI are discussed below.

The detector used by SOFI is a Rockwell Hg:Cd:Te 1024x1024 Hawaii array manufactured by Rockwell Scientific, with 18.5 micron pixels. The array is read out in four quadrants. The average quantum efficiency is 65 % (Fig. 2.4). The dark current is very low, 20 e⁻/hour, and the readout noise with the IRACE controller in DCR (Double Correlated Read) mode is 12 e⁻. In NDR (Non Destructive Read), values as low as ~3 e⁻ have been reached with integration of one minute. About 0.1% of the pixels are bad. An updated mask of bad pixels is available in the SOFI web page.

The gain of the array ~5.4 e/ADU. The well depth of the array is around 170,000 electrons (32,000 ADU). Although the array non-linearity is limited to less than 1.5 % for a signal up to 10,000 ADU (Fig. 2.5), we recommend that observers keep the exposure short enough so that the *background* does not exceed 6,000 ADU. This is due to the bias of the array, which has a complicated dependence on the flux when the flux is above 6,000 ADU.

NOTA BENE: A star of K~10 mag produces 9,000-10,000 ADU at the central pixel under seeing ~0.6-0.7 arcsec and average conditions, in Large Field imaging mode. This is the case with many of the photometric standards and since these values are close to the non-linearity limit, the user may consider defocusing the telescope for standard star observations or switching to one of the imaging modes with finer pixel scale for observations that can not

The DCS has to cope with backgrounds that range from a fraction of a ADU/sec/pix, as seen in the spectroscopic modes to six hundred ADU/sec/pix, as seen during summer with the K_s filter and the

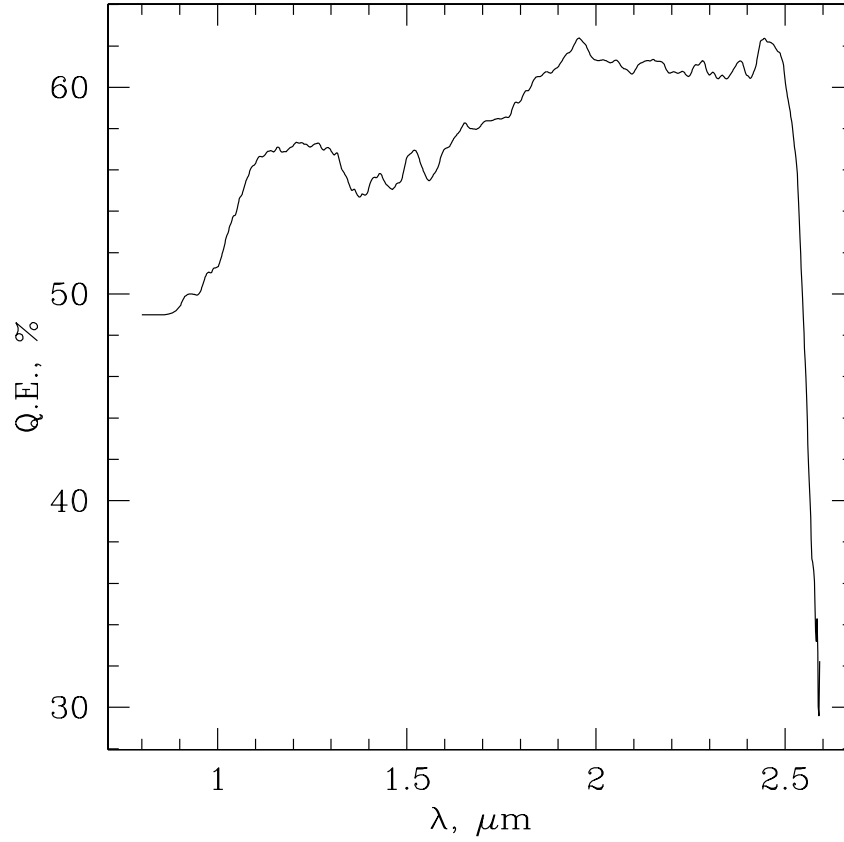


Figure 2.4: Quantum Efficiency of the SOFI detector at $T = 78$ K. The peak Q.E. is at $1.970 \mu\text{m}$, and the long wavelength cut-off is at $2.579 \mu\text{m}$.

large field objective. The high count rate from the background, particularly at K_s , limits a single integration to about 10 seconds, and it can drop down even to 6 seconds, depending on the humidity, clouds cover, etc. Thus, in order to accumulate sufficient photons without saturating the detector, the observer should acquire many individual integrations and then average them on the fly into a final frame. The control system allows to do exactly that.

In this context, we define DIT, the *Detector Integration Time*, as the amount of time during which the signal is integrated onto the detector diodes, and NDI as the number of detector integrations that are obtained and **averaged** together. These averaged frames make up the raw data, and in normal co-add mode, they are the smallest block of data presented to the user. So that the total exposure time of a single image (= single raw data file, available to the user) is $\text{NDI} \times \text{DIT}$. Note that the RTD (Real Time Display) can show either a single DIT or the averaged $\text{NDI} \times \text{DIT}$.

Additionally, some of the science templates allow the user to obtain several consecutive exposures of $\text{NDI} \times \text{DIT}$, as defined by NINT and NJIT parameters.

NOTA BENE: The counts in a raw data file always correspond to DIT seconds. However, a single raw data file represents total integration time of $\text{NDI} \times \text{DIT}$, because the counts in the file are the average of NDI sub-integrations, each of DITseconds!

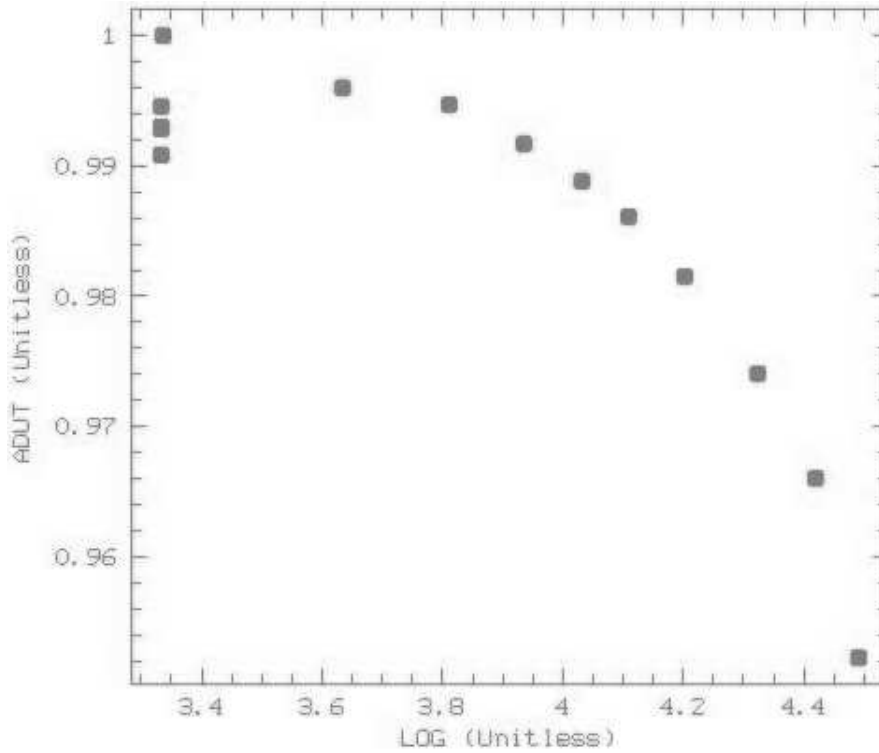


Figure 2.5: Deviation from the linearity of the SOFI detector: measured fraction of the expected level for linear response versus the measured illumination level (in log of ADUs) for May 10, 2012.

2.5.1 Readout Modes

Unlike optical CCDs, the charges in individual pixels of infrared detector arrays are not shifted from pixel to pixel during the read-out process. Instead, each pixel is independently read and each column is individually reset. This enables one to develop several readout methods. For SOFI, two readout methods are available: DCR (Double Correlated Read) and NDR (Non-Destructive Read). DCR incurs less overheads and is suitable for modes where the dominant source of errors is the Poisson noise of the sky emission: imaging and low resolution spectroscopy. NDR has a lower read-out noise and is suitable for modes for which the read-out noise is comparable to the other sources of error: low and medium resolution spectroscopy.

- **Double Correlated Read:** Here the voltage is sampled twice, once at the beginning of the integration and a second time at the end of the integration. This method is called Double Correlated Read or DCR for short. It is commonly used in high background situations where integrations are forced to be short. In principle this method is susceptible to $1/f$ noise; however, the Poisson noise from the background will be, by far, the most important noise source. The minimum integration time for SOFI in this mode is 1.182 seconds.
- **Non-Destructive Read:** In Non-Destructive Read or NDR the array is sampled several times after the reset. We will call the number of samples **NSAMP**. In other words, the array is read **NSAMP** times during each individual DIT. The flux in each pixel is computed by fitting a linear function to the voltage as a function of time. The fitted slope – equivalent to the photon rate – is then multiplied by the integration time. There are several ways the signal can be read in the NDR mode. The one used in SOFI concentrates the read-out at the beginning and at the end of the integration, in this way the noise is minimised. This mode is called Fowler sampling. Unfortunately, with the number of readings the glowing produced by the heat from

the shift-registers at the border of the array increases. For $\text{NSAMP} = 60$ the photon noise of the glowing starts to compete with the read-out noise and above 60 becomes the dominant source of noise. For this reason the **NSAMP must not exceed 60**. To further minimise the noise, keeping **NSAMP** small, the analogue signal can be sampled several times for each reading, this is done with the parameter **NSAMPIX**. The read-out noise is reduced approximately by the square root of the number of samples **NSAMP**. This method is also less susceptible to 50 Hz pickup noise. The drawback of the NDR is that it takes slightly longer to process the data although the disadvantage becomes smaller as the exposure time increases. The shortest integration time in this mode is given by $1.64 \times \text{NSAMP}$. NDR is only available for the spectroscopic modes. For most applications we recommend $\text{NSAMP} = 30$ and $\text{NSAMPIX} = 4$.

There is a noiseless component of the signal introduced by the reset which is somewhat unstable a short time after the reset. Thus, there is a small delay of 100 ms between the reset and the first read. This component decays with time, for this reason dark frames taken with long DIT have negative counts.

The Hawaii array, like the NICMOS III array is read out simultaneously in four quadrants. This leads to a characteristic jump in the count level at rows 1 and 513. This jump depends on the DIT and on the incident flux on the array. The jump is stable with time, and it disappears after the sky subtraction.

2.5.2 Features on the Detector

There is a number of features of the detector. The most prominent ones are a few circular depressions in the quantum efficiency with radius of ~ 20 -40 pixels. Their presence is not a source of concern because they flat-field out. In addition, there are some variable features, i.e. dark regions caused by dust particles that have fallen inside the instrument. Rotation of the instrument can remove them.

Occasionally, the detector array shows elevated noise pattern. The effect is probably caused by interference and despite the replacement of most electronic components, it remained unreproducible. Usually, it is confined to a single quadrant. The noise amplitude is only a few ADU and it adds negligible error in imaging or low-resolution spectroscopy mode. However, in high-resolution spectroscopy mode, especially when observing faint targets, it is advisable to select nodding and jittering parameters placing the target in the unaffected quadrants.

2.5.3 Windowed Reading

The SOFI detector can be windowed. Each window is defined by the starting pixel coordinates and the size of the windowed region. The entire array is still read out by the IRACE controller; however, only the windowed section is transferred to the workstation. This leads to a slight overhead decrease. For position angle $\text{PA} = 0^\circ$ starting coordinates of (1,1) correspond to the North-East corner and they increase toward South-West, respectively.

To reduce further the overhead, one has to resort to hardware windowing. This is available only for the Fast Photometry / Burst observations. In this case only part of the array is read and transferred and the significant reduction of the overheads comes from reading a smaller fraction of the array. To take advantage of having four Analog-to-Digital Converters, the windowed part of the array is placed at the centre of the detector (512,512) so the four ADC's can be used simultaneously. The window size determines both the minimum DIT and the readout time.

2.6 Calibration Unit

The SOFI calibration unit is located inside the telescope adaptor. It contains a halogen lamp for internal spectroscopic flat fielding, Xenon and Neon lamps for wavelength calibration. Line identifications are given in the Appendix A.

The halogen lamp is used to flat field spectroscopic data. These flats are called Nasmyth flats. Note that the Nasmyth flats introduce an extra spatial response by the lamp, that has to be removed by additional night sky spectroscopic flat.

The halogen lamp is **not suitable for flat fielding of imaging data**. Use the dome flat field lamp for this (see Sec. 3.5.2 for details).

2.7 Instrument Performance and the Exposure Time Calculator

Table 2.5 we list the approximate zero-point of the broad band filters with the large field objective as measured over many nights starting from July 1998 to May 2005. Limits are similar for the small field and spectroscopic objectives. For more recent detection limits in the spectroscopic modes and for detection limits with the refer to the SOFI web page.

Filter	ZP	Average Background		Detection Limit	
		Mag.	sq. arc second	Point Source	Extended Source
<i>Z</i>	22.6
<i>J</i>	23.2	15.5-16.1	...	22.7	22.1
<i>J_s</i>	23.1
<i>H</i>	23.0	13.4-14.7	...	21.8	21.4
<i>K_s</i>	22.4	12.8-13.3	...	20.8	20.4

Table 2.5: Measured SOFI performance for the broad band filters: 1 hour exposure.

The point source detection limits are based on the following assumptions:

- Signal to noise ratio of 5, computed over 21 pixels,
- Pixel scale: 0.288 arc seconds/pixel.
- 1 hour exposure made up of 60 one minute integrations.
- Seeing 0.75 arc seconds.
- Airmass 1.2
- Backgrounds of 16.0, 16.0, 14.2 and 13.0 in *J*, *J_s*, *H* and *K_s*, respectively.

The extended source detection limits were made for an aperture with a diameter of 3 arc seconds and for a S/N of 4.

The values of Table 2.5 can be re-scaled to different S/N ratios, fluxes (*F*) and integration times (*t*) keeping in mind that for background limited performances $S/N \propto F \times \sqrt{t}$.

Background limited performances are reached when $ADU \gg 30$, assuming read-out noise $12 e^-$ in the double correlated read-out mode, and a gain of about $5.4 e^-/ADU$.

NOTA BENE: The detection limits given above vary with the background level which is strongly sensitive to the humidity and and air temperature. The background can easily change within 50%.

2.8 Instrumental Overheads

The fraction of time spent not collecting photons is defined as the instrumental overhead. A good conservative rule of thumbs states that it is typically about 30% of the integration time.

There are several sources of telescope and instrument overhead:

- acquisition overhead which depends on the instrument mode, the selected acquisition template, the brightness of the target, the angular distance from the previous target; the simple imaging preset takes on average 3 min and all other presets - 6 min; the more challenging spectroscopic acquisitions of faint targets can require up-to 5-10 minutes
- the time necessary to offset the telescope between different jitter positions and any remaining time to complete the last DIT that was started during the telescope offset; recall that the sequencer is not interrupted unless the DIT is changed; in case of guiding there is an extra overhead necessary to reacquire the guiding star; therefore, the intervals between the offsets have to be as long as possible, to minimise these losses and this can be achieved with increasing NDIT (i.e. $\text{NDIT} \times \text{DIT} \sim 2\text{-}3$ min for imaging, in case of excellent sky conditions), but not too long as to compromise the sky subtraction; for practical purposes the average telescope movement overhead is estimate to be 15 sec
- there is a 0.1 second delay between the reset of the array and the first read at every DIT and 1.7 sec to read the array; this implies that to minimise the overhead the user should try to use as long DIT as possible (but not too long as to saturate the detector or to work in the non-linear regime above 10000 ADU); therefore, each DIT “costs” 1.8 sec of overhead, if the entire detector is being read
- the movement of the telescope rotator (that is used to change the orientation of SOFI) is relatively slow but it can not be quantified because it is never known in advance what would be the preceding rotator position and how much it will have to move, so this overhead is simply included in the average acquisition time; in cases of templates that rotate SOFI between exposures it is assumed that the angle between the individual exposures will be always small, i.e. less than 15 deg, which takes less than ~ 10 sec; typically this adds negligible extra time to the time that the telescope takes to move and so it is ignored.

It is difficult to give a precise estimate of the instrumental overhead, but here is a common example. A one hour exposure made up of 60 one minute, unguided exposures (with typical for H , K or K_s DIT = 10 and NDIT = 6) with telescope offsets for jittering in between will take a total of ~ 80 minutes. Guided exposures will take 35 to 50% overheads.

Chapter 3

Observing in the IR

3.1 The IR Sky

Observing in the IR is more complex than observing in the optical. The difference arises from a higher and more variable background, by stronger atmospheric absorption and telluric emission throughout the 1 to 2.5 micron wavelength region.

Short-ward of 2.3 microns, the background is dominated by non-thermal emission, principally by aurora, OH and O₂ emission lines. The vibrationally excited OH lines are highly variable on a time scale of a few minutes. Pronounced diurnal variations also occur. The lines are strongest just after sunset and weakest a few hours after midnight. A complete description and an atlas of the sky emission lines can be found in the paper by Rousselot et al. (2000, A&A 354, 1134).

Long-ward of 2.3 microns, the background is dominated by thermal emission from both the telescope and the sky, and is principally a function of the temperature. The background in K_s can vary by a factor of two between the winter and summer months but is more stable than the J or H band background on minute-long time-scale. It also depends on the cleanliness of the primary mirror. Imaging in broadband K_s and the wide field objective can result in backgrounds of 600-700 ADU/sec, depending strongly on the temperature and humidity.

The IR window between 1 and 2.5 microns contains many absorption features that are primarily due to water vapour and carbon dioxide in the atmosphere. These features are time varying and they depend non-linearly with airmass. The atmosphere between the J and H bands and between the H and K bands is almost completely opaque. The atmospheric transmittance between 1 and 2.5 microns as seen by SOFI is plotted in Appendix B. As the amount of water vapour varies so will the amount of absorption. The edges of the atmospheric windows are highly variable which is important for the stability of the photometry in J and K_s filters.

These difficulties have led to the development of specific observing techniques for the IR. These techniques are encapsulated in the templates that are used to control SOFI and the telescope. In this section we link common observational scenarios with specific templates and we give some concrete examples. In later sections we discuss some finer points. The templates are described in detail in The SOFI TSF Parameters Reference Guide.

3.2 Imaging

It is not unusual for the objects of interest to be hundreds or even thousands of times fainter than the sky. Under these conditions it has become standard practice to observe the source (together with the inevitable sky) and subtract from it an estimate of the sky. Since the sky emission is generally

variable, the only way to obtain good sky cancellation is to do this frequently. The frequency depends on the wavelength of observation (and respectively on the nature of the sky background emission) and on meteorological conditions. Ideally, one would like to estimate the sky more quickly than the time scale of the sky variations. While this could be done quickly with the traditional single- and especially double-channel photometers, the overhead in observing with array detectors and the necessity of integrating sufficient photons to achieve background limited performance are such that the frequency is of the order of once per minute. This sky subtraction technique has the additional advantage that it automatically removes offsets due to fixed electronic patterns (bias) and dark current.

NOTA BENE: The *sky* and the *object+sky* have to be sampled equally; integrating more on the *object+sky* than on the *sky* will not improve the overall signal-to-noise ratio because the noise will be dominated by the *sky*.

There are two standard techniques to estimate the sky. The first is appropriate for angularly large objects or crowded fields and the second is appropriate for angularly small (in comparison with the field of view) objects or uncrowded fields.

3.2.1 Selecting the best DIT and NDIT

Selecting the best DIT and NDIT is a complex optimisation problem and it depends on the nature of the program: targets, required signal-to-noise, frequency of sky sampling, etc. Therefore, it is hard to give general suggestions and the users should exercise their judgement and discuss their choices with the telescope operator who is supporting their run.

The first constraint is to keep the signal from the target on the linear part of the detector array dynamic range, which is below 10,000 ADU. The minimum DIT of 1.182 sec allows to observe without problem stars of ~ 10 mag, under average seeing/humidity conditions. Keeping in the linear regime any bright stars that may have randomly fallen into the field is also desirable because these stars will become a source of horizontal lines due to the cross-talk. The cross-talk is correctable well only if the sources that cause it are not saturated.

Unfortunately, it is not always possible to use small DIT, because the smaller DIT increases greatly the overhead – up to 50%, in the cases of 1.2-2 sec DIT. For comparison, observations with DIT=8 sec and NDIT=8 or 64 sec of integration in total, require about 80 sec so the overhead is 25%. The user may choose to accept the large overheads or the presence of increased “lines” in the image, that are hard (but not completely impossible) to correct.

The sky background is another factor that has to be accounted for when selecting a DIT. It is the strongest in K_s band when it amounts to 400-800 ADU per second, depending strongly on the humidity. The sky background can easily saturate the array by itself if the user selects a big DIT of 15-30 sec or more, depending on the filter.

Furthermore, this background is not constant. It varies on a time scale of a 1-3 minutes becoming a source of systematic uncertainties. To account for them the user must monitor these variations on the same time scale. This is done by alternatively observing the target and a clear sky field next to the target (the next sections contain some useful tips how to optimise these sampling in case of different types of targets). The frequency of the sky sampling is determined by the product DIT \times NDIT, plus the overhead. The user should try to keep this time close to 1-2 min in average sky conditions or ~ 3 min in exceptionally good conditions. To verify the choice of sky sampling the observer should frequently subtract sequential images from one another and monitor how large the average residual is. Ideally, it should be smaller than or comparable to the expected Poisson noise but this is rarely the case. Usually, a few tens or hundreds of ADU are considered acceptable by most users.

Finally, the total integration time is accumulated by obtaining a certain number of images, usually

specified by the total number of exposures and/or the number of jittered images at each position, depending on which template is used. In case when relatively large integrations are necessary, it is simply a matter of increasing the number of exposures. However, in the cases when the total required time can be accumulated in less than 5-7 exposures it might become difficult to create a good sky for the sky subtraction, especially if the field is crowded because the sky image may contain residuals from the stellar images that will produce “holes” in the sky-subtracted data. This situation will require to increase artificially the number of exposures to 5 or 7. It might be possible to compensate this increase by decreasing correspondingly NDIT to keep the total integration time constant. Still, there will be some increase in the overheads.

Summarising, under average conditions, for faint targets, one can safely use DIT=30, 15 and 10 sec for J , H and K_s filters, respectively. The narrow band filters can tolerate DIT of 60-120 sec. The array level will be dominated entirely by the sky. Brighter targets require to reduce these times down to 12, 8 and 6 sec for the broad band filters, for 12-15 mag stars, and all the way down to the minimum DIT of 1.182 sec for ~ 10 mag stars. The users may even have to consider splitting their observations into “shallow” and “deep” sequences, optimised for different magnitude ranges.

These issues are discussed again in Sec. 3.6.

3.2.2 Small Objects or Uncrowded Fields

If the object of interest is an uncrowded stellar field, it is not necessary to take separate sky observations. In this case one can dither within the field and use the object frames to create sky frames. As a rule, the offsets should be greater than 10 arc-seconds, and if very deep exposures are required, the offset vector should not be replicated. For example, an offsetting scheme that is based on a rectangular grid of points will in a deep exposure show faint negative images arranged symmetrically around each real image. As a last resort, these can be minimised with appropriate masking during the data reduction.

Similar to the uncrowded fields, small objects can be observed in the same manner. In this case the user should make sure that the offsets are at least 2-3 times bigger than the size of the object.

A pseudo-random offsetting scheme (see sec. 3.6) is used in the template **SOFI_img_obs_AutoJitter**. The offsets are restricted to be within a square box centred on the object. The dimension of the box is defined by the parameter “Jitter Box Width.” This and other parameters of the template are listed in Table 3.1.

If the value of the “Return to Origin” parameter is true (T), the telescope, at the end of the template, moves back to where it was at the start of the template. In general, auto-guiding is not required because of the frequent telescope offsets.

For a more detailed discussion about guiding options and the algorithm used to compute the offsets, please refer to section 3.6.

If you wish to enter the offsets individually, use the template **SOFI_img_obs_Jitter**. This template is discussed fully in Appendix C.

3.2.3 Large Objects or Crowded Fields.

For objects larger than $\sim 40\%$ of the field or for very crowded fields, it is necessary to image the sky and object separately. Unfortunately, it is common that the sky frames will contain other objects, and it is not uncommon that one of these objects will be in the same region of the array as the science object. To avoid this it is standard to obtain several sky images on different locations, usually in the context of object-sky pairs. This technique assumes that the sky fields are sufficiently uncrowded

Parameter signature	Value
Exposure Name	Hubble_Deep_Field
DIT	10
NDIT	6
Number of exposures	60
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Filter wheel 1	K_s
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined Offset? (T/F)	F
Jitter Box Width (arcsec)	40
Return to Origin? (T/F)	T

Table 3.1: Parameters of the SOFI_img_obs_AutoJitter template with commonly used values.

that on any given position of the array most of the images will have sky, and only a minority will have objects. Clearly, minimum three sky images are necessary for this technique but the experience shows that a reasonable minimum number of the sky images (and respectively, the object-sky pairs) is 5-7, to ensure a good removal of the objects from the sky frames. Note that this may lead to an extra overhead because in some cases the NDIT has to be reduced artificially (contrary to the optimisation strategy discussed in Sec. 2.8) to a number below the optimal, just to split the total integration into 5-7 images, adding an extra overhead for the telescope offsets.

This technique is encapsulated in the template called **SOFI_img_obs_AutoJitterOffset**. The parameters of this template are listed in Table 3.2. While preparing such observations, please keep in mind that the total duration of an imaging OB *in Service Mode* can not exceed 1 hour (in Visitor Mode the duration of an OB is not constrained). The example corresponds to the maximum duration acceptable.

Parameter signature	Value
Exposure Name	NGC6118
DIT	10
NDIT	6
Number of exposures	36
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Filter wheel 1	K_s
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined Offset? (T/F)	F
Jitter Box Width (arcsec)	20
Sky Offset Throw (arcsec)	600
Rotate Pupil? (T/F)	T

Table 3.2: The parameters of the template SOFI_img_obs_AutoJitterOffset with typical values.

In this template, one observes alternatively the object and sky. There are 36 exposures in total,

18 on the object and 18 on the sky. The sky positions are randomly chosen to lie on a circle that is centred on the object. The diameter of this circle is set by the parameter “**Sky Offset Throw.**” Each of the six sky positions is different. Each of the six object positions are also different – they are randomly distributed within a square box that is centred on the object. The box size is set by the parameter “**Jitter Box Width.**” Auto-guiding is set by the parameter “**Combined offset.**” If selected, guiding is done when the telescope is pointing to the object field. As the NTT tracks very well and as one usually spends no more than a few minutes on a single position, it is not recommended.

At the end of the template the telescope returns to the original position, making it especially convenient for obtaining a sequence of images in different filters on the same field.

The parameter “**Rotate Pupil**” is an option which allows you to rotate the instrument in between the sky and object positions so that the one obtains better sky cancellation. For imaging in filters with high backgrounds, that is the K_s and narrow band filters with central wavelengths greater than 2.2 microns, we recommend that you set this option to T (true).

For a more detailed discussion about guiding options, the algorithm used to compute the offsets, and the reasoning behind the option to rotate the pupil, please refer to Sec. 3.6.

If you wish to enter the offsets manually, use the template **SOFI_img_obs_JitterOffset**. If one wishes to use a more complex pattern that does not involve observing the object and the sky alternatively, use either **SOFI_img_obs_Jitter** or **SOFI_img_obs_GenericImaging**. These templates¹ are discussed fully in The SOFI TSF Parameters Reference Guide.

NOTA BENE: The offset to a “clear” *sky* field introduces extra 130-150% overhead, because we have to spend the same amount of time on the *sky* as on the *object+sky*, plus the usual 30-50% overheads for readout, telescope offsets, etc. Therefore, it is recommended – if possible – to take advantage of the SOFI Large Field mode and to observe such targets in the same way as compact targets (see Sec. 3.2.2).

3.2.4 Maps of Large Fields

To cover a large area of the sky with a map the template **SOFI_img_obs_AutoJitterArray** is available. This template allows to define an array of positions (NEXPO) through a list of offsets in RA and DEC and to randomly jitter (number of jitter positions defined by NJIT) around each of the offset positions. The NJIT images around each offset are completed before moving to the next offset position.

However, in some cases it is preferable first to obtain a single image at each offset position, then to “shift” the entire offset pattern by a small jitter and to repeat it. This can be done with the template **SOFI_img_obs_AutoJitterArray_1** that has identical parameters, takes the same number of images on the same locations as **SOFI_img_obs_AutoJitterArray** but differs only in the order in which these images are taken.

The template **SOFI_img_obs_AutoJitterArray** may be preferable if it is more important to make sure that at least some of the offset positions are observed with the required depth, while the template **SOFI_img_obs_AutoJitterArray_1** may be preferable if covering the entire mapped area (albeit at a shallower depth) has higher priority than the depth of the mosaic.

The parameters of both templates are listed in Table 3.3.

¹ As of 10-05-2012 these two template must never be used with guiding option set to “S”! See the template description for details.

Parameter signature	Value
Exposure Name	SOFI_Map
DIT	10
NDIT	6
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
NJITT	6
NEXPO	9
Jitter Box Width (arcsec)	40
Filter wheel 1	K_s
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined Offset? (T/F)	F
Return to Origin? (T/F)	T
RA Offset List (arcsec)	0 450 450 0 -450 -450 0 450 450
DEC Offset List (arcsec)	0 0 0 -450 0 0 -450 0 0

Table 3.3: Parameters of the SOFI_img_obs_AutoJitterArray template with commonly used values.

3.2.5 Imaging of Moderately Large Object

This section describes imaging of moderately large objects, i.e. objects comparable to the SOFI field of view. We discuss this case after the mapping of large areas because it uses the same template **SOFI_img_obs_AutoJitterArray_1** as the area mapping.

There are two distinct cases:

(i) Observations of “semi-extended” objects as large of a quarter to half of the SOFI field of view; these may be moderately distant galaxies, small Milky Way clusters or LMC/SMC clusters. They are too big for simple jittering mode observations but still leave room for more efficient observing strategy where the observer does not need to sample clear sky 50% of the observing time. Instead, the user can observe the sky simultaneously with the target, adopting clever offset pattern that would move the target in the centres of the four quadrants, or in the centres of the two halves of the array. In the latter case, a suitable rotation offset may be necessary to apply during the acquisition, so the side of the array is aligned with the major axis of the object.

A typical example of a 4-point observation is shown in Fig. 3.1. The figure shows the target – a round object with diameter 90 arcsec – as it will appear on the RTD and on the SOFI images. The best choice is to use the template **SOFI_img_obs_AutoJitterArray_1**. As it is described in Section 3.2.4 this template allows the user to take a sequence of jittered images around user-defined positions: first, all positions are imaged ones, then a random jitter is added to the entire pattern and it is repeated as many times as necessary. Note that the offsets are executed **before** the images are taken, and the offsets are defined along RA and DEC, in arcseconds. On the figure, the images in the sequence are numbered to show the order in which the target will move during the observations. The zero number is the location of the target acquisition and it will be discussed further.

The parameters of both templates are listed in Table 3.4. Note that for 4-point observation the NEXPO parameter **must** be equal to four! Therefore, the total integration per position is controlled by the time spent on each image (= DIT×NDIT, usually 1-3 minutes) and by the number of jittered images taken at each individual position NJITT. In this example, the total integration time is: NJIT×NEXP×NDIT×DIT = 12×4×10×6=2880 sec.

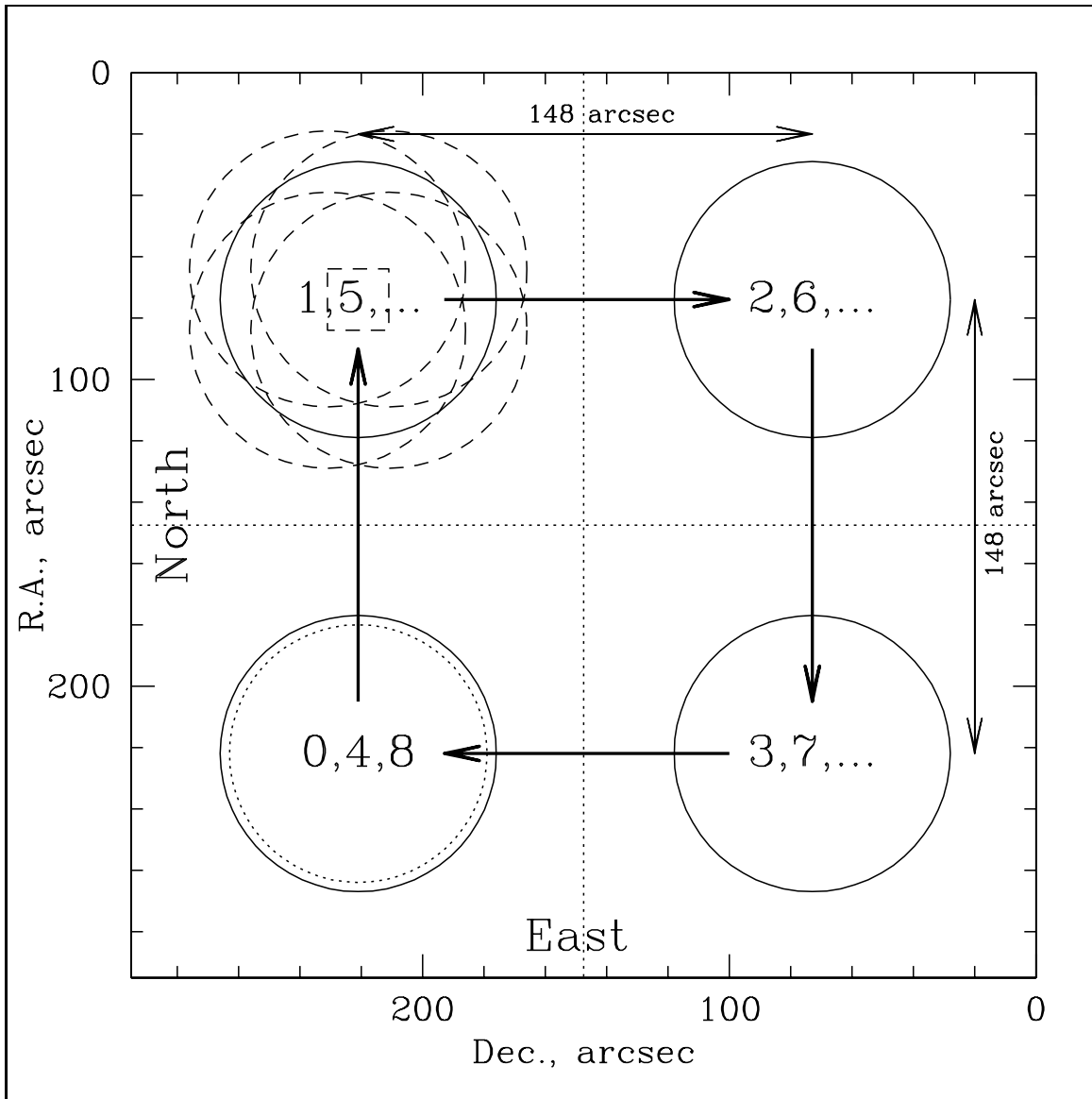


Figure 3.1: Example of 4-point observation scheme of a “semi-extended” objects without extra overhead for observation of clear sky. Remember, than although the target is moving on the RTD, it is really the telescope that is moving!

In general, the template allows to define manually the offsets for each individual image. However, this is inconvenient if the number of images is large. To simplify the matter, we suggest to the user to define a “closed” loop, where after four images the telescope points back to the original position. This means that the sum of the offsets along the RA is zero, same as the sum of the offsets along the DEC.

Fortunately, the autojitter feature of the template ensures that the the next batches of four images will not be taken at exactly the same positions as the previous ones. In the example the jitter box is 20 arcsec, and for simplicity it is shown only on the first of the four positions. Four dashed circles show where the target would be if the random jitter offsets move the target centre to the edges of the jitter box. This has to be remembered well, because too large a jitter box may move the target out of the field of view or it may lead to overlapping of the target edges in sequential images. The overlap is better to be avoided because it may affect the sky subtraction.

The user should try to keep the target as close to the centre as possible without allowing the overlap

Parameter signature	Value
Exposure Name	SOFL4-point
DIT	6
NDIT	10
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
NJITT	12
NEXPO	4
Jitter Box Width (arcsec)	20
Filter wheel 1	K_s
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined Offset? (T/F)	F
Return to Origin? (T/F)	T
RA Offset List (arcsec)	148 0 -148 0
DEC Offset List (arcsec)	0 148 1 -148

Table 3.4: Parameters of the SOFLimg_obs_AutoJitterArray_1 template with commonly used values for 4-point observation of a “semi-extended” object.

mentioned above. In fact, this example allows some modification: the offsets can be reduced from 148 to 110-120 arcsec and the jitter from 20 to 10 arcsec. The small jitter is acceptable because during every cycle we take four images at different array locations reducing the effect from the bad pixels (the major reason for jittering in this case). Note that this is usually not the case with **SOFLimg_obs_AutoJitter** and template observations of uncrowded targets where the jitter is needed to produce a sky image for the sky subtraction.

There is one more – albeit simple – complication in this observing strategy: the offsets are executed before the images are taken, so the acquisition must place the target at the location of the fourth image in the cycle. This is the position shown in the figure with a dashed circle and marked with zero. Finally, the thick arrows show the apparent movement of the target on the RTD from one position to another. North and East are also marked.

This example can easily be generalised - one can observe a slightly small target at five positions, including one at the centre of the array. In this case NEXPO=5. On the contrary, an edge-on galaxy with a major axis of 3-4 arcmin a minor axis of 1-2 arcmin must be observed at only two array positions, so NEXPO=2. In this case the user should define a suitable rotator offset to align the major axis with the array edge.

(ii) Observations of objects that need 2-3 fields to cover the entire target with additional images to sample the sky.

In principle, the template **SOFLimg_obs_AutoJitterOffset** can accommodate such observations if they are carried out in a sequence of one field at a time. However, this means 100% overhead for the sky sampling. Instead, the user can reduce the overhead by alternating between one sky position and multiple target positions. A typical case is observing a large 6-10 arcmin galaxy in a sequence: sky - target field 1 - target field 2 - sky - target field 1 - target field 2 - sky - ... This is a 3-point pattern, two of the three pointings are science targets and one is sky. Naturally, this scheme is more efficient than taking a separate sky for each target, but the sky is monitored less frequently.

Alternatively, one can alternate between two or more different sky positions, to improve the sky

subtraction: sky1 - target field 1 - target field 2 - sky2 - target field 2 - target field 1 - sky1 - ... The last pattern is in effect a 6-point sequence. An example for this observation scheme is given in Fig. 3.2 and the offsets are listed in Table 3.5. Note that the rotation offset has to be defined in the acquisition template. Also, the acquisition has to put the telescope at the location marked on the figure with zero. However, the first image in the sequence is actually a sky, on the location marked with number one. The parameter NEXPO has to be set to six. The sequence in which the images are taken is shown in the figure on the right. The target is marked with a thick dark line, and it is covered by two overlapping fields. The jitter boxes around each position are not shown. One 6-point cycle will obtain two images of 1 min at each target field, 12 cycles will accumulate 24 min of integration on each of the two galaxy fields and 12 min on each of the two sky fields. The total observation will require about 1.5hr, including the overheads.

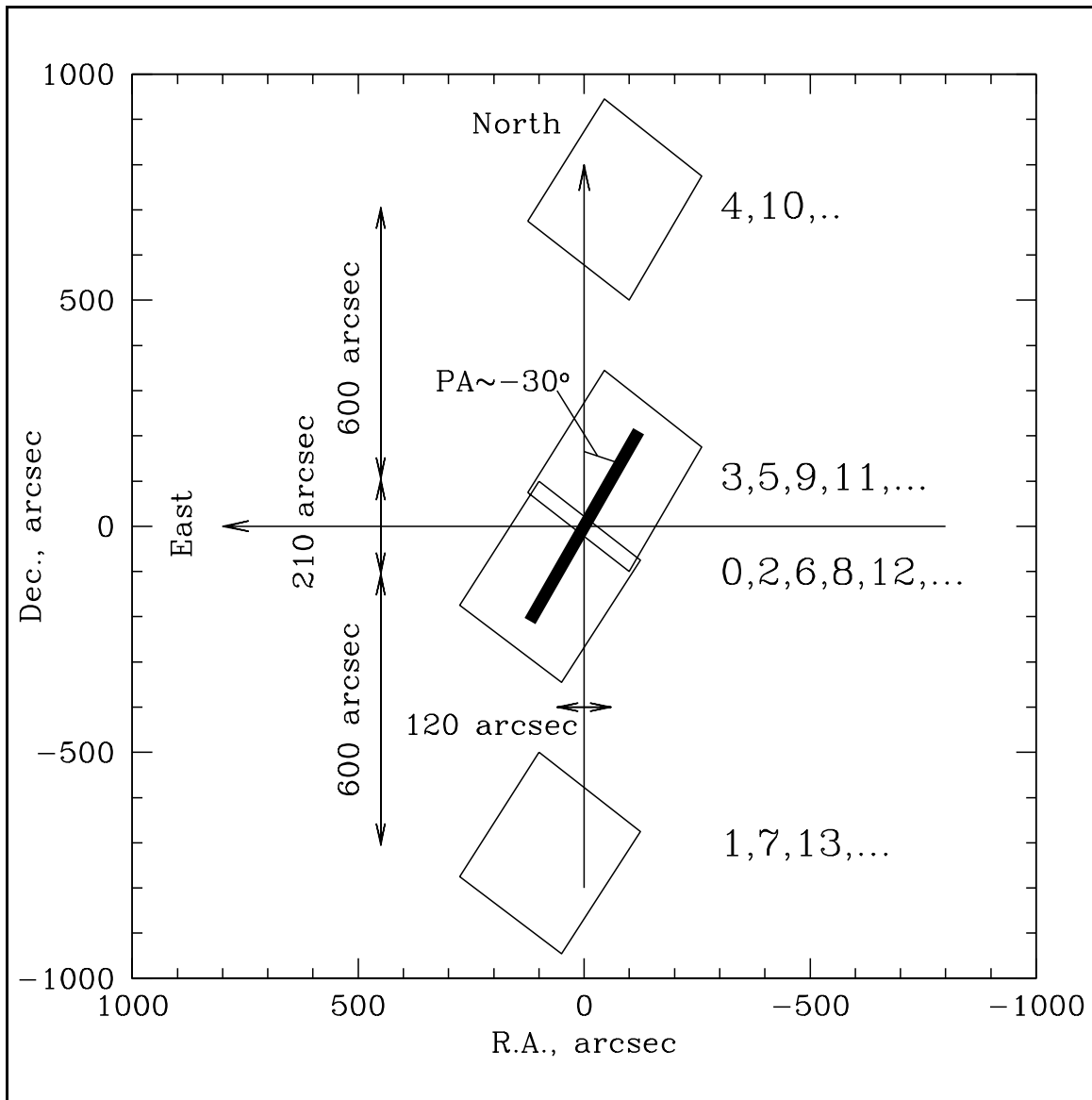


Figure 3.2: Example of a 6-point observation of “semi-extended” objects alternating between two target fields and one of two different skies.

Nota Bene: The reduction of data obtained with the patterns described in this subsection is usually more complicated and requires special attention and more manual intervention than the “classical” observations! In addition, the sky sampling is rarefied. However, in most of the cases these is a small price to pay for 20-30% more data. The users must read carefully Section 5.2.4 where the sky

Parameter signature	Value
Exposure Name	SOFL_6-point
DIT	6
NDIT	10
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
NJITT	12
NEXPO	6
Jitter Box Width (arcsec)	10
Filter wheel 1	K_s
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined Offset? (T/F)	F
Return to Origin? (T/F)	T
RA Offset List (arcsec)	0 0 -120 0 0 120
DEC Offset List (arcsec)	-600 600 210 600 -600 -210

Table 3.5: Parameters of the SOFLing_obs_AutoJitterArray_1 template with commonly used values for 4-point observation of a “semi-extended” object.

subtraction is described.

3.2.6 Faint Objects Around Bright Objects

The typical scenario that fits this situation is when one observes a faint companion of a bright star. As the luck will have it, the larger the contrast between the two objects, the more appealing the science case for such observations.

The bright star dictates a choice of strategy with minimal DIT and numerous frequent jitter offset – to minimise the artefacts rising from the electronic cross talk effect that can not be corrected satisfactorily in case if the star is heavily saturated. Such a strategy comes with a significant disadvantage – with the minimum DIT = 1.182 sec the overheads can easily reach 200%. Furthermore, the internal reflections inside the telescope and the instrument will cause optical ghosts from the bright source that are not negligible, and therefore, will also require the minimum possible DIT to minimise the ghosts.

On the contrary, the faint star calls for just the opposite strategy – use the longest DIT that does not saturate the detector array just from the sky (which often can dominate the signal even at the peak pixel of the faint star’s image), and move the telescope as little as possible to minimise the overheads.

The minimum DIT in combination with the simple jittering and offset templates can certainly help to minimise the crosstalk and the ghosts but they will not move them across the field of view. In other words, the affected areas will be the same on all jittered (or offset) images, with respect to the bright source. This may significantly affect the scientific output of the observations, i.e. the depth of the image will vary, there will be cross-talk produced star-like artefacts that will move together with the bright star, making any common proper motion studies unreliable.

One solution to this problems is to brake the observations into a few series taken at different position angles, i.e. three series of images divided by 120 deg. While this will certainly help to minimise the

effects, an alternative is to acquire all the images of the object at different angles. True, it will add an extra overhead for moving the rotator between the exposures but this extra overhead is partially compensated by having only one acquisition as opposed to the few extra acquisitions that will be associated with the different series of images.

Some extra taxes to be paid when using this strategy are:

- complication in the data reduction which will include not just the usual image alignment but derotation as well
- degradation of the image quality toward the edges of the field of view because the field distortions differ in the different parts of the array
- depending on the used rotation procedure, the corners of the field of view may be lost, so the final image will be a circle with 4.92 arcmin diameter (in case of Large Field mode) and not a 4.92×4.92 arcmin square.

Given these limitations, the templates that rotate the instrument for every exposure are suitable only for imaging of faint targets within ~ 1 arcmin from a bright source.

The templates that perform the observations described here are called `SOFI_img_obs_AutoJitterRot` and `SOFI_img_obs_JitterRot`.

3.2.7 Objects with Fast Variability

The typical observational problems addressed here are occultations, transits and fast variables. In the “normal” case the minimal DIT is 1.182 sec and it is determined by the fraction of the array being read, i.e. the detector windowing. The “normal” templates do allow to window the detector down but this does not reduce the minimal DIT because the windowing is applied to the data by the Sofi observing software **after the entire array was read**. However, the detector can be read only partially if the so called “hardware windowing” is applied. At the moment this is available to the user only in **FastPhot** mode.

In addition to the windowing, to speed up the photometry the new mode allows to save time by storing only the individual array readouts instead of their difference and by using the fast `RdRstRd` (Read-Reset-Read) readout mode. The readout time is in order of 0.002-0.008 seconds, i.e. between every integration a few milliseconds are spent to read the detector and to transfer the image.

The hardware window is centred (at 512,512) at the detector to take advantage of the simultaneous work of the four ADCs (analog-to-digital converters). The window can be rectangular. Due to the detector readout pattern, windowing along the rows and columns have different effect on the DIT reductions. If you need a rectangular window, it is better to extend it along the rows and to have the narrow side along the columns.

There is no telescope movement during a time series, so it is best to make sure the object is placed on an array location clear of bad pixels. The pixel scale is $0.288 \text{ arcsec px}^{-1}$ so windows smaller than $\sim 4 \times 4 \text{ px}^2$ are likely to lead to aperture losses and will not allow to measure reliably the sky variations. The preliminary tests indicated that the monitoring of the sky is crucial for the accuracy of the photometry. Furthermore, having a reference source in the field (so it is observed simultaneously with the target) helps to improve the photometric accuracy further. It is best if the reference source is of comparable brightness or brighter than the target.

It is advisable to split long time series into smaller bits, typically 5-10 min. The minimum separation between the different series (without moving the telescope or a filter wheel) is 8 sec.

The data product is a 3-dimensional cube where every slice along the third axis is essentially 1-dimensional image or detector read. The first two dimensions have the meaning of the usual rows and columns. The file header contains the UT at the start of the sequence and the UT at the end (as time of file creation). The ESO time distribution system has absolute accuracy better than tenth of a sec and even much better relative accuracy. This will be tested and more detail will be made available in the future.

To take full advantage of this more, the targets have to be sufficiently bright, to avoid the necessity to rebin the data later on.

The standard calibration plan for FastPhot includes darks and flats (normally dome flats, sky flats if possible) taken with the same detector configuration. This mode is not supported by the SofI pipeline.

The time to execute a series of 1000 images, using the minimal DIT for a given window size is given in Table 3.6.

dx	dy	min DIT	T _{exec}
px	px	sec	sec
1024	1024	2.1033	2105
512	1024	1.2645	2104
1024	512	1.0519	1053
512	512	0.6324	633
128	512	0.3179	318
512	128	0.1584	158
128	128	0.0798	79
64	128	0.0667	66
128	64	0.0401	40
64	64	0.0336	33
32	64	0.0303	30
64	32	0.0170	17
32	32	0.0153	16
16	32	0.0145	16
32	16	0.0079	17
16	16	0.0075	33

Table 3.6: Minimal DIT and time to complete a series of N_{DIT} = 1000 for different window sizes. The readout mode was RdRstRd.

Currently, the mode is offered only in Visitor Mode because it typically require real time decisions by the user. Members of the community interested in using the new mode are strongly encouraged to contact the SofI Instrument Scientist in advance for further details.

Nota Bene The FastPhot mode is not fully commissioned! More details will be made available on the SofI web page.

3.2.8 Standard Stars

The IR window between 1 and 2.5 microns contains several large absorption features that are primarily due to water vapour and carbon dioxide in the atmosphere. The edges of the atmospheric windows are highly variable. Unfortunately, the edges of some IR filters, particularly *J* and *K_s*, are defined by these absorption features rather than the transmission curves of the filters themselves. Thus, when the column density of water vapour is variable, accurate photometry can be difficult to achieve. On good nights (generally when the humidity is low and it is cold) it has been possible to

achieve better than 1% photometry; however, on most nights this should be considered as the best limit and the typical accuracy is 3-5%. Good planning of the observation and careful data reduction has allowed some users to reach with SOFI relative photometry of $\sim 0.3\%$!

To get good photometry you should choose standard stars that are as close as possible to your program objects and you should observe them before and after you observe your program objects. The classical, optical method of determining extinction coefficients by observing standards over a wide range of airmasses works less well in the infra-red as the extinction coefficients are generally not a linear function of airmass and may vary with the conditions, i.e. with the humidity.

During the night you should observe at least three different standards and you should observe at least one standard every two hours. On nights where the humidity (which is a rough measure of the water vapour column) is varying considerably, let's say by 40% in one hour, you will need to observe standards more frequently.

When observing standards, two or more images (typically five: one in the centre of the array and four in the centres of each quadrant) of the standard star are obtained with a telescope shift(s) in between. In this way one image can be used as the reference sky of the other(s). If several images are obtained, the uniformity of the flat-field illumination can be checked.

Several standard star lists are currently available, but it should be noted that each list was created using detectors and filters that differ from the detectors and filters used with SOFI. Note also that most standards were observed with single channel photometers with very wide apertures. Thus close companion stars were probably included. The standards listed in Carter and Meadows (1995, MNRAS, 276, 734) have proved to be very useful, although they may prove to be too bright for SOFI, and to require defocusing the telescope. The more recent NICMOS standards by Person et al. (1998, AJ, 116, 2475) are more suitable as one does not need to defocus the telescope, except for seeing better than 0.6 arcsec. They are the standards of choice for most observers. Predefined OBs are available at the telescope for most of these stars. A list of references to these and other standard star lists as well as useful papers regarding the transformations between one photometric system and another are available from the SOFI web page.

It is important that you take care where the standard star lies on the array. Use the RTD (Real Time Display) to find areas that are clear of bad pixels. Furthermore, check if other objects are in the field. The offset between exposures should not be such that other objects in the field interfere with flux of the standard when one frame is subtracted from another. Alternatively, the standard can be observed in several different parts of the array. A pattern of five exposures with the standard observed once in the centre of the array and once in each of the quadrants is one example. Table 3.7 demonstrates the usage of the **SOFI_img_cal_StandardStar** template.

In this template the offsets are along array columns and rows and the units are arc-seconds, unlike other templates where the offsets are given along the RA and Dec. The offsets are relative. This template is very similar to the template **SOFI_img_obs_Jitter**. Indeed, one can use the **SOFI_img_obs_Jitter** template to do the same observation; however, we strongly encourage observers to use the **SOFI_img_cal_StandardStar** template as it is both easier to use and it tags the resulting images as standard star frames for archiving purposes and for easier identification by the zero point calculation tool available at the telescope.

In general is not required for standard star observations, and it is useful to set the "Return to Origin" parameter to true.

Parameter signature	Value
Exposure Name	sj9104
DIT	5
NDIT	10
Number of exposures	5
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Filter wheel 1	J
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined Offset? (T/F)	F
Return to Origin? (T/F)	T
X offset list (arcsec)	0 75 -150 0 150
Y offset list (arcsec)	0 75 0 -150 0

Table 3.7: Parameters of the standard star template with commonly used values.

3.3 Polarimetry

In Polarimetry, a Wollaston prism and a mask wheel consisting of alternating opaque and transmitting strips are inserted into the beam. The widths of the transmitting sections are about 40 arc seconds, whereas the widths of the opaque sections are slightly larger. Thus, to cover the whole field one needs to take three separate images shifted by 30 arc seconds from each other. This technique is embodied in the template **SOFI_img_obs_Polarimetry**. The parameters for this template are listed in Table 3.8.

Parameter signature	Value
Exposure Name	Planetary_Disk
DIT	1.8
NDIT	10
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of exposures	3
Filter wheel 1	J
Filter wheel 2	open
Combined Offset? (T/F)	F
Return to Origin? (T/F)	T
X offset list (arcsec)	0 0 0
Y offset list (arcsec)	-64 32 32
Rotator Offset	0

Table 3.8: Parameters of the polarimetric template with commonly used values.

This example takes three exposures in the *J* band. By default, the polarimetric mode uses the large field objective. Offsets are relative to array rows and columns, are independent of the rotator angle and are in arc-seconds. In this example, three exposures with offsets along the Y-direction are taken. In this way, an entire field can be covered. Observers may wish to include additional exposures with

offsets in the X-direction to help with sky subtraction.

To measure the Stokes parameters and hence the degree and direction of linear polarisation, one needs at least one additional set of observations with a position angle different from the first. The usual practice is to take two sets of observations with a rotational offset of 45 degrees. This can be done with the ‘‘Rotator Offset’’ parameter. The approximate rotation centres for the imaging modes are listed in Table 3.9.

Objective	Rotation Centre (x,y)
Large Field	517, 504
Small Field	501, 505
Spectroscopic Objective	537, 502

Table 3.9: Approximate mechanical rotation centres for the imaging modes. Note that these values change after each instrument intervention.

There are two ways to do this. In the first method, one includes two observation templates for each polarimetric observation, one with the ‘‘Rotator Offset’’ angle set to zero and the second with the ‘‘Rotator Offset’’ angle set to 45. Alternatively, one can set the rotator angle through the acquisition template (See Appendix C) and keep the ‘‘Rotator Offset’’ angle zero in both cases. Since the rotator axis does not correspond to the centre of the array, the later method is preferable when one is trying to determine the linear polarisation of a single object. For those who wish to map the polarisation over a large field, either method can be used.

3.4 Spectroscopy

3.4.1 Small Objects and Uncrowded Fields

In spectroscopy, like imaging, accurate sky cancellation is important. If the object is small enough, the object can be observed at different slit positions. Sky cancellation is then achieved by subtracting one frame from another. This is a very efficient method because it is not necessary to spend extra time integrating off the target. This technique is embodied in the template **SOFI_spec_obs_AutoNodOnSlit**. Typical parameters for this template are listed in Table 3.10.

In this template the telescope nods the object between two positions along the slit that are ‘‘Nod Throw Along Slit’’ arc seconds apart. For convenience we will call one of these positions, position A and the other position B. At position A the object is observed for 9 minutes and the observer will receive 3 frames (specified by NINT), each the average of 3 (specified by NDIT) exposures of 60 seconds (specified by DIT). The telescope then moves to position B, where it again integrates for nine minutes producing three frames, each the average of three 60 second exposures. This completes one cycle. The number of A-B pairs – also called cycles – is defined by the parameter ‘‘Number of AB or BA cycles.’’ In this particular case the number of cycles is 4, so the total exposure time is $60^{sec} \times 3 \times 3 \times 2 \times 4 = 72$ minutes. If the number of cycles is greater than 1, the telescopes moves in the following pattern: A, B, B, A, A, B, B, A, etc. Between consecutive A-A or B-B positions, there is no telescope offset, but one can include random offsets between consecutive A-B or B-A positions. In other words, every A-B or B-A cycle can be shifted with respect to the previous. The size of the region along the slit in which the random offsets are placed, is defined by the parameter ‘‘Jitter Box Width.’’ In the example given here the 1st and 2nd A positions will differ. Additionally the 1st and 3rd B positions will also differ.

In principle, the nodding frequency should be similar to that used during direct imaging. However, the desire to reach background limited performance forces longer exposures in the spectroscopic modes.

Parameter signature	Value
Exposure Name	Hubble_Deep_Field
DIT	60
NDIT	3
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Spectro Mode	LONG_SLIT_RED
Which Slit	long_slit_1
Combined Offset? (T/F)	T
Jitter Box Width (arcsec)	10
Return to Origin? (T/F)	T
Nod Throw Along Slit (arcsec)	60
NINT	3
Number of Cycles	4

Table 3.10: Parameters of the template `SOFI_spec_obs_AutoNodOnSlit` with commonly used values.

As a general rule, one should limit the time spent at any one position to less than 15 minutes.

Although guiding is optional, we recommend that you guide. To do so, set the parameter “`Combined offset`” to true.

The Non Destructive Read-Out is recommended for faint objects that require long DIT. This is done with the templates `SOFI_spec_obs_AutoNodNonDestr` and `SOFI_spec_obs_GenSpecNonDestr`. In addition to the parameters of Table 3.10 these two templates include those showed in Table 3.11. The array is read within each DIT a number of times equal to `NSAMP`, and for each read-out the signal is sampled `NSAMPPIX` times. The minimum DIT that can be used for a given `NSAMP` is: $\text{NSAMP} \times 1.64 \text{ sec}$.

Parameter signature	Value
NSAMP	30
NSAMPPIX	4

Table 3.11: Parameters in the templates with Non Destructive Read-Out and commonly used values.

3.4.2 Extended Objects and Crowded Fields

For more complex nodding patterns or for observations of either extended objects or crowded fields, please use the template `SOFI_spec_obs_GenericSpectro`. This template is discussed in detail in Appendix C. It allows to place the slit on a user-defined sequence of positions, alternating between the *object+sky* and the *sky*.

Similar to the imaging of extended objects, an effort should be made to take advantage of the large slit length of 4.92 arcmin, to obtain observations with the simple nodding procedure described above, because the integrating of the clear sky increases the total time spent on the object by a factor of two, for the same integration time on the target.

NOTA BENE: Frames taken with the NDR mode must be flat-fielded using dome flats taken with the same mode. For this reason a specific template exists: `SOFI_spec_cal_DomeFlatNonDestr`.

3.4.3 Telluric Standards and Flux Calibration

The issue of flux calibrating IR spectra is still very much an issue of some debate. This section is written not as the definitive method to flux calibrate IR spectra, but rather it is written to illustrate the problems involved, and various methods for their solution.

As the IR window is dominated by time varying atmospheric absorption features that depend non-linearly with airmass, **it is necessary to divide object spectra with the spectrum of what we shall call a telluric standard.** The standard should be at an airmass that is as close as possible to that of the science target (usually $\Delta \text{secz} \leq 0.05-0.1$), it should be observed immediately before or after the target, and it should be observed with the same instrument set up (grism, slit, objective, etc.). From a technical point of view, the telluric standards are observed in exactly the same way as the science targets. Since they are usually bright stars with no nearby companions of comparable magnitude, the template `SOFI_spec_obs_AutoNodOnSlit` is perfectly suited for this.

It is very important that object and standard star spectra are accurately aligned. A small misalignment will result in poor cancellation of the atmospheric absorption features. Misalignment could be caused by slit misalignment in between the observations, by instrument flexure or, when the seeing is smaller than the slit width, by inaccurate centring of objects on the slit. In general, this problem can be remedied later during the data reduction by re-aligning of the spectra: first by using the telluric features to determine the offset in wavelength direction, and second - by geometric transformation of the image in spatial direction. However, these are complicated and time-consuming steps, that involve heavy modification of the data prior to any scientific measurements, and it is preferable to avoid them, if possible.

The division by the telluric standard a number artificial “emission” features introduces into the target spectrum, because of the intrinsic absorption features of the telluric standard itself. By choosing carefully the spectral type of the telluric standards one can minimise the number of such features but there are no stars with completely featureless IR spectra and therefore some such “emissions” will always remain. To remove them, **after the division by the telluric standard, one should multiply by the absolute spectral energy distribution of the telluric standard.** However, this information is usually not available, so one has to model the spectral energy distribution of the standard. For hot stars, a black body fit to the star may be appropriate, but for stars of spectral type later than B, a more accurate description of the energy distribution may be required.

Another possibility is to use a solar analog for telluric standard, taking advantage of the fact that the true intrinsic IR solar spectrum of the Sun is available. It was corrected by observing the Sun at different zenith angles and extrapolating to $\text{secz}=0$ – an extremely time-consuming technique that is not reasonable to the usually faint science targets. The IR Solar spectrum is available from the National Solar Observatory (“An atlas of the solar spectrum in the infrared from 1850 to 9000 cm^{-1} ” Livingston W. & Wallace L. N.S.O., Technical Report #91-001, July 1991). A more detailed description on how to use solar analogues as telluric standards and an IRAF-based tool can be obtained from Maiolino, Rieke & Rieke (1996, AJ, 111, 537). The tool is available at the Sofi web page: http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/sofi_scripts/ in a tarball that also contains related material, including lists of stars with reliable spectral class determinations that can be used as telluric standards.

The final step is to flux calibrate the object spectrum. The flux calibration is done by **integrating the object spectrum over one or more broadband filter band-passes and scaling the result with the respective broadband magnitudes.** This method typically gives accuracy of order of 5-10%. It is very different to what is done to flux calibrate spectra in the optical. In the optical, atmospheric absorption is a relatively smooth function of wavelength, so it is sufficient to divide extracted spectra by a smoothed standard star spectrum and then multiply by the smoothed absolute energy distribution of the standard. In the IR, accurate spectroscopic standards do not

exist. For an example of IR flux calibration in the more general case of extended objects see Ivanov et al. (2000, ApJ, 545, 190).

The choice of which standard to use depends on which part of the spectrum is of interest. All stars have absorption lines, so the idea is to choose a star which does not have strong features near the wavelength of interest. Hot stars provide relatively featureless spectra; however, they have strong hydrogen absorption lines, so they should not be used as flux calibrators if the region around the hydrogen lines is of interest. Later type stars such as G stars have weaker hydrogen lines, but are contaminated by multiple weak absorption lines. These stars have the additional advantage of being very numerous. Stars of later type should not be used as these stars contain numerous weak lines throughout their spectra. However, such stars have very weak hydrogen absorption and may be usefully employed to determine the strength of hydrogen absorption in the hot stars. In general, IR standards are significantly brighter than optical standards, nevertheless, stars should be fainter than seventh magnitude.

A moderately-high resolution IR spectral library for stars with well-known parameters, including spectral types, is available from Ivanov et al. (2004, ApJS, 151, 387).

A list of O, F and G type stars selected from the HIPPARCOS catalogue, with magnitudes appropriate for SOFI spectroscopy is available at the telescope and on the SOFI web page, together with a web-based tool for automated search for stars with given spectral type and brightness near the science target: <http://www.eso.org/sci/facilities/paranal/sciops/catsearch.html>

A similar tool *stdsopMain*, with extra capability to automatically create OBs for observing of telluric standards, is available on the SOFI instrument workstation but **we recommend that the users create their own OBs with the P2 tool to retain consistent program ID.**

Nota Bene: Until recently IR spectrophotometric standards were not available. However, Cohen et al. (2003, AJ, 125, 2645) published “supertemplates” for a series of 33 stars, that span wavelength range from the UV to mid-IR (40 μ m). Tests show excellent agreement between the synthetic photometry derived from these spectra and the observed magnitudes of the stars, making this spectral library useful for flux calibration of IR spectra.

3.5 Calibration Frames

3.5.1 Darks (Biases)

The concept of bias frames with IR arrays does not have the same meaning as CCD bias frames. With IR arrays, a zero second exposure is not possible. It is better to think of all exposures without direct illumination as dark frames. Note that in this context the biases are an integral and indistinguishable part of the darks.

Dark frames are taken with the template **SOFI_img_cal_Darks**. The parameters for this template are listed in Table 3.12.

In this template, one can enter a list of DITs and NDITs. If the number of exposures is greater than the number of elements in either of these lists, the list is repeated until the correct number of exposures have been completed.

The structure seen in the darks is quite complicated. In general, it is not a linear function of time. The signal is made of several components: shading, a component which depends on the DIT and on the incident flux; heat from the readout amplifiers, commonly referred to as amplifier glow; and classical dark current from the random generation of electron/hole pairs. The heat status from the readout amplifiers is a function of the number of reads only, whereas the dark current is a linear

Parameter signature	Value
Exposure Name	Dark_Frames
DIT LIST	10 10 20 20 10 10
NDIT LIST	6 6 3 6 3 6
Number of exposures	12
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1

Table 3.12: Parameters in the SOFI_img_cal_Darks template with commonly used values.

function of time.

Other than depending on the DIT, the shading pattern seen in darks frames also depends on the incident flux. Thus, subtracting a dark frame from a science frame with the same DIT will not remove the shading pattern perfectly. This issue is discussed in more detail in the next section.

The experience has shown that the best way to remove the dark current is by taking data in exactly the same way and with exactly the same exposure time as the science observations. This is routinely done for the removal of the sky emission, so the removal of the dark is naturally removed when the sky is removed. The only case when the removal of the dark remains an issue is the flat fields (see the next section for discussion).

Darks appear to be stable over the period of a typical observing run. Thus, it is sufficient to take darks only once during the run, or more often – as a verification of the detector array status.

Nota bene: Darks are not part of the calibration plan, for the reasons described above. If you still need darks, prepare your own calibration OBs.

3.5.2 Flat Fields

As in the case of visible observations, one must correct for differences in pixel sensitivities. The method for creating flat fields is identical in both the spectroscopic and imaging modes. They are created by exposing alternatively an illuminated and an unilluminated dome panel. The flat-field images are constructed from the difference of the two, normalised to one. This technique is especially important when working beyond 2.3 microns, where it removes the thermal component of the signal, which does not depend on the intensity of the flat-field lamp. Spectroscopic flats must be taken with the same slit as the observations.

Imaging flat-fields can also be obtained from the twilight sky or from the observations themselves. However, dome flats represent better the low frequency sensitivity variations of the array.

Twilight spectroscopic flats cannot be used to flat field your data, although they can be used to determine the slit transmission function for spectroscopy of extended sources. The slit transmission function describes the transmission variation along the slit. It is usually affected by magnetised or sticky dust particles that attach to the slit edges. These particles can make the slit effectively narrower. Normally, steps are taken to clean the slit during interventions so the slit transfer function is highly uniform. This effect can be neglected for spectroscopy of point sources.

Spectroscopic flat fields can be taken with the halogen lamp in the calibration unit (Nasmyth flat). Again, one takes an image with the lamp on and off. The flat is the difference between the two, normalised to unity. The optical path length between the Nasmyth flat field lamp and the front window of SOFI is considerably smaller than the path length between the incandescent lamp on the

floor of the dome and the front window of SOFI. Thus, one will see less of the strong atmospheric absorption features in the Nasmyth flats.

There are four templates to create dome flat fields, two for imaging **SOFI_img_cal_DomeFlats** and **SOFI_img_cal_SpecialDomeFlats**, and two for spectroscopy **SOFI_spec_cal_DomeFlats** and **SOFI_spec_cal_DomeFlatNonDestr**. Please note that the last template must be used **ONLY** to flat frames taken with the NDR Mode. Table 3.13 gives an example for imaging dome flats.

Parameter signature	Value
Exposure Name	Imaging.Dome.Flats
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of exposures	1
DIT LIST	3
NDIT LIST	60
Filter wheel 1	J
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING

Table 3.13: Parameters in the SOFI_img_cal_DomeFlats template with commonly used values.

The template actually takes four exposures in total: the first with the lamp off, then two exposures with the lamp on and then a final exposure with the lamp off. The intensity of the dome flat field lamp is controlled by a panel on the *tcs* machine. The instrument operator will adjust the voltage until a level of a few thousand (usually 3000-4000 ADU) counts is reached when the lamp is on.

In this template, one can enter a list of DITs, and NDITs. If the number of exposures is greater than the number of elements in either of these lists the list is repeated until the the correct number of exposures have been completed. **Note that each “exposure” generates a sequence of 4 frames**, as explained before. If the parameter "Number of Exposures" is >1, a corresponding sequence of frames will be generated. The template **SOFI_spec_cal_DomeFlatNonDestr** contains in addition the parameters of Table 3.11.

It was pointed out earlier that shading, one of the components that make up a dark frame, is a function of the incident flux. This means that the method of creating dome flats described above does not remove perfectly the shading pattern. The effect for most observations is small, and manifests itself as a discontinuity of a few percent in the ZP across the centre of the array. For most programs, this is not a major problem. Nevertheless, we have developed a technique to remove this residual, and this technique is embodied in the template **SOFI_img_cal_SpecialDomeFlats**. This template, in addition to the four frames taken with the template **SOFI_spec_cal_DomeFlats**, takes frames with the mask partially obscuring the array. Therefore, each exposure generates a sequence of 8 frames. In addition to them, 1-2 or more images are taken at the beginning in order to adjust and check the illumination level (set to 4000-6000 ADU), so the total number of the obtained images by this template may vary. Note that only the last 8 images are actually the flat field. Half of this sequence is shown in Fig. 3.3. If the parameter "Number of Exposures" is >1, a corresponding sequence of frames will be generated (so generally, use 1). These frames are used to estimate the shading. Further details are given in the section describing data reduction.

Flat Field frames for Large Field imaging with broad band filters are taken by the observatory staff as a part of the calibration plan in the morning after the observations, in the filters that were used during the night. Narrow band flats are responsibility of the users. It is also a responsibility of the users to verify the quality of the flats (and any other day-time calibrations).

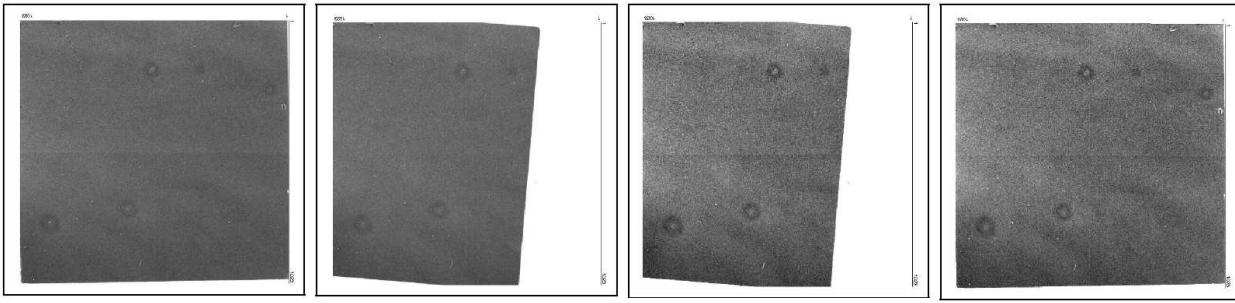


Figure 3.3: Examples of Special Dome Flat images. From left to right: lamp off, lamp off with mask, lamp on with mask, lamp on.

The stability of the flats was studied over a period of about a year, and it appears the degradation over time is minimal as can be seen in Fig. 3.4. Therefore, it's not necessary to take flats every day, especially if high precision photometry is not required to achieve the scientific goals of your program. However, it is advisable to take at least two sets of flats to check their consistency.

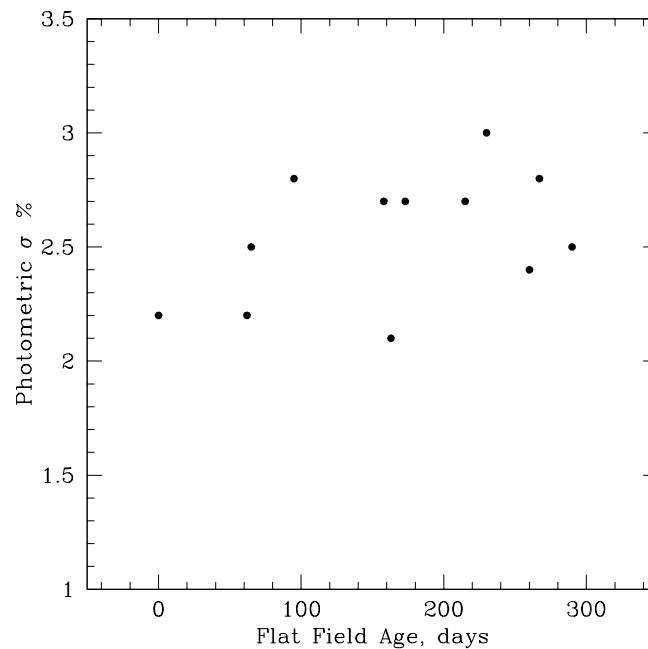


Figure 3.4: Stability of the Special Dome Flat images: accuracy of the photometry as a function of the flat field “age”.

3.5.3 Illumination Corrections

It is often the case that the flat field generated from the dome, the twilight sky or the night sky itself does not represent completely the low frequency sensitivity variations of the array. In other words, the dome flat field contains a residual large scale variation that is usually due to the fact that the illumination of the screen can never represent accurately the illumination of the sky. In addition, the illumination of the screen can change with time due to the variations in the lamps used to illuminate it, and due to variation of the screen and dome temperature (for $\lambda \geq 2.3$ micron).

This problem can be solved by observing a bright star (preferably a standard) in a rectangular grid

pattern of 9-16 positions across the array. After sky subtraction and flat fielding with the uncorrected flat, the intensities (not the magnitudes!) are fitted with a low order polynomial. Next, the surface generated by this fit is normalised to one. Then, the flat is multiplied by this surface to create a corrected flat field.

Alternatively, one can think of this surface as a second flat field, and instead of correcting the dome flat, the user can divide the data by the regular dome flat and then by this surface. This surface is sometimes called an illumination correction. Without this correction the low frequency variability from edge to edge of the array is between 1 and 3 % depending on the filter, the objective and the type of flat used. After the illumination correction the low frequency variability improves to better than 1%.

Illumination Correction surfaces can be downloaded from the SOFI Web page. The observatory staff produces them typically once a month. Templates for observing a standard star in an array of 16 positions are available in the impex directory of SOFI on wsofi (in the calib directory: **SOFI_Illum_Correction_J**, **SOFI_Illum_Correction_H**, etc.).

NOTA BENE: The illumination correction frame refers **only** to the flat field that was used to create it, and it must not be used with other flat fields! The reason for this is that the illumination correction removes the effect of different dome screen illumination and in general every time the user obtains dome flats the intensity of the lamps is slightly different.

3.5.4 Arcs

Wavelength calibration can be done with the Xenon and Neon lamps in the adaptor. The Xenon lamp has a better distribution of lines and is by itself sufficient to calibrate the wavelength scale of data taken with the low resolution grisms. Both the Neon and Xenon lamps should be used to calibrate the medium resolution grisms.

Arcs are taken with the template **SOFI_spec_cal_Arcs**. The parameters in this template are displayed in Table 3.14. In this template, one can enter a list of DITs, NDITs, grisms, slits and lamps. If the number of exposures is greater than the number of elements in either of these lists, the list is repeated until the correct number of exposures have been completed. The choices for the lamps are: Xe - Xenon, Ne - Neon, B - Both and N - None. The choices for the slit are: 0.6, 1 and 2. The choices for the spectral mode are: B and R for low resolution grisms, H and K for the high resolution grism. Although, the template enables one to do this, it is best to keep things simple, i.e. just take arcs with one slit and if another slit is required, run the template again with that slit.

The Xenon lamp warms up quickly and produces high background level in the arcs, notably when calibrating the red grism. To minimise the background, the user may want to split the time when the lamp is on giving the lamp some time to cool. To do this, just decrease the number of exposures and repeat the OBs with a few minutes delay, as many times as the NEXP was reduced.

Nota Bene: Remember to take images with the lamps off, to remove the pattern due to dark/bias, scattered light and thermal background. This makes it easier to identify the weaker emission lines.

In general you should avoid taking arcs during the night since the high illumination levels of the strongest emission lines in the arc images may cause persistence problems. It is therefore advisable in such cases that you use the arcs in the morning after your run to do the wavelength calibration and then use the atmosphere to check the zero point of the calibration during the night. Both OH emission and sharp atmospheric absorption lines can be used. The P1 branch lines of OH are very good for wavelength calibration. Avoid the blended Q branch lines. Sample spectra for Xenon, Neon, and the main OH emission lines are displayed in appendix A.

Predefined OBs for wavelength calibration with the correct exposure times are available at the tele-

Parameter signature	Value
Exposure Name	Arc_Calibration
DIT LIST	2 2 2 2
NDIT LIST	5 5 5 5
Number of exposures	4
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Spectro Mode	B R
Slit List	1
Lamp List	N Xe N Ne

Table 3.14: Parameters in the SOFI_spec_cal_Arcs template with commonly used values.

scope.

3.6 Finer Points

Some of the information in this section is available through the manual but here it is concentrated for the convenience of the user.

3.6.1 Choosing DIT, NDIT and NINT

The appropriate value of DIT depends on the intensity of the source and the background. First and foremost, the DIT must be kept short enough that moderately bright objects of scientific interest do not saturate. On the other hand, in order to maximise the S/N ratio, one would like to work in background limited instrument performance, i.e. the sky intensity in a single frame should be sufficiently high that the sky shot noise will dominate over the detector read-out noise. In the K_s band, where the sky is bright, this can be reached with DITs of one second. Through the spectroscopic modes, where the background is lower, longer DITs are necessary. For most observations this will mean a DIT between 10 and 60 seconds.

In addition, shorter DIT means larger overheads because every DIT is associated with certain fixed overheads (i.e. for reset and read-out), and the shorter the DIT, the more DITs it will take to complete the desired total integration time.

Once the DIT is chosen, the NDIT can be determined from the desired offsetting frequency. During good weather conditions one can stay on the object as long as ~ 2 min for broad band imaging (~ 5 min for narrow band imaging and ~ 15 min for spectroscopy) before switching to the sky, and in this case $\text{NDIT} = 2 \text{ min} / \text{DIT}$ is appropriate. Note that in case of short DITs the NDIT may have to be reduced further by a factor of 1.2-1.3 to take into account the overheads. If conditions are not so good, e.g. the sky intensity is fluctuating rapidly, or if one simply does not wish to integrate that long, then a lower value of NDIT can be used. As a general rule, the switching frequency should be higher for larger fields of view or longer wavelengths.

Another consideration that has to be taken into account in choosing the NDIT is the need to have a sufficiently high number of jittered images (1) to create a good sky image without stellar residuals, and (2) to remove well the array cosmetic defects and the cosmic rays (albeit there are few cosmic rays because of the short integrations for individual images in the IR). Three jittered images is the lowest meaningful number to ensure the array cosmetics removes well, but having 5-7 or more images

is preferable, especially if they have to be used for creating a *sky* frame.

An additional restriction on the offsetting frequency applies to imaging observations. Unlike equatorial telescopes, the pupil plane of the NTT rotates relative to the image plane and, for reasons which may be due to internal reflections within the objectives, this causes pupil ghosts (images of the telescope support structure near the primary mirror) to appear in sky subtracted images. **In other words, the traces of the spider supporting the secondary mirror do not fully subtract on images taken more than 5-10 min apart.** The effect is greatest when the parallactic angle changes quickest and this occurs when the telescope is near the meridian. The effect is worst for images taken in K_s and the large field objective. To minimise this effect, the beam switching frequency should be such that the parallactic angle does not change by more than 0.5 degrees between exposures. For “jittered” images the residual effect is at or below the shot noise level when applying standard reduction techniques which involve running sky measurements over ≈ 10 frames. For applications requiring large object-sky offsets, it is advisable to use a strategy which automatically equalises the parallactic angle of the two. This strategy can be employed in two SOFI templates: **SOFI_img_obs_AutoJitterOffset** and **SOFI_img_obs_JitterOffset** – just set the parameter “rotate pupil” to “true”. It is possible to reduce the effect of the ghosts during the data reduction by running sky measurements over only 2-4 frames but this makes the data reduction somewhat more complicated.

Finally, in the spectroscopic mode only, one can accumulate a number of NINT exposures in one nodding position (A or B) before moving the telescope to the other side of the nod. This feature allows to minimise the overheads for moving the telescope between the two nodding positions. These overheads are bigger than the ones related to the small jitter offsets.

3.6.2 Autoguiding

The NTT tracks very well over most of the sky. For IR imaging and some spectroscopic observations, the exposures are short enough that there is no need to guide. Guiding is necessary only for spectroscopy of faint targets, with exposures of 10-15 min or longer.

Nota Bene: If you choose to guide you should be aware that the guiding mode is currently fixed to *box2star* (as opposed to *star2box*) and this option has to be chosen in templates that allow to specify the guiding mode.

3.6.3 SOFI Observing Modes

For SOFI, there are several standard observing modes. For example, there are modes for wide field imaging and long slit spectroscopy. The observing mode defines which optical elements are inserted into the beam. As an example, in the large field imaging mode, the grism wheel is set to open, the mask wheel is set to the large field mask, and the objective wheel is set to the wide field objective. In addition, the collimator – used for focusing the instrument – is moved to a predefined position, depending on the mode and the filter. The standard modes available with SOFI are listed in the following Table 3.15.

3.6.4 Template Parameters - Signatures and Keywords

Template parameters can be described either by their signature, which are aimed at being self-explanatory, or by their keyword, which are more meaningful to the OS (Observing Software). When a template is added within P2, the signature form of the parameter is displayed. When the OB is passed to the BOB (Broker of Observing Blocks) for execution, signatures are translated into keywords. Although users will mostly be faced with signatures, they will later see the corresponding

Observing Mode	Description
LARGE_FIELD_IMAGING	Imaging with the large field objective
LARGE_FIELD_SO_IMAGING	Imaging with the large field objective
SMALL_FIELD_IMAGING	Imaging with the small field objective; decommissioned
POLARIMETRY	Polarimetric mode
LONG_SLIT_RED	Long slit spectroscopy with the red grism
LONG_SLIT_BLUE	Long slit spectroscopy with the blue grism
LONG_SLIT_K	2.0-2.3 micron, long slit spectroscopy with the medium resolution grism
LONG_SLIT_H	1.5-1.8 micron, long slit spectroscopy with the medium resolution grism
DARK	Dark exposures
IMAGING_LF_FE	Imaging with the large field objective and the focal elongator; decommissioned
UNDEFINED	Undefined Mode

Table 3.15: The currently supported SOFI observing modes.

keywords, either when BOB is executing the template, or in the FITS image headers where some of these keywords are stored.

As an example, setting the filter on filter wheel 1 is done with the signature "Filter wheel 1", whereas the corresponding keyword is `INS.FILT1.ID`.

3.6.5 File Naming and Exposure Number

The parameter "Exposure Name" determines the name the output files will be given on the instrument workstation. The raw data become available to the user under different names, consisting of an abbreviation of the instrument and a time stamp indicating the time the file was created, i.e. `SOFI.2005-02-03T10:17:22.354.fits`, where the time is UT. However, the exposure name is stored in the fits header and some users may find it convenient to use it in order to rename the files to somewhat more intuitive names for easier data reduction.

The exposure names consist of a base name, defined by the parameter "Exposure Name," plus an extension of the type `_0003.fits`. The extension is automatically set. If the base file name is found on disk, the sequential number is automatically incremented. For example, if there exists a file "fileName_0024.fits" on disk, the next file generated by the template is "fileName_0025.fits". If no file with the same base name is found on disk, the sequential number is automatically set to 0001 (see for an example Table 3.16).

NOTA BENE: Do not use spaces, slashes or other non-alphanumeric characters in file names and target names, otherwise the SOFI OS will report an error. However, the underscore symbol ("_") is acceptable.

Signature	Keyword	Default	Description	Input Value
Exposure Name	DET EXP NAME	–	Exposure Base filename	ngc1068
Number of Exposures	SEQ NEXPO	–	Number of exposures for the template	10

Table 3.16: File naming signatures and keywords. The first filename will be `ngc1068_0001.fits` if no file with the "ngc1068" string was found on disk, and the last file will be `ngc1068_0010.fits`.

3.6.6 Detector Window

In all templates, **except acquisition templates**, it is possible to window the detector. The window is described with the coordinates of the origin and the window size (Table 3.17). Note that the entire array is read, the windowing parameters only determine which part of the image will be stored in the raw data file. Therefore, the windowing does not decrease appreciably the overheads.

Signature	Keyword	Default Value	Explanation
Number of rows	DET WIN NX	1024	Number of rows
Number of columns	DET WIN NX	1024	Number of columns
First column of window	DET WIN STARTX	1	First column of window
First row of window	DET WIN STARTY	1	First row of window

Table 3.17: Detector Window Signatures and Keywords.

Chapter 4

Phase 2 Preparation and Observing with SOFI

4.1 General Issues

Observations at the NTT are carried out in “classical” mode, i.e. **there is no support astronomer on the site**. The observers are expected to arrive a few days before their run begins, and to use this time to familiarise with the instrument. At all times there are day- and night-time telescope operators who provide technical support. Astronomical support can be obtained by e-mail contact with the instrument scientist (names and contact details are listed in the instrument web pages), well in advance before the observation started because immediate response cannot be guaranteed. Limited help on short-time scale may be available from the support astronomers at Paranal via phone but it is also not guaranteed. **Clearly, the early preparation for the observations is the key for success!** Some further hints:

- The visitors will be given a visitor account to use the computers in the visitors’ office; there are possibilities to link a laptop to the local network (WiFi and Ethernet).
- The visitors are encouraged to call the La Silla safety number officer at 4444 in case of medical or other emergency, and to warn a member of the observatory staff for any outdoor activity; remember that **the closest hospital is a few hours away!**
- The visitors are expected to contact their support telescope operator upon arriving.
- The visitors are expected to read carefully the La Silla web page to acquaintance themselves with the rules of the operation of the observatory and the instrument they are going to use.
- The visitors should not edit by hand an OB or write a script to produce them, as it may cause software crashes.

Detailed visitor guidelines are available at the La Silla web page:

<http://www.eso.org/sci/facilities/lasilla/sciops/observing/VA.GeneralInfo.html>

Starting P104 the new p2 tool, see:

www.eso.org/p2ls

must be used to prepare the observations.

4.2 The VLT environment: P2, BOB, OS, TCS, DCS, ICS

Observations are described by observing blocks (OBs) – essentially scripts for driving the telescope and the instruments. The OBs consist of one or more templates – in effect, commands in the scripts. The templates have a number of parameters, much like the functions in computer programming.

A typical science OB is made up of four components: a target, an acquisition template, an observation description, and a set of constraints. The observation description itself is made up of one or more observation and calibration templates. The constraints describe the acceptable conditions under which an observation can be carried out in service mode, and therefore the constraints section can be ignored by the visitor mode observers.

Templates are the simplest unit of an observation. They are split into three categories: acquisition, observing and calibration. The templates are described fully in Appendix C.

The P2 tool is a online web interface based on java that enables observers to create lists of targets and observing descriptions and then associate them with an acquisition template to form an OB. The day-time telescope operator will start you up on P2, and will help you to optimise the observations but **we highly recommend the users to check the following link:**

<http://www.eso.org/sci/observing/phase2.html> The OBs can also be prepared (and actually should be) before coming to the Observatory and can be checked on site before the beginning of the observations.

Once the OBs are created, they can be *checked in*: once this is done, they will be visible in the visitor Observing Tool or vOT under the PI account. From here, OBs can be sent to the execution sequence from which they can be retrieved by BOB, or they can be exported as ASCII files. The exported format of the files allows to ftp them from one machine to another.

NOTA BENE: Under no circumstances edit the OBs by hand with a text editor or with scripts! The software is sensitive to the format of the OB and extra (or not enough) spaces/tabs often cause crashes!

When observing, vOT is used to select OBs for execution. Observing blocks are selected by sending them to the execution sequence; OBs will be executed sequentially. If you want a different OB than the one on top of the sequence to be executed, you will have to move it up manually.

Note that once an OB is transferred to BOB, any changes that the observer make to the OB in the P2 have no effect on the observations, unless the OB is loaded into BOB again.

In general, putting the OBs in the execution sequence is all the observer needs to do. Some observer intervention may be required during the execution of an observing block. For example, the user may need to select the correct object to be put in the slit or to indicate an offset during the acquisition. Everything else is either automatic or it is done by the support staff of the observatory. The users are advised to use the available time for inspection of the raw data.

BOB is a very versatile tool. It can be used to display the contents of an OB, it can skip templates within an observing block, it can be used to pause at a template, and it even allows to edit the parameters of an OB (but only in the templates that have not been started yet). Thanks to this facility it is possible to fine tune the execution of an OB in real time. However, with this freedom comes the danger to make a mistake, especially during a long and stressful observing night. Therefore, we recommend to make changes within the P2, check in the OB, and then to re-load the edited OB into

BOB. This provides an extra layer of verification because the P2 has some built-in tools for verification of the OBs. Please note that if you wish to edit an already checked-in OB, you will need to *check it out* first.

BOB runs a Tcl script which sends commands to the OS (Observing Software) which then sends commands to the TCS (Telescope Control System), the DCS (Detector Control System) and the ICS (Instrument Control System). The status of the instrument is displayed in the OS GUI (Graphical User interface) and the obtained raw images are displayed on the RTD (Real Time Display). The results are also sent to the archive machine and the off-line workstation *wg5off*. This will be discussed in more detail later.

When creating OBs, it is useful to keep them as simple as possible. Do not create complex OBs which switch from one mode to another. OBs with multiple science templates with different filters, imaging the same target are acceptable. Likewise it is possible to create spectroscopic OBs with different grisms in different templates, but with the same slit. This is because the different slits project on slightly different places on the detector.

4.3 Arriving at the Telescope.

The NTT is run by a section of the joint control room (the former library) located in the building just below the hotel. Consult the map on the La Silla web-page to find out where it is located. Please be aware that you will be sharing the control room with the observers at two other telescopes and some understanding and patience may be required.

When you arrive at the control building, you will be confronted by a vast number of terminals. Most of these are of no concern to you. You will be primarily interested in the terminals at the right end of the room NTT section. Starting from right to left there is an empty table with a network connection that can be used to plug-in a laptop. The first desktop computer next to it is *wg5off*, where you can run *IRAF*, *MIDAS*, *eclipse*, *IDL*, web-browsers and other useful programs. This is a dual monitor machine. From here you can access and examine your data. Note that the system clocks of these computers are set to the UT.

To the left of the *wg5off* is the terminal *wh5dhs* where you run the VOT. Immediately to the left of *wg5dhs* is a sequence of three vertically placed dual screen terminals for the EFOSC, and SOFI workstations (in this order). The bottom screens show their OSs, BOBs, etc., while the upper screens show their RTDs. For convenience we can run some applications on the screens of the instruments that are not in use at the moment. To the left of the SOFI workstation are the computers of the telescope. The user will have no interaction with them.

When you arrive at the telescope, the entire system will have been started for you by the daytime operator. During the entire night you will be accompanied a telescope/instrument operator (TIO). It is important to emphasise that the telescope/instrument operator will be the one controlling BOB and SOFI. You as an observer will be selecting OBs for execution via vOT and inspecting the incoming data.

NOTA BENE: The TIO is responsible for the operation and safety of the telescope and instrument; he has the authority to refuse an action that would jeopardise this safety, this include closing the telescope in case of dangerous weather conditions.

On the terminal, *wg5off*, do the following:

- The machine will be started by the support staff, together with all the necessary applications, so there is no need to log in.
- In a workspace of your choice you can run IRAF, MIDAS, PYTHON, IDL, and other reduction software. Your data is stored as FITS files in the directory */data/raw/YYYY-MM-DD/*, where *YYYY-MM-DD* is the date at the beginning of the observing night. We suggest that you do all your reductions in this directory.
- This machine (*wg5off*) also runs the important *dataSubscriber* program, responsible for transferring the data from the instrument work station to the archive and to the off-line machine. Please make sure NOT to stop it!
- There is a number of quick-look data reductions tools, available on *wg5off*: for imaging and spectroscopic observations, for photometric calibration. They are described in the instrument web page.

On the terminal *wg5dhs*, do the following:

- If required, login to this terminal as *visitor*. The day-time telescope operator will provide you with the password.
- In a workspace of your choice, click the left mouse button on the background window free area of the screen to start a terminal and from the prompt type “vOT” to start the vOT tool. Next, you will need to login with your username and password. They are the same as the ones you were given to login in the web-letter service, where you had found the notification for the approval of your program. If you do not have the username and the password, ask the day-time telescope operator or contact the User Support Department (usd-help@eso.org).

You are now ready to observe, and even to obtain some calibrations before the beginning of your observations. You can use them to reduce the data during the night.

4.3.1 Image Analysis

The optics of the NTT are actively controlled, and to get the best out of the NTT, frequent corrections to the primary are necessary. The telescope aberrations are determined by examining images of a guide star recorded through a Hartmann mask. The process is called image analysis. Its outcome is used to compute corrections to the optical configuration, which are applied as deformations to the primary mirror and as displacements of the secondary along three axes. Preferably, the image analysis should be done with the telescope pointing in the direction of your next target.

As most exposures in the IR are relatively short, the image analysis is done off line. Time spent doing image analysis is not spent on your source, thus image analysis should be used wisely. It usually takes 10 to 20 minutes to perform an image analysis and verify the result. Once done the telescope will be at the correct focus. It is no longer necessary to check the focus after image analysis.

The image analysis is done by the TIO and the only action a visitor can take is to request an image analysis of the image if the quality is degrading. How frequently one does an image analysis depends on the required image quality. Here are some hints on how frequently it should be done.

i) If image quality needs to be better than 1.0", image analysis should be done if the altitude changes by 20 degrees (sometimes 30 degrees; trust the advice of the telescope operator). For long integrations, where the source is changing altitude, one should perform the image analysis at an altitude that is 20 degrees higher if the source is rising or 20 degrees lower if the source is setting. The analysis should be repeated once the altitude of the object differs from the altitude where the previous image analysis was done by more than 20 degrees. In this way the object is never observed more than 10 degrees from the last image analysis.

ii) If the image quality is not critical, (this is usually the case for spectroscopic observations) then perform the image analysis once at the beginning of the night. It should be done with the telescope at an altitude where the bulk of the observations will be done.

iii) Good focus is very important. Poor focus can result in distorted images if any aberrations are present. In fact, distorted images are a good sign that either the telescope needs to be refocused or another image analysis needs to be done. The focus depends linearly with the temperature of the four serurriers. The coefficient is 0.079 mm of M2 movement per degree. In good seeing, a displacement in M2 of 0.03 mm induces a noticeable degradation in the image quality. So if the temperature has changed by 0.5 degrees or more since the last focus, one should refocus. On La Silla, the temperature changes quickest during the first hour of the night.

iv) During long integrations, greater than five minutes, it is possible to run the NTT in "closed loop" mode, that is, the optics are adjusted during the observations. For SOFI, integrations are usually shorter, and this is usually not needed.

NOTA BENE: The best way to find out if it is time for another image analysis is to monitor the image quality of the incoming raw data and to request an image analysis if you notice degradation! As a general guideline, however, you should trust the TiO's experience.

4.3.2 Focusing

The telescope and the instrument are focused separately. Nominally, the instrument is focused automatically to a pre-defined collimator position depending on the instrument configuration, and this process is fully transparent to the user. However, recently the collimator has shown tendency to get stuck, and currently it is set at a fixed position, optimal for wide field imaging and spectroscopy.

The telescope focusing can be done with a focus pyramid that produces five separate images. There is a specific template to do the focusing. Note that the image analysis normally corrects for the telescope focus (so the focus pyramid is used mostly for health check), and the only change of the focus that actually has to be done manually is the correction for the temperature of the primary mirror. All this is done by the TIO. The role of the visitor is to monitor the quality of the incoming raw data for any degradation.

4.4 The SOFI OS GUI Panel

The SOFI OS GUI panel (Figure 4.1) is used by the instrument operator to show the status of the instrument. It can also be used to setup the instrument, the detector, execute exposures and even point the telescope. A detailed description of this panel is given in the OS manual (VLT-MAN-ESO-14100-1510). In addition to the OS GUI

panel, there are various panels related to ICS, DCS and IRACE showing the status of the instrument and the detector. These panels are not of direct interest to the observer.

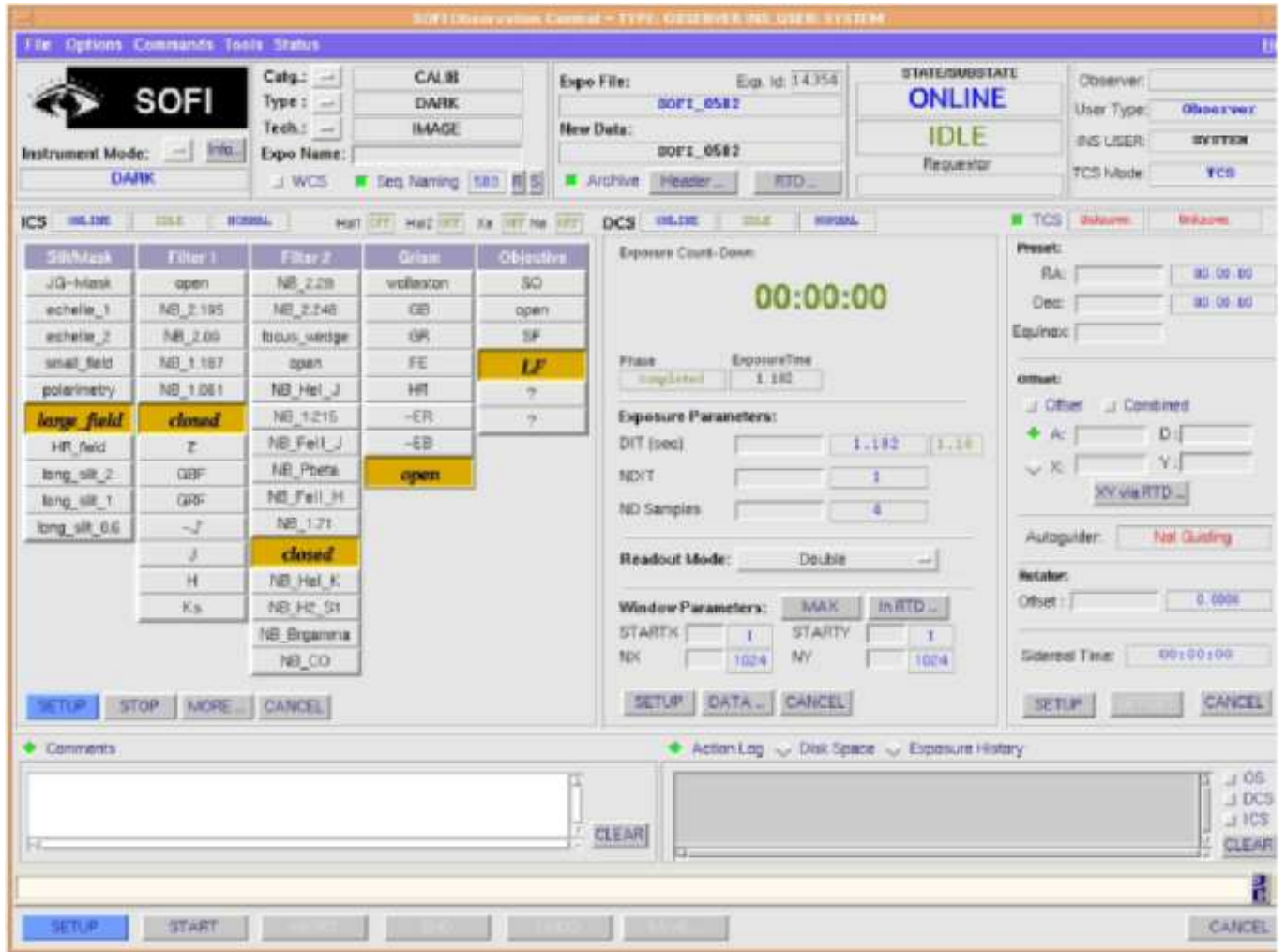


Figure 4.1: OS of SOFI

We do not allow observers to use this panel to start exposures, as it is far more efficient and safer to observe with OBs running on BOB. Furthermore, data archiving, a fundamental requirement at the NTT and the VLT, is only possible if frames are accurately classified because the images taken via the OS panel have incomplete headers. The OS GUI is used only for trouble shooting.

4.5 RTD

The SOFI RTD (Real Time Display; Figure 4.2) is a versatile tool of frequent use. As the array is continuously read out, images are continuously displayed on the RTD. During acquisition or when the instrument is idle, the data is still taken (with the last used DIT and NDIT) but it is not written to disk.

A very useful feature of the RTD is the ability to store a frame into memory buffer

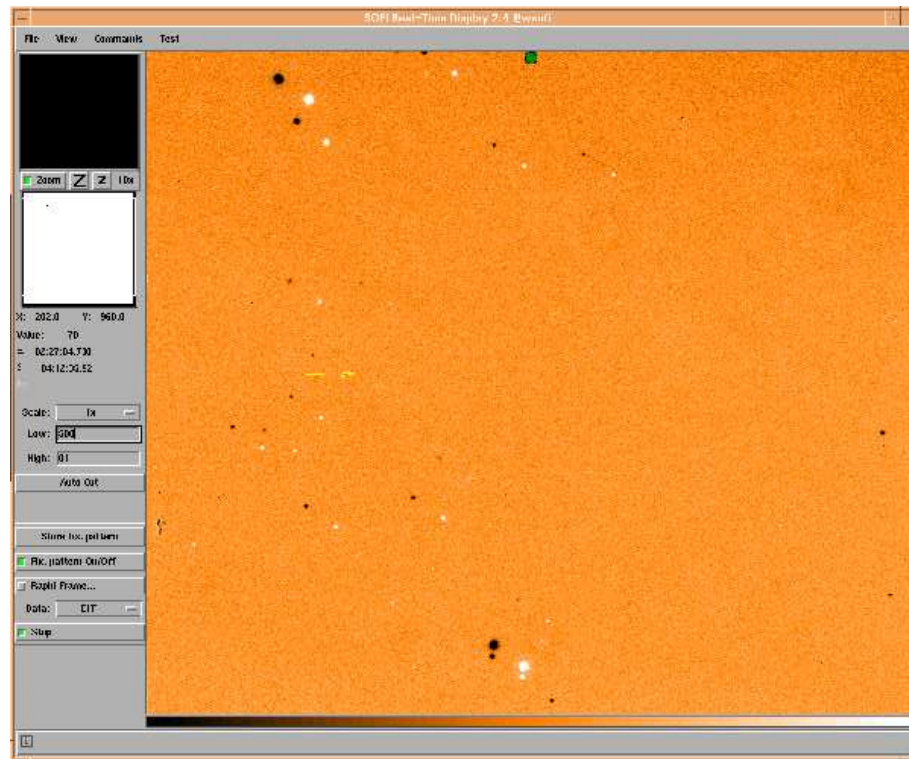


Figure 4.2: Real Time Display, with the “Store Fixed Pattern” option enabled. Note the offset between the negative and the positive star pattern. This offset is necessary to avoid self-cancellation of the objects in the field. The green dot on the top of the image flashes every time when the RTD updates the image. The yellow arrow is a remnant from a previous acquisition – it indicates the direction and the size of the offset that was carried out.

and to subtract that frame from incoming frames. In effect, this feature allows a simple sky subtraction. This is essential for recognising faint objects on backgrounds with significant structure. The stored frame is called a “fixed pattern” and to activate this feature, first click on the button labeled Store Fixed Pattern on the RTD screen and the current image will be stored. Next, click on ON/OFF.... to subtract the stored images from the following ones. You’ll probably have to click on the button labeled Auto Cuts or edit the cut levels directly to see the image. Clicking on Store Fixed Pattern again will load a new fixed pattern.

Naturally, it is best to store a fixed pattern before the telescope is centred on your science target. Acquisition templates “Move to ...” take an image away from the target, at an offset position defined by the user in the template. To expedite the process, one can even store a fixed pattern during the movement of the telescope (first make sure the DIT and NDIT have been updated to the values in the acquisition template!) but this may lead to dark/negative lines instead of dots on the subtracted image.

The RTD can display either each individual frame or it can display the average of NDIT frames (INT mode). The option is selected from the button labeled DIT or INT, respectively. The RTD can also display the image zoomed or reduced by an integer factor: clicking on the button to the right of label “Scale” will rise a pop-up menu with the available options.

The menus on the top left hand corner of the RTD can be used to activate sub-windows with various tools, which the observer may find useful:

- The pick object sub-window can determine the centroid and FWHM of selected objects, in pixels. It is useful for monitoring of the image quality.
- The statistics sub-window can be used to display the statistics of a region that can be defined by the cursor. It is useful for checking if the counts in the core of an image are below the non-linearity limit of 10,000 ADU.
- The cuts sub-window can be used to plot a trace that can be defined by the cursor, which is especially useful to check the count levels in spectra and the general background level. Zooming within this sub-window is possible with the left mouse button, the right mouse button un-zooms.

Last but not least, the RTD display is also used with the acquisition templates to place objects of interest into slits in spectroscopic modes, into the clear portion of the polarimetric mask in polarimetric mode, or to a specific position on the array in imaging mode. The interaction between the RTD and the acquisition templates is done by the TIO, after the user specifies his/her requirements.

4.6 The Data Flow Path

The raw data passes from the detector to the IRACE controller where pre-processing of the data occurs. Exactly what the pre-processing does depends on the readout method. For example, in non-destructive readout mode the pre-processing involves fitting a line through the individual readouts and multiplying the slope by the DIT.

Next, the data is written to the disk of the instrument workstation *wsofi* as FITS files, and it is displayed on the RTD.

Finally, copies of the data are sent to the archive machine *wg5dhs*, and to the off-line computer *wg5off*. During this transfer the file names are changed from “exposure name plus number” to “instrument name plus time stamp”. The archival machine is used to produce the data package the user is given at noon after the last night of the observing run. The offline machine is where the user can access the data, in a directory */data/raw/YYYY-MM-DD/*.

The data transfer processes are transparent to the user.

4.6.1 Sofi Pipeline

A VLT compliant SOFI pipeline is available at the NTT. It carries out on-line data reduction and derives zero points automatically, in nearly real-time. The observer can find the products on *wg5off* in the */data/reduced/YYYY-MM-DD/* sub-directory. The products of the on-line pipeline are not science grade because they were produced with out of date calibrations. They are meant only for quick and rough data quality assessment. A proper reduction can be carried out with the pipeline recipes which are available via Gasgano and esorex on the off-line machine *wg5off*, in interactive regimen. The user is also provided with number of tools for data reduction (i.e. MIDAS and IRAF scripts). They are described in the Sofi web page at:

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/index.html>
(see in particular the *Useful Scripts* link).

The pipeline recipes carry out the following reduction steps:

- **Imaging:** cross talk removal, sky subtraction, flat fielding, alignment and combination of images specified in a list provided through Gasgano or esorex; observations with large offsets at clear *sky* fields are not reduced.
- **Spectroscopy:** sky subtraction, flat fielding, alignment and combination of the 2-dimensional spectra, subtraction and wavelength calibration of an 1-dimensional spectrum of the brightest object. Naturally, this recipe works best for single bright objects, the fainter the object, the higher the chance of failure.
- **Zero Points:** cross talk removal, sky subtraction, flat fielding, simple aperture photometry. Note that the tool is sensitive to the position of the star on the array and if it is not centred within ~ 20 arcsec it may fail to recognise the standard.

The pipeline can be downloaded from: <http://www.eso.org/sci/software/pipelines/>

The most recent instructions how to use it are available on the SOFI web page:

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/SofIPipeline.html>

4.6.2 The Archive

All data taken at the NTT are archived into the ESO archive. For more details see:

<http://archive.eso.org/>

4.6.3 The Calibration Plan

The calibration plan aims to ensure that all data stored in the archive can be fully calibrated, regardless of the specific needs of the program. Indeed, it is easy to imagine a situation that the data may be used later – after their proprietary period expires – from the archive for the purpose of another object. The calibration plan ensures that any limitations in that respect will not come from the lack of appropriate calibrations.

Therefore, daily calibrations are taken as a matter of routine, both by the observatory support staff and in Visitor mode runs by the observers wishing to calibrate their own data. The standard procedure is to run every morning after the end of the night time observations a special tool *calobBuilt* that inspects all files taken during the night and determines what calibrations are necessary. Furthermore, this tool prepares a calibration OB that can be loaded directly into BOB and executed.

However, Visitor mode observers are encouraged not to rely on our calibration plan to provide them with these data, but to take their own calibrations in the afternoons before their observations start. This has a number of advantages: they would do their own quality control of the calibrations, and they will have at their disposal appropriate calibrations if they wish to reduce the data during the night.

The calibration plan includes:

- **Imaging:** special dome flats for all filters and imaging modes used.
- **Spectroscopy:** arcs and flats for each spectroscopic mode and slit combination; the flats are accompanied by images taken with lamps off.
- **Polarimetry:** flats.

Nota bene: Darks are not a part of the calibration plan.

4.7 At the End of the Night and at the End of Your Run

The raw data are immediately available in the ESO archive. However, if you want to save data, i.e. reduced one, we suggest you save a copy of your data at the end of each night. Do not wait until the end of your run. It is common to generate several Gb of data per night with SOFI! Leaving the backup until the end of the run may mean that you could miss your plane! The raw data is stored in the directories `/data/raw/YYYY-MM-DD/` on the machine `wg5off`, the reduced on `/data/reduced/YYYY-MM-DD/`.

At the end of your run you should fill out the *End of Mission Report*. This can be done via any web browser. You will find the form on the La Silla home page:

<http://www.eso.org/sci/facilities/lasilla/sciops/EoM/>

Chapter 5

Data Reduction

This chapter does not aim to be the ultimate guide to reducing SOFI data but it does outline the general steps, it does provides useful tips and in effect it can be used as a mini-cook book. It represents the experience the NTT team has gained in reducing data from SOFI. This section of the manual will evolve with time and we are very keen to hear of any suggestions that people may have.

The first and the last principle of data reduction is to look at the data. The data do contain all the answers you seek. We urge the users to experiment with different techniques, combination, and rejection algorithms, and to look at the quality of the final data product. We also warn the users, that pipelines, as good as they might be, are dangerous tools because they separate the observer from the data. Therefore, the products from any automatic data reduction tools must be treated with caution. First and foremost, one must understand the data: the specific instrumental effects, the stability of the weather conditions during the observations, etc.

5.1 Basic Concepts

Here we list some basic concepts for the SOFI data:

- There are two read out methods available with SOFI: NDR (non-destructive) and DCR (double-correlated). Data taken in one readout method should NOT be used to calibrate data taken in another (an exception are the spectroscopic arcs).
- The exposure time keyword represents the total integration time for a single image, i.e. DITxNDIT. As the image you receive is the average of NDIT exposures of DIT seconds, the correct number to use for flux calibration is DIT and not the value given in the exposure time keyword.

5.2 Imaging

5.2.1 Inter-quadrant Row Cross Talk

This is the first step of the data reduction but some specific cases described below require post-processing of the final image to remove cross talk residuals.

A bright source imaged on the array produces a “ghost” that affects all the lines where the source is and all the corresponding lines in the other half of the detector. For

instance, if the bright source is on row 300, the cross-talk affects the row $300 + 512 = 812$ and vice versa, if the source is on row 750, the other affected row will be $750 - 512 = 238$. The cross talk is seen either on high S/N data or if the bright source reaches near saturation level.

Though the effect is not completely understood it is well described and can be easily corrected. The intensity of the ghost is in fact 1.4×10^{-5} times the integrated flux of the line. An IRAF script for removing the cross talk is available from the SOFI web-page: http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/sofi_scripts/. The script has to be given input and output lists of images.

The script can remove the cross talk only if the source that causes it is not saturated! In case of saturation, the counts of the source will be decreased by the non-linearity and the subtraction of the ghost will be incomplete. However, even in those cases the users can improve the data. The cross talk manifests itself on the final reduced image as: (i) residual of the ghost, and (ii) residual of the line leaks, i.e. horizontal lines with increased counts. While there is nothing one can do for the former effect, short of masking, the latter one can be removed. This can be achieved in a few steps:

(i) block-average the image along the horizontal lines; note that the IRAF task *blkavg* is not well suited for this purpose because it does not allow to use median averaging or to impose rejection and lower threshold limits. Instead, it is better to use *imcombine* feeding it with a list that contains each row of the image as individual entry: `final_image[1,*]`, `final_image[2,*]`, `final_image[3,*]`, etc. The use of a median averaging, and a suitable rejection algorithm with low value of the upper threshold is necessary to exclude from the averaging any sources. The result from this step will be a 1-dimensional vector-image.

(ii) expand the vector-image to the full size of the final image by replicating it; the IRAF task *blkrep* is well suited for this purpose. The result from this step is a 2-dimensional image with a size matching the size of the final image.

(iii) subtract the image produced in the previous step from the final image; ideally, this should remove the residuals from the horizontal lines. However, one should exercise caution because any imperfect flat fielding may affect the procedure, leaving artificial gradients across the image.

As explained above, in many cases the cross talk effect can not be removed completely, but sometimes the field of view can be rotated so that the ghost image can be relocated to an area of the field away from the target. It is enough to relocate the bright source on a different row or column than the target, and on a different row than the one corresponding row in the upper/lower half of the array.

5.2.2 Masking the Bad Pixels

The SOFI detector array suffers only minor cosmetic defects – it has $\sim 0.1\%$ or ~ 1000 bad pixels. Some bad pixels masks are available in the SOFI web page:

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/bad.pix.html>

We consider the following types of bad pixels:

- Dead pixels : they are computed from the histograms of flats taken during SOFI Calibration Plan; the distribution of the pixels values is fitted with a Gaussian (mean, sigma), and we select the dead pixels as the ones whose values are 4σ below the mean value.

- Hot pixels : they are computed from the histograms of darks taken during SOFI Calibration Plan; the distribution of the pixels values is fitted with a Gaussian (mean, sigma), and we select the hot pixels as the ones whose values are 4σ above the mean value.
- Noisy Pixels : they are computed from the histograms of the sigma image obtained when averaging the flats and the darks of the previous 2 steps.
- Frame Pixels : the mask also marks those pixels as bad that are masked by the edges of the focal plane field mask.

NOTA BENE: The masking must occur before the images are shifted and aligned!

It is up to the user how to exclude the masked pixels. Perhaps, the simplest way is to use two masks. The first one will have zero values at the bad pixels and one at good ones. Every raw data frame has to be multiplied by this mask, ensuring that all bad pixels will have the same value. The second mask will have zero values at the good pixels and a large negative value, i.e. $-1e6$ at the bad values. This mask will have to be added to the row frames after the multiplication. It will set the bad pixels to a value outside the “good” pixel values for SOFI.

Of course, the fractional offset during the image alignment may move some of this artificial negative flux into the acceptable values but the large value makes sure that this will happen only rarely.

The final step of the application of the bad pixel mask is to set an appropriate lower threshold during the combination of the aligned images.

Indeed, some data reduction tools like IRAF have the capability to incorporate a bad pixel mask and the shifting of the images into the combination step. However, the *imcombine* task can carry out only integer shifts and this will degrade the spatial resolution of the image. This is the reason we recommend a 2-step shifting – first with fractional offsets, and then with integer offsets (see Section 5.2.6). The only possibility to apply the bad pixel mask in this case is before any of the offsets, as described here.

5.2.3 Subtracting the Dark/Bias Frame

First and foremost, a SOFI dark frame changes with the DIT, so if a dark frame is subtracted from other frames, the DIT must be identical.

Secondly, the underlying bias pattern to any image is a function of the flux incident on the array, so it is not always useful to subtract separately the dark frame. In fact, the bias and the dark contributions to the total signal are almost impossible to separate without special efforts. In cases where the mean count level is low, for example narrow band images, the dark frame is a good approximation to the underlying bias pattern. In cases where the mean flux level is high, for example broad band images, the dark frame is a poor approximation to the underlying bias pattern, so subtracting it serves no real purpose.

Therefore, we suggest to avoid separate dark subtraction. Instead, the user can subtract the darks together with the sky, as described in the next section. The advantage of this approach is that the overall illumination is nearly constant, so the varying bias level problem mentioned above will be avoided.

The only case when the dark/bias subtraction is still necessary is the reduction of the flat field calibration images. See the description of the special procedure for that in

Section 5.2.5.

5.2.4 Sky Subtraction

This is a critical operation, which depends on how many frames are available to create a *sky* image. The experience shows that at least 5-7 images are necessary to achieve acceptable sky subtraction, and this number depends strongly on the nature of the sky field: if the field is crowded with stars (i.e. the Galactic plane), this may not be enough. In such cases we recommend to make sure that at least 12-15 images at different locations are available. In general, the more frames you have, the better the results (but up to a point, because of such effects as pupil rotation and increasing overheads).

There are two basic strategies for sky subtraction:

(i) creating a single *sky* image, and then subtracting it from each individual frame. Here is a step-by-step description of the procedure to remove the sky from a stack of n jittered frames, where n is large.

(a) For each image, take the 10 images that were taken closest in time. These images will be used to estimate the sky. In case of jittered observations, when the target field is on all images, these 10 images will indeed include the target. In case of observations alternating between the *object* and the *sky* the list must include only the 10 closest *sky* images.

(b) Scale these 10 images to a common median. To first approximation, this is the *sky* image. However, it is possible to improve it. One possibility is simply to reject the most deviating values, i.e. rejection of the highest and the lowest 1 or 2 values, 3σ rejection, etc. We strongly encourage the user to experiment here. Another option is a more sophisticated 2-iteration rejection. After making a zero-order *sky* as described above, return to the list of the images you used to create this *sky*. From each of them subtract the *sky*. The result will be images with nearly zero background and only the stars above the background. Study the distribution of pixel values with image histograms and determine the border value between the background and the stars. Usually, the limit is ~ 3 times the background standard deviation. Next, mask out the values above this limit. The easiest way to do it is just to replace the values above the chosen limit by a high or low value (i.e. $1e5$ to $-1e5$) so when you combine the “masked” images to produce a second iteration *sky*, these pixels will remain out of the range of the “good” values used in the combination. The “good” values are usually defined as parameters of the combining task. Be careful before modifying the files with replacement. It is better to copy them into dummy files first!

(c) Scale the *sky* image to have the same median as the image from which it will be subtracted and then carry out the subtraction.

This method reduces the negative traces of stellar images on the sky-subtracted data.

The number of 10 images used to create the sky is somewhat arbitrary. It can be 5-7-... or even all available images.

(ii) subtracting from each image the average of the previous and the next image. In the more general case, the average (or the median) of the previous 2-3-... and the following 2-3-... images is subtracted. Strictly speaking, this is a specific case of the method described above. However, our experience shows that in most cases using only the two images nearest in time is sufficient.

The advantage of this technique is that it takes into account sky variation on time scale

a few times $\text{NDIT} \times \text{DIT}$. In the specific case of SOFI it also removes very well the pupil ghosts (usually observed as diagonal “stripes” in the South-East corner of the array) discussed earlier. Both these advantages are pronounced better if the number of images we use to subtract the sky is kept small - ideally, just the preceding and the succeeding images. While the reason for this is obvious for sky variations, it is not the case with the pupil ghost until one remembers that the pupil rotates and the ghost moves across the field of view. Therefore, it is best to sky-subtract using images that have the ghost nearly at the same position as on the target image.

The problem with this method is that some negative residuals from stellar images will remain on the sky-subtracted images. Indeed, with only the preceding and the succeeding images to subtract the sky, it is impossible to remove the “negative stars”. This difference with the previous method makes us to consider it separately. Note that the “negative stars” can be dealt with later, on the stage of combining the shifted object images. Given enough object images and sufficiently low crowding, this technique works well.

This method is applicable for observations taken in both simple jittering mode, and with offsets on a clear *sky* region.

A special attention is required to observations obtained following the schemes described in Section 3.2.5. As pointed there, the sky sampling in those observing strategies is somewhat rarefied and the user should check the residuals from the sky subtraction.

In case of observations in which the target has been kept on the array during every image, the user should mask the target out before constructing a sky image. This can be done in many ways, perhaps the most straightforward solution is to follow the same procedure as for masking out the bad pixels. In the case of 4-point observations described in Section 3.2.5 you have to create four pairs of masks: one with zeros in the target quadrant, and one everywhere else, and another, with a large negative value, i.e. $-1e6$, at the target quadrant, and zero everywhere else. Then every image has to be multiplied by that of the first masks that matches the location of the cluster, and summed with that of the second masks that again, matches the target quadrant. Remember to produce those manipulations on a copy of the original data. Finally, a suitable low threshold in *imcombine* will make sure that the target has been masked out from the sky.

An alternative to this masking is to use a “floating” sky, constructed from the six (in case of the 4-point scheme) images immediately preceding and succeeding every image. Those images have the target located on different locations than on the sky-subtracted image. As usual, the “global” average is preferable if the sky background is stable, and the “floating” sky subtraction is preferable if the sky varies.

NOTA BENE: The sky subtraction of every observation has to be considered individually! The recipes provided here refer to the most general cases. Every observation is different and covering all possibilities is out of the scope of this manual.

We provide some tools that help to carry out the sky subtraction. Check the *Data Reduction* section of the SOFI web page, and if necessary, ask the telescope operator what options are available to receive help—usually they are limited to calling Paranal on the phone or to sending an e-mail to the SofI Instrument Scientist (the contact details are available on the SofI web page).

5.2.5 Flat Fields and Illumination Corrections

Creating a good flat field for SOFI data is difficult. The simplest way of creating flat fields is to subtract an image of the dome flat field screen with the dome lamp off from an image with the dome lamp on. This flat field has two short-comings. First, the shade pattern of the array is a function of the overall flux, so the shade pattern in the image with the lamp on is different from that in the image with the lamp off. Thus, the difference of the two will contain a residual shade pattern. Second, the illumination of the dome panel is slightly different from that of the sky. Furthermore, it changes with time because of the ageing of the lamps used to illuminate the screen, and in the case of observations with $\lambda \geq 2.3$ micron, due to the variations of the dome temperature. Both these effects are at the 1-2% level and both can be removed.

The residual shade pattern can be removed if the `SOFI_img_cal_SpecialDomeFlat` template was used (see sec. 3.5.2). This template takes the usual sequence of images with the dome lamp on and off and, in addition, it takes the same set of images with the mask wheel vignetting the array (Figure 3.3). The vignetted part of the array is relatively free of scattered light, so it can be used to estimate and remove the shade pattern. However, this estimate of the shade pattern is valid for the frame that is vignetted and it is not valid for un-vignetted frame. The difference can be estimated from a region that is common to both the vignetted and un-vignetted frames. The method is outlined below.

- (i) Collapse the following image sections along detector rows to form one dimensional images: columns 50 to 150 in the vignetted image, columns 500 to 600 in vignetted image and columns 500 to 600 in the un-vignetted image. For reference, call these images A, B and C, respectively.
- (ii) The difference between the shade pattern in the vignetted and un-vignetted images is C-B.
- (iii) The shade pattern in the un-vignetted frame is then $A + (C - B)$ and this one dimensional image should be subtracted from the un-vignetted frame.
- (iv) Steps (i) to (iii) should be done for the images with the lamp on and the images with the lamp off. Once the shading pattern is removed for both, one subtracts the image with lamp off from the image with the lamp on to form the flat field.

Small IRAF and MIDAS procedures performing these steps are available at the telescope and from the SOFI web page in the *Data Reduction* section, named `special_flat.cl` and `midas_specialflat.prg`, respectively.

The illumination correction removes the difference between the illumination pattern of the dome flat screen and the sky. It is derived from a grid of 9-16 observations of a star (preferably a standard) across the field of view. The illumination correction is created by fitting a plane to the fluxes (not magnitudes!) of the star, after it has been flat fielded. Therefore, each illumination correction surface refers to a particular special dome flat and it can not be used with another special dome flat!

Special Dome Flats and Illumination Correction Surfaces for all broad band filters for Large Field mode are prepared by the observatory staff monthly. The latest versions can be downloaded from the Web Page of SOFI:

<ftp://ftp.eso.org/web/sci/facilities/lasilla/sofi/fits/Archive/>

and fits images for a number of 16×16 grids can be found at:

ftp://ftp.eso.org/web/sci/facilities/lasilla/sofi/fits/Archive_Illum_Corr/

An IRAF script for producing Illumination Correction Surfaces – `illumination.cl` – is offered to the users. It can be downloaded from the SOFI web page:

http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/sofi_scripts/

Twilight sky flats and flats created from the observations themselves can also be used to flat field the data. However, these frames suffer from the same problems as the uncorrected dome flats and, in the case of the shade pattern, it is not clear how they can be corrected.

Last but not least - the Special Dome Flats and the Illumination Correction Surfaces are mode and filter specific, i.e. they have to be prepared separately for each instrument mode and filter combination.

5.2.6 Image Alignment and Combination

Most data reduction packages have tools for alignment of images and for their combination. For example in IRAF one can use *imexam*, *imalign* and/or *imcombine*. We refer the user to the corresponding user manuals.

Note that in case of IRAF, the task *imcombine* can not handle fractional offsets. The workaround is to use a 2-step shifting procedure. First, shift the images by the fractional part only with the task *imshift*. If you intend to do photometry on the final image, set the interpolation parameter of *imshift* to linear in order to conserve the flux during the shifting. In addition, set the boundary value that will be given to pixels with no value to constant, and set that constant to a large negative value, i.e. $-1e6$. Later, during the combination these pixels can be excluded with an appropriate lower cut-off threshold. If you use any of the outer options for the boundary pixels, you will end up “creating” data.

The second step of the shift can be carried out with *imcombine*, using the offset facility of this task. The final image will have the largest possible dimensions, and the signal-to-noise ratio will most likely degrade toward the edges because of the smaller and smaller number of images that cover the outer regions. Remember to set the lower threshold of *imcombine* appropriately, to exclude both the marked bad pixels and the pixels that were given bogus values at the edges of the shifted images.

5.3 Long Slit Spectroscopy

5.3.1 Inter-quadrant Row Cross Talk

This step is a must for reducing spectroscopic observations because unlike the imaging, the spectra are always nearly aligned along the array columns, so a ghost of the upper half of the spectrum always appears on the top of the lower one, and vice versa, unlike the imaging when the observers may be fortunate not to have bright star ghosts on top of their targets (or they can even rotate the instrument to avoid it). The effect is particularly dangerous when one looks for emission lines because their ghosts are emission line like and can be confusing.

More details on the cross talk can be found in Section 5.2.1.

5.3.2 Sky Subtraction

In most cases, the user simply has to subtract one image from another. This should remove most of the night sky emission. The subtraction of the average of every preceding and succeeding image works particularly well. The residuals can be removed during the extraction of a 1-dimensional spectrum.

The tool for sky subtracting of imaging data works for the spectroscopy as well.

5.3.3 Flat Fields

Spectroscopic flats are taken with the dome flat field screen. Like imaging flats, one subtracts an image with the dome lamp off from an image with the dome lamp on. These flats also suffer from the shortcoming that the shade pattern is not perfectly removed.

However, the slits do not cover the entire chip; there is a region of approximately 50 pixels wide which is free of direct illumination and it may be possible to use this region to estimate the residual shade pattern and to correct the flat.

5.3.4 Removing Slit Curvature

Spectra taken with all SOFI grisms show slit curvature. For the red and blue grisms, the curvature amounts to a few pixels from the middle of the slit to the edge, and is well fitted with a quadratic. For point source observations, removing the slit curvature is an unnecessary step. However, extended objects may require to correct it in order to avoid degradation of the spectral resolution.

The correction is done by doing a 2-dimensional wavelength solution to arc spectra and both MIDAS and IRAF have tasks to do this. For the red and blue grisms a quadratic in the slit direction and a cubic or quadratic in the wavelength direction is adequate. There are no cross terms. An excellent description of the method and the IRAF tools to do it can be found in the cookbook *2D-Frutti Reductions with IRAF* by Mario Hamuy and Lisa Wells. The relevant section is cited verbatim as Appendix A in *A User's Guide to Reducing Slit Spectra with IRAF* by Phil Massey, Frank Valdes and Jeannette Barnes. These documents are available on the WWW but for user's convenience we have placed copies of them at:

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/doc/index.html>

For the medium resolution grism, the slit curvature is larger. As yet, we have not tried to do a 2 dimensional wavelength solution to arc spectra taken with this grism.

5.3.5 Arcs

For the red and blue low resolution grisms, it is sufficient to use the Xenon lamp for wavelength calibration. For the medium resolution grism, both the Xenon and Neon lamps should be used.

The main Xenon and Neon lines are identified in Appendix A. For the medium resolution and red grisms, a cubic fit to calibrate the dispersion is adequate. For the blue grism, a quadratic fit is better.

5.3.6 Removing of the Atmospheric Absorption Features and Flux Calibration

The IR spectra are dominated by atmospheric absorption features. To remove these features it is customary to observe a star with a featureless spectrum at a similar airmass to that of the target. This star is called telluric standard. The reduction involves the following steps:

- the spectrum of the target has to be divided by the spectrum of the telluric standard. This will remove the atmospheric absorption because – presumably – it affects both spectra in the same way.
- the science spectrum has to be multiplied by the true intrinsic spectrum of the standard. This will remove the artificial “emission” features introduced to the target spectrum by the intrinsic absorption in the spectrum of the standard.

Removing the atmospheric features is a critical operation. In general, the spectra of the atmospheric standard and the spectra of the science target will not have exactly the same wavelength scale. They may differ by as much as half a pixel. This could be caused by internal flexure within the instrument or by the science object and the standard being on different parts of the slit. Thus, before dividing the standard into the science target, one should use the atmospheric absorption features in both to realign the spectra. The IRAF task *telluric* may be useful for this.

The true intrinsic spectrum of the standard is usually not known with any great precision. As discussed in Section 3.4.3, usually the telluric standards are selected to be wither solar analog (because a good Solar spectrum is available) or early type stars (because their spectra are easier to model). The IR Solar spectrum is available from the National Solar Observatory (“An atlas of the solar spectrum in the infrared from 1850 to 9000 cm⁻¹” Livingston W. & Wallace L. N.S.O., Technical Report #91-001, July 1991). A more detailed description how to use solar analogues as telluric standards and an IRAF-based tool can be obtained from Maiolino, Rieke & Rieke (1996, AJ, 111, 537). An associated IRAF tool is available on the Sofi web-site:

http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/sofi_scripts/

Recently, the fist library of flux-calibrated near-infrared spectra became available (Rayner et al. 2009, ApJS, 185, 289; Cushing et al. ApJ, 2005, 623, 1115), and it can be used to remove the artificial “emissions” for telluric of spectral types other than the Solar. For an example on how to treat early spectral type telluric standards, and for an empirical library of spectra of early type stars refer to Hanson et al. (1996, ApJS, 107, 281).

Good empirical IR spectrophotometric standards are not available. The only remaining way to do flux calibration is to observe the science target in one or more broad band filters and to scale the spectrum so that it agrees with the broadband fluxes. This method can achieve an absolute calibration of 5-10%. For an example of IR flux calibration in the more general case of extended objects see Ivanov et al. (2000, ApJ, 545, 190). See also the discussion on the flux calibration in Section 3.4.3.

5.3.7 Alignment and Combination

There are two general approaches to combine spectra:

- (i) to combine 2-dimensional images after an appropriate geometric correction that will align the wavelength axis and straighten the spectra

(ii) to combine 1-dimensional spectra (usually after wavelength calibration) that were extracted from the 2-dimensional images.

Every one of these methods has certain advantages and disadvantages, and the user will have to decide what is the best way to reduce the data. Obviously, the first method is easier for reducing extended objects because it will become trivial to extract spectra at different locations of these objects. However, the geometric correction involves heavy modification of the data, and some degradation of the resolution may be possible. The second method is easier to apply on point source spectra, and it has an added advantage of giving a straightforward estimate of the observational error - this is just the standard deviation obtained during the combination.

5.4 Polarimetry

As in the case of imaging, you should start to use a bad pixel mask. Some bad pixel masks can be retrieved from:

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/badpix.html>

Then, you have to apply the flat field to correct for the pixel-to-pixel sensitivity variations across the chip. The flat field exposures should be performed using the same filter as the object exposures and the Wollaston prism in the optical path. Since the whole instrument rotates in order to provide a certain orientation of the Wollaston prism, separate exposures at different orientations are not required. We recommend to take dome flats, because when a sky map is used, some ghosts of the object on consecutive images may result in a wrong flat-field, especially at the position(s) of the object. Exposures with the lamp OFF and with the lamp ON must be combined to obtain the final flat field:

$$FF = (lamp_ON_a - lamp_OFF_a) + (lamp_ON_b - lamp_OFF_b)$$

The final flat field is polarised, i.e. the median value of the upper field is different from the median value of the lower field. However, this problem can be overcome by an independent normalisation of the flat field in each field. A further artificial polarisation is introduced due to the deviation of the transmission ratio of the Wollaston prism (ideal: 50%:50% = intensity upper:lower beam). Three alternatives for the correction of the wrong transmission ratio are proposed (to see them in detail, please read the technical report available at

http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/Georges_Sofi_Polarimetry.pdf

- Use the transmission ratios between the two beams of the Wollaston: of $C_{K_s} = 0.968$ and $C_J = 0.954$, respectively for K_s and J . Ratios for the other filters have not been determined as of yet.
- Since the sky is also affected by this effect, a good measure for the intensity ratio of both beams can be derived from the median of the lower/upper image. This method may not be applicable if the polarised object covers a large fraction of the lower/upper image (because in this case the median value may be influenced by a possibly polarised object) and moreover, this method only works properly if the background radiation of the sky is unpolarised so the images are not affected by scattered moonlight.

- Fitting a $\cos(2\theta)$ function to the intensity of the object on the image as a function of the orientation (θ) of the Wollaston prism (see Ageorges 2000;
<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/polarimetry/>).

The next step is the sky subtraction: since the intensity of the background radiation is an additive component to the intensity of the object, the subtraction of the sky will also result in the removal of the contribution of the sky polarisation to the net polarisation of the object. Therefore, the sky has to be estimated and subtracted independently for both the upper and the lower field.

Finally, the remaining instrumental polarisation is expected to be mainly caused by the reflection on the tertiary mirror M3. It should therefore depend on the altitude angle of the telescope. However, the remaining instrument polarisation was found to be $< 0.3\%$ for the K_s band, and $< 0.4\%$ for the J band (see Appendix B of the technical report for more details). Since the (statistical) error of the results from which these limits have been derived is in the same order of magnitude, a possible contribution by the mirror M3 could not be extracted.

If you plan to do polarimetric observations with SOFI we recommend to read carefully the three following reports:

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/report.ps>

http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tecdoc/tech_rep_polarimetry.ps

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/polarimetry/>

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/reduction/polarimetry/sofipolmode.html>

Please note that this manual gives only a short description, consult the polarimetry reports from the SOFI web page for more details.

Appendix A

Calibration Arcs

The adapter contains both Xenon and Neon lamps. The Xenon lamp produces an even spread of lines for both the red and blue grisms. It is well suited for wavelength calibration. Fig. A.1 show the main Xenon lines for the blue grism. There are two electronic ghosts, caused by the very bright lines near one micron, between 1.35 and 1.4 microns. Fig. A.2 shows the main Xenon lines for the red grism. The continuum in the red is thermal emission from the lamp.

Files with the Xe and Ne line wavelengths are available from the SofI web page:

<http://www.eso.org/sci/facilities/lasilla/instruments/sofi/tools/Neon-XenonLines.html>

Both the Neon and Xenon lamps should be used to calibrate the medium resolution grism. Figures A.3, A.4, A.5 and A.6 show the main Xenon and Neon lines.

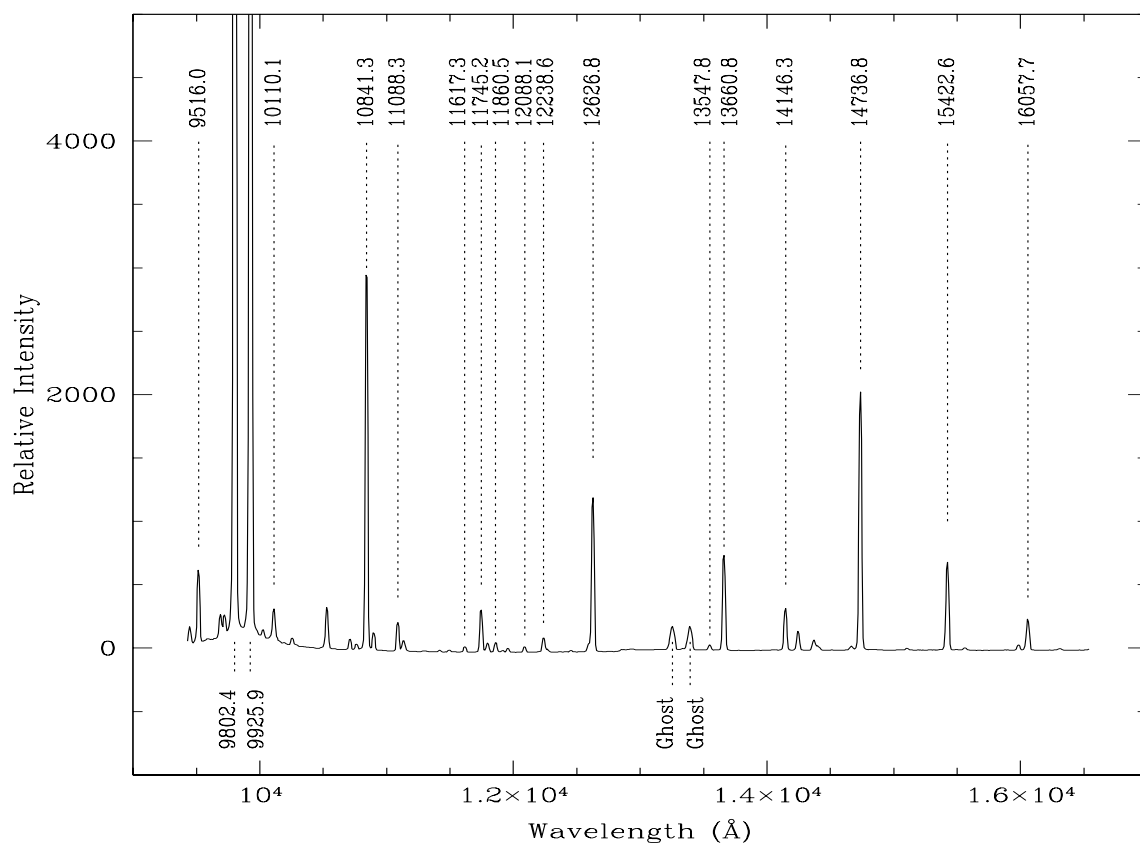


Figure A.1: A Xenon arc spectrum taken with the blue grism. The main lines are marked.

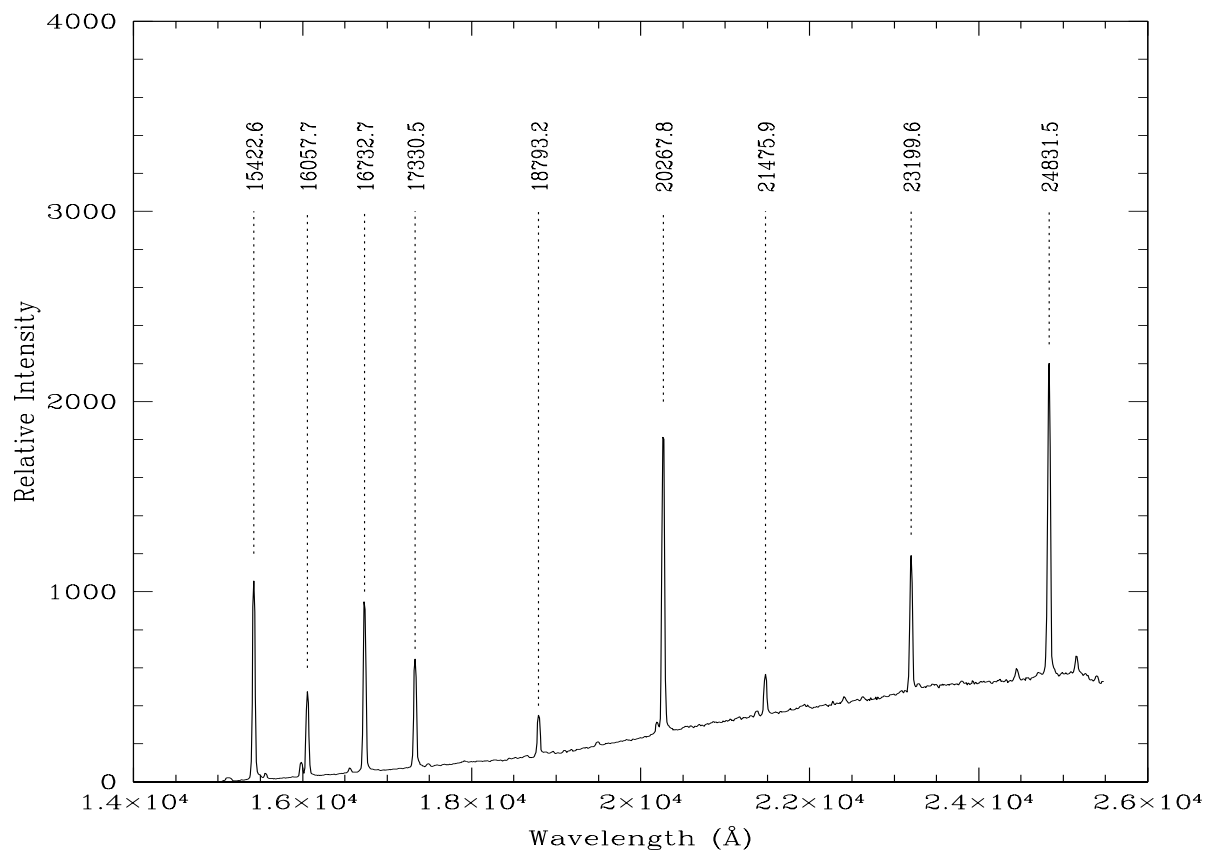


Figure A.2: A Xenon arc spectrum taken with the red grism. The main lines are marked.

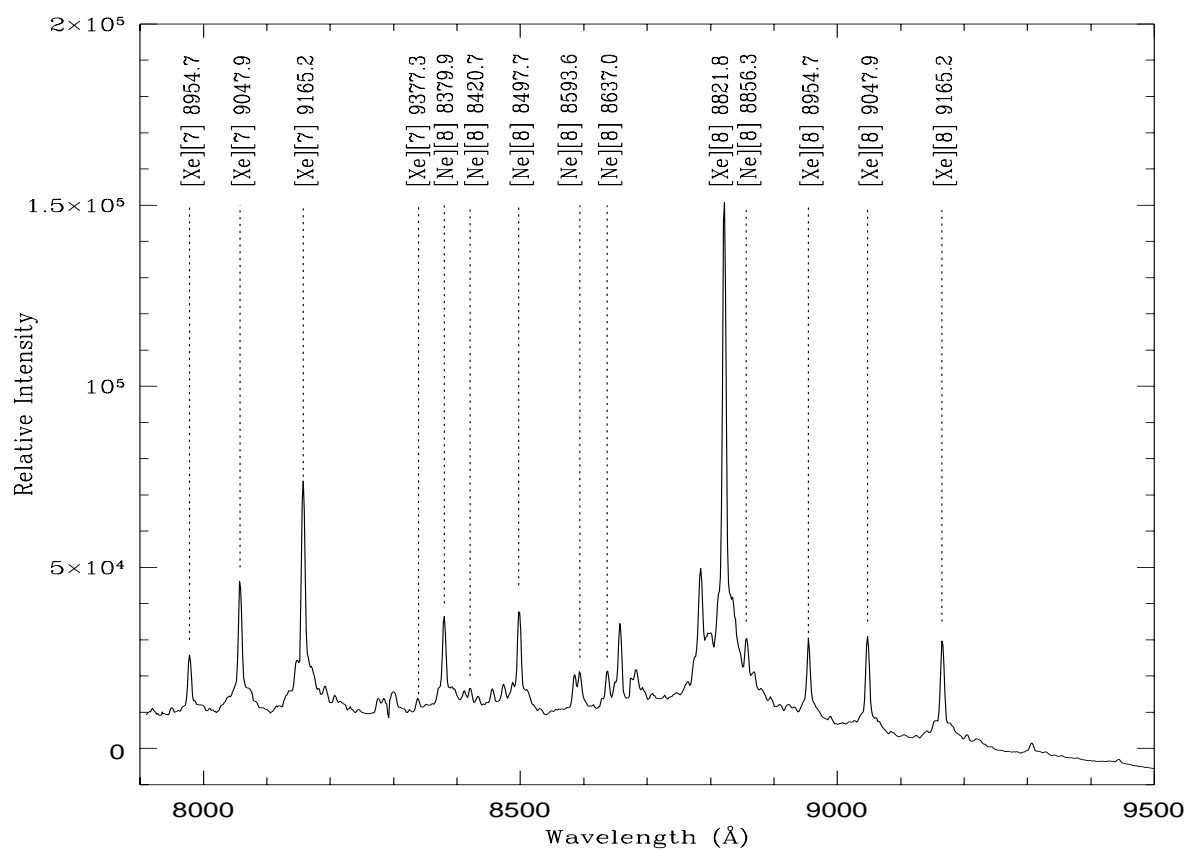


Figure A.3: A Xenon and Neon arc spectrum taken with the medium resolution grism at the Z atmospheric window. The main lines are marked.

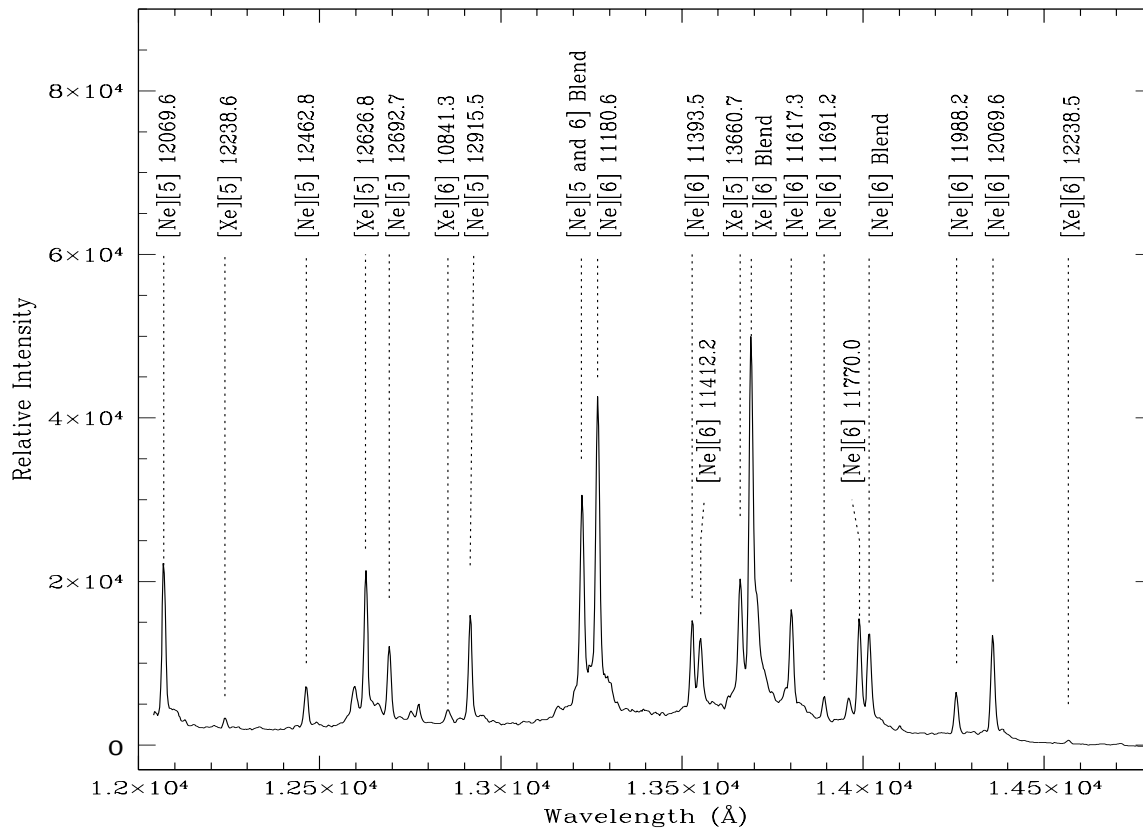


Figure A.4: A Xenon and Neon arc spectrum taken with the medium resolution grism at the J atmospheric window. The main lines are marked.

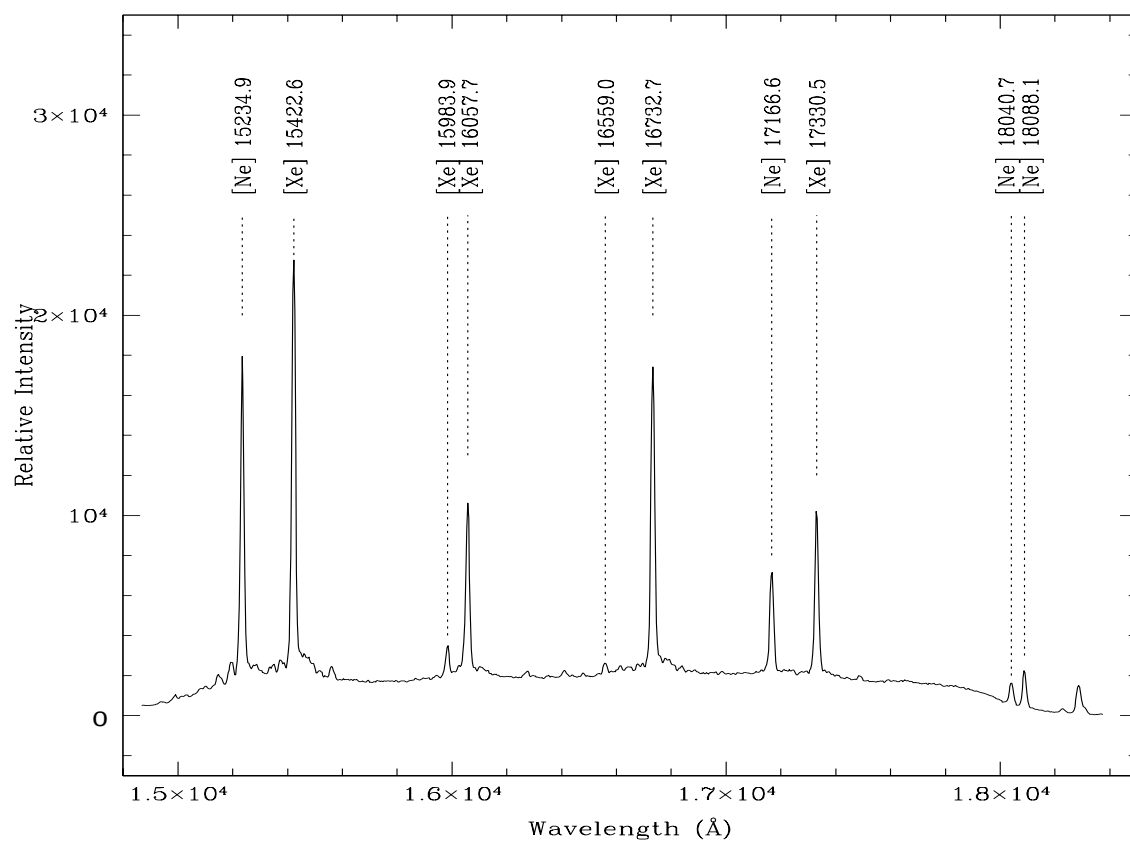


Figure A.5: A Xenon and Neon arc spectrum taken with the medium resolution grism at the H atmospheric window. The main lines are marked.

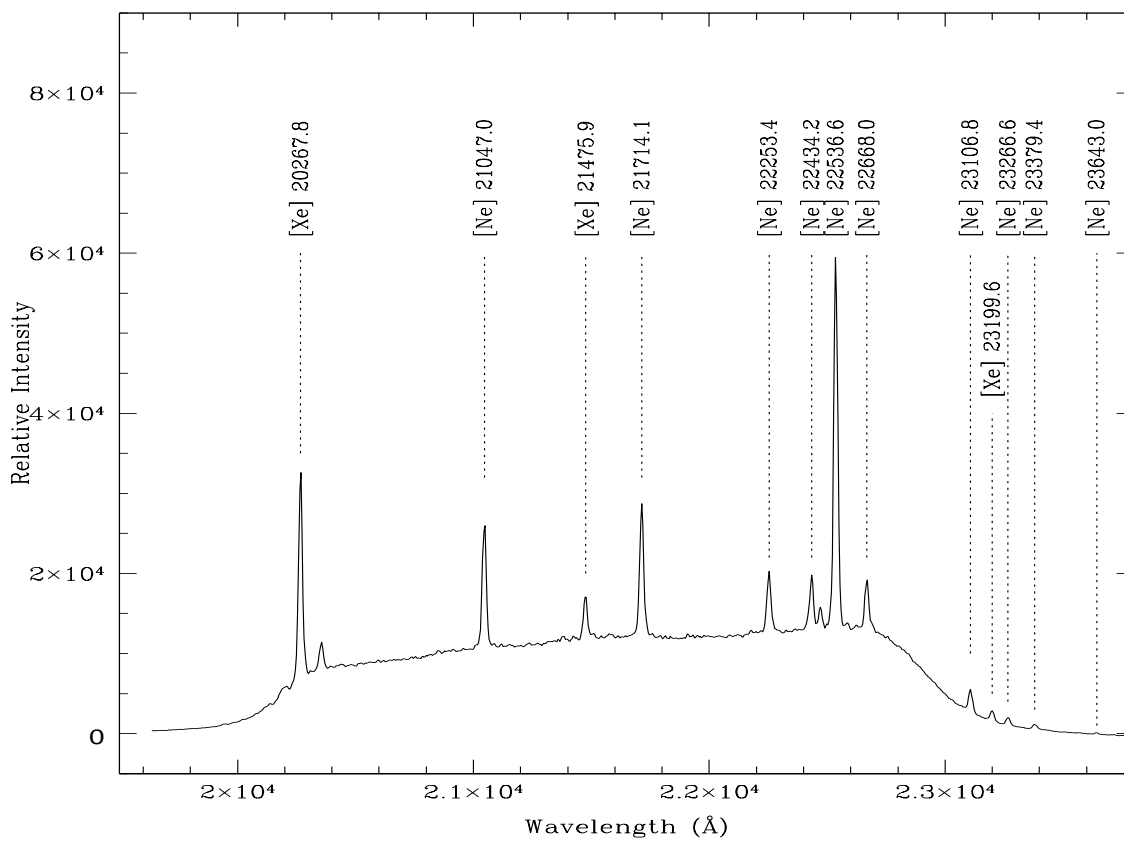
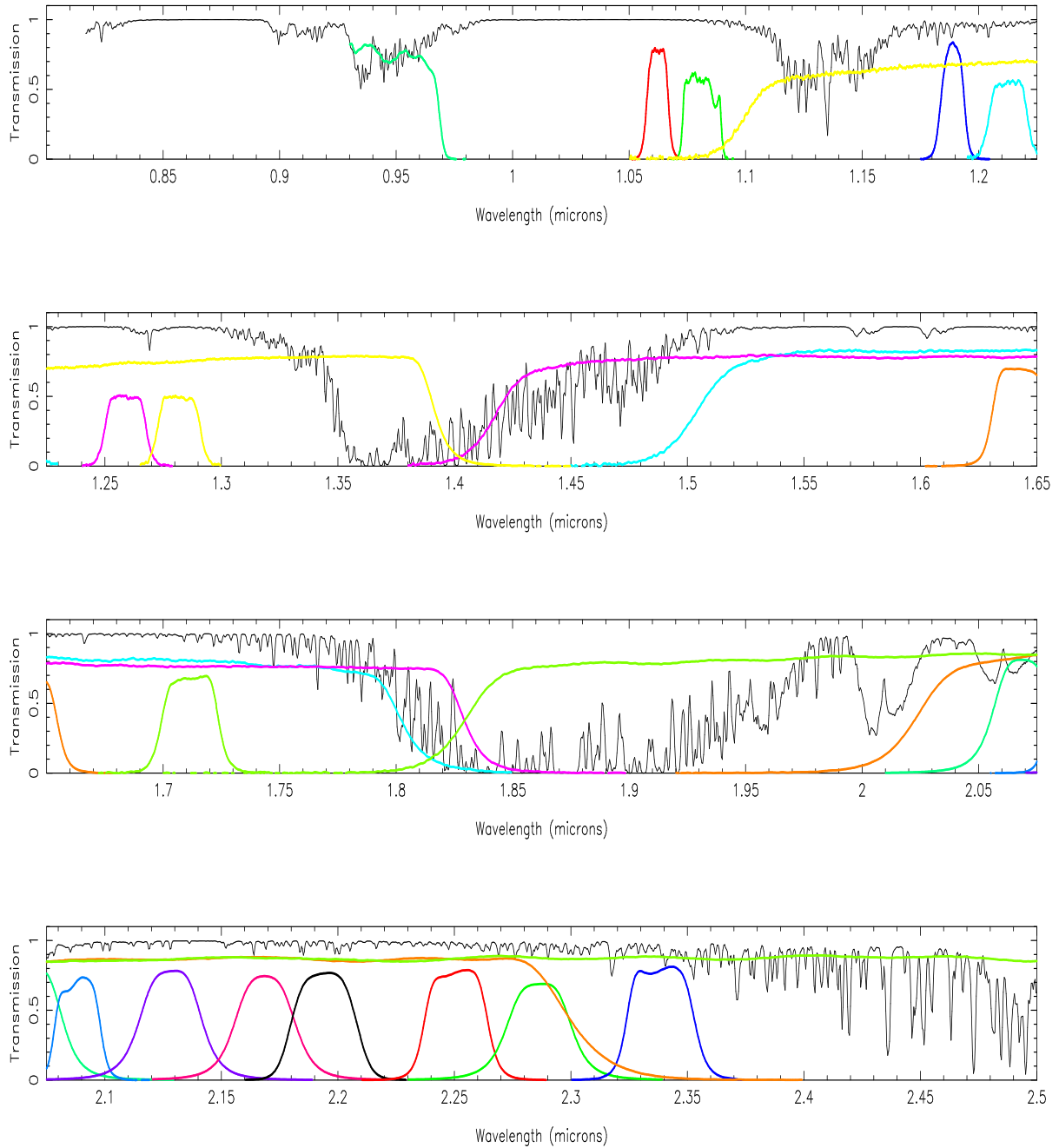


Figure A.6: A Xenon and Neon arc spectrum taken with the medium resolution grism at the K atmospheric window. The main lines are marked.

Appendix B

Atmospheric Absorption

The atmospheric transmission in the 0.8 to 2.5 micron region is plotted as a function of wavelength (Fig. B.1), together with pass-bands for the SOFI filters.

Figure B.1: The atmospheric transmission at a resolution of 8\AA .

Appendix C

SOFI Templates: A Reference Guide

C.1 General Points

This section provides a summary of the currently released SOFI templates. It assumes that the reader is already familiar with the capabilities of the instrument. Information about SOFI may be found on-line at the NTT Web site (<http://www.ls.eso.org/>). Note that undocumented and unreleased, prototype templates may be available at the telescope. Contact the SOFI support astronomer for further details.

SOFI has three types of templates: observation templates (OT) for science observations, Calibration templates (CT) for calibration exposures, and Acquisition templates (AT) for target acquisition. Tables C.1, C.2, C.3, and C.4 associates templates with observational scenarios.

Not all parameters have default values. In those cases where there is no default value, the value is “NODEFAULT”.

IMPORTANT! For those parameters that require a string to be entered, for example the target name or the exposure name, spaces and special characters must NOT be used, only letters, numbers and “_” (underscore symbol) are allowed.

The SOFI templates are maintained by the SOFI instrument scientist and by the Science Operations team.

The following sections will review the SOFI templates in order of the instrument modes: imaging, spectroscopy, polarimetry.

This Appendix supersedes the SOFI Template Manual (ESO document LSO-MAN-ESO-40100-0007).

Type of Acquisition	Template(s) to use
Simple telescope preset (point and shoot)	SOFI_img_acq_Preset
Preset telescope and move an object onto a pixel	SOFI_img_acq_MoveToPixel
Preset telescope and center an object in a slit	SOFI_img_acq_MoveToSlit
Preset telescope and position an object for polarimetry	SOFI_img_acq_Polarimetry

Table C.1: Short guide for acquisition templates

Type of Imaging	Template(s) to use
Imaging of uncrowded fields or point-like objects	SOFI_img_obs_AutoJitter or SOFI_img_obs_Jitter
Imaging of faint objects around bright objects	SOFI_img_obs_AutoJitterRot or SOFI_img_obs_JitterRot
Imaging of crowded fields or extended objects	SOFI_img_obs_AutoJitterOffset or SOFI_img_obs_JitterOffset
Map of extended fields	SOFI_img_obs_AutoJitterArray or SOFI_img_obs_AutoJitterArray_1
Imaging requiring complex telescope offsets and/or guiding options	SOFI_img_obs_GenericImaging
Fast Photometry	SOFI_img_obs_FastPhotJitt
Imaging Polarimetry	SOFI_img_obs_Polarimetry

Table C.2: Short guide for imaging and polarimetry templates

Type of Spectroscopy	Template(s) to use
Spectroscopy of point-like or moderately extended objects	SOFI_spec_obs_AutoNodOnSlit
As above but in Non Destructive read-out mode	SOFI_spec_obs_AutoNodNonDestr
Spectroscopy of extended objects (i.e. wider than ~ 2 arc-minutes), or complex sequences of slit positions	SOFI_spec_obs_GenericSpectro
As above but in Non Destructive read-out mode	SOFI_spec_obs_GenSpecNonDestr

Table C.3: Short guide for spectroscopic templates

Type of calibration	Template(s) to use
Darks	SOFI_img_cal_Darks
Imaging Dome Flat Fields	SOFI_ima_cal_DomeFlats
Special Imaging Dome Flats	SOFI_ima_cal_SpecialDomeFlats
Standard Star (imaging)	SOFI_img_cal_StandardStar
Polarimetric Dome Flat Fields	SOFI_img_cal_PolarimDomeFlats
Arcs (spectroscopy)	SOFI_spec_cal_Arcs
Spectroscopic Dome Flats (for different readout modes)	SOFI_spec_cal_DomeFlats or SOFI_spec_cal_DomeFlatsNonDestr
Spectroscopic Adapter Flats	SOFI_spec_cal_Flats

Table C.4: Short guide for calibration templates

C.2 SOFI Imaging Templates

C.2.1 SOFI Imaging Acquisition Templates

Acquisition Templates determine how a target field is acquired by the telescope. In case of SOFI some instrument parameters can be set up in the Acquisition Templates, such as the instrument mode. There are two acquisition template for simple imaging:

SOFI_img_acq_Preset

This acquisition template (Table C.5) does a simple telescope preset (i.e. pointing) in case of imaging observations, to the coordinates of the Target associated with the Observing Block. The Differential tracking rates can be set using TEL.TARG.ADDVELALPHA and TEL.TARG.ADDVELDELTA. There is no operator intervention, no acquisition image is saved on disk.

To save time, the instrument and detector parameters – most importantly the instrument mode and the filters – can be set to values that will be used in the first Observation Template attached with this Observing Block. An exception is the NDIT which should be set to 1, to save time, because even after the telescope points at the target, the next template will not start until the NDITxDIT integrations are completed.

It is the best template to use when pointing is not critical. The NTT typically points with an error of a few arc second. Most observers find this template sufficient for the purposes of their programs.

SOFI_img_acq_MoveToPixel

This acquisition template (Table C.6) has two parts. First, it presets the telescope to the coordinates of the Target associated with the Observation Block. Next, it requests the operator to interact with the RTD so that the user-specified object is moved to a user-specified location on the array.

In order for objects to be clearly seen, the telescope can do a small offset after the preset. The size and direction of the offset are template parameters (TEL.TARG.OFFSETALPHA and TEL.TARG.OFFSETDELTA). The operator is then prompted to store a fixed pattern that is subsequently subtracted from incoming images via the real-time display (RTD) features. After storing the fixed pattern, the telescope returns to the original position. At this point in time the operator is offered the possibility to change DIT and NDIT. This is

Parameter signature	Header Keyword	Value	Description
DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Add Velocity Alpha	TEL.TARG.ADDVELALPHA	0.0	Additional tracking velocity in RA (arcsec/sec)
Add Delta Velocity	TEL.TARG.ADDVELDELTA	0.0	Additional tracking velocity in DEC (arcsec/sec)
Rotation Offset on Sky	TEL.ROT.OFFANGLE	0.0	Rotation offset on the sky (degrees)

Table C.5: SOFI_img_acq_Preset.

useful if the target can not be identify securely because of the low signal-to-noise. If the operator changes the values, the telescope offsets again and the operator is required to store another fixed pattern before the telescope returns to the nominal position. This loop can continue until the operator has identified the target.

To move the target from one position of the array to another, the operator simply draws an arrow on the screen with the left hand-side button of the mouse. At this point a window, which lists the pixel co-ordinates at the start and the end of the arrow, will appear. The operator can accept the offsets, cancel, or edit the co-ordinates directly. If the offsets are accepted, the telescope offsets by the desired amount. Finally, the operator is given the possibility to redraw the arrow for refining the position of the target if necessary. Once the operator is satisfied, the template finishes.

Guiding is possible (defaulted to the Telescope Control System's option *Box-To-Star*; parameter value B) only if SEQ.COMBINED.OFFSET is set to T and if guiding is on before the start of the template. The experience has shown that the guiding is necessary if the observer plans to stay at one pointing more than 5-10 min. This is rarely the case in imaging mode but can often happen in spectroscopic observations.

The SEQ.PRESET parameter allows the user to use the fine-tunning pointing functionality of the template without presetting the telescope. On other words, the user can improve the position of the target on the array, without full telescope preset. This option is particularly useful if acquisition was aborted for some reason or if multiple OBs to observe the same target are executed one after another.

The parameter SEQ.SAVE allows to save an acquisition image at the prize of some additional overhead for file transfer and saving.

The interactive pop-up windows are usually displayed before new images have arrived on the RTD. Therefore, operators are strongly advised to carefully check that a new image has arrived before clicking on these windows (e.g. for storing a fixed pattern, for changing the DIT/NDIT configuration). The arriving of a new image on the RTD is marked by a flashing green dot in the middle of the upper part of the RTD window.

Parameter signature	Header Keyword	Value	Description
DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Alpha Offset (arcsec)	TEL.TARG.OFFSETALPHA	0.0	Alpha Offset for sky subtraction (arcsec)
Delta Offset (arcsec)	TEL.TARG.OFFSETDELTA	0.0	Delta Offset for sky subtraction (arcsec)
Add Velocity Alpha	TEL.TARG.ADDVELALPHA	0.0	Additional tracking velocity in RA (arcsec/sec)
Add Delta Velocity	TEL.TARG.ADDVELDELTA	0.0	Additional tracking velocity in DEC (arcsec/sec)
Rotation Offset on Sky	TEL.ROT.OFFANGLE	0.0	Rotation offset on the sky (degrees)
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	False	T - guiding ON, F - guiding OFF
Preset Telescope ? (F/T)	SEQ.PRESET	True	T - full preset, F - fine-tuning of the pointing
Save Image ? (F/T)	SEQ.SAVE	False	T - preserve the acq. image, F - no

Table C.6: SOFI_img_acq_MoveToPixel.

C.2.2 SOFI Imaging Science Templates

SOFI_img_obs_AutoJitter

This observation template (Table C.7) offsets the telescope between exposures according to a random pattern of offsets automatically generated within the template. It is ideal for long integrations on “empty” fields or fields containing isolated point sources, and does not require a long list of offsets to be defined by the observer.

The offsets are distributed randomly within a box whose size is defined by SEQ.JITTER.WIDTH (in arc seconds), with the condition that the distance between any two points in a series of ten values is greater than certain minimum specified within the template. This is intentionally done to ensure that the 5 frames before and after any frame are not too close spatially and can be safely used for creating a sky-image without large residuals for the sky subtraction.

By default, there is no telescope offset before the first exposure. If SEQ.RETURN is set to True (T), the telescope slews back to its original position at the end of the template; if not, the telescope is not moved. This feature is useful if more exposure of the same field are taken (i.e. with different filters) without running the corresponding acquisition template. In this case the parameters should be set to True.

If SEQ.COMBINED.OFFSET is set to T and if guiding was started before the start of the template, guiding is *Box-To-Star*. Considering that in imaging mode the time between the offsets is limited typically to 1-3 min by the sky background variation, the guiding is rarely used in imaging mode.

The value of SEQ.JITTER.WIDTH corresponds to the full width of the box in which the offsets are generated. The choice of box size is a compromise between constraints:

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Jitter Box Width (arcsec)	SEQ.JITTER.WIDTH	40	Jitter box size
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True

Table C.7: SOFI_img_obs_AutoJitter.

(1) Too wide a box may lead to insufficient image overlap. The signal-to-noise ratio will degrade toward the edged of the final combined image.

(2) Too small a value may lead to poor sky subtraction near extended objects, including bright stars with noticeable wings. Setting this parameter to zero is equivalent to staring the observations.

Consider the following example (Table C.8): the template will produce 6 fits files. Each image is jittered regarding to the previous one, and the jitter offsets are chosen randomly inside a box of 20 arcsec (Jitter Box Width) around the central position. Each of them corresponds to the average of 6 exposures (NDIT) of 10 sec (DIT). At the end of the exposures, the telescope moves back to the preset position (Return to origin = True).

SOFI_img_obs_AutoJitterRot

This observation template (Table C.9, C.10) is identical with the SOFI_img_obs_AutoJitter with the additional capability to move the rotator by a given offset. The offset is relative with respect to the current position. The user is strongly advised to check if the rotator has sufficient range to complete the observations and to monitor the rotator angle variation during the execution of the template. The offset parameter is a list, the default value is list of one zero, the valid range is from -360 to 360 degr.

SOFI_img_obs_AutoJtterOffset

This observation template (Table C.11) allows the user to move the telescope alternatively between the object and a nearby patches of sky. However, when pointing to the object, the position of the telescope is still randomly distributed within a box of size `SEQ.JITTER.WIDTH` (in arcsec).

The offset skies are at a constant distance (defined by `SEQ.SKYTHROW`, in arcsec) from the original telescope position, but at an angle randomly distributed between 0 and 360 degrees (i.e. the offset skies are distributed on a circle surrounding the initial telescope position).

Parameter signature	Value
Exposure Name	NGC6118
DIT	10
NDIT	6
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures?	6
Filter wheel 1	K_S
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined offset ? (T/F)	F
Jitter Box Width (arcsec)	20
Return to origin ? (T/F)	T

Table C.8: SOFI_img_obs_AutoJitter. Example.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Jitter Box Width (arcsec)	SEQ.JITTER.WIDTH	40	Jitter box size
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
Rotation offsets list (degr) on Sky	SEQ.OFFANGLE	0	List of rotator offsets on the sky

Table C.9: SOFLimg_obs_AutoJitterRot.

Parameter signature	Value
Exposure Name	NGC6118
DIT	10
NDIT	6
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures?	6
Filter wheel 1	K_S
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined offset ? (T/F)	F
Jitter Box Width (arcsec)	20
Return to origin ? (T/F)	T
Rotation offsets list (degr) on Sky	30

Table C.10: SOFI_img_obs_AutoJitterRot. Example.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Jitter Box Width (arcsec)	SEQ.JITTER.WIDTH	40	Jitter box size
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
Sky Offset Throw (arcsec)	SEQ.SKYTHROW	300	Radius of the sky offsets region (arcsec)
Pupil Rotation	SEQ.ROTPUPIL	T	T - rotate, N - not rotate

Table C.11: SOFI_img_obs_AutoJtterOffset.

Parameter signature	Value
Exposure Name	NGC6118
DIT	6
NDIT	10
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures?	10
Filter wheel 1	K_S
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined offset ? (T/F)	F
Jitter Box Width (arcsec)	20
Return to origin ? (T/F)	T
Sky Offset Throw (arcsec)	600.
Rotate Pupil ?	T

Table C.12: SOFI_img_obs_AutoJitterOffset. Example.

By default, there is no telescope offset before the first exposure.

If `SEQ.RETURN` is set to True (T), the telescope slews back to its original position at the end of the template. If not, the telescope is not moved.

It is assumed that odd numbered exposures are on the object, and are consequently identified as SCIENCE frames. Even numbered exposures are assumed to be sky frames, and are consequently identified as OTHER.

If the number of exposures is even, $(\text{SEQ.NEXPO})/2$ pairs of object/sky frames are produced. If the number of exposures is odd, then an extra frame is taken, and $(\text{SEQ.NEXPO} + 1)/2$ pairs of object/sky frames are taken. Users are encouraged to give the parameter `SEQ.NEXPO` an even value.

Guiding is only possible for the object frames, if `SEQ.COMBINED.OFFSET` is set to True (T; must be done by the observer during the preparation of the Observing Block) and guiding is started before the template starts (done by the operator after warning from the observer. By default, there is no guiding for the sky frames.

This template can also be used to observe small extended objects by setting the throw small enough so that the object is always within the field of view. In such a case, however, `SEQ.ROTPUPIL` should be set to False (F).

Using this template with a jitter box size parameter `SEQ.JITTER.WIDTH` set to zero is equivalent to staring on the object.

For the pupil rotation see the discussion in Sec. 3.6.

Consider the following example (Table C.12): the sky is taken at 600 arcsec of the object, and a small jitter of 20 arcsec around the target is done. The number of exposures corresponds to the TOTAL number of exposures, that means SCIENCE exposures + SKY exposures. So here, you will have 5 images of your objects, and 5 images of the sky. Each image is the average of 10 (NDIT) exposures of 6 sec (DIT). No guiding is used.

SOFI_img_obs_Jitter

This observation template (Table C.13) allows the user to offset the telescope between exposures according to a list of predefined offsets (`SEQ.OFFSETALPHA.LIST` and `SEQ.OFFSETDELTA.LIST`), giving more freedom and flexibility. Complicated patterns can be used, pointing the telescope exactly on

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
RA offsets list (arcsec)	SEQ.OFFSETALPHA.LIST	NODEFAULT	list of offsets along RA
Dec offsets list (arcsec)	SEQ.OFFSETDELTA.LIST	NODEFAULT	list of offsets along Dec

Table C.13: SOFI_img_obs_Jitter.

areas of interest. However, this comes at a price of more difficult preparation of the observation. The offsets in the list are **relative** to the previous position, they are in RA and DEC, and they are in arcsec. Note that the telescope offset are executed **before** the corresponding image is taken. In other words, if you want the first image to be at the position you have acquired the target, the first offset has to be (0,0).

If the number of the images is longer than the number of the offsets, the software loops back to the beginning of the list and executes again the first offset. If the number of the images is shorter than the number of the offsets, the software executes only as many offsets as necessary, starting from the beginning of the list.

The RA and DEC lists can have different length but this can be extremely confusing. It is good practice to use either lists with a single value or lists of equal length.

If `SEQ.COMBINED.OFFSET` is set to True (must be done by the observer during the preparation of the Observing Block) and if guiding was started before the start of the template (will be done by the operator after a warning from the observer), the guiding is *Box-To-Star*. Given the short time spent on a single position – typically 2-3 min, determined by the sky level variations – guiding is usually rarely used.

If `SEQ.RETURN` is set to True (T), the telescope slews back to its original position at the end of the template. If not, the telescope is not moved.

Consider the following example (Table C.14): it gives to you 5 images at 5 different positions. Each image is the average of 3 (NDIT) exposures of 20 sec (DIT). The first image is taken at the position at which the target was acquired.

The following example (Table C.15) will produce 10 images, repeating the sequence of 5 offsets. Note that the sums of the offsets along RA and along Dec are non-zero, so the second offset pattern will be shifted relative to the first one by 2.5 arc to the East and 2.5 arcsec to the North.

Parameter signature	Value
Exposure Name	NGC6118
DIT	20
NDIT	3
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures	5
Filter wheel 1	J
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined offset ? (F/T)	T
Return to origin ? (T/F)	T
RA offset list (arcsec)	0 -50 0 100 0
DEC offset list (arcsec)	0 -50 100 0 -100

Table C.14: SOFI_img_obs_Jitter. Example.

Parameter signature	Value
Number of Exposures?	10
RA offset list (arcsec)	0. 5. -10. 2.5 5.
DEC offset list (arcsec)	0. 5. -10. 2.5 5.

Table C.15: SOFI_img_obs_Jitter. Offset Example.

This example illustrates an interesting usage of this template – to define a “closed loop” offset sequence that will bring the telescope at the original position (or close to the original position, i.e. 2-5 arcsec away, to imitate jittering, minimizing the effects of the array cosmetic defects). Therefore, these template can execute sequences of the type: target1-sky-target2-target1-sky-target2-..., saving from the time spent on the sky. Unlike some other templates such as the SOFI_img_obs_AutoJitterArray, the freedom to choose the offsets means that the targets and the sky can form an irregular pattern.

SOFI_img_obs_JitterRot

This observation template (Table C.16, C.17) is identical with the SOFI_img_obs_Jitter with the additional capability to move the rotator by a given offset. The offset is relative with respect to the current position. The user is strongly advised to check if the rotator has sufficient range to complete the observations and to monitor the rotator angle variation during the execution of the template. The offset parameter is a list, the default value is list of one zero, the valid range is from -360 to 360 degr.

SOFI_img_obs_JitterOffset

This observation template (Table C.18) is very similar to SOFI_img_obs_AutoJitterOffset. The only difference is that the user fully specifies the offsets in `SEQ.OFFSETALPHA.LIST` and `SEQ.OFFSETDELTA.LIST`, instead of executing randomly distributed offsets.

It is assumed that offsets are given so that the telescope alternates between object (odd numbered frames) and sky. The object frames are identified as `SCIENCE`, and the sky frames are identified as `OTHER`. The telescope offsets are executed before the taking the images.

Guiding is only possible for the object frames, if `SEQ.COMBINED.OFFSET` is set to True (T; must be done by the observed during the preparation of the Observing Block) and guiding is started before the template starts (done by the operator after warning from the observer).

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
RA offsets list (arcsec)	SEQ.OFFSETALPHA.LIST	NODEFAULT	list of offsets along RA
Dec offsets list (arcsec)	SEQ.OFFSETDELTA.LIST	NODEFAULT	list of offsets along Dec
Rotation offsets list (degr) on Sky	SEQ.OFFANGLE	0	List of rotator offsets on the sky

Table C.16: SOFI_img_obs_JitterRot.

Parameter signature	Value
Exposure Name	NGC6118
DIT	20
NDIT	3
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures	5
Filter wheel 1	J
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined offset ? (F/T)	T
Return to origin ? (T/F)	T
RA offset list (arcsec)	0 -50 0 100 0
DEC offset list (arcsec)	0 -50 100 0 -100
Rotation offsets list (degr) on Sky	30

Table C.17: SOFI_img_obs_JitterRot. Example.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
RA offsets list (arcsec)	SEQ.OFFSETALPHA.LIST	NODEFAULT	list of offsets along RA
Dec offsets list (arcsec)	SEQ.OFFSETDELTA.LIST	NODEFAULT	list of offsets along Dec
Pupil Rotation	SEQ.ROTPUPIL	T	T - rotate, N - not rotate

Table C.18: SOFI_img_obs_JitterOffset

Compare the following example (Table C.19) with the example for SOFI_img_obs_AutoJitterOffset template: the images are now taken in a user-defined sequence instead of in a random jittering pattern. The first offset (0,0) is executed before taking the first image, the second offset (5,5) is executed before taking the second image and so on. The number of exposures corresponds to the TOTAL number of exposures. Each image is the average of 10 (NDIT) exposures of 6 sec (DIT). No guiding is used.

SOFI_img_obs_AutoJitterArray and SOFI_img_obs_AutoJitterArray_1

These templates (Table C.20) are very similar to SOFI_img_obs_Jitter. They define a sequence of offsets in RA and DEC but in addition to that they allow to jitter randomly around each of the defined offset (array) positions. They are useful to build maps of large areas of the sky.

SEQ.OFFSETALPHA.LIST and SEQ.OFFSETDELTA.LIST are used to define the map (array) on the sky. The offsets between different positions are relative and they do not have to be equidistant (i.e. irregularly shaped maps/arrays are acceptable). If you have less offsets than number of exposures – defined by the parameter SEQ.NEXPO, – BOB will start again with the first offset. It repeats simply the list of offsets you enter until it has the same number of exposures. It is possible to randomly jitter around each position of the map within a box whose size is set by SEQ.JITTER.WIDTH. The number of jittered frames at each offset (array) position is SEQ.NJITT and the telescope moves to the next map (array) position only after they are completed. The total number of frames acquired by this template will be the product of the number of array positions by the number of jitters around each one of them, i.e. SEQ.NEXPO x SEQ.NJITT.

A modified version of this template is called SOFI_img_obs_AutoJitterArray_1. It has identical parameter set, and it takes the same number of images as the original template but in a different order: exposures are taken at the positions of the array as defined by the offsets lists, then the random offset is applied and the entire array pattern is executed again. This template is preferable to the original one because it gives better sky subtraction. Also, if the weather conditions are poor (i.e. the

Parameter signature	Value
Exposure Name	NGC6118
DIT	6
NDIT	10
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures?	5
Filter wheel 1	Ks
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined offset ? (T/F)	F
Return to origin ? (T/F)	T
RA offset list (arcsec)	0. 5. -10. 2.5 5.
DEC offset list (arcsec)	0. 5. -10. 2.5 5.
Rotate Pupil ?	T

Table C.19: SOFI_img_obs_JitterOffset. Example.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
NJITT	SEQ.NJITT	1	Number of Jittered exposures around each Array position
Jitter Box Width (arcsec)	SEQ.JITTER.WIDTH	40	Jitter box size
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
RA offsets list (arcsec)	SEQ.OFFSETALPHA.LIST	NODEFAULT	list of offsets along RA
Dec offsets list (arcsec)	SEQ.OFFSETDELTA.LIST	NODEFAULT	list of offsets along Dec

Table C.20: SOFI_img_obs_AutoJitterArray and SOFI_img_obs_AutoJitterArray_1

Parameter signature	Value
Exposure Name	NGC6118
DIT	6
NDIT	10
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures?	4
NJITT	3
Jitter box width (arcsec)	20
Filter wheel 1	Ks
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined offset? (T/F)	F
Return to origin? (T/F)	T
RA offset list (arcsec)	50 0 50 0
Dec offset list (arcsec)	50 50 0 50

Table C.21: SOFI_img_obs_AutoJitterArray and SOFI_img_obs_AutoJitterArray_1. Example.

sky background variations are strong, there are thin clouds, etc.), the entire mapped area is imaged under relatively more similar conditions in comparison with the SOFI_img_obs_AutoJitterArray template. On the other hand, if the conditions grow even worse and make it impossible to complete the mapping, the observer may be left with a shallow image over the entire area as opposed to a deep image of part of the mapped area. The observer must decide which strategy is better suited to the goals of the program before deciding which template to use.

Consider the following example (Table C.21): you will have 4 main positions defined in the RA and DEC offset lists. Around each of the offset position, you will have 3 images with jitter (inside a box of 20 arcsec defined by the parameter `Jitter Box Width`). In total this means 3×4 images = 12 images. Each of them is the average of 10 exposures (NDIT) of 6 sec (DIT), so the total exposure time is $10 \times 6 \times 3 \times 4 = 12$ min (without the overheads).

The template SOFI_img_obs_AutoJitterArray will take a sequence of 3 jittered images around the first offset position, then move to the next offset position and take 3 more jittered images and so on until the entire mapped region is covered. The template SOFI_img_obs_AutoJitterArray_1 will instead take one image at each 4 offset positions, then it will repeat the entire offset pattern NJITT=3 times, adding before each cycle a small jitter offset.

SOFI_img_obs_GenericImaging

This template (Table C.22) has some similarity with SOFI_img_obs_Jitter, but it allows the user to do any sequence of telescope offsets whether they be combined offsets (offsetting the telescope and the guide probe; necessary in case of observations with guiding) or non-combined offsets (offsetting the telescope alone; when guiding is not used, i.e. the telescope pointing changes on average at least once every 5-10 minutes). However, with this flexibility comes complexity. This is a complicated template and it is meant to be used in situations in which none of the other templates is suitable.

Telescope offsets and guiding options are defined as lists in `SEQ.OFFSETALPHA.LIST`, `SEQ.OFFSETDELTA.LIST` and `SEQ.GUIDING.LIST`. The offsets are relative to the previous position, they are in RA and DEC and they are in arcsec.

There are three guiding options: N, B, and S. The option N must be used only for observations without guiding, i.e. the offsets are non-combined. This works well for short total integrations (up to 5-10 min) at each individual position. The template will fail if this option is selected, AND the

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
List if NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
RA offsets list (arcsec)	SEQ.OFFSETALPHA.LIST	NODEFAULT	list of offsets along RA
Dec offsets list (arcsec)	SEQ.OFFSETDELTA.LIST	NODEFAULT	list of offsets along Dec
Obs Type (O or S)	SEQ.OBSTYPE.LIST	NODEFAULT	O - object, S - sky
Guiding (N B S)	SEQ.GUIDING.LIST	NODEFAULT	N - no guiding, B - Box-To-Star, S - Star-To-Box DON'T USE

Table C.22: SOFI_img_obs_GenericImaging

guiding is started. The option B stands for *Box-To-Star*, and it must be used with guiding, i.e. the offsets are combined with with appropriate counter-movement of the guiding arm to maintain the pointing. This is necessary for long total integrations (more than 5-10 min) at each individual position. The option S stands for *Star-To-Box*. It has been disabled (i.e., the template fails if the guiding option is set to S) but it was kept in the software for historic reasons.

Note: With large combined offsets, the guide probe may not be able to follow the one guide star. In such a case, the guiding system will automatically find another star, and resume guiding.

If any of the entries in `SEQ.GUIDING.LIST` are different from N, and if guiding is off at the beginning of the template, a pop up window prompts the operator to start guiding.

A list is also provided for NDIT (`SEQ.NDIT.LIST`), but not for DIT.

The observation type can be defined for each image, and is entered as a list in `SEQ.OBSTYPE.LIST`. O stands for Object and S stands for Sky. This is to make it easier for the user to classify the images during the data reduction.

The total number of exposures is defined in `SEQ.NEXPO`. This number can differ from the number of elements in the aforementioned lists. Lists do not need to have the same length. If the number of exposures is larger than the number of elements in a list, the list is restarted from the beginning as many times as needed until the correct number of frames have been acquired. The lists can have any length; however having lists of different lengths can become extremely confusing. It is good practice to use either lists of one value, and/or lists of equal length.

At the end of the template, the telescope is returned to the original position, if `SEQ.RETURN` is set to True (T). If not, the telescope is not moved at the end of the template.

Consider the following example (Table C.23). It contains 7 offsets with the option to guide *Box-To-Star* (the only other option is *No Guiding*; *Star-To-Box* has been disabled). The object is always

Parameter signature	Value
Exposure Name	NGC6118
DIT	20
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures?	7
List of NDIT	3 2 3 2 3 2 3
Filter wheel 1	Ks
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Return to origin? (T/F)	T
RA offset list (arcsec)	0 -150 140 150 -140 -250 250
DEC offset list (arcsec)	0 -150 140 150 -140 -250 250
Obs Type (O or S)	O S O S O S O
Guiding (N B S)	B

Table C.23: SOFI_img_obs_GenericImaging. Example.

observed with a NDIT equal to 3, whereas the sky is always observed with NDIT equal to 2. At the end there are 4 images for the object, each of them is the average of 3 (NDIT) expositions of 20 sec (DIT). And 3 images of the sky, each of them is the average of 2 (NDIT) expositions of 20 sec (DIT).

SOFI_img_obs_FastPhotJitt

This observation template (Table C.24) takes a fast photometry sequence. It has built capability to apply a random jitter within a predefined box on the sky and to repeat the fast photometry sequence. It offsets the telescope between bursts according to a random pattern of offsets automatically generated within the template. However, the implication when using this template is that the detector will be heavily windowed to minimize the minimal DIT . Thus, usually the jitter box size will be extremely small and placing the target on a detector location clear of bad pixels may be a better strategy. Therefore, it is recommended to set the jitter box size to zero.

The offsets are distributed randomly within a box whose size defined by `SEQ.JITTER.WIDTH` (in arc seconds), with the condition that the distance between any two points in a series of ten values is greater than certain minimum specified withing the template. This is intentionally done to ensure that the 5 frames before and after any frame are not too close spatially and can be safely used for creating a sky-image without large residuals for the sky subtraction.

By default, there is no telescope offset before the first exposure. If `SEQ.RETURN` is set to True (T), the telescope slews back to its original position at the end of the template; if not, the telescope is not moved. This feature is useful if more exposure of the same field are taken (i.e. with different filters) without running the corresponding acquisition template. In this case the parameters should be set to True.

To minimize the overheads, this template is used without telescope guiding and it does not verify if the telescope is guiding or not.

If the `DET.BURST.MODE` parameter is set to True, then the individual reads are stored as opposed to False when the reconstructed images (i.e. the difference between the final and the initial detector reads) are stored.

The `DET.STORE.CUBE` determines how to store the data: if it is set to True, the observing software generates a single data cube in which every layer is a 1-dimensional image. The cube contain NDIT layers

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	0.05	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	3	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	100	Number of columns in the window
Number of rows	DET.WIN.NY	100	Number of rows in the window
First column of window	DET.WIN.STARTX	463	First column of window
First row of window	DET.WIN.STARTY	463	First row of window
Turn on Burst mode ? (T: on, F: off)	DET.BURST.MODE	T	T - store detector reads F - store reconstructed images
Store DITs in a cube (T) on in individual files (F)	DET.STORE.CUBE	T	T - store data in 3-D fits cube F - store data in individual fits files
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Jitter Box Width (arcsec)	SEQ.JITTER.WIDTH	5	Jitter box size
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Date of event YYMMDD, set to 0 for no event	SEQ.EVENT.DATE	0	Event date for timed observations
Time of event HHMMSS, set to 0 for no event	SEQ.EVENT.TIME	0	Event time for timed observations

Table C.24: SOFI_{img_obs_FastPhotJitt}.

if `DET.BURST.MODE` is set to False and $2 \times \text{NDIT}$ if `DET.BURST.MODE` is set to True. If `DET.STORE.CUBE` is set to False, then `NDIT` or $2 \times \text{NDIT}$ individual fits files are generated, respectively. The fastest photometry and the most convenient data handling is achieved with both `DET.BURST.MODE` and `DET.STORE.CUBE` set to True.

The `DET.WIN.STARTX` and `DET.WIN.STARTY` are ignored because by default the usable window is centered on the detector to speed up the readout – this allows to use simultaneously the four analog-to-digital converters.

Finally, for the purpose of occultation observations, the template can perform timed observations that will start at a moment defined by the parameters `SEQ.EVENT.DATE` and `SEQ.EVENT.TIME`, in UT. To forfeit this option, set both of them to zero and the observations will start at the moment the template is executed by the operator.

Nota bene This is a preliminary version of the template and it is likely to evolve with time. Currently, it is only available the telescope. Consult your support astronomer or the Sofi Instrument Scientist for details.

Consider the following example (Table C.25): the template will produce 6 data cube files. Each image is located on the same place (i.e. no jittering) and it consists of 2000 reads (because `NDIT` = 1000 and `DET.BURST.MODE` is set to True), each with `DIT` = 0.1 sec. The detector is windowed down to $100 \times 100 \text{ px}^2$. The sequence will begin immediately after the start of the template and the set up of the instrument (as opposed as at predefined time). At the end of the exposures, the telescope moves back to the preset position (`Return to origin` = True).

Parameter signature	Value
Exposure Name	Star
DIT	0.1
NDIT	1000
Number of columns	100
Number of rows	100
First column of window	463
First row of window	463
Turn on Burst mode ? (T: on, F: off)	T
Store DITs in a cube (T) on in individual files (F)	T
Filter wheel 1	K_S
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Jitter Box Width (arcsec)	0
Return to origin ? (T/F)	T
Number of Exposures?	1
Date of event YYMMDD, set to 0 for no event	0
Time of event HHMMSS, set to 0 for no event	0

Table C.25: SOFI_img_obs_FastPhotJitt. Example.

C.2.3 SOFI Imaging Calibration Templates

SOFI_img_cal_Darks

This calibration template (Table C.26) produces dark frames, both for imaging and spectroscopy. The instrument is set to the DARK mode, where both filter wheels and the grism wheel are in the closed position.

The number of frames is defined in `SEQ.NEXPO`. Dark frames can be taken with different DIT and NDIT values, defined as lists in `SEQ.DIT.LIST` and `SEQ.NDIT.LIST`.

Note that `SEQ.DIT.LIST` and `SEQ.NDIT.LIST` do not need to have the same length. If the number of exposures is larger than the number of elements in a list, then the list is restarted from the beginning, as many times as needed until the correct number of frames have been acquired.

This template should be used to take dark frames at the end of the night for all the detector integration times that have been used during the night, though it is not strictly necessary since the dark signal will be removed along with the subtraction of the sky background.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
List of DITs	SEQ.DIT.LIST	NODEFAULT	Detector Integration Time; individual exposure (sec)
List if NDIT	SEQ.NDIT.LIST	NODEFAULT	Number of DITs averaged into an individual image

Table C.26: SOFI_img_cal_Darks

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
List of DITs	SEQ.DIT.LIST	NODEFAULT	Detector Integration Time; individual exposure (sec)
List if NDIT	SEQ.NDIT.LIST	NODEFAULT	Number of DITs averaged into an individual image
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode

Table C.27: SOFI_img_cal_DomeFlats and SOFI_img_cal_SpecialDomeFlats.

SOFI_img_cal_DomeFlats and SOFI_img_cal_SpecialDomeFlats

Imaging dome flats are taken with these calibration template (Table C.27). For each element in `SEQ.DIT.LIST` and `SEQ.NDIT.LIST`, four images are taken: one with the dome lamp off, two with the dome lamp on and a fourth with the dome lamp off. The intensity of the dome lamp is controlled manually.

`SOFI_img_cal_SpecialDomeFlats` takes an additional set of images with the aperture wheel partially masking the array. These images are used to estimate the bias pattern. The reasons for using this template rather than the first are explained in the *SOFI Users manual*. The intensity of the dome lamp is controlled semi-automatically: the template maintains a database and once it is populated for a given configuration, the lamp intensities are taken from there.

The parameter sets for both templates are identical.

SOFI_img_cal_StandardSatr

This calibration template (Table C.28) is intended for standard star observations in one filter. If a standard has to be observed in a few different filters, this template has to be called as many times as necessary, with the appropriate filter defined in each template call.

The telescope can be offset to position the object at different positions on the array. **Note** that unlike the other templates where the offsets are defined along the RA and Dec, here they are defined along detector rows and columns, and are entered into `SEQ.OFFSETX.LIST` and `SEQ.OFFSETY.LIST`. Offsets are in arcsec.

Guiding is possible (defaulted to *Box-To-Star*) if `SEQ.COMB.OFFSET` is set to True (T) and if guiding was started before the start of the template.

If `SEQ.RETURN` is set to True (T), the telescope slews back to its original position at the end of the template. If not, the telescope is not moved.

Although it is possible to observe standard stars with any of the imaging science templates, using this one is recommended because it sets up properly the header keywords necessary for the pipeline data reduction of the standards.

In the next example (Table C.29) `SOFI_img_cal_StandardStar` is used to obtain 5 images: the first on the acquisition position and the others defined by the list of offsets. At the end you have an “x”

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Instrument Mode	INS.IMODE	NODEFAULT	Instrument Mode
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
X offsets list (arcsec)	SEQ.OFFSETX.LIST	NODEFAULT	list of offsets along rows
Y offsets list (arcsec)	SEQ.OFFSETY.LIST	NODEFAULT	list of offsets along columns

Table C.28: SOFLimg_cal_StandardSatr.

pattern.

Parameter signature	Value
Exposure Name	S9104
DIT	2
NDIT	5
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures?	5
Filter Wheel 1	J
Filter wheel 2	open
Instrument mode	Large_field
Combined offset ? (F/T)	T
Return to origin ? (T/F)	T
X offset (arcsec)	0 45 -90 0 90
Y offset (arcsec)	0 45 0 -90 0

Table C.29: SOFI_ima_obs_StandardStar. Example.

C.3 SOFI Polarimetric Template:

C.3.1 SOFI Polarimetric Acquisition Template

SOFI_img_acq_Polarimetry

This acquisition template (Table C.30) is very similar to the SOFI_img_acq_MoveToPixel template. The polarimetric mask is displayed on the RTD and is superimposed on the image of the field. Then, the operator is prompted to define an offset by drawing an arrow on the RTD holding the left mouse button. The offset fine-positions the object into the transparent region of the mask.

Since the POLARIMETRY mode uses the large field objective, the instrument mode is not a parameter of the template.

C.3.2 SOFI Polarimetric Science Template

SOFI_img_obs_Polarimetry

This observation template (Table C.31) is used for polarimetry. It is quite simple, as it works as an imaging template.

The number of exposures must correspond to the number of the offsets, specified in `SEQ.OFFSETX.LIST` and `SEQ.OFFSETY.LIST`. **Note** that unlike most of the imaging templates, here the offsets are defined along detector rows and columns in arcsec, so that the users can move the object easily along the strips of the polarimetric mask.

Most importantly, the user have to decide which angle are necessary: one should observe the same object at several angles to determine the Stokes parameters (cf. section 2.4). All frames can be given an orientation relative to the previous position angle, by setting `SEQ.ROT.OFFANGLE` parameter (in degrees). When the template starts, the instrument is rotated on the sky by `SEQ.ROT.OFFANGLE`, and remains at this position until the end of the template. After the last exposure, the instrument is rotated back to the original position.

With this scheme, it is possible for the user to sample the object and the sky as desired for one rotator position, and then restart the template with another orientation on the sky for another series of exposures. At least two different orientations, separated by 45 degrees, are required for

Parameter signature	Header Keyword	Value	Description
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Alpha Offset (arcsec)	TEL.TARG.OFFSETALPHA	0.0	Alpha Offset for sky subtraction (arcsec)
Delta Offset (arcsec)	TEL.TARG.OFFSETDELTA	0.0	Delta Offset for sky subtraction (arcsec)
Add Velocity Alpha	TEL.TARG.ADDVELALPHA	0.0	Additional tracking velocity in RA (arcsec/sec)
Add Delta Velocity	TEL.TARG.ADDVELDELTA	0.0	Additional tracking velocity in DEC (arcsec/sec)
Rotation Offset on Sky	TEL.ROT.OFFANGLE	0.0	Rotation offset on the sky (degrees)
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	False	T - guiding ON, F - guiding OFF
Preset Telescope ? (F/T)	SEQ.PRESET	True	T - full preset, F - fine-tuning of the pointing
Save Image ? (F/T)	SEQ.SAVE	False	T - preserve the acq. image, F - no

Table C.30: SOFI_img_acq_Polarimetry.

computing the Stokes parameters. This implies that **the template must be called at least twice within an Observation Block**, or with two different Observation Blocks, both with two different rotator positions. The most likely situation will be to set `SEQ.ROT.OFFANGLE` to 0 degrees in the first template, and then to 45 degrees in the second template.

Guiding is possible (defaulted to *Box-To-Star*) if `SEQ.COMBINED.OFFSET` is set to T and if guiding was started before the start of the template. If `SEQ.RETURN` is set to True (T), the telescope slews back to its original position at the end of the template. If not, the telescope is not moved.

Note that the instrument mode is not a parameter of this template because all polarimetric observations are taken in Large Field mode.

This example (Table C.32) gives the typical values of the parameters. Note that to determine the Stokes parameters, the template will have to be executed twice with a rotator offset of 45 degree.

C.3.3 SOFI Polarimetric Calibration Template

SOFI_img_cal_PolarimDomeFlats

This template (Table C.33) takes polarimetric dome flats for a list of rotation angles.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures	SEQ.NEXPO	1	Number of exposures in the sequence
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	F	T - guiding ON, F - guiding OFF
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
X offsets list (arcsec)	SEQ.OFFSETX.LIST	NODEFAULT	list of offsets along rows
Y offsets list (arcsec)	SEQ.OFFSETY.LIST	NODEFAULT	list of offsets along columns
Rotator Offset ?	SEQ.ROT.OFFANGLE	45	Rotator offset (degrees)

Table C.31: SOFI_img_obs_Polarimetry.

Parameter signature	Value
Exposure Name	NGC6118
DIT	20
NDIT	3
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposures?	5
Filter wheel 1	Ks
Filter wheel 2	open
Instrument Mode	LARGE_FIELD_IMAGING
Combined offset ? (F/T)	T
Return to origin ? (T/F)	T
X offset list (arcsec)	0 -50 0 100 0
Y offset list (arcsec)	0 -50 100 0 -100
Rotator offset?	0

Table C.32: SOFI_img_obs_Polarimetry. Example.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DIT (individual exposure)	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
NDIT (number of DIT)	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures ?	SEQ.NEXPO	1	Number of exposures in the sequence
List of Rotator angles	SEQ.ROT.OFFANGLE	NODEFAULT	Rotator offset list (degrees)
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True

Table C.33: SOFI_img_cal_PolarimDomeFlats.

C.3.4 SOFI Spectroscopic Templates

C.3.5 SOFI Spectroscopic Acquisition Templates

SOFI_img_acq_MoveToSlit

This acquisition template (Table C.34) is very similar to the SOFI_img_acq_MoveToPixel template. The selected slit is drawn on the real-time display (RTD) superimposed on the image of the field. In most cases, operators will use the option to move the selected object to the center of the slit. Since each slit has a different position on the detector, the slit name (i.e. the 0.6, 1 or 2 arcsec slits) is explicitly stated in the template via the INS.WHICHSLIT parameter.

This acquisition has two parts. First, it presets the telescope to the coordinates of the Target associated with the Observation Block. Next, it takes and displays an image on the RTD. Finally, it requests the operator to interact with the RTD so that the user-specified object is moved to either a center of the slit, or to a user-specified location on the slit.

To make the target clearly visible, the telescope can do a small offset after the preset to obtain a “sky” image. The size and direction of the offset are template parameters (TEL.TARG.OFFSETALPHA and TEL.TARG.OFFSETDELTA). The operator is then prompted to store a fixed pattern that is subsequently subtracted from the incoming images via the real-time display (RTD) features. After storing the fixed pattern, the telescope returns to the original position. At this point in time the operator is offered the possibility to change DIT and NDIT. This is useful if the target can not be identify securely because of the low signal-to-noise. If the operator changes the values, the telescope offsets again and the operator is required to store another fixed pattern before the telescope returns to the nominal position. This loop can continue until the operator has identified the target.

To move the target from one position of the array to another, the operator has two options:

(1) If the target must be centered on the slit, one can use the build in feature by clicking on the “center” button on the RTD dialog box. The software will draw an arrow from the location of the selected object to the center of the slit and the operator will be asked to confirm the offset. This option is faster and it works well if the science observation is done via simple nodding along the slit.

(2) The operator can simply draw an arrow on the screen holding the left hand-side button of the mouse. At this point a window, which lists the pixel co-ordinates at the start and the end of the arrow, will appear. The operator can accept the offsets, cancel, or edit the co-ordinates manually. If the offsets are accepted, the telescope offsets by the desired amount. Finally, the operator is given the possibility to redraw the arrow for refining the position of the target if necessary. Once the operator is satisfied, the template finishes.

Nota Bene: There is no way of checking if the slit that is chosen in the Acquisition Template and the slit that will be used subsequently in the Observation Description Templates are the same. It is, therefore, of utmost importance that the astronomer ensures that they are the same.

Since the spectroscopic modes use the large field objective by default, the instrument mode is not a parameter of the template.

Guiding is possible (defaulted to the Telescope Control System's option *Box-To-Star*) **only** if `SEQ.COMBINED.OFFSET` is set to T (this must be done by the observer during the preparation of the Observing Block) and if guiding is on before the start of the science templates (this will be done by the operator after warning by the observer, before starting the execution of the Observing block). The experience has shown that the guiding is necessary if the observer plans to stay at one pointing more than 15 min. This is rarely the case in imaging mode but can often happen during spectroscopic observations.

The `SEQ.PRESET` parameter allows the observer to use the fine-tuning pointing functionality of the template without presetting the telescope. In other words, the user can change the position of the target on the array, without full telescope preset. This option is particularly useful if acquisition was aborted for some reason or if multiple OBs to observe the same target are executed one after another. In this case, it is advisable to re-acquire the target with a false preset every 1-2 hours to ensure it has not drifted away from the slit.

The parameter `SEQ.SAVE` allows to save an acquisition image at the cost of small additional overhead for file transfer and saving.

The interactive pop-up windows are usually displayed before new images have arrived on the RTD. Therefore, operators are strongly advised to carefully check that a new image has arrived before clicking on these windows (e.g. for storing a fixed pattern, for changing the DIT/NDIT). The arriving of a new image on the RTD is marked by a flashing green dot in the middle of the upper part of the RTD window.

C.3.6 SOFI Spectroscopic Science Templates

`SOFI_spec_obs_AutoNodOnSlit` and `SOFI_spec_obs_AutoNodNonDestr`

These observation templates (Table C.35 and C.36) nod the telescope between two positions (A and B) along the slit. A cycle is a pair of two observations: AB or BA. Cycles are repeated in ABBA sequences. E.g. 3 cycles corresponds to an ABBAAB sequence, 4 cycles correspond to an ABBAABBA sequence, etc. Each observation consist of one or more exposures, as defined in the template by the `NINT` parameter. The total number of frames corresponds to $\text{SEQ.NABCYCLES} \times \text{SEQ.NINT} \times 2$.

The mean amplitude of the nod is defined by `SEQ.NODTHROW` (in arcsec). The first exposure (A) is done after offsetting the object along the slit by $-(\text{SEQ.NODTHROW})/2$ arcsec, which can be either negative or positive. The second exposure (B) is therefore $+(\text{SEQ.NODTHROW})/2$ arcsec from the initial position along the slit.

In addition to nodding, random (typically small – a few tens of arcsec or smaller) offsets can be added in the middle of a cycle to locate the spectra at different positions on the array, improving the bad pixel removal and the flat fielding. A sequence of 4 cycles with jittering will result in the

Parameter signature	Header Keyword	Value	Description
DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Filter wheel 1	INS.FILT1.ID	NODEFAULT	Filter wheel 1 position
Filter wheel 2	INS.FILT2.ID	NODEFAULT	Filter wheel 2 position
Which Slit ?	INS.WHICHSLIT	NODEFAULT	Slit
Alpha Offset (arcsec)	TEL.TARG.OFFSETALPHA	0.0	Alpha Offset for sky subtraction (arcsec)
Delta Offset (arcsec)	TEL.TARG.OFFSETDELTA	0.0	Delta Offset for sky subtraction (arcsec)
Add Velocity Alpha	TEL.TARG.ADDVELALPHA	0.0	Additional tracking velocity in RA (arcsec/sec)
Add Delta Velocity	TEL.TARG.ADDVELDELTA	0.0	Additional tracking velocity in DEC (arcsec/sec)
Rotation Offset on Sky	TEL.ROT.OFFANGLE	0.0	Rotation offset on the sky (degrees)
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	False	T - guiding ON, F - guiding OFF
Preset Telescope ? (F/T)	SEQ.PRESET	True	T - full preset, F - fine-tuning of the pointing
Save Image ? (F/T)	SEQ.SAVE	False	T - preserve the acq. image, F - no

Table C.34: SOFI_{img_acq_MoveToSlit}.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Spectro Mode	INS.SMODE	NODEFAULT	Spectroscopic Mode
Which Slit ?	INS.WHICHSLIT	NODEFAULT	Which Slit
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	False	T - guiding ON, F - guiding OFF
Jitter Box Width (arcsec)	SEQ.JITTER.WIDTH	40	Jitter box size
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
Nod Throw (arcsec)	SEQ.NODTHROW	60	Nod throw (arcsec)
NINT	SEQ.NINT	1	Number of exposures in each A or B position
Number of AB or BA cycles ?	SEQ.NABCYCLES	1	Number of AB cycles

Table C.35: SOFI_{spec_obs_AutoNodOnSlit}.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
NSAMP	DET.NSAMP	4	Number of Samples
NSAMPPIX	DET.NSAMPPIX	4	Sample Number per Reading
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Spectro Mode	INS.SMODE	NODEFAULT	Spectroscopic Mode
Which Slit ?	INS.WHICHSLIT	NODEFAULT	Which Slit
Combined offset ? (F/T)	SEQ.COMBINED.OFFSET	False	T - guiding ON, F - guiding OFF
Jitter Box Width (arcsec)	SEQ.JITTER.WIDTH	40	Jitter box size
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
Nod Throw (arcsec)	SEQ.NODTHROW	60	Nod throw (arcsec)
NINT	SEQ.NINT	1	Number of exposures in each A or B position
Number of AB or BA cycles ?	SEQ.NABCYCLES	1	Number of AB cycles

Table C.36: SOFI_spec_obs_AutoNodNonDestr.

following sequence:

$$A(B+\epsilon_1)(B+\epsilon_1)(A+\epsilon_2)(A+\epsilon_2)(B+\epsilon_3)(B+\epsilon_3)(A+\epsilon_4)$$

where ϵ_N are the random jittering offsets. They are generated within the interval defined by `SEQ.JITTER.WIDTH` (in arcsec). If `SEQ.JITTER.WIDTH` is set to zero, then the template will just nod between A and B. It is recommended to use a jitter box width smaller or equal to the nodding length to avoid any chance of overlapping of the spectra on sequential images.

If `SEQ.COMBINED.OFFSET` is set to T and if guiding was started before the start of the template the telescope will guide with *Box-To-Star*. Guiding is recommended for observations longer than 15 min but close to the meridian or the zenith it may be required even for shorter sequences.

The templates have been optimized to minimize overheads. However, after a telescope offset, one or more DITs are skipped to make sure the telescope has stabilized after the move. Otherwise, when the telescope is not moved between frames (i.e. when `NINT`>1, or at the beginning of a new cycle) all DIT-s are kept.

At the end of the templates, the telescope returns to the original position if `SEQ.RETURN` is true (T). If not, the telescope is not moved at the end of the template.

SOFI_spec_obs_AutoNodOnSlit and SOFI_spec_obs_AutoNodNonDestr differ only in the control over the readout mode. The second template allows to define `NSAMP` and `NSAMPPIX` (as explained in the section 3.4.1). Otherwise, they are identical. The only two additional parameters with respect to the `SOFI_spec_obs_AutoNodOnSlit` are `DET.NSAMP` and `DET.NSAMPPIX`, respectively the number of samples and the sample number per reading. During the integration the signal is sampled a number of times defined by `DET.NSAMP`. The larger the number of samples the lower the read out noise. However `DET.NSAMP` must not exceed 60, above this value in fact the glowing of the

Parameter signature	Value
Exposure Name	NGC6118
DIT	40
NDIT	2
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Spectro Mode	LONG_SLIT_BLUE
Which slit	long_slit_1
Combined offset ? (F/T)	T
Jitter Box Width (arcsec)	20
Return to origin ? (T/F)	T
Nod Throw (arcsec)	100.
NINT	3
Number of AB or BA cycles	3

Table C.37: SOFI_spec_obs_AutoNodOnSlit. Example.

Parameter signature	Value
Exposure Name	NGC6118
DIT	50
NDIT	2
NSAMP	30
NSAMPPIX	4
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Spectro Mode	LONG_SLIT_BLUE
Which slit	long_slit_1
Combined offset ? (F/T)	T
Jitter Box Width (arcsec)	20
Return to origin ? (T/F)	T
Nod Throw (arcsec)	100.
NINT	3
Number of AB or BA cycles	3

Table C.38: SOFI_spec_obs_AutoNodNonDestr. Example.

shift registers become the dominant source of noise. For long integrations `DET.NSAMP` between 20 and 40 is advisable. A way to further low the read out noise keeping `DET.NSAMP` small enough is to sample the video signal more than once. The number of sampling of the video signal is `DET.NSAMPIX`, 4 is a good number for this parameter. For long integrations DIT should not exceed 300 sec.

The following examples (Table C.37 and C.38) show typical parameters values. According to the NINT and number of cycles they produce $3 \times 3 \times 2 = 18$ files, each of them corresponding to the average of the DIT x NDIT.

SOFI_spec_obs_GenericSpectro and SOFI_spec_obs_GenSpecNonDestr

These observation templates (Table C.39 and C.40) are for spectroscopy and they have the flexibility to do any sequence of telescope offsets whether they be combined offsets (offset the telescope and guide probe) or non-combined offsets (offset the telescope alone). `SOFI_spec_obs_GenSpecNonDestr` is identical to `SOFI_spec_obs_GenericSpectro` except for the non-destructive readout more.

With this flexibility of these templates comes complexity. They are meant to be used in situations

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures ?	SEQ.NEXPO	1	Number of exposures
Spectro Mode	INS.SMODE	NODEFAULT	Spectroscopic Mode
Which Slit ?	INS WHICHSLIT	NODEFAULT	Which Slit
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
X offset list (arcsec)	SEQ.OFFSETX.LIST	NODEFAULT	List of offsets in X (arcsec)
Y offset list (arcsec)	SEQ.OFFSETY.LIST	NODEFAULT	List of offsets in Y (arcsec)
Observation Type (S or O)	SEQ.OBSTYPE.LIST	NODEFAULT	Observation type: O - object S - sky
Guiding (N B S)	SEQ GUIDING LIST	NODEFAULT	N - no guiding, B - Box-To-Star, S - Star-To-Box DON'T USE

Table C.39: SOFI_spec_obs_GenericSpectro.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
DET.DIT	DET.DIT	NODEFAULT	Detector Integration Time; individual exposure (sec)
DET.NDIT	DET.NDIT	NODEFAULT	Number of DITs averaged into an individual image
NSAMP	DET.NSAMP	4	Number of Samples
NSAMPPIX	DET.NSAMPIX	4	Sample Number per Reading
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures ?	SEQ.NEXPO	1	Number of exposures
Spectro Mode	INS.SMODE	NODEFAULT	Spectroscopic Mode
Which Slit ?	INS WHICHSLIT	NODEFAULT	Which Slit
Return to Origin ? (T/F)	SEQ.RETURN	T	Returns the telescope to the original pointing if True
X offset list (arcsec)	SEQ.OFFSETX.LIST	NODEFAULT	List of offsets in X (arcsec)
Y offset list (arcsec)	SEQ.OFFSETY.LIST	NODEFAULT	List of offsets in Y (arcsec)
Observation Type (S or O)	SEQ.OBSTYPE.LIST	NODEFAULT	Observation type: O - object S - sky
Guiding (N B S)	SEQ GUIDING LIST	NODEFAULT	N - no guiding, B - Box-To-Star, S - Star-To-Box DON'T USE

Table C.40: SOFI_spec_obs_GenSpecNonDestr.

that cannot be accommodated by the nodding templates described above, i.e when the target is too extended and to obtain sky spectrum one has to move further than the slit length allows.

Telescope offsets and guiding options are defined as lists in `SEQ.OFFSETX.LIST`, `SEQ.OFFSETY.LIST` and `SEQ.GUIDING.LIST`. Telescope offsets are relative to the previous position, they are defined **along detector lines (X) and columns (Y)**, and they are in arcsec. **The slit is oriented along the N-S direction, i.e. along the X axis** (N to the left-negative and S to the right-positive) and E-W along the Y axis (E down-negative and W up-positive) for Rotation Offset on Sky `TEL.ROT.OFFANGLEPA=0 deg`. **This parameter is equivalent to the “classical” Position Angle PA on the sky**, i.e. if you want to align the slit along the major axis of a galaxy with `PA=32 deg`, you have to set `TEL.ROT.OFFANGLEPA=32 deg`.

There are three guiding options: N, B, and S. The option N must be used only for observations without guiding, i.e. the offsets are non-combined. This works well for short total integrations (up to 5-10 min) at each individual position. The template will fail if this option is selected, AND the guiding is started. The option B stands for *Box-To-Star*, and it must be used with guiding, i.e. the offsets are combined with with appropriate counter-movement of the guiding arm to maintain the pointing. This is necessary for long total integrations (more than 5-10 min) at each individual position. The option S stands for *Star-To-Box*. It has been disabled (i.e., the template fails if the guiding option is set to S) but it was kept in the software for historic reasons.

Note: With large combined offsets, the guide probe may not be able to follow the one guide star. In such a case, the guiding system will automatically find another star, and resume guiding.

The observation type can be defined for each image, and entered as a list in `SEQ.OBSTYPE.LIST`. O stands for Object and S stands for Sky.

In addition to observations that require large sky offsets, these templates are very useful for slit scanning across an object (i.e. simulation of 3D spectroscopy) by defining a list of offsets in the y direction (perpendicular to the slit direction). A third case when these templates can be used is when the user wants to obtain spectra at more than two positions across the slit - then the offsets will be only in the x direction.

These templates cannot by nature be easily optimized. After each exposure, one or more DITs are skipped.

At the end of the templates, the telescope is returned to the original position if `SEQ.RETURN` is set to true (T). If not, the telescope is not moved.

The lists can have any length; however having lists of different lengths can become extremely confusing. It is good practice to use either lists of one value, or lists of equal length.

The following examples (Table C.41 and C.42) demonstrate a typical usage of the generic templates. Six images will be obtained: 1- the object at the acquisition position; 2- the object in an offset position (-15 arcsec); 3- the sky, the object is off the slit; 4- the object at another offset position; 5- the sky again; and 6- the object at the original position. A Fig. with the object positions is also shown.

Each of the obtained 6 files is the average of 2 (NDIT) exposures of 40 sec (DIT). The guiding option asks to put the box to the star.

C.3.7 SOFI Spectroscopic Calibration Templates

`SOFI_img_cal.Darks`

The template (Table C.26) for imaging darks can be used for obtaining spectroscopic darks, if necessary. However, notice that in practice this is never necessary because: (i) in case of observations the

Parameter signature	Value
Exposure Name	NGC6118
DIT	40
NDIT	2
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of Exposure	6
Spectro Mode	LONG_SLIT_BLUE
Which slit	long_slit_1
Return to origin ? (T/F)	T
X offset list (arcsec)	0 -60 0 120 0 -60
Y offset list (arcsec)	0 0 60 -60 60 -60
Obs Type (O or S)	O O S O S O
Guiding (N B S)	B

Table C.41: SOFI_spec_obs_GenericSpectro. Example.

Parameter signature	Value
Exposure Name	NGC6118
DIT	40
NDIT	2
NSAMP	4
NSAMPPIX	10
Number of columns	1024
Number of rows	1024
First column of window	1
First row of window	1
Number of exposures	6
Spectro Mode	LONG_SLIT_BLUE
Which slit	long_slit_1
Return to origin ? (T/F)	T
X offset list (arcsec)	0 -60 0 120 0 -60
Y offset list (arcsec)	0 0 60 -60 60 -60
Obs Type (O or S)	O O S O S O
Guiding (N B S)	B

Table C.42: SOFI_spec_obs_GenSpecNonDestr. Example.

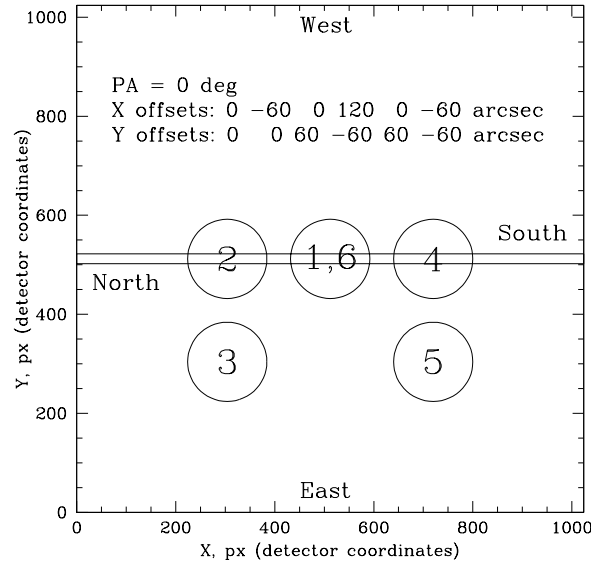


Figure C.1: Positions of the offsets listed in Tables C.42 and C.41 on the detector.

dark (including the “bias”) is subtracted together with the sky, and (ii) in case of the calibrations the respective templates take frames with lamps off to be used as “darks”.

SOFI_spec_cal_Arcs

This calibration template (Table C.43) takes arc spectra with the calibration unit. The calibration mirror is automatically inserted at the beginning of the template, and is automatically removed at the end.

The number of arc frames is defined by `SEQ.NEXPO`. Arc frames can be taken with different DIT and NDIT values, defined as lists with `SEQ.DIT.LIST` and `SEQ.NDIT.LIST`. `SEQ.SPECTROMODELIST` defines the spectroscopic mode, entered as a list with either B, R, Z, J, H, K or NB_1.061 standing for the available spectroscopic modes. In addition, different slits can be called sequentially, with `SEQ.SLIT.LIST`. The allowed values are 1, 2 or 0.6, for `long_slit_1`, `long_slit_2` and `long_slit_0.6`, respectively. Finally, the arc lamps are entered into `SEQ.LAMP.LIST`. Valid values are N, Xe, Ne or B. N stands for None (this is why there is no need to take separate dark frames; but make sure to obtain at least one image with no lamp at each configuration!), Xe stands for Xenon, Ne stands for Neon, and B stands for Both, i.e. Xenon and Neon lamps simultaneously.

Nota Bene: Taking images with no lamp (the option N) has to be specified explicitly in the `SEQ.LAMP.LIST`! This is responsibility of the user. Remember, that these images are necessary to remove the bias/dark/scattered light contribution.

Note that there is no need to take the arc calibration with non-destructive readout because the DIT can always be increased to ensure that the arc spectrum has a sufficient signal-to-noise ratio. There is also no reason to match the readout mode of the arc calibration with the readout mode of the data.

SOFI_spec_cal_DomeFlats and SOFI_spec_cal_NonDestrDomeFlats

Spectroscopic dome flats are taken with these two calibration templates (Table C.44 and C.45). The only difference between the two templates is the readout mode - it has to match the readout mode used to obtain the scientific observations.

For each element in `SEQ.DIT.LIST` and `SEQ.NDIT.LIST`, the templates take four images: the first with the dome lamp off, the next two with the dome lamp on, and the fourth one with the dome

lamp off. The intensity of the dome lamp is controlled manually.

Note that in case of `SOFI_spec_cal_NonDestrDomeFlats` it is not crucial that `NSAMP` and `NSAMPPIX` are the same for the science frames and calibrations. In fact, this is sometimes impossible because the `DIT` of the scientific exposure is rather long, often in order of 5-10 min. Good values for the parameters of the `NonDestrDomeFlat` templates are: `DIT=10`, `NSAMP=6`, `NSAMPPIX=4`.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures ?	SEQ.NEXPO	1	Number of exposures
List of DITs	SEQ.DIT.LIST	NODEFAULT	Detector Integration Time; individual exposure (sec)
List of NDIT	SEQ.NDIT.LIST	NODEFAULT	Number of DITs averaged into an individual image
Spectral Mode List	SEQ.SPECTROMODELIST	NODEFAULT	List of Spectroscopic Modes: B R Z J NB_1.061 H K
Slit List	SEQ.SLIT.LIST	NODEFAULT	Slit list: 0.6, 1, 2
Spectral Lamp List	SEQ.LAMP.LIST	N	Lamp list: N - none, B - both Xe - Xenon, Ne - Neon

Table C.43: SOFI_spec_cal_Arcs.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
Number of exposures ?	SEQ.NEXPO	1	Number of exposures
List of DITs	SEQ.DIT.LIST	NODEFAULT	Detector Integration Time; individual exposure (sec)
List of NDIT	SEQ.NDIT.LIST	NODEFAULT	Number of DITs averaged into an individual image
Spectral Mode List	SEQ.SPECTROMODELIST	NODEFAULT	List of Spectroscopic Modes: B R Z J NB_1.061 H K
Slit List	SEQ.SLIT.LIST	NODEFAULT	Slit list: 0.6, 1, 2

Table C.44: SOFI_spec_cal_DomeFlats.

Parameter signature	Header Keyword	Value	Description
Exposure Name	DET.EXP.NAME	SOFI	File name prefix
Number of columns	DET.WIN.NX	1024	Number of columns in the window
Number of rows	DET.WIN.NY	1024	Number of rows in the window
First column of window	DET.WIN.STARTX	1	First column of window
First row of window	DET.WIN.STARTY	1	First row of window
NSAMP	DET.NSAMP	4	Number of Samples
NSAMPPIX	DET.NSAMPPIX	4	Sample Number per Reading
Number of exposures ?	SEQ.NEXPO	1	Number of exposures
List of DITs	SEQ.DIT.LIST	NODEFAULT	Detector Integration Time; individual exposure (sec)
List of NDIT	SEQ.NDIT.LIST	NODEFAULT	Number of DITs averaged into an individual image
Spectral Mode List	SEQ.SPECTROMODELIST	NODEFAULT	List of Spectroscopic Modes: B R Z J NB_1.061 H K
Slit List	SEQ.SLIT.LIST	NODEFAULT	Slit list: 0.6, 1, 2

Table C.45: SOFI_spec_cal_NonDestrDomeFlats.

Appendix D

Frame Types

There are several basis frame types, which are identified by three keywords in the FITS header. The full list is given in Table D.1.

Keyword	Value	Type
DPR CATG	SCIENCE	Exposure on target (object)
DPR CATG	OTHER	Exposure off target (sky)
DPR CATG	CALIB	Calibration frame
DPR TECH	IMAGE	Exposure in imaging mode
DPR TECH	SPECTRUM	Exposure in spectroscopic mode
DPR TECH	POLARIMETRY	Exposure in polarimetric mode
DPR TYPE	LAMP	Arc spectrum (comparison lamp)
DPR TYPE	DARK	Dark frame
DPR TYPE	STD	Standard Star
DPR TYPE	FLAT	Flat Field Frame

Table D.1: FITS header keywords defining the content of the images.

Appendix E

Photometric Standards

The majority of users prefer to calibrate their data with the NICMOS photometric standard stars (Persson et al. 1998 A.J. 116, 2475). Pre-prepared OBs for observations of most of the standard stars from this list are available at the telescope. For convenience, these standards that are observable from the Southern hemisphere are listed in Tables E.1 and E.2. Finding charts with 2x2 arcmin field of view are shown in Figs. E.1, E.2, E.3, E.4 and E.5. North is up and East is to the left. OBs for all Person standards that can be observed from La Silla (including the extremely red ones) are available in a special queue, accessible via the Observing Tool.

Note that under seeing conditions better than ~ 0.6 arcsec the peak of the images for stars with ~ 10 mag, at any broad band filter, approach the non-linearity limit of the detector (10000 ADU), even for the minimum DIT. This is the case of many red standards, listed in Table E.2. Some defocusing via the secondary mirror M3 may be necessary – typically, shifting the M2 by 0.1 mm in either direction is sufficient. We caution against larger defocusing that would turn the stars into doughnuts because this increases the relative contribution of the sky, degrading the photometric accuracy.

No	HST	RA (J2000)	DEC	J	$\sigma(J)$	H	$\sigma(H)$	K	$\sigma(K)$	K_S	$\sigma(K_S)$
9101	P525-E	00:24:28.3	+07:49:02	11.622	0.005	11.298	0.005	11.223	0.008	11.223	0.005
9103	S294-D	00:33:15.2	-39:24:10	10.932	0.006	10.657	0.004	10.596	0.005	10.594	0.004
9104	S754-C	01:03:15.8	-04:20:44	11.045	0.005	10.750	0.005	10.693	0.010	10.695	0.005
9105	P530-D	02:33:32.1	+06:25:38	11.309	0.010	10.975	0.006	10.897	0.006	10.910	0.005
9106	S301-D	03:26:53.9	-39:50:38	12.153	0.007	11.842	0.005	11.772	0.010	11.788	0.006
9108	P533-D	03:41:02.4	+06:56:13	11.737	0.009	11.431	0.006	11.337	0.008	11.336	0.005
9109	S055-D	04:18:18.9	-69:27:35	11.552	0.002	11.326	0.002	11.255	0.027	11.269	0.002
9111	S361-D	04:49:54.6	-35:11:17	11.246	0.006	11.031	0.006	10.992	0.033	10.980	0.006
9113	S252-D	05:10:25.6	-44:52:46	11.059	0.005	10.776	0.005	10.708	0.034	10.713	0.005
9115	S363-D	05:36:44.8	-34:46:39	12.069	0.007	11.874	0.005	11.826	0.007	11.831	0.005
9116	S840-F	05:42:32.1	+00:09:04	11.426	0.009	11.148	0.009	11.077	0.014	11.058	0.008
9118	S842-E	06:22:43.7	-00:36:30	11.723	0.011	11.357	0.009	11.264	0.016	11.261	0.010
9119	S121-E	06:29:29.4	-59:39:31	12.114	0.006	11.838	0.005	11.765	0.009	11.781	0.005
9121	S255-S	06:42:36.5	-45:09:12	12.719	0.004	11.434	0.004	—	—	11.372	0.004
9125	S005-D	07:19:38.6	-84:35:06	10.885	0.007	10.598	0.006	10.514	0.013	10.522	0.008
9129	S209-D	08:01:15.4	-50:19:33	10.914	0.007	10.585	0.006	10.487	0.021	10.496	0.009
9132	S312-T	08:25:36.1	-39:05:59	11.949	0.006	11.669	0.005	11.608	0.004	11.609	0.004
9133	S495-E	08:27:12.5	-25:08:01	11.521	0.007	11.048	0.008	10.965	0.016	10.960	0.010
9134	P545-C	08:29:25.1	+05:56:08	11.881	0.007	11.624	0.005	11.575	0.005	11.596	0.006
9135	S705-D	08:36:12.5	-10:13:39	12.362	0.010	12.098	0.011	—	—	12.040	0.014
9136	S165-E	08:54:21.7	-54:48:08	12.489	0.008	12.214	0.008	12.138	0.018	12.142	0.011
9137	S372-S	09:15:50.5	-36:32:34	11.153	0.007	10.891	0.007	10.830	0.019	10.836	0.010
9138	S852-C	09:41:35.8	+00:33:12	11.354	0.006	11.041	0.006	10.981	0.015	10.982	0.008
9140	S262-E	09:45:42.8	-45:49:40	11.409	0.011	11.085	0.008	—	—	11.022	0.012
9141	S708-D	09:48:56.4	-10:30:32	11.081	0.008	10.775	0.008	10.715	0.035	10.718	0.010
9143	P550-C	10:33:51.8	+04:49:05	12.344	0.007	12.121	0.005	12.067	0.006	12.081	0.005
9144	S264-D	10:47:24.1	-44:34:05	11.642	0.009	11.335	0.008	11.263	0.018	11.280	0.010
9146	S217-D	12:01:45.2	-50:03:10	12.323	0.007	11.002	0.005	10.931	0.003	10.936	0.004
9147	S064-F	12:03:30.2	-69:04:56	12.111	0.007	11.803	0.007	11.722	0.013	11.724	0.007
9150	S791-C	13:17:29.6	-05:32:37	11.661	0.008	11.310	0.007	11.250	0.014	11.267	0.008
9153	P499-E	14:07:33.9	+12:23:51	11.947	0.008	11.605	0.008	11.560	0.013	11.540	0.008
9154	S008-D	14:23:45.5	-84:09:58	11.232	0.007	10.990	0.007	10.904	0.009	10.915	0.008
9155	S867-V	14:40:58.0	-00:27:47	12.045	0.008	11.701	0.005	11.622	0.005	11.633	0.005
9157	S273-E	14:56:51.9	-44:49:14	11.341	0.007	10.924	0.005	10.851	0.004	10.849	0.004
9160	S870-T	15:39:03.5	+00:14:54	10.914	0.008	10.701	0.008	10.649	0.010	10.659	0.009
9164	P565-C	16:26:42.7	+05:52:20	12.180	0.007	11.895	0.006	11.842	0.007	11.844	0.006
9170	S875-C	17:27:22.2	-00:19:25	11.132	0.005	10.835	0.005	10.739	0.006	10.744	0.005
9172	S279-F	17:48:22.6	-45:25:45	12.477	0.009	12.118	0.006	12.026	0.006	12.031	0.006
9173	S024-D	18:18:46.2	-80:06:58	11.039	0.007	10.778	0.007	10.693	0.009	10.711	0.008
9175	S071-D	18:28:08.9	-69:26:03	12.252	0.006	11.916	0.007	11.834	0.011	11.839	0.007
9178	S808-C	19:01:55.4	-04:29:12	10.966	0.007	10.658	0.008	10.566	0.014	10.575	0.008
9181	S234-E	20:31:20.4	-49:38:58	12.464	0.011	12.127	0.008	12.095	0.007	12.070	0.007
9182	S813-D	20:41:05.1	-05:03:43	11.479	0.005	11.142	0.005	11.082	0.010	11.085	0.005
9183	P576-F	20:52:47.3	+06:40:05	12.247	0.004	11.940	0.004	11.873	0.007	11.880	0.005
9185	S889-E	22:02:05.7	-01:06:02	12.021	0.005	11.662	0.004	11.586	0.012	11.585	0.005
9186	S893-D	23:18:10.0	+00:32:56	11.403	0.009	11.120	0.006	11.045	0.006	11.055	0.006
9187	S677-D	23:23:34.4	-15:21:07	11.857	0.003	11.596	0.003	11.538	0.009	11.542	0.003

Table E.1: Infrared photometric standard stars (Persson et al. 1998, AJ, 116, 2475).

Name	RA (J2000)	DEC	J	$\sigma(J)$	H	$\sigma(H)$	K	$\sigma(K)$	K_S	$\sigma(K_S)$
BRI 0021	00 24 24.6	-01 58 22	11.835	0.008	11.086	0.007	10.552	0.010	10.561	0.008
T832-38078	03 04 02.0	+00 45 52	11.833	0.010	11.248	0.008	10.890	0.013	10.913	0.008
LHS 191	04 26 20.1	+03 36 40	11.621	0.013	11.058	0.012	10.667	0.020	10.717	0.016
IRAS 537 W	05 40 10.5	-07 27 38	12.974	0.009	11.032	0.009	9.981	0.013	10.013	0.011
IRAS 537 S	05 40 15.4	-07 28 46	13.855	0.012	12.087	0.009	10.972	0.014	11.035	0.012
LHS 2026	08 32 30.5	-01 34 37	12.066	0.006	11.497	0.005	11.129	0.007	11.156	0.008
LHS 2397a	11 21 49.2	-13 13 10	11.897	0.008	11.190	0.007	10.691	0.008	10.709	0.010
cskd-8	12 31 13.3	-63 40 21	11.741	0.010	9.913	0.010	9.151	0.011	9.137	0.013
cskd-9	12 31 16.7	-63 40 11	11.372	0.010	9.788	0.010	9.161	0.011	9.154	0.013
cskd-12	12 31 30.0	-63 42 41	11.585	0.007	9.506	0.005	8.617	0.007	8.608	0.007
cskd-15a	12 31 44.5	-63 49 08	10.879	0.010	9.286	0.007	8.625	0.010	8.617	0.011
cskd-16	12 31 57.8	-63 42 21	11.986	0.011	9.982	0.008	9.137	0.011	9.124	0.010
cskd-18	12 31 59.2	-63 41 42	9.365	0.010	9.058	0.008	8.999	0.011	8.992	0.010
cskd-20	12 32 04.0	-63 43 46	9.903	0.011	8.590	0.008	8.089	0.011	8.077	0.010
cskd-21	12 32 10.9	-63 43 16	10.382	0.006	9.918	0.006	9.647	0.007	9.658	0.007
cskd-34	12 31 23.6	-63 46 45	12.494	0.015	10.537	0.014	9.677	0.016	9.668	0.019
cskd-37	12 31 30.1	-63 47 12	12.489	0.015	10.415	0.014	9.528	0.016	9.513	0.019
cske-23	12 31 56.0	-63 37 43	10.624	0.008	8.762	0.006	7.917	0.008	7.889	0.008
cskf-12	12 31 30.1	-63 51 03	9.671	0.008	8.858	0.007	8.561	0.009	8.561	0.008
cskf-13a	12 31 39.5	-63 51 03	10.631	0.010	9.861	0.009	9.486	0.011	9.506	0.011
cskf-14a	12 31 45.9	-63 49 36	10.095	0.008	8.851	0.007	8.387	0.009	8.361	0.008
T868-53850	15 00 26.4	-00 39 29	11.589	0.008	10.993	0.008	10.633	0.009	10.657	0.012
T868-110639	15 10 17.0	-02 41 05	12.612	0.009	11.865	0.010	11.336	0.009	11.353	0.011
L134	15 53 38.4	-04 39 04	11.408	0.009	9.627	0.010	8.875	0.011	8.898	0.013
Oph N9	16 27 13.3	-24 41 34	17.549	0.022	12.389	0.011	9.527	0.014	9.620	0.018
L547	18 51 15.6	-04 16 02	11.872	0.010	9.831	0.011	8.888	0.014	8.870	0.016
BRI 2202	22 05 36.0	-11 04 27	11.652	0.010	11.083	0.009	10.721	0.013	10.736	0.011

Table E.2: Additional red infrared photometric standard stars (Persson et al. 1998, AJ, 116, 2475).

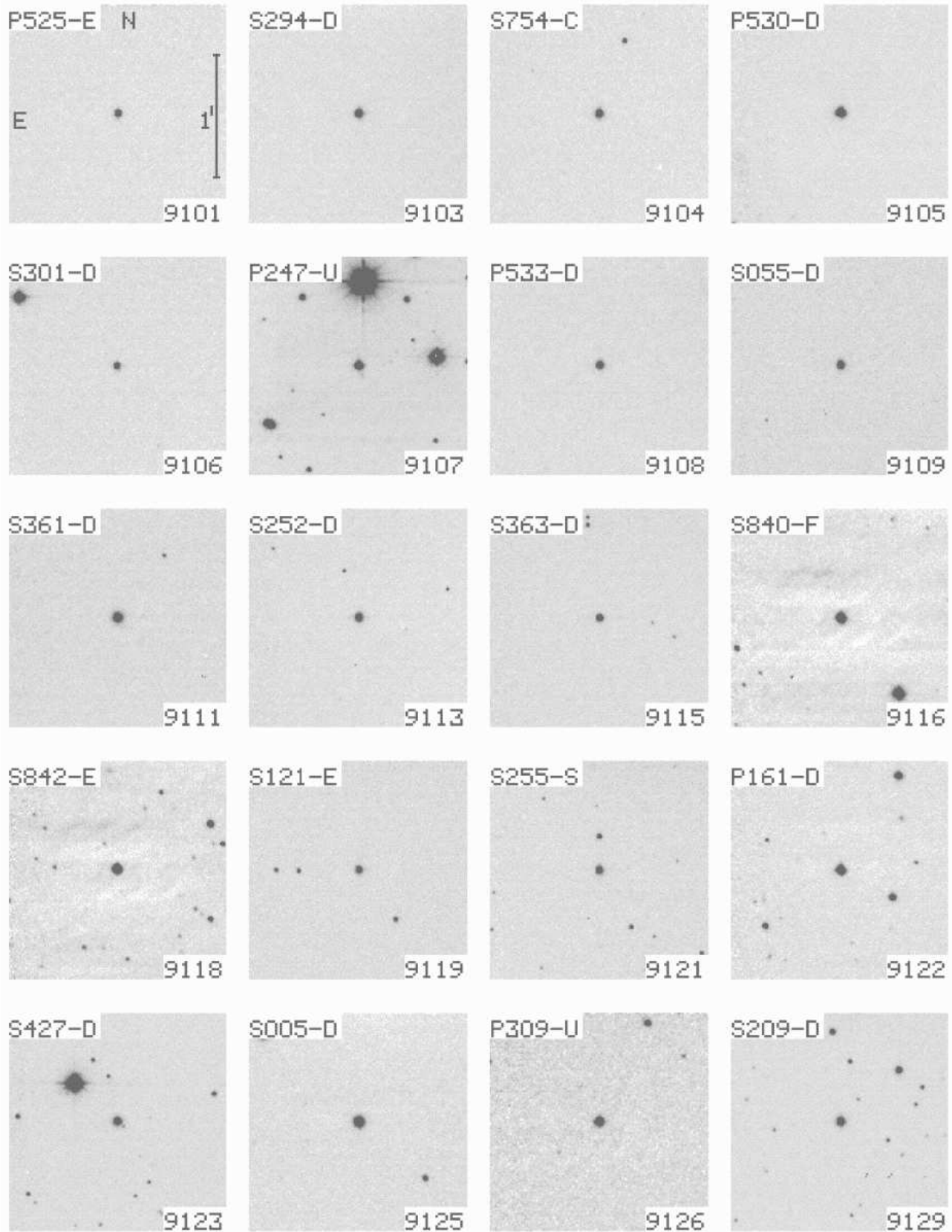


Figure E.1: Finding charts for the photometric standards of Persson et al. (1998). I.

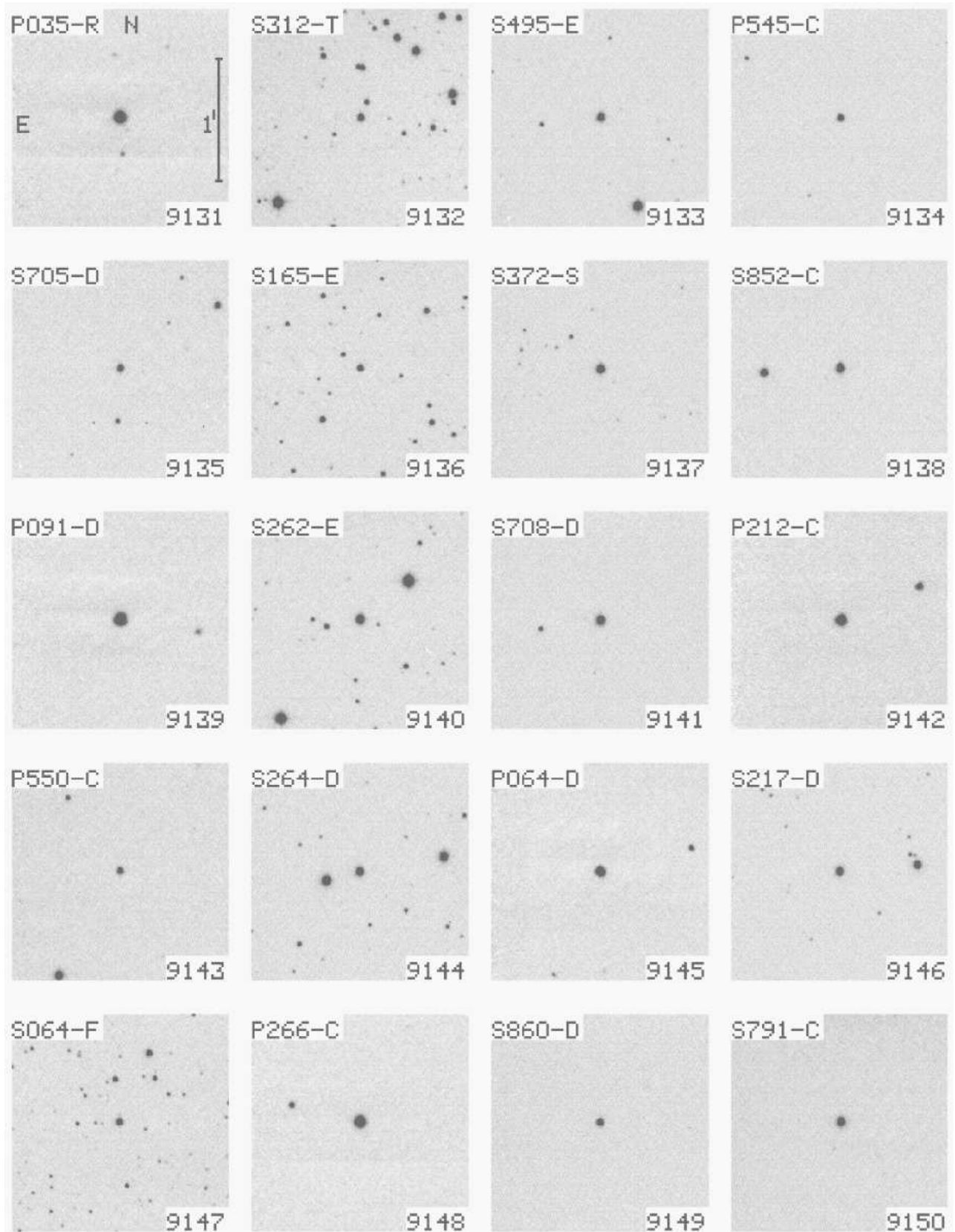


Figure E.2: Finding charts for the photometric standards of Persson et al. (1998). II.

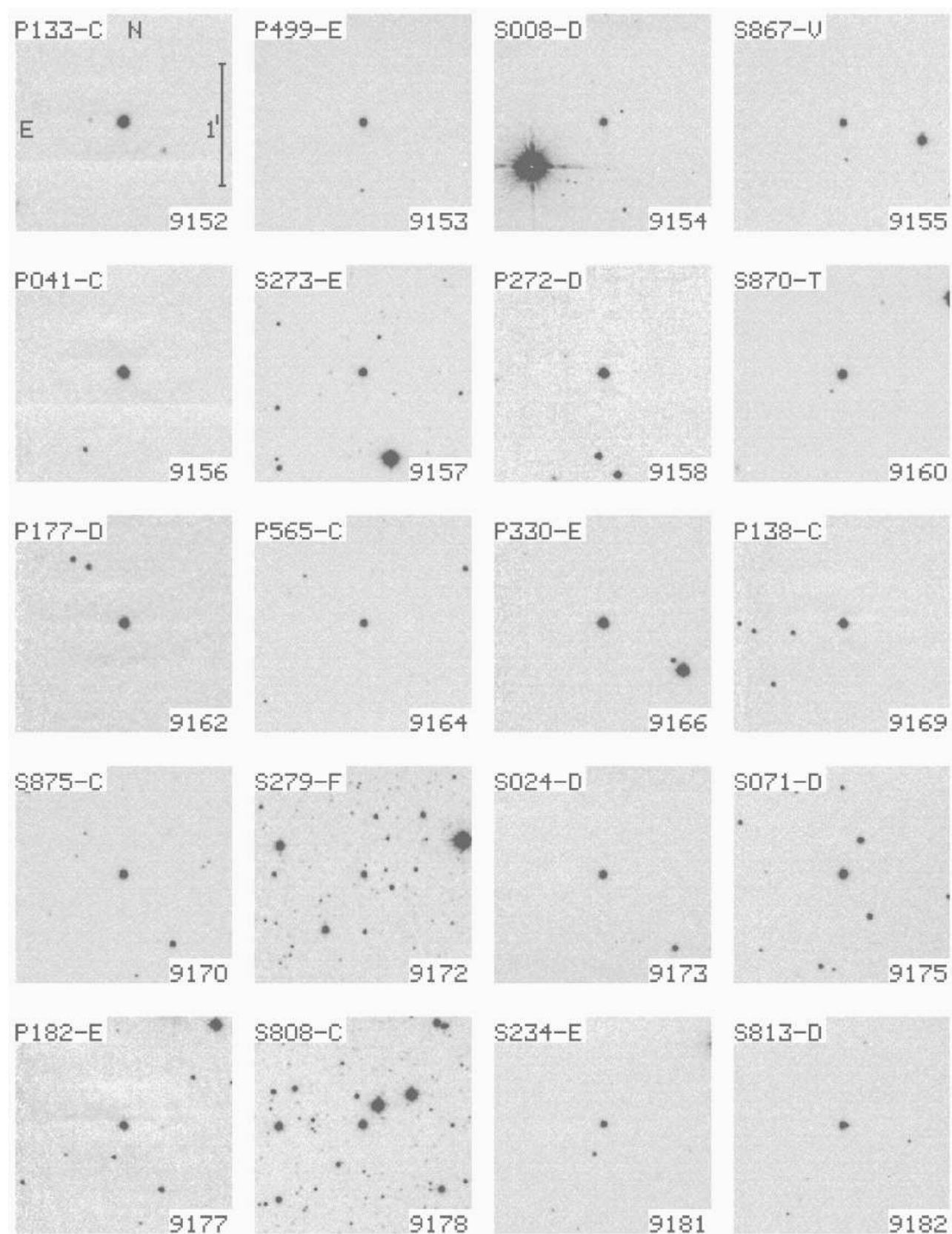


Figure E.3: Finding charts for the photometric standards of Persson et al. (1998). III.

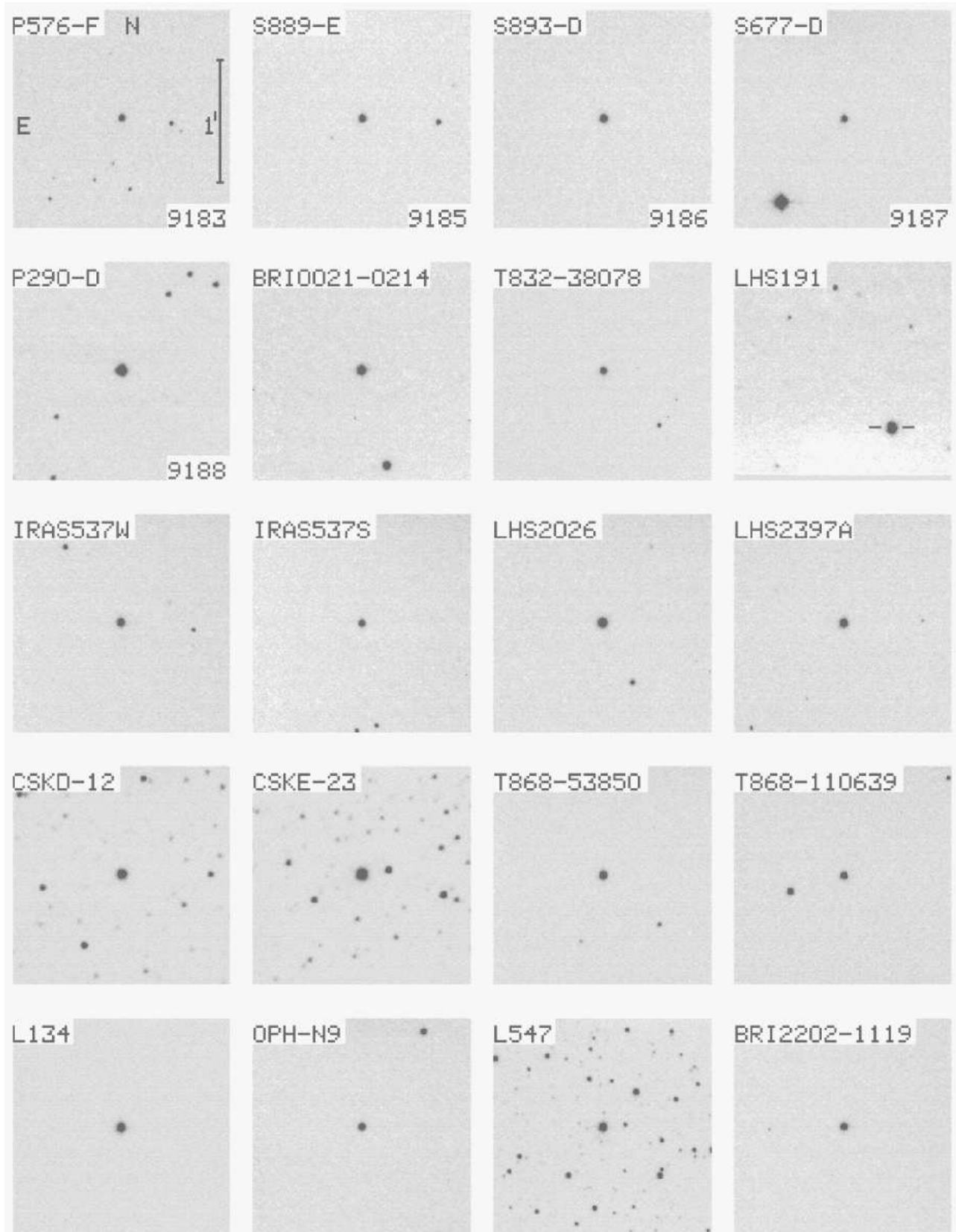


Figure E.4: Finding charts for the photometric standards of Persson et al. (1998). IV.

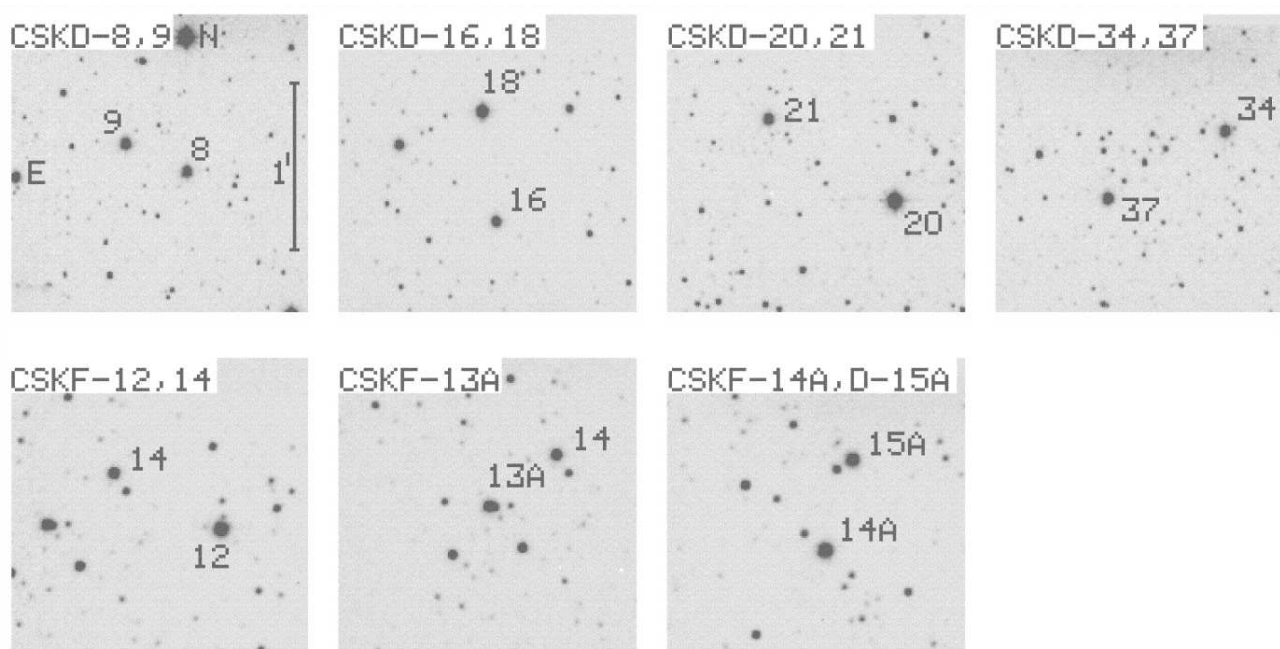


Figure E.5: Finding charts for the photometric standards of Persson et al. (1998). V.

___oOo___