



European Southern Observatory

Organisation Européenne pour des Recherches Astronomiques dans l'Hémisphère Austral  
Europäische Organisation für astronomische Forschung in der südlichen Hemisphäre

Paranal Observatory

Very Large Telescope

Visible Multi-Object Spectrograph

**VIMOS**

**User Manual**

Doc. No. VLT-MAN-ESO-14610-3509

Issue: 78.1

Date: June 6, 2006

**Change Record**

| Issue/Rev. | date          | sections affected                       | Reason/Remarks                                     |
|------------|---------------|-----------------------------------------|----------------------------------------------------|
| draft      | Oct 24, 2002  | all                                     | Draft delivered by VIRMOS Consortium               |
| 1.0        | Feb 11, 2003  | all                                     | first release for P71 (phase II) and P72 (Phase I) |
| 1.1        | July 14 2003  | all                                     | release for P72 (phase II) and P73 (Phase I)       |
| 1.2        | Jan 12, 2004  | all                                     | release for P73 (phase II) and P74 (Phase I)       |
| 1.3        | Jun 14, 2004  | all                                     | release for P74 (phase II) and P75 (Phase I)       |
| 5.0        | Nov 26, 2004  | all                                     | release for P75 (phase II) and P76 (Phase I)       |
| 76.0       | Mar 21, 2005  | 2.7; 3.5.7; 3.6.6;<br>3.7.7             | updated release for P76 (Phase I)                  |
| 76.1       | June 30, 2005 | 2.1;2.2;2.3;2.6;3.4;<br>3.5;3.6;3.7;4.6 | release for P76 (Phase II)                         |
| 77.0       | Aug 31, 2005  | 3.6.1, 3.7.2                            | release for P77 (Phase I)                          |
| 77.2       | Dec 01, 2005  | Appendix                                | release for P77 (Phase II)                         |
| 78.0       | Mar 01, 2006  | none                                    | release for P78 (Phase I)                          |
| 78.1       | Jun 06, 2006  | all                                     | release for P78 (Phase II)                         |

Editors: G. Marconi ESO, Paranal Science Operations; e-mail: gmarconi@eso.org  
S. Bagnulo ESO, Paranal Science Operations; e-mail: sbagnulo@eso.org  
M. Rejkuba ESO, User Support Department; e-mail: mrejkuba@eso.org

# Contents

|          |                                               |           |
|----------|-----------------------------------------------|-----------|
| <b>1</b> | <b>Introduction</b>                           | <b>1</b>  |
| 1.1      | Scope                                         | 1         |
| 1.2      | More Information on VIMOS                     | 1         |
| 1.3      | Contact Information                           | 1         |
| 1.4      | Period of Validity of this Manual             | 2         |
| 1.5      | Version notes, change track                   | 2         |
| 1.5.1    | Version 78.1                                  | 2         |
| 1.6      | Acknowledgements                              | 2         |
| <b>2</b> | <b>Instrument Characteristics</b>             | <b>3</b>  |
| 2.1      | Overview                                      | 3         |
| 2.1.1    | Instrument modes                              | 3         |
| 2.1.2    | Field of view orientation                     | 5         |
| 2.2      | Imaging                                       | 5         |
| 2.2.1    | Image Quality                                 | 6         |
| 2.2.2    | Imaging Filters                               | 6         |
| 2.3      | Multi-Object Spectroscopy (MOS)               | 6         |
| 2.3.1    | Spectral range                                | 7         |
| 2.3.2    | Spectral resolution                           | 7         |
| 2.3.3    | Multiplex characteristics                     | 7         |
| 2.4      | Integral Field Unit (IFU)                     | 7         |
| 2.4.1    | Fiber arrangement and spatial characteristics | 7         |
| 2.4.2    | Spectral range                                | 10        |
| 2.4.3    | Spectral resolution                           | 10        |
| 2.4.4    | Multiplex characteristics                     | 10        |
| 2.5      | Detectors                                     | 10        |
| 2.5.1    | Detector read-out modes                       | 11        |
| 2.5.2    | Controllers                                   | 11        |
| 2.5.3    | Read-out noise and conversion factor          | 11        |
| 2.5.4    | Fringing                                      | 11        |
| 2.5.5    | Shutters                                      | 12        |
| 2.5.6    | More information                              | 12        |
| 2.6      | Calibration Units                             | 12        |
| 2.7      | Flexures                                      | 13        |
| <b>3</b> | <b>Observing with VIMOS</b>                   | <b>15</b> |
| 3.1      | Interaction with the telescope                | 15        |
| 3.1.1    | Guide stars and Guidecam software tool        | 15        |
| 3.1.2    | Guide Stars: Policy                           | 16        |
| 3.1.3    | Cumulative offsets                            | 17        |
| 3.1.4    | Telescope and Instrument Focus                | 17        |
| 3.2      | Atmospheric Dispersion                        | 17        |
| 3.3      | Fringing and Flexures                         | 17        |
| 3.4      | Imaging (IMG)                                 | 18        |
| 3.4.1    | Pre-Imaging                                   | 18        |
| 3.4.2    | Science imaging                               | 19        |

|          |                                                                            |           |
|----------|----------------------------------------------------------------------------|-----------|
| 3.4.3    | IMG mode: calibration plan . . . . .                                       | 19        |
| 3.4.4    | IMG mode: pipeline . . . . .                                               | 19        |
| 3.4.5    | IMG mode: Exposure Time Calculator . . . . .                               | 19        |
| 3.5      | Multi-Object Spectroscopy (MOS) . . . . .                                  | 20        |
| 3.5.1    | Pre-imaging . . . . .                                                      | 20        |
| 3.5.2    | Rotator on sky . . . . .                                                   | 20        |
| 3.5.3    | VIMOS Mask Preparation Software (VMMPS) . . . . .                          | 20        |
| 3.5.4    | Finding charts . . . . .                                                   | 20        |
| 3.5.5    | Visitors using MOS . . . . .                                               | 20        |
| 3.5.6    | MOS mode: grism+filter setup . . . . .                                     | 22        |
| 3.5.7    | MOS mode: attached night time calibrations . . . . .                       | 22        |
| 3.5.8    | MOS mode: calibration plan . . . . .                                       | 22        |
| 3.5.9    | MOS mode: pipeline . . . . .                                               | 22        |
| 3.5.10   | MOS mode: Exposure Time Calculator . . . . .                               | 23        |
| 3.6      | Integral Field Unit (IFU) . . . . .                                        | 23        |
| 3.6.1    | Acquisition accuracy . . . . .                                             | 23        |
| 3.6.2    | Dithering/jitter pattern . . . . .                                         | 23        |
| 3.6.3    | IFU mode: grism + filter setup . . . . .                                   | 23        |
| 3.6.4    | IFU mode: attached night time calibrations . . . . .                       | 23        |
| 3.6.5    | IFU mode: calibration plan . . . . .                                       | 23        |
| 3.6.6    | IFU mode: pipeline . . . . .                                               | 23        |
| 3.6.7    | IFU mode: Exposure Time Calculator . . . . .                               | 24        |
| 3.7      | VIMOS overheads . . . . .                                                  | 24        |
| <b>4</b> | <b>VIMOS templates</b> . . . . .                                           | <b>25</b> |
| 4.1      | Introduction . . . . .                                                     | 25        |
| 4.2      | Template names . . . . .                                                   | 25        |
| 4.3      | Things to know - Reminder . . . . .                                        | 25        |
| 4.4      | List of templates . . . . .                                                | 26        |
| 4.5      | Acquisition templates . . . . .                                            | 26        |
| 4.5.1    | VIMOS_img_acq_Preset . . . . .                                             | 26        |
| 4.5.2    | VIMOS_mos_acq_Mask . . . . .                                               | 28        |
| 4.5.3    | VIMOS_mos_acq_Standard . . . . .                                           | 29        |
| 4.5.4    | VIMOS_ifu_acq_Preset . . . . .                                             | 30        |
| 4.6      | Science templates . . . . .                                                | 31        |
| 4.6.1    | VIMOS_img_obs_Offset . . . . .                                             | 31        |
| 4.6.2    | VIMOS_mos_obs_Offset . . . . .                                             | 32        |
| 4.6.3    | VIMOS_ifu_obs_Offset . . . . .                                             | 33        |
| 4.7      | Calibration Templates . . . . .                                            | 34        |
| 4.7.1    | Nighttime Imaging Calibrations . . . . .                                   | 34        |
| 4.7.1.1  | VIMOS_img_cal_Photom . . . . .                                             | 34        |
| 4.7.1.2  | VIMOS_img_cal_Astrom . . . . .                                             | 35        |
| 4.7.2    | Nighttime MOS Calibrations . . . . .                                       | 35        |
| 4.7.2.1  | VIMOS_mos_cal_NightCalib . . . . .                                         | 35        |
| 4.7.2.2  | VIMOS_mos_cal_Standard . . . . .                                           | 35        |
| 4.7.3    | Nighttime IFU Calibrations . . . . .                                       | 36        |
| 4.7.3.1  | VIMOS_ifu_cal_Specphot . . . . .                                           | 36        |
| 4.7.3.2  | VIMOS_ifu_cal_NightCalib. . . . .                                          | 36        |
| <b>A</b> | <b>VIMOS Filters</b> . . . . .                                             | <b>39</b> |
| <b>B</b> | <b>VIMOS Grisms</b> . . . . .                                              | <b>44</b> |
| <b>C</b> | <b>Template Signature Files for Observation and Acquisition.</b> . . . . . | <b>49</b> |
| C.1      | Acquisition Template . . . . .                                             | 49        |
| C.2      | Observation Template. . . . .                                              | 49        |
| C.3      | Template Signature Files for Calibrations . . . . .                        | 49        |

## **D Abbreviations and Acronyms**

**53**

# List of Figures

|     |                                                                              |    |
|-----|------------------------------------------------------------------------------|----|
| 2.1 | VIMOS Opto-Mechanical layout . . . . .                                       | 4  |
| 2.2 | Image orientation . . . . .                                                  | 5  |
| 2.3 | IFU head - View from the telescope and detail of the fiber modules . . . . . | 9  |
| 2.4 | IFU mode: fibers numbering on the masks . . . . .                            | 9  |
| 2.5 | QE curve of one VIMOS CCD . . . . .                                          | 11 |
| 3.1 | Guidecam Graphical User Interface . . . . .                                  | 16 |
| 3.2 | Atmospheric Dispersion: Example . . . . .                                    | 18 |
| 3.3 | VMMPS Graphical User Interface . . . . .                                     | 21 |
| 4.1 | P2PP: example of a pre-imaging OB . . . . .                                  | 27 |
| 4.2 | P2PP: Example of a MOS OB . . . . .                                          | 28 |
| 4.3 | P2PP: VIMOS ADP browser window. . . . .                                      | 29 |
| A.1 | filter transmission curves . . . . .                                         | 40 |
| A.2 | filter transmission curves . . . . .                                         | 41 |
| A.3 | Order Sorting Filter transmission curves . . . . .                           | 42 |
| A.4 | Order Sorting Filter transmission curves . . . . .                           | 43 |
| B.1 | Efficiency curve of the low resolution grisms. . . . .                       | 45 |
| B.2 | Efficiency curve of the medium resolution grism. . . . .                     | 46 |
| B.3 | Efficiency curve of the high resolution grisms. . . . .                      | 47 |
| B.4 | Efficiency curve of the high resolution grisms. . . . .                      | 48 |

# List of Tables

|      |                                                                                                                      |    |
|------|----------------------------------------------------------------------------------------------------------------------|----|
| 2.1  | VIMOS opto-mechanical characteristics . . . . .                                                                      | 4  |
| 2.2  | VIMOS Imaging characteristics and performance . . . . .                                                              | 6  |
| 2.3  | Standard grism / filter combinations . . . . .                                                                       | 7  |
| 2.4  | MOS: spectral resolution and spectral ranges observed using the mask for spectrophotometric standard stars . . . . . | 8  |
| 2.5  | Spectral coverage $\Delta\lambda$ for three positions of the slit in field of view (MOS mode) . . . . .              | 8  |
| 2.6  | IFU Modes . . . . .                                                                                                  | 8  |
| 2.7  | Spectral characteristics in IFU mode . . . . .                                                                       | 10 |
| 2.8  | IFU: spectral coverage in the various quadrants . . . . .                                                            | 10 |
| 2.9  | Basic characteristics of the VIMOS CCDs . . . . .                                                                    | 11 |
| 2.10 | Detector Readout Modes . . . . .                                                                                     | 11 |
| 2.11 | Detector read-out noise and conversion factors . . . . .                                                             | 12 |
| 3.1  | VIMOS Calibration Plan: IMG Mode . . . . .                                                                           | 19 |
| 3.2  | VIMOS Calibration Plan: MOS Mode . . . . .                                                                           | 22 |
| 3.3  | VIMOS Calibration Plan: IFU Mode . . . . .                                                                           | 24 |
| 4.1  | Science and Nighttime Calibration Templates . . . . .                                                                | 26 |
| A.1  | Characteristics of the VIMOS filters . . . . .                                                                       | 39 |





# Chapter 1

## Introduction

### 1.1 Scope

The VIMOS User's Manual is intended to provide information on the following topics:

- Overall description of the instrument (Sect. 2)
- Observing with VIMOS (Sect. 3)
- Observing and calibration templates (Sect. 4)

### 1.2 More Information on VIMOS

All VIMOS Manuals are available on the VIMOS instrument WEB pages together with the most updated informations on the instrument:

<http://www.eso.org/instruments/vimos/>

Information and software tools for the preparation of service and visitor mode observations with VIMOS are available at:

<http://www.eso.org/observing/p2pp/ServiceMode.html>

<http://www.eso.org/observing/p2pp/VIMOS/VIMOS-P2PP.html>

In particular, information concerning mask preparation with VMPS and guide star selection with Guidecam is given in the following additional manuals:

VMPS  
Guidecam

Visiting astronomers will find further instructions on the Paranal Science Operations WEB page and the Paranal Observatory home page:

<http://www.eso.org/paranal/>

<http://www.eso.org/paranal/sciops/>

### 1.3 Contact Information

In case of specific questions related to Service Mode observations and proposal preparation please contact the ESO User Support Departement:

[usd-help@eso.org](mailto:usd-help@eso.org)

For visitor mode observations please contact the Paranal Science Operations Team:

`paranal@eso.org`

## 1.4 Period of Validity of this Manual

This manual is valid for Phase II preparation of Period 78.

## 1.5 Version notes, change track

### 1.5.1 Version 78.1

Version 78.1 contains more detailed and accurate information about the spectral range that can be observed in spectroscopic modes (MOS and IFU) with the various grisms. Table 2.11 has been updated. Hyperlinks have been added to the text. The various sections of the manual have been re-arranged.

## 1.6 Acknowledgements

The following institutes have participated to the development of VIMOS:

- Laboratoire d'Astrophysique de Marseille, Marseille, France
- Observatoire de Haute-Provence, Saint-Michel l'Observatoire, France
- Laboratoire d'Astrophysique de l'Observatoire Midi-Pyrénées, Toulouse, France
- Osservatorio Astronomico di Capodimonte, Napoli, Italy
- Istituto di Fisica Cosmica e Tecnologie Relative, Milano, Italy
- Osservatorio Astronomico di Brera, Milano, Italy
- Istituto di Radioastronomia, Bologna, Italy
- Osservatorio Astronomico di Bologna, Bologna, Italy

The first edition of this User Manual was delivered by the VIRMOS Consortium (P.I. O. Le Fèvre).

# Chapter 2

## Instrument Characteristics

### 2.1 Overview

VIMOS is the Visible Multi-Object Spectrograph for the Very Large Telescope (VLT) of the European Southern Observatory (ESO). This instrument has been built to provide the ESO community with a wide field spectrograph with high throughput and multiplex, dedicated to deep surveys.

VIMOS operates in the 0.37–1 micron domain in 3 main observing modes: direct imaging (IMG), multi-slit spectroscopy (MOS), and integral field spectroscopy (IFU). Because of the large field of view, VIMOS optical path is split into four channels. The four quadrants of the instrument are operated in parallel (i.e., it is not possible to mix modes or filters or grisms between quadrants).

The field of view of the four channels is  $4 \times 7' \times 8'$  in IMG and MOS and up to  $54'' \times 54''$  in IFU mode in low resolution spectroscopy. It provides an image scale of  $0.205''/\text{pixel}$  in MOS and  $0.67''/\text{resolution element}$  in IFU mode. Each of the 4 channels has a  $2048 \times 4096$  pixels EEV CCD, with the 4k pixels being used along the dispersion to maximize the spectral coverage. VIMOS is installed on the Nasmyth B focus of UT3 (Melipal). The VIMOS mechanical and optical characteristics are summarized in Table 2.1.

A schematic opto-mechanical layout is shown in Fig. 2.1.

- The focal plane is divided in 4 identical channels and is equipped with: the mask exchange units (MEU) which host up to 15 pre-punched masks for MOS spectroscopy per quadrant; the IFU fiber-head, and the IFU masks;
- The folding mirror section hosts the four folding mirrors and the flexure compensators.
- The Filter/Camera Section: the beam section hosts optical analyzers (filters and/or grisms), the camera, and the exposure shutter in front of the CCD.
- FEU: the filter exchange units (one for each channel) are located below the main body of the instrument below the Filter/Camera Section; each unit hosts a juke-box of 10 filters permanently mounted.
- GEU: the grism exchange units are located below the Filter/Camera Section; they host a carousel of six grisms permanently mounted.

#### 2.1.1 Instrument modes

VIMOS offers the following modes for observations:

- Imaging (IMG) (see Sect. 2.2)
- Multi-Object Spectroscopy (MOS) (see Sect. 2.3)
- Integral Field Spectroscopy (IFU) (see Sect. 2.4)

Table 2.1: VIMOS opto-mechanical characteristics

|                        |                                                                                                                                                                |
|------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Instrument Location    | Nasmyth B VLT-UT3 (Melipal)                                                                                                                                    |
| Opto-mechanical layout | 4 beams, each a complete focal-reducer, F/1.88 output                                                                                                          |
| Wavelength coverage    | 0.37 to 1 microns                                                                                                                                              |
| Detectors              | $4 \times 2048 \times 4096$ , $15 \mu\text{m}/\text{pixel}$                                                                                                    |
| Spectral Resolution    | $R = 180$ to $2500$ (1 arcsec slit)                                                                                                                            |
| Filters                | 10 per channel, <i>U'BVRIz</i> , OS-red, OS-blue, GG435, GG475, 170mm diameter                                                                                 |
| Grism                  | 6 per channel, LR red, LR blue, MR, HR blue, HR orange, HR red <sup>[1]</sup>                                                                                  |
| Flexures               | Passive compensation, motion of $\pm 2$ pixels over $360^\circ$ rotation                                                                                       |
| Masks                  | 15 masks simultaneously loaded in instrument at maximum (per channel) <sup>[2]</sup> .<br>Slits of any position and shape, width $> 0.6''$ , length $< 30''$ . |
| Multiplex              | 840 simultaneous slits, $10''$ long at $R = 200$<br>210 simultaneous slits, $10''$ long at $R = 2000\text{-}2500$                                              |
| Integral Field Unit    | $54'' \times 54''$ field, 6400 fibers with $0.67''$ sampling<br>$27'' \times 27''$ field, 6400 fibers with $0.33''$ sampling                                   |

<sup>[1]</sup> From October 2005, a set of four new and more efficient HR red VPHG gratings is installed and available in the 4 channels.

<sup>[2]</sup> Some restrictions apply in operations; currently the number of masks that can be mounted at the instrument cabinet is 8. Please refer to the VIMOS web pages for updates.

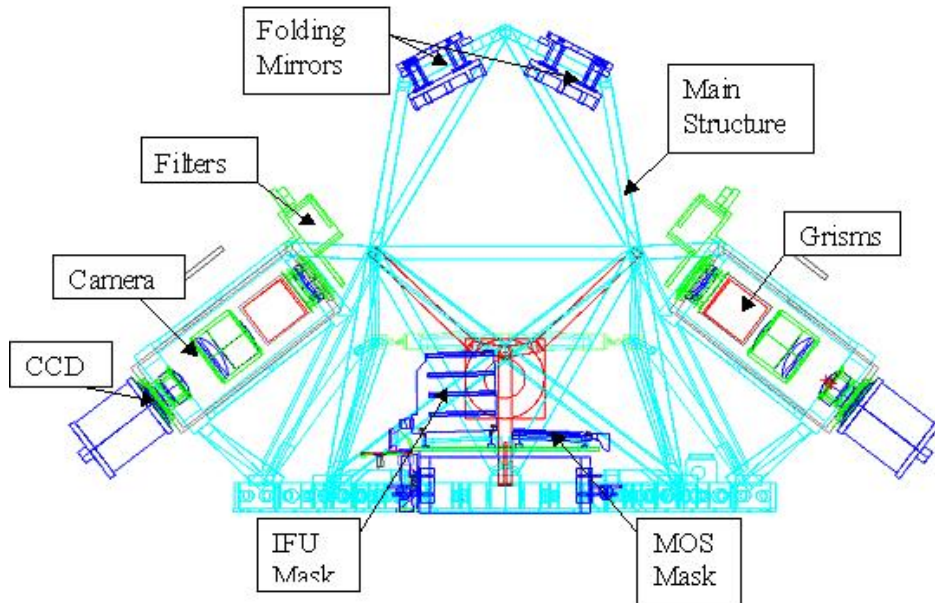


Figure 2.1: VIMOS Opto-Mechanical layout

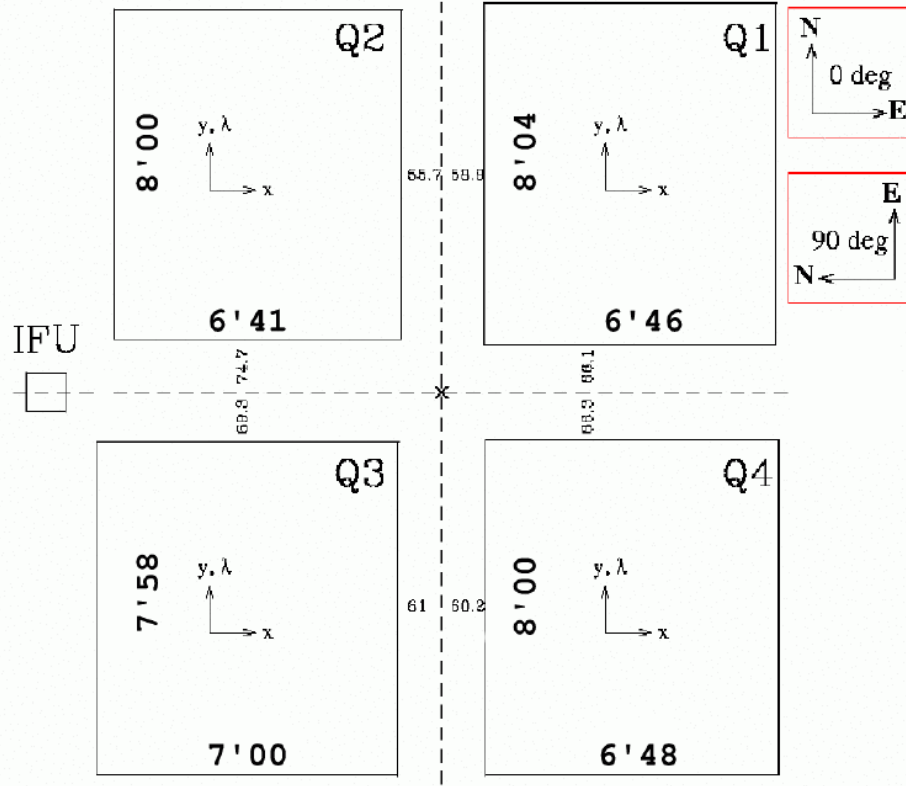


Figure 2.2: Image orientation. Layout of the imaging field of view. Indicated are the orientation on sky for different rotator angles, the image and wavelength axes, the location of the IFU head.

IMG mode uses broadband filters  $U'BVRIZ$ . MOS and IFU modes share the set of grisms and order sorting filters given in Table 2.3. In MOS mode, a set of four laser punched masks made of INVAR is inserted in the focal plane. In IFU mode, fibers are arranged along pseudo-slits in the so-called “IFU masks”. IFU masks may be inserted in the instrument focal plane (in the same position as occupied by the INVAR masks in the MOS mode). The spectral characteristics of the two spectroscopic modes are therefore similar, but not identical.

### 2.1.2 Field of view orientation

Figure 2.2 shows the orientation of the field of view of the instrument. At  $0^\circ$  rotator angle the  $x$  axis of the CCDs is oriented E and the  $y$  axis is oriented N. Increasing rotator angle moves the compass counterclockwise, e.g., at  $90^\circ$  rotator angle the  $x$  axis is oriented S and the  $y$  axis oriented E. The wavelength axis in spectroscopy is the  $y$  axis, with wavelength increasing with  $y$ .

## 2.2 Imaging

Direct imaging is used to produce wide field images (VIMOS is the largest imager on the VLT) using broad band filters. This mode has two functions: produce science images, and produce the mandatory pre-images needed to prepare masks (see Sect. 3).

The gap between the quadrants is approximately  $2'$ . The pointing coordinates correspond approximately to the center of the gap.

Table 2.2: VIMOS Imaging characteristics and performance

|                               |                                                                                                                                      |
|-------------------------------|--------------------------------------------------------------------------------------------------------------------------------------|
| Field of view                 | 4×7' × 8' separated by 2' gaps<br>illuminated CCD area: 2048 × 2350 pix <sup>2</sup><br>slightly vignetted on the corners outside 8' |
| Pixel scale                   | 0''.205/pixel                                                                                                                        |
| Image Quality                 | 80 % encircled energy in 0''.4 within the whole field                                                                                |
| Filters                       | <i>U'</i> , Bessel <i>BVRI</i> , Gunn <i>z</i> (see section A)                                                                       |
| F-ratio                       | 1.88                                                                                                                                 |
| Limiting Magnitudes in 1h exp | V=27.7 I=26.6 5 $\sigma$ detection for a point source<br>V=26.5 I=25.4 5 $\sigma$ detection in 3'' for a galaxy                      |

Table 2.2 summarizes the main characteristics of the imaging mode.

### 2.2.1 Image Quality

For the most updated information please have a look at the Quality Control and Data Flow Operations [WEB pages](#)

The focal plane of VIMOS is not strictly flat because of the complexity of the optical elements involved. This results in slight focus variations across the field, which translates in a variable PSF from the center to the periphery. The sky concentration effects are in any case small and negligible for flat-field and photometric calibrations.

### 2.2.2 Imaging Filters

VIMOS provides positions for 10 broadband filters in any of the four channels. Presently installed filters inside the instruments are *U'*, Bessel *BVRI*, Gunn *z*, and some order separation filters. Only the *U'BVRIz* filters are available for imaging programmes. The special *U'* is not a standard Bessel filter (see Figures A.1). With the Bessel *I* and Gunn *z* filters the CCDs show fringes of approximately 7%.

The complete list of filters together with the transmission curves are presented in Appendix A (see Table A.1, and Figs. A.1 to A.4).

ESO will not accept requests for user's supplied filters for VIMOS.

## 2.3 Multi-Object Spectroscopy (MOS)

A variety of grisms with different wavelength ranges and dispersions is available in MOS mode. Order separation filters are available to eliminate the overlap, *for a given slit*, between orders 1 (the order used for science) and order 2 (always present on the detector). This is done by restricting the wavelength range of a particular grism to less than one octave in wavelength. The 0th order contamination is present only in multi-layered mode, i.e., when observing with the LR and the MR grisms.

Table 2.3 indicates the recommended and allowed grism + filter combinations. **Any other combinations of grisms and filters should be requested with a waiver.** Defining in the OBs combinations not allowed and for which a waiver was not requested will result in these OBs being rejected by the User Support Department.

ESO will not accept requests for user's supplied grisms or filters for VIMOS.

Transmission curves are available in Appendix A for the filters, and Appendix B for the grisms.

Table 2.4 shows the spectral range covered in MOS modes by the various grisms, as obtained by the automatic pipeline when reducing the spectrophotometric standard star. This spectral range is valid for slits close to the center of the field.

Table 2.3: Standard grism / filter combinations for MOS and IFU modes

| Grism     | Filter  |                              |
|-----------|---------|------------------------------|
|           | Default | Allowed                      |
| LR blue   | OS-blue | Free                         |
| LR red    | OS-red  | Free                         |
| MR        | GG475   | OS-red, OS-blue <sup>1</sup> |
| HR blue   | Free    | –                            |
| HR red    | GG475   | –                            |
| HR orange | GG435   | –                            |

<sup>1</sup> The grism-filter configuration MR+OS-blue is not supported by the pipeline

### 2.3.1 Spectral range

In MOS mode, slits can be positioned at any position in the imaging field of view. With high-resolution grisms, the observed spectral range changes slightly according to the position of the slit. Table 2.5 shows the spectral ranges covered by the detector for each grism (and default order sorting filter) for the slits located at three different positions in the field of view. With low resolution grisms, the spectral coverage is independent of the position of the slit.

### 2.3.2 Spectral resolution

The spectral resolution obtained with a 1" slit width with the various grisms is given in Table 2.4.

### 2.3.3 Multiplex characteristics

A particular feature of the MOS mode is the possibility of stacking multiple layers of slits along the dispersion direction. This is only possible with the low resolution grisms (up to four slits along the dispersion direction) and the intermediate resolution grisms (up to two slits), for which the spectra are short enough (in number of pixels) to allow stacking. With the high-resolution grisms, only one slit can be placed along the dispersion direction.

However, please note that in the multi-layer mode there is some order overlap: the second order of one slit may overlap with the first order of the adjacent slit. There is also 0<sup>th</sup> order contamination between different layers of slits.

For faint objects, order overlap is usually not a problem: only the sky will contribute significantly to the overlap, with a contribution that cannot be distinguished from the first order. This however constrains the positioning of the slits in the focal plane: stacked slits along the dispersion direction have to be strictly parallel and of the same length to allow the second orders to overlap nicely with the first order of the adjacent slits.

With the low-resolution grisms, the use of *no filter* instead of the default order sorting filters OS-blue and OS-red will result in longer spectra and second order overlap. It is **strongly recommended to reduce the multiplex/number of slits along the dispersion direction** to avoid/reduce multiple overlaps. The multiplex can be reduced by increasing the distance between slits.

## 2.4 Integral Field Unit (IFU)

### 2.4.1 Fiber arrangement and spatial characteristics

The Integral Field Unit (IFU) is a dedicated opto-mechanical system aimed at producing spectra of a contiguous area of up to  $54 \times 54$  arcsec<sup>2</sup>. The sky image is projected onto a micro-lens array with a choice of two spatial samplings (magnifications): 0.33" or 0.67" per resolution element (fibre). Each 2D input micro-lens is coupled to an optical fiber. The output of the fibers is re-arranged on a linear set of micro-lenses to produce an entrance

Table 2.4: MOS: spectral resolution and spectral ranges observed using the mask for spectrophotometric standard stars

| Grism               | Default Order sorting filter | Spectral range          | Number of pixels | Spectral Res. (1" slit) | Dispersion (Å/ pixel) | Spectral multiplex |
|---------------------|------------------------------|-------------------------|------------------|-------------------------|-----------------------|--------------------|
| LR blue             | OS-Blue                      | 370-670 nm              | 550              | 180                     | 5.3                   | 4                  |
| LR red              | OS-Red                       | 550-950 nm              | 550              | 210                     | 7.3                   | 4                  |
| MR                  | GG475                        | 480-1000 nm             | 2000             | 580                     | 2.5                   | 2                  |
| MR                  | OS red                       | 550-970 nm              | 2000             | 580                     | 2.5                   | 2                  |
| HR blue             | Free                         | 415-620 nm <sup>1</sup> | 4096             | 2050                    | 0.51                  | 1                  |
| HR orange           | GG435                        | 520-760 nm <sup>2</sup> | 4096             | 2150                    | 0.6                   | 1                  |
| HR red <sup>4</sup> | GG475                        | 630-870 nm <sup>3</sup> | 4096             | 2500                    | 0.6                   | 1                  |

<sup>1</sup> With HR blue grisms the spectral range in quadrant 1 (for a standard star) is 400-620 nm

<sup>2</sup> With HR orange grisms the spectral range in quadrant 1 (for a standard star) is 505-745 nm

<sup>3</sup> With HR red grisms the spectral range in quadrant 1 (for a standard star) is 650-865 nm

<sup>4</sup> New HR\_red VPHG grisms have been installed from October 5, 2005.

Table 2.5: Spectral coverage  $\Delta\lambda$  for three positions of the slit in field of view (MOS mode)

| Grism     | Order sorting filter | $\Delta\lambda$ | $\Delta\lambda$ | $\Delta\lambda$ |
|-----------|----------------------|-----------------|-----------------|-----------------|
|           |                      | CCD center      | CCD top         | CCD bottom      |
| LR red    | OS-red               | 550 - 950 nm    | 550 - 950 nm    | 550 - 950 nm    |
| LR blue   | OS-blue              | 360 - 670 nm    | 360 - 670 nm    | 360 - 670 nm    |
| MR        | GG475                | 500 - 1000 nm   | 500 - 1000 nm   | 500 - 950 nm    |
| HR blue   | Free                 | 415 - 620 nm    | 468 - 674 nm    | 370 - 566 nm    |
| HR red    | GG475                | 630 - 870 nm    | 690 - 930 nm    | 570 - 805 nm    |
| HR orange | GG435                | 520 - 760 nm    | 577 - 817 nm    | 463 - 703 nm    |

slit to the spectrograph, at the mask focal plane location. Each IFU mask (one for each channel of VIMOS) has a set of 4 fibre-slits (or pseudo-slits), for a total of 6400 fibers.

The special arrangement of fibers for VIMOS IFU head is showed in Fig. 2.3. The fibers numbering has been arranged following the scheme: mask - line - module - fiber (i.e. 4 - B - 5 - 67). The arrangement of the fibers into the 4 pseudo-slits mask (one for each channel) is reported in Fig. 2.4. The detailed IFU table is distributed to the users together with the data release.

A masking-shutter (IFU shutter) is provided in front of the input micro-lens array; this allows the exclusive use of the central 1/4<sup>th</sup> of the field when IFU is used with intermediate and high spectral resolution grisms. In this configuration, only 1 fiber-slit per quadrant (mask) is used, for a total of 1600 fibers. The templates for IFU mode offer the option to insert a shutter to limit the field to the central part also for low resolution grisms. When small field can be accommodated, this avoids second order contamination by the sky or background objects from the 3 pseudo-slits corresponding to the outer part of the field of view.

Table 2.6 shows the spatial characteristics of the IFU mode.

Table 2.6: Integral Field Unit: spatial characteristics

| Grism                          | Shutter | Spatial sampling (arcsec/fibre) | Field of view arcsec <sup>2</sup> | Number of fibres | Pseudo-slits per quadrant |
|--------------------------------|---------|---------------------------------|-----------------------------------|------------------|---------------------------|
| LR blue, LR red                | OFF     | 0.67"/fibre                     | 54" x 54"                         | 80 x 80          | 4                         |
| LR blue, LR red                | ON      | 0.67"/fibre                     | 27" x 27"                         | 40 x 40          | 1                         |
| LR blue, LR red                | OFF     | 0.33"/fibre                     | 27" x 27"                         | 80 x 80          | 4                         |
| LR blue, LR red                | ON      | 0.33"/fibre                     | 13" x 13"                         | 40 x 40          | 1                         |
| MR, HR blue, HR orange, HR red | ON      | 0.67"/fibre                     | 27" x 27"                         | 40 x 40          | 1                         |
| MR, HR blue, HR orange, HR red | ON      | 0.33"/fibre                     | 13" x 13"                         | 40 x 40          | 1                         |



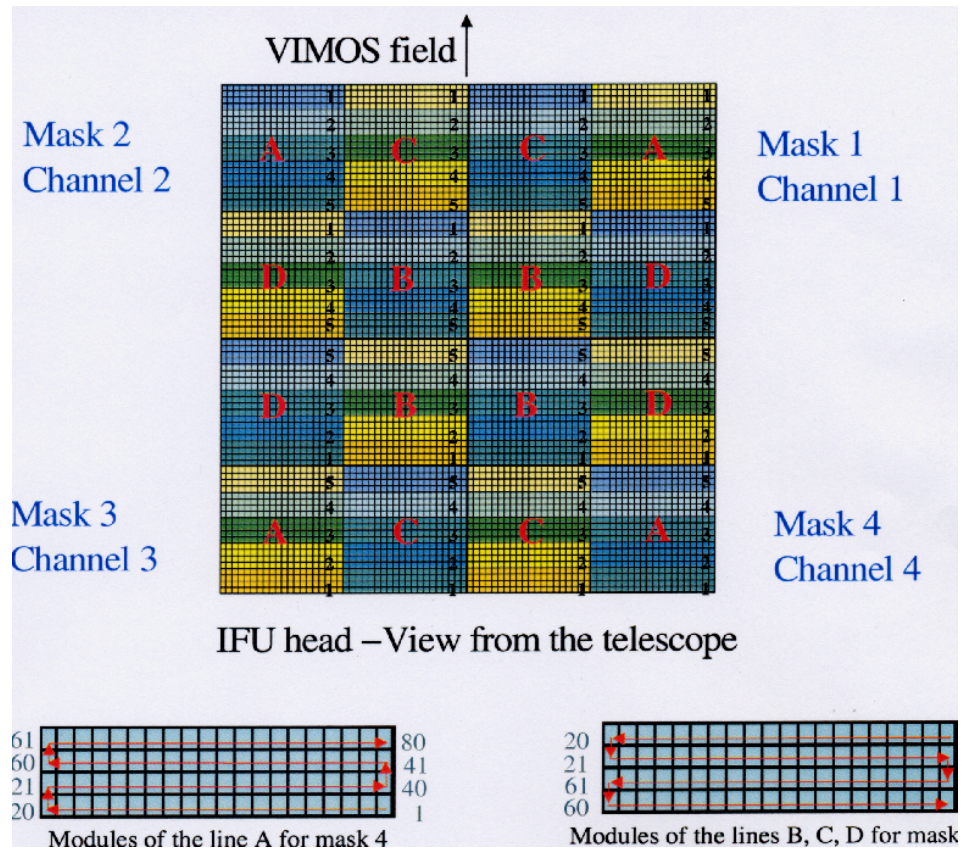


Figure 2.3: IFU head - View from the telescope and detail of the fiber modules

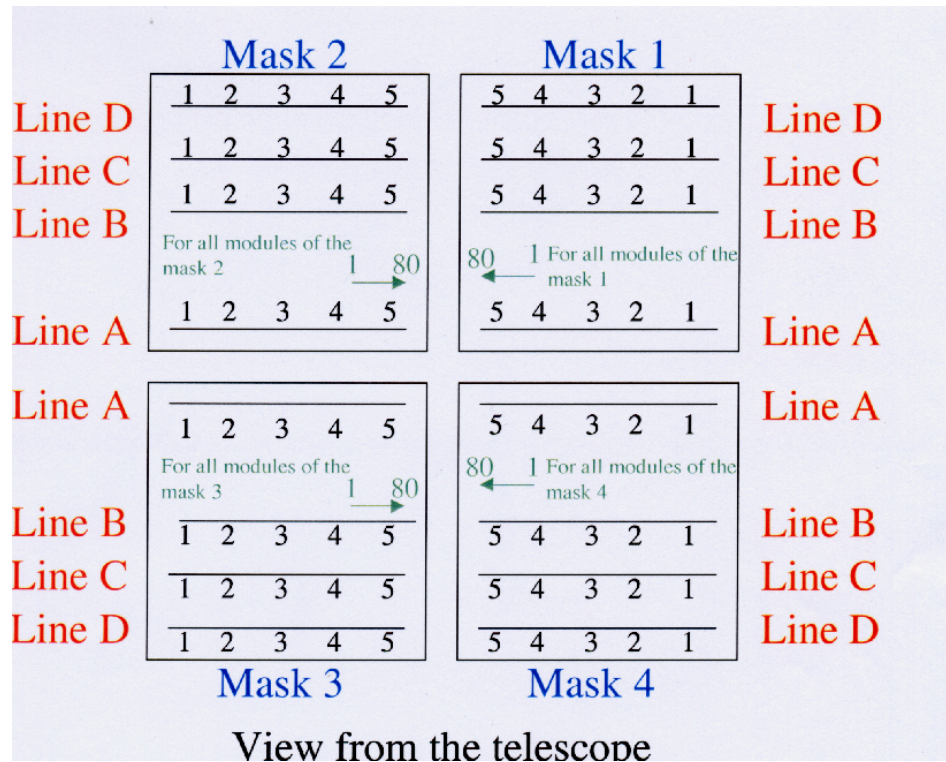


Figure 2.4: IFU mode: fibers numbering on the masks

Table 2.7: Spectral characteristics in IFU mode

| Grism     | Default Order<br>sorting filter | Spectral range<br>common to all quadrants | Spectral Res.<br>(1" slit) | Dispersion<br>(Å / pixel) | Spectral<br>multiplex |
|-----------|---------------------------------|-------------------------------------------|----------------------------|---------------------------|-----------------------|
| LR blue   | OS-Blue                         | 400-670 nm                                | 220                        | 5.3                       | 4                     |
| LR red    | OS-Red                          | 590-915 nm                                | 260                        | 7.3                       | 4                     |
| MR        | GG475                           | 490-1015 nm                               | 720                        | 2.5                       | 1                     |
| HR blue   | None                            | 415-620 nm                                | 2550                       | 0.51                      | 1                     |
| HR orange | GG435                           | 525-740 nm                                | 2650                       | 0.6                       | 1                     |
| HR red    | GG475                           | 645-870 nm                                | 3100                       | 0.6                       | 1                     |

**Note:** new HR red VPHG grisms have been installed for October 5, 2005.

Table 2.8: IFU: spectral coverage in the various quadrants

| Grism             | Q1            | Q2            | Q3            | Q4            |
|-------------------|---------------|---------------|---------------|---------------|
| LR blue + OS blue | 390 - 670 nm  | 395 - 675 nm  | 400 - 680 nm  | 400 - 680 nm  |
| LR red + OS red   | 570 - 915 nm  | 580 - 930 nm  | 590 - 930 nm  | 585 - 930 nm  |
| MR + GG475        | 485 - 1015 nm | 490 - 1020 nm | 490 - 1020 nm | 490 - 1015 nm |
| HR blue + free    | 415 - 620 nm  | 415 - 620 nm  | 415 - 620 nm  | 415 - 620 nm  |
| HR orange + GG435 | 525 - 745 nm  | 525 - 745 nm  | 525 - 745 nm  | 525 - 740 nm  |
| HR red + GG475    | 635 - 860 nm  | 635 - 865 nm  | 635 - 860 nm  | 645 - 875 nm  |

## 2.4.2 Spectral range

Table 2.7 shows the spectral range covered in IFU modes by the various grisms, as obtained by the automatic pipeline when reducing the spectrophotometric standard star.

With low and intermediate resolution grisms the spectral coverage is independent of the position of the IFU pseudo-slits. With the high-resolution grisms, the spectral coverage is determined also by the (fixed) positions of the pseudo-slits in the focal plane, but it varies slightly from quadrant to quadrant. The spectral coverage for each quadrant derived from the automatic reduction of standard star observations with VIMOS pipeline is listed in the Table 2.8.

## 2.4.3 Spectral resolution

In IFU mode, spectral resolution is fixed by the combination grism + fiber size (the fiber FWHM is about 3.2 pixels, corresponding to 0.66"). The spectral resolution is about 1.25 times the spectral resolution obtained in MOS mode with a 1" slit width. The spectral resolution for various grism+filter combinations is given in Table 2.7.

## 2.4.4 Multiplex characteristics

The same multiplex capability described in MOS mode is also used in IFU mode. With the low-resolution grisms there are 4 pseudo slits of fibers along the dispersion direction. With intermediate and high-resolution grisms only the central pseudo-slit is used (in practice by masking the entrance field of the IFU so that the other pseudo-slits do not receive sky light).

## 2.5 Detectors

The detector mosaic of VIMOS consists of four 2k×4k EEV44 CCDs, thinned and back side illuminated.

Table 2.9 indicates the basic characteristics of the VIMOS CCDs. Figure 2.5 shows the quantum efficiency curve of one of the CCDs.

Table 2.9: Basic characteristics of the VIMOS CCDs

| Parameter                                                         | Value     |
|-------------------------------------------------------------------|-----------|
| Number of pixels                                                  | 2048×4096 |
| Pixel size ( $\mu\text{m}$ )                                      | 15        |
| Dark current at $-120^\circ\text{C}$ ( $e^-/\text{px}/\text{h}$ ) | <7        |
| Deviation from Linearity (Full well; %)                           | < 0.5     |
| Charge transfer efficiency                                        | >0.999999 |

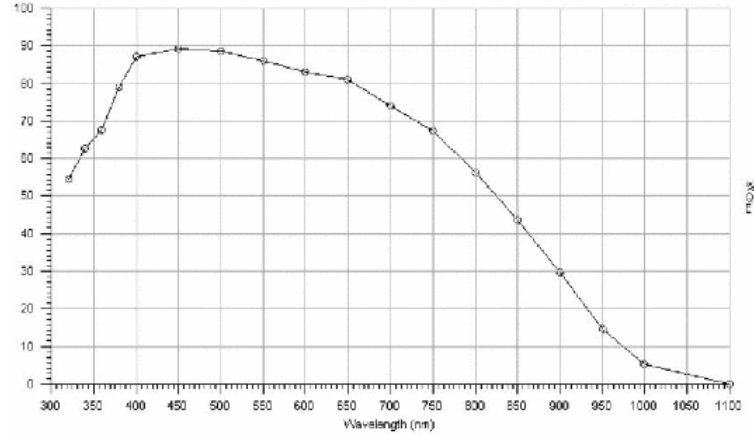


Figure 2.5: QE curve of one VIMOS CCD

### 2.5.1 Detector read-out modes

The CCDs have two read-out ports but only one is used. This increases the read-out time but as a benefit the gain, offset and read-out noise are constant along the image. The default readout modes are indicated in Table 2.10.

Table 2.10: Detector Readout Modes and Characteristics

| VIMOS Mode   | Size    | Speed                  | Binning | Gain | Port   | Windowing     | Readout time |
|--------------|---------|------------------------|---------|------|--------|---------------|--------------|
| Imaging      | 4×2k×2k | 225kpx s <sup>-1</sup> | 1×1     | Low  | 1-port | Not supported | 30s + 15s    |
| Spectroscopy | 4×2k×4k | 225kpx s <sup>-1</sup> | 1×1     | High | 1-port | Not supported | 45s + 20s    |

### 2.5.2 Controllers

The 4 CCDs are controlled by two FIERA controllers. The detectors of channel 2 and 3 are controlled by FIERA A (VIMOSA), the detectors of channel 1 and 4 by FIERA B (VIMOSB).

### 2.5.3 Read-out noise and conversion factor

Table 2.11 indicates the readout noise and gain values for the 4 chips. Note that these values are subject to change and are indicated here only for reference. Exact values have to be found in the image headers.

RON and conversion factor monitoring is available through the [QC1 WEB pages](#)

### 2.5.4 Fringing

With the  $z$ -Gunn and Bessel  $I$  filters the fringing amplitude is around 7% in imaging and in spectroscopy the amplitude of fringing is of the order of 7-8% at maximum.

Table 2.11: Detector read-out noise and conversion factors. For reference only, subject to change

| Chip  | Read-out Mode                                           | RON [e <sup>-</sup> rms] | K [e <sup>-</sup> /adu] |
|-------|---------------------------------------------------------|--------------------------|-------------------------|
| CCD 1 | 225 kpx s <sup>-1</sup> , low gain (IMG mode)           | 5.3                      | 1.70                    |
| CCD 2 | 225 kpx s <sup>-1</sup> , low gain                      | 5.2                      | 1.86                    |
| CCD 3 | 225 kpx s <sup>-1</sup> , low gain                      | 5.3                      | 1.95                    |
| CCD 4 | 225 kpx s <sup>-1</sup> , low gain                      | 4.9                      | 1.80                    |
| CCD 1 | 225 kpx s <sup>-1</sup> , high gain (MOS and IFU modes) | 4.3                      | 0.51                    |
| CCD 2 | 225 kpx s <sup>-1</sup> , high gain                     | 4.6                      | 0.56                    |
| CCD 3 | 225 kpx s <sup>-1</sup> , high gain                     | 4.7                      | 0.57                    |
| CCD 4 | 225 kpx s <sup>-1</sup> , high gain                     | 4.3                      | 0.52                    |

Despite this relatively strong fringing, it is important to note that the structure is quite stable and dithered images taken in a sequence can easily remove the structure (see Sect. 3.3).

### 2.5.5 Shutters

The shutter (one for CCD camera) guarantees uniform illumination of the CCD to the 1 % level or better for exposure times as short as 1 sec (the shortest possible exposure time).

### 2.5.6 More information

More (updated) information on detector characteristics cosmetic and performance on sky at:

<http://www.eso.org/projects/odt/>  
<http://www.eso.org/observing/dfo/quality/>

## 2.6 Calibration Units

VIMOS has three calibration units illuminating a calibration screen at the back of the Nasmyth shutter (closing the Nasmyth tunnel).

Each unit has the following lamps:

- Flatfield lamps:
  - QTH10 halogen lamp for imaging screen flats in *VRIz*, and MOS spectroscopic screen flats in LR red and blue
  - QTH50 halogen lamp for imaging screen flats in *U'B*, MOS screen flats in MR, HR red and blue, and all IFU screen flats
- Arc lamps:
  - He
  - Ar
  - Ne

The updated atlas of the spectral lines for the different grisms could be find at

<http://www.eso.org/instruments/vimos/inst/atlas/index.html>

## 2.7 Flexures

The passive flexure compensation of VIMOS is optimized to reach a reasonable compromise in every position of the Nasmyth rotator. An astatic compensator system is installed in quadrants 2, 3, and 4. Image motion due to instrument flexures under gravity are currently within  $\pm 2$  pixels for the four channels.

See Sect. 3.3 for recommendations on how to deal with fringing and flexures.





## Chapter 3

# Observing with VIMOS

All observations with VIMOS are done via “observing blocks” (OBs). OBs consist of one acquisition template including target information followed by “observing templates” related to the observing mode. The users will fill out the parameter fields (“keywords”) of the templates (e.g. grisms, filters, etc). All the preparations are done with the Phase II proposal preparation tool [p2pp](#).

In addition, VIMOS masks have to be prepared with the VIMOS mask manufacturing preparation software [VVMPS](#).

The strategy behind observing blocks and templates is to prepare the observations well in advance to minimize any interactive steps during the observations (optimization and service mode compatibility). The execution of the OBs will be mostly automatic and the execution will be done by telescope and instrument operators and staff astronomers. Direct interaction at execution time is only required for the target acquisition and to assess the quality of the data.

The preparation of visitor and service mode observations requires special care, and the relevant documentation should be consulted. See Sect. 1 for instructions on how to retrieve the corresponding information.

### 3.1 Interaction with the telescope

#### 3.1.1 Guide stars and Guidacam software tool

For telescope guiding and active optics correction it is necessary to find a guide star of magnitude in the range 11-13.5 within the Nasmyth field of view. The guide star is picked up by the Guide Probe, which is attached to the telescope *adapter*. (Note that the adapter is independent of the telescope *rotator* to which the instrument is attached.)

In IMG and MOS mode, it is not always possible to find a guide star such that the Telescope Guide Probe is not vignetting the VIMOS field of view. Therefore, users are asked to select the guide star for their observations. The choice of the guide star is **mandatory** for pre-imaging, for imaging, and for MOS programs. In fact, the guide star for MOS is selected at the time of pre-imaging: guide star coordinates and position of guide probe (whether positive or negative) is stored in the pre-image headers, and subsequently propagated by VVMPS (see Sect. 3.5.3), the VIMOS Mask Preparation Software, into the Aperture Definition in Pixel (ADP) files. At the time of the MOS acquisition, this information is extracted from the ADP files and the same guide star as used for pre-imaging will be re-used for the MOS observations. Only in some exceptional cases, the telescope operator may have to select another guide star: see Sect. 3.1.2. Users need not select the guiding star for observations carried out in IFU mode.

The position and the amount of the vignetting introduced in the field of view by the guide probe can be controlled with the help of the [Guidacam](#) software, that allows one to visualize the geometry of the VIMOS field of view and of the Telescope Guide Probe, and allows users to select the guide star most suitable to their observation.

Figure 3.1 shows the Guidacam Graphical User Interface displaying the geometry of the VIMOS field of view

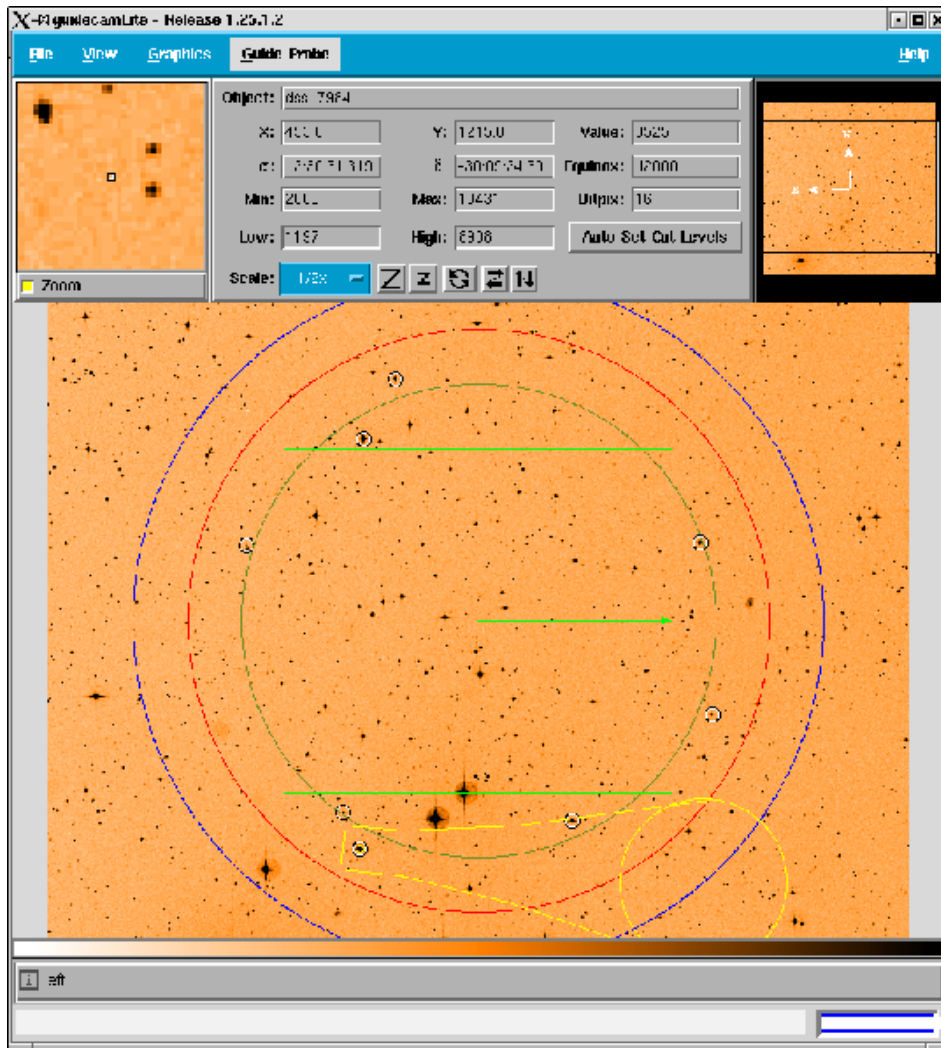


Figure 3.1: Guidecam Graphical User Interface. Superimposed to the field of view are the VIMOS field of view layout, Guide Probe drawing and USNO selected stars in the field.

and the Guide Probe at the Nasmyth focal plane of the telescope.

### 3.1.2 Guide Stars: Policy

Starting from P72, ESO applies the following policy regarding guide stars:

- Selection of a Guide Star for VIMOS observations is mandatory for MOS and IMG mode, to assess and control the amount and position of vignetting that may be introduced by the Guide Probe; see Phase 2 Webpages for detailed informations.
- If no guide star is provided by the user (i.e. the option **CATALOGUE** of the **Get Guide Star From P2PP** parameter of the **VIMOSifu.acq.Preset** template is selected), the guide star will be selected by the Telescope Operator at execution time.
- If the user provides a guide star (by selecting the option **SETUPFILE** of the **Get Guide Star From P2PP** parameter of the **VIMOSimg.acq.Preset** template), efforts will be made to use this guide star for imaging, pre-imaging and for the subsequent MOS observations, although this cannot be entirely guaranteed. Observing conditions (seeing, Moon, etc.) may indeed force to select a different guide star than the one selected by the user. In addition, it may happen that a guide star suitable for preimaging observations will not be usable for the spectroscopic follow-up if e.g. the seeing has changed.



- It is requested that the users selecting a guide star for their observations use Guidecam to generate the finding chart to be delivered as part of their Phase II material.

### 3.1.3 Cumulative offsets

For small telescope offsets, the telescope may keep the same guide star. However, observations involving cumulative offsets greater than about  $30''$  often result in significant vignetting of the guide probe, or even in the loss of the guiding star. For these reasons, all dither positions for offset patterns used in the OBs must fit inside a  $30'' \times 30''$  size box. If larger offsets are necessary, users need to prepare separate OBs and use appropriate guide star for each of them.

### 3.1.4 Telescope and Instrument Focus

The telescope focus is automatically set by the active optics system. No intervention is required by the observer. Defocussing of the telescope is not possible during the observations. The instrument focus is corrected automatically for the different thickness of the various filters, and for temperature variations (autofocus).

## 3.2 Atmospheric Dispersion

VIMOS is not equipped with Atmospheric Dispersion Corrector in none of its modes. The impact of atmospheric dispersion will be mostly noticed in MOS mode. In imaging, some image elongation may be noticed at high airmasses in the U or B band. In IFU, the effect will mostly be noticed with the LR.blue grism and with the  $0.33''/\text{fiber}$  sampling.

In MOS mode, the effects can lead to severe slit losses if some precautions are not taken. To minimize atmospheric dispersion effects at high airmass, and assuming that the observations are performed close to meridian, the slits need to be oriented N-S, forcing the rotator angle to be set at  $90^\circ$ . This rotator angle is the default one to be used for pre-imaging, and the same rotator angle will be used during spectroscopic follow-up. The VIMOS operational scenario foresees to carry out MOS observations as close as possible to Meridian, within  $\pm 2$  hrs. Although this represents a major operational constraint, this is the only safe mode to ensure that slit losses are minimized in a systematic way for the blue setups, provided that slits are oriented N-S.

More information and simulations can be found on the [VIMOS web page](#). Figure 3.2 shows a worst case scenario of MOS observations with LR.blue (spectral coverage 370 – 670 nm) for a field with declination  $+20^\circ$ .

## 3.3 Fringing and Flexures

The amount of flexures reported in Sect. 2.7 has a number of operational consequences:

- The accuracy of the slit positions determined from pre-imaging depends on the rotator position. This is for the time being one of the reasons why slit widths narrower than  $0''.6$  are not allowed.
- Accuracy of sky subtraction / flat fielding: due to fringing at red wavelengths, accuracy of sky subtraction / flat fielding may be limited by flexures. With dithering / jittering sequences, median running sky frames can be built in imaging or spectroscopy allowing accurate sky subtraction and/or flat fielding of the data
- In spectroscopy, arcs and flats are taken during daytime and at constant rotator angle. Depending on the variation in rotator angle between daytime calibrations and night time observations, some flat field residuals and wavelength calibration offsets may be introduced due to the flexures. In order to minimize/eliminate these effects, night time calibrations can be taken with an adequate template (see Sect. 4). **Attached flat and arcs calibrations are mandatory for all IFU observations. Attached arcs are mandatory for MOS science observations with grisms LR blue and HR blues, and highly recommended for MOS science observations with grisms MR, HR red, HR orange, and HR red.** Execution time is charged to the user program.

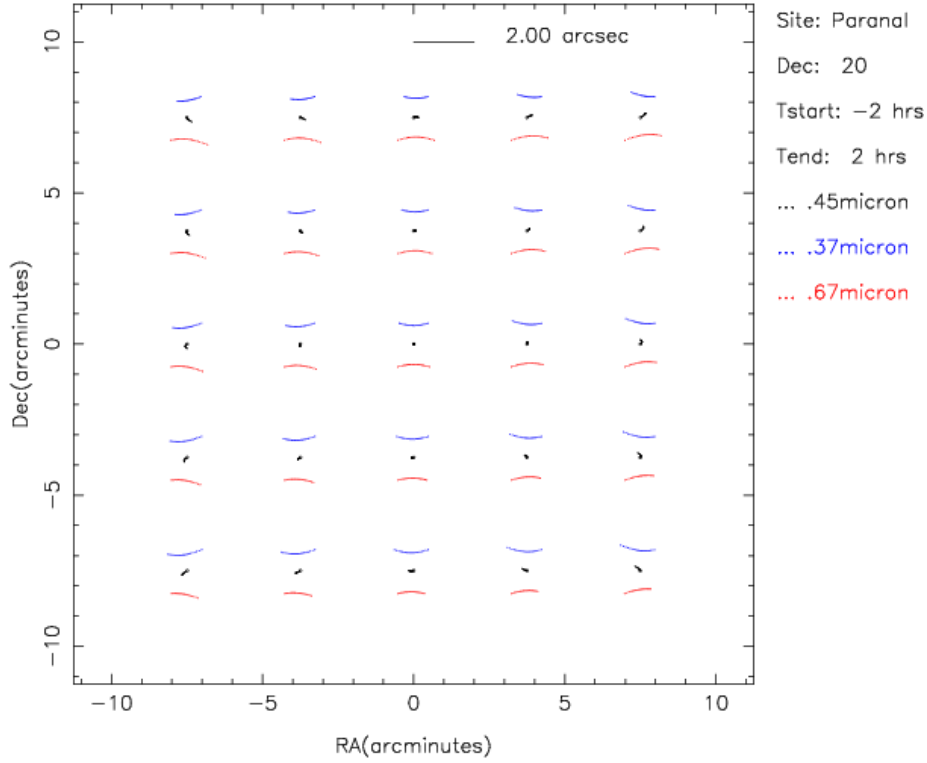


Figure 3.2: Example of atmospheric dispersion at Paranal for a field at declination  $+20^\circ$ , observed with the LR\_blue grism from  $\pm 2$  hrs from Meridian. The blue and red trails show the evolution of the blue and red wavelengths during the exposure. Atmospheric spectra start with N-E orientation and move to N-W through a N-S orientation while crossing the Meridian. The maximum amplitude of the dispersion is  $\sim 2''$  at  $\pm 2$  hrs from Meridian. Slits clearly have to be oriented N-S as a default orientation to minimize the slit losses, assuming that all MOS observations will be carried out within  $\pm 2$  hrs from Meridian.

### 3.4 Imaging (IMG)

Imaging mode is used for both pre-imaging runs and science imaging runs. In both cases, users are requested to use Guidecam for the selection of guide stars. See Sect. 3.1.1 concerning Guidecam and Sect. 3.1.2 concerning the policy regarding Guide Stars.

#### 3.4.1 Pre-Imaging

VIMOS pre-Imaging is mandatory prior to MOS follow-up. Pre-imaging is carried out in service mode only. The following instructions for preparation of pre-imaging OBs should be strictly followed.

**Program ID** It is mandatory for Phase I to request a specific pre-imaging run in the proposal form, and for Phase II to prepare the pre-imaging OBs by using the corresponding Program ID.

**Filter used** The correction of the distortions is derived from *R* filter observations only. To minimize the slit losses due the distortion of the instrument, the use of the *R* filter for preimaging is mandatory. Pre-imaging with filters other than *R* may be requested with a waiver.

**Observation Category: PRE-IMAGE** It is mandatory to set the Observation Category to PRE-IMAGE in the observation template VIMOS\_img\_obs.Offset (see Sect. 4). Failure to do so will result in delays to process and deliver the pre-imaging data, and may result in MOS observations not being carried out.

Table 3.1: VIMOS Calibration Plan: IMG Mode

| Calibration                  | Number      | Frequency | When?         | Products          | Accuracy | Comment                |
|------------------------------|-------------|-----------|---------------|-------------------|----------|------------------------|
| Bias                         | 5           | Daily     | Day time      | Master Bias       | RON      |                        |
| Darks                        | 3           | Monthly   | Day time      | Dark current      |          |                        |
| Sky Flats                    | 5           | 3-7 days  | Twilight time | Master Flat       | 2%       | <i>U'BVRIZ</i> filters |
| Photom. Std.                 | 4 quadrants | Daily     | Night time    | Zero Point        | 5%       | <i>U'BVRI</i> filters  |
| Photom. Std.<br>high airmass | 1 quadrant  | Weekly    | Night time    | Extinction coeff. | 10%      | <i>U'BVRI</i> filters  |
| Astrometric fields           | 4           | 6 months  | Night time    | Distortion map    | 1 pixel  | <i>R</i> filter        |

**Rotator on sky** The rotator on sky, otherwise called position angle on sky, cannot be freely selected by the user. This is in order to minimize the slit losses due to atmospheric dispersion (Sect. 3.2). The default value assigned to this parameter in P2PP should therefore not be modified. The default value is 90°, and corresponds to an orientation of the field with E to the top, N to the left (see Fig. 2.2). In *special cases* it may be possible to deviate from these rules *if* there is a valid scientific reason; in that case a waiver must be specifically requested to ESO before you submit your Phase II information.

### 3.4.2 Science imaging

**Observation Category: SCIENCE** It is mandatory to set the Observation Category to SCIENCE in the observation template VIMOS\_img\_obs\_Offset (see Sect. 4), otherwise data will be treated as pre-imaging and will not be processed as regular imaging data.

**Rotator on sky** Users are encouraged, but not forced, to use the default value of the rotator angle, i.e., 90°.

### 3.4.3 IMG mode: calibration plan

Table 3.1 indicates the nature, frequency and accuracy of VIMOS calibrations expected to be taken as part of the calibration plan. Only calibration frames required for the reduction of science data are reported. More calibrations will be taken by the operation staff for the purpose of maintaining and monitoring the instrument configuration and performance, e.g. mask to CCD transformations etc. These calibrations are not reported here. Calibration monitoring, including photometric zero points, can be found on the pages of the Quality Control group:

<http://www.eso.org/observing/dfo/quality/>

### 3.4.4 IMG mode: pipeline

See <http://www.eso.org/qc/pipeline-status.html> for the information concerning the latest status of the ESO pipeline.

### 3.4.5 IMG mode: Exposure Time Calculator

The approximate limiting magnitudes in direct imaging mode can be obtained from the exposure time calculator available on the ESO WEB pages <http://www.eso.org/observing/etc/>.

## 3.5 Multi-Object Spectroscopy (MOS)

### 3.5.1 Pre-imaging

Pre-Imaging is mandatory for all MOS observations. Pre-Imaging is carried out in service mode (see Sect. 3.4.1).

### 3.5.2 Rotator on sky

Because of the atmospheric dispersion effects (Sect. 3.2), the orientation on sky for MOS and pre-imaging observations is not a parameter set by the user, but is forced by the system to be 90°. All service and visitor mode users who have particular constraints for the orientation of their fields should request a waiver.

### 3.5.3 VIMOS Mask Preparation Software (VMMPS)

The use of VMMPS (Sect. 3.5.3) is compulsory for the creation of the so-called ADP (Aperture Definition in Pixels) files which are subsequently loaded into p2pp when preparing Observation Blocks. ADP files will be used in Paranal to manufacture the masks.

Here are some general features of VMMPS, for a quick reference guide. Consult the [VIMOS web page](#) for instructions to retrieve the package and the documentation.

- VMMPS allows to define a mask in two situations:
  - From a contributed catalogue, consisting of list of targets with accurate *relative* astrometry (better than 0''/3 rms). Input catalog is a list of RA-DEC coordinates. This catalogue needs to have, with the same *relative* astrometric accuracy, coordinates of bright objects / stars in the field of view. These objects / stars, after detection on the pre-images, will be used to register the target coordinates to the VIMOS coordinate reference frame. Some of these bright stars will be used as reference for acquisition.
  - From a catalog of target positions (in pixel coordinates) **extracted from a VIMOS image**.
- VMMPS maximizes the number of slits that can be assigned to the targets in the input catalogue, in one or more masks, depending on user's choice.
- Upon realization of a first mask definition, VMMPS allows manual interaction to unselect targets, define curved slits, tilted slits, etc. However, note that manual modification of the mask definition should be done only by experienced users and with care.
- VMMPS allows to define some objects as compulsory
- VMMPS requires that some bright stars are defined as reference targets for acquisition purposes.

Figure 3.3 shows the VMMPS Graphical User Interface.

### 3.5.4 Finding charts

Finding charts for service mode OBs need to be prepared with VMMPS and show all targets with slits overplotted.

### 3.5.5 Visitors using MOS

Visitors using the MOS mode are requested to prepare Observation Blocks ahead of their scheduled visitor run, and upon completion of their pre-imaging observations (like in service mode). This is to guarantee that visitors will have a sufficient knowledge of the VIMOS Mask Preparation Software (VMMPS, Sect. 3.5.3) prior to their visit at Paranal, and to allow mask manufacturing ahead of time. Limited additional support

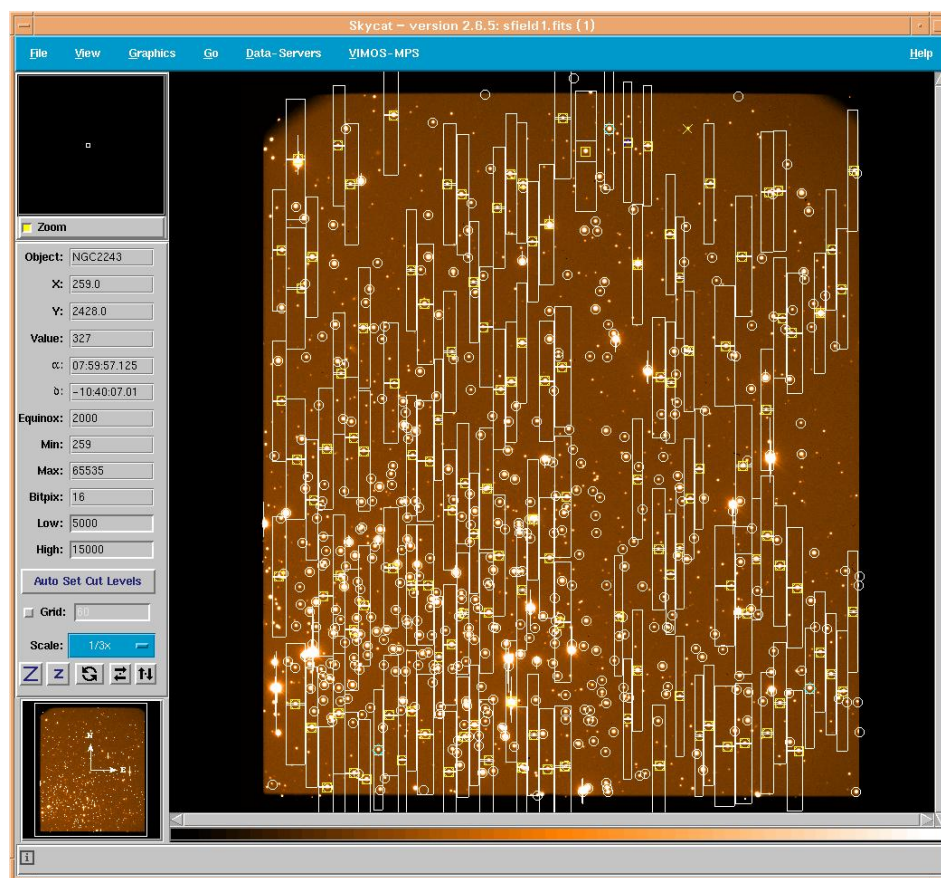


Figure 3.3: VMMPs Graphical User Interface: Example of ADP file superimposed to field of view

Table 3.2: VIMOS Calibration Plan: MOS Mode

| Calibration             | Number      | Frequency                                 | When?      | Products          | Accuracy <sup>1</sup> | Comment                      |
|-------------------------|-------------|-------------------------------------------|------------|-------------------|-----------------------|------------------------------|
| Bias                    | 5           | Daily                                     | Day time   | Master Bias       | RON                   |                              |
| Darks                   | 3           | Monthly                                   | Day time   | Dark current      |                       |                              |
| Screen Flats            | 5           | Daily                                     | Day time   | MOS Flat          | 5%                    | All used masks               |
| Screen Arcs             | 1           | Daily                                     | Day time   | Dispersion coeff. | 0.3 pix               | All used masks               |
| Screen Arcs             | 1           | Attached to science OB with HR or LR blue |            | Dispersion coeff. | 0.3 pix               | All used masks               |
| Specphot. Standard Star | 4 quadrants | Daily                                     | Night time | Response function | 10%                   | All grism+filter setups used |

<sup>1</sup> Internal accuracy, not considering instrument flexures (see Sect. 2.7).

**IMPORTANT NOTE** The execution time for all attached night time calibrations (screenflats and arcs) is charged to the user.

is guaranteed for mask preparation on the mountain, up to 3 masks per night, to be prepared at least 48 hrs in advance. The instrument has a capacity of 15 masks (per quadrant) loaded simultaneously (one slot is normally occupied by a mask dedicated to observations of spectrophotometric standard stars). However, some restrictions may apply for technical reasons in the number of masks which can be hosted simultaneously into the cabinets. Until P77, only 8 masks can be simultaneously accommodated in the instrument cabinet (1 of which is the mandatory spectroscopic mask), and the same is currently foreseen for P78. *Please refer to the Call for Proposals and to the instrument WEB pages for the updated status of the instrument situation.* Note that masks cannot be exchanged at night.

### 3.5.6 MOS mode: grism+filter setup

In MOS mode only one filter+grism combination per OB is permitted (the only filter exchange allowed is the one between acquisition and science template). Users who want to observe the same targets with different filter+grism combinations are requested to submit separate OBs for different filter+grism combinations and to consider the respective overheads.

### 3.5.7 MOS mode: attached night time calibrations

To correct for instrument flexures, the use of attached nighttime arcs (see Sect. 3.5.8), is MANDATORY for all MOS OBs using LR blue and HR blue grisms, and strongly recommended for all setups. Users may also attach screenflats. The execution time for these nighttime calibrations is always subtracted from the total allocated time.

### 3.5.8 MOS mode: calibration plan

Table 3.2 indicates the nature, frequency and accuracy of VIMOS calibrations expected to be taken as part of the calibration plan. Only calibration frames required for the reduction of science data are reported. More calibrations will be taken by the operation staff for the purpose of maintaining and monitoring the instrument configuration and performance, e.g. mask to CCD transformations etc. These calibrations are not reported here. Calibration monitoring data can be found on the pages of the Quality Control group:

<http://www.eso.org/observing/dfo/quality/>

### 3.5.9 MOS mode: pipeline

See <http://www.eso.org/qc/pipeline-status.html> for the information concerning the latest status of the ESO pipeline.

### 3.5.10 MOS mode: Exposure Time Calculator

The approximate limiting magnitudes in MOS mode can be obtained from the exposure time calculator available on the ESO WEB pages <http://www.eso.org/observing/etc/>.

## 3.6 Integral Field Unit (IFU)

### 3.6.1 Acquisition accuracy

In IFU mode acquisition is “blind”, i.e., limited by the accuracy of the telescope pointing ( $\sim 3''$  accuracy).

### 3.6.2 Dithering/jitter pattern

With the large field of view of the VIMOS-IFU, techniques equivalent to “shift-and-add” classically used in imaging, are a possible mode of operations. Dithering/jitter pattern is a list of differential pointing positions around a central position. Several short images (with exposure times set to have sky background limited exposures) are taken at each position of a dithering pattern, which are combined to reconstruct an image with the total required exposure time. This allows to eliminate the sky contribution, cosmic rays, the CCD fringing (particularly strong in the red), cosmetic defects on the detectors, and to produce very accurate flat fielding for correction of pixel to pixel variations. The dithering pattern depends on the mean size of the observed objects: it should be larger than the mean size of the objects for the method to work efficiently.

### 3.6.3 IFU mode: grism + filter setup

In IFU mode, only one filter+grism combination per OB is permitted. Users who want to observe the same targets with different filter+grism combinations are requested to submit separate OBs for different filter+grism combinations and to consider the respective overheads.

### 3.6.4 IFU mode: attached night time calibrations

To better take into account the effects of instrument flexure the use of attached nighttime calibrations, arcs and flats (see Sect. 3.6.5), is MANDATORY for all IFU OBs. This template will provide 1 or 3 flat-field images (for HR grism: 3 screenflats are mandatory, for the other cases, user can choose the number) as well as 1 arc, taken at the same rotator position of the science target. The execution time for these nighttime calibrations is always subtracted from the total allocated time.

### 3.6.5 IFU mode: calibration plan

Table 3.3 indicates the nature, frequency and accuracy of VIMOS calibrations expected to be taken as part of the calibration plan. Only calibration frames required for the reduction of science data are reported. More calibrations will be taken by the operation staff for the purpose of maintaining and monitoring the instrument configuration and performance, e.g. mask to CCD transformations etc. These calibrations are not reported here. Calibration monitoring data can be found on the pages of the Quality Control group:

<http://www.eso.org/observing/dfo/quality/>

### 3.6.6 IFU mode: pipeline

See <http://www.eso.org/qc/pipeline-status.html> for the information concerning the latest status of the ESO pipeline.



Table 3.3: VIMOS Calibration Plan: IFU Mode

| Calibration                | Number              | Frequency                   | When?      | Products              | Accuracy <sup>1</sup> | Comment                      |
|----------------------------|---------------------|-----------------------------|------------|-----------------------|-----------------------|------------------------------|
| Bias                       | 5                   | Daily                       | Day time   | Master Bias           | RON                   |                              |
| Darks                      | 3                   | Monthly                     | Day time   | Dark current          |                       |                              |
| Screen Flats               | 3                   | Daily                       | Day time   | IFU Flat in HR Orange |                       | Monitoring                   |
| Screen Arc                 | 1                   | Daily                       | Day time   | IFU Arc in HR Orange  |                       | Monitoring                   |
| Screen Flats               | 1 to 3 <sup>2</sup> | Attached to all science OBs |            | IFU Flat              | 5 %                   |                              |
| Screen Arcs                | 1                   | Attached to all science OBs |            | Dispersion coeff.     | 0.3 pix               |                              |
| Specphot.<br>Standard Star | 4 quadrants         | Daily                       | Night time | Response function     | 10%                   | All grism+filter setups used |

<sup>1</sup> Internal accuracy, not considering instrument flexures (see Sect. 2.7).

<sup>2</sup> With grism HR blue, 3 screen flats are always mandatory.

**IMPORTANT NOTE** The execution time for all attached night time calibrations (screenflats and arcs) is charged to the user.

### 3.6.7 IFU mode: Exposure Time Calculator

The approximate limiting magnitudes in IFU mode can be obtained from the exposure time calculator available on the ESO WEB pages <http://www.eso.org/observing/etc/>.

## 3.7 VIMOS overheads

The current overall efficiency of VIMOS is around 65-70 % depending which mode is used. The way in which some functions are operated may change at any time. Consequently, we cannot give here an estimate of the overheads that is reliable on the long term. The user is referred to the current Call for Proposals and to the VIMOS web pages for the most updated estimate of the overheads:

<http://www.eso.org/instruments/vimos/overview.html>



# Chapter 4

## VIMOS templates

### 4.1 Introduction

This section describes the VIMOS templates. Templates are to be used during phase II for service and visitor programmes.

Templates are characterized by signatures to be filled in by the users with P2PP. Templates are grouped together to build an Observation Block (OB), which typically consists of one acquisition template and one or more observing or calibration templates. At the telescope, OBs are passed to the instrument and executed according to the user defined parameter values.

Unless specified otherwise, the information provided in this part applies to both service and visitor modes.

### 4.2 Template names

VIMOS templates are divided in groups according to the instrument mode and to the observations to be performed. This is reflected in the name of the various templates.

The template naming scheme is the following:

$$\text{VIMOS\_}<mode>\_<type>\_<description>.\text{tsf}$$

where

|                    |                                                                                                                                                                                  |
|--------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <i>mode</i>        | is the mode of the instrument (either : <i>img</i> , <i>mos</i> , <i>ifu</i> )                                                                                                   |
| <i>type</i>        | is the type of the template (either: <i>acq</i> , <i>obs</i> or <i>cal</i> )                                                                                                     |
| <i>description</i> | is a string identifying the purpose of the template<br>( <i>Offset</i> for an observation, <i>Preset</i> for a preset,<br><i>NightCal</i> for attached night calibrations, etc.) |

### 4.3 Things to know - Reminder

- See Sect. 2.1.2 for a description of the orientation on sky of VIMOS.
- The default orientation on sky for MOS and pre-imaging templates is 90° (see Sect. 3.2). Any departure from this value should be requested with a waiver.
- Offset convention. The templates make extensive use of telescope offsets, which are entered manually as lists. The convention is that offsets are relative. In some templates, the offsets can be defined in detector

coordinates (e.g. X, Y), or in RA / DEC. All offsets are to be defined in arcseconds. The offsets refer to the target. Examples are provided further down in the description of the templates.

- Only some combinations of filters and grisms are allowed. See Sect. ?? and Table 2.3.

## 4.4 List of templates

Table 4.1 lists the VIMOS templates.

Table 4.1: Science and Nighttime Calibration Templates

| Type                  | Name                     |
|-----------------------|--------------------------|
| Acquisition templates | VIMOS_img_acq_Preset     |
|                       | VIMOS_mos_acq_Mask       |
|                       | VIMOS_mos_acq_Standard   |
|                       | VIMOS_ifu_acq_Preset     |
| Science templates     | VIMOS_img_obs_Offset     |
|                       | VIMOS_mos_obs_Offset     |
|                       | VIMOS_ifu_obs_Offset     |
| Calibration templates | VIMOS_img_cal_Photom     |
|                       | VIMOS_img_cal_Astrom     |
|                       | VIMOS_mos_cal_NightCalib |
|                       | VIMOS_mos_cal_Standard   |
|                       | VIMOS_ifu_cal_NightCalib |
|                       | VIMOS_ifu_cal_Specphot   |

## 4.5 Acquisition templates

### 4.5.1 VIMOS\_img\_acq\_Preset

This template is used for acquisition of a field in imaging mode, for both imaging and pre-imaging OBs.

Figure 4.1 shows the P2PP window with an example of a pre-imaging Observing Block (OB).

The coordinates of the target have to be entered in the *Target* window at the bottom of the P2PP window.

The parameters for this template are the following:

- **Filter.** This filter should be the same as for the subsequent imaging template. The available filters are: U, B, V, R, I, z.
- **Rotator on sky.** This parameter defines the orientation of the instrument on sky. **This parameter should always be set to '90' for pre-imaging templates, unless a waiver is requested.** The orientation corresponding to 90 degrees is such that the slits will be oriented N-S for the spectroscopic follow-up.
- **Alpha offset** and **Delta offset.** These are optional additional offsets in alpha and delta added to the coordinates of the object specified in the target panel at the bottom of the main p2pp window. The coordinates at the center of the instrument (i.e. the central point between the four quadrants) are the coordinates specified in the Target Package plus the offsets here specified.
- **Get Guide Star From.** Two options are offered for this parameter:
  - **CATALOGUE.** This is the default option for which the guide star will be searched through a catalog and will be selected at the telescope by the Telescope Operator. If this option is selected the subsequent fields (**Guide Star RA / DEC**) will be ignored.

ObsBlock: Target\_preimaging: VIMOS

File Edit Synthesise FindingCharts

Name: target\_preimaging

Status: partially defined

\* Exposition Time: 00:00:00.000

User Priority: 1

OD Name: target\_preimaging

User Comments:

| VIMOS_img_offset        | 1         | VIMOS_img_offset                       | 1          |
|-------------------------|-----------|----------------------------------------|------------|
| Alpha offset            | 0         | Exposure time (seconds)                | 0          |
| Delta offset            | 0         | Observation Category                   | PHOTOMETRY |
| Rotator velocity        | 0         | Number of Exposures per Telescope Off. | 0          |
| Additional Velocity RA  | 0         | Number of Telescope Offsets            | 0          |
| Additional Velocity DEC | 0         | Offset Coordinates                     | 0          |
| Get Guide Star From     | SETUPFILE | List of offsets (arcsec) in X or RA    | 0          |
| Guide Star RA           | 0         | List of offsets (arcsec) in Y or DEC   | 0          |
| Guide Star DEC          | 0         | Return to Original (Y/N)               | Y          |
| Filter                  | R         | Filter                                 | R          |

Target Constraint Set Time intervals Scheduling Time Calibration Requirements

Name: target

Class: Irregular

Right Ascension: 20:00:00.000

Declination: -20:00:00.000

Equinox: 2000

Epoch: 2000.0

Proper motion RA: 0.0

Proper motion DEC: 0.0

Diff RA: 0.0

Diff DEC: 0.0

Figure 4.1: P2PP: example of a pre-imaging OB

- **SETUPFILE**. If this option is selected, the guide star that will be used is the one for which the coordinates are defined in the (Guide Star RA / DEC) fields. A dedicated tool (Guidecam) is delivered by ESO to help the users to determine the best possible guide star. This option is mandatory for service mode imaging and pre-imaging OBs.
- **Guide Star RA and Guide Star DEC**. Coordinates of the guide star when the **SETUPFILE** option of the (Get Guide Star From) field is selected.
- **Guide Star PPOS**. Position of the Guide Probe when the **SETUPFILE** option of the (Get Guide Star From) field is selected.
- **Optional additional velocity (moving targets)** in alpha and delta (**Additional Velocity RA** and **Additional Velocity DEC**). This velocity should be expressed in arcsec/sec. These parameters are mostly intended for imaging or possibly IFU programmes of moving targets. If no additional velocity is needed, these parameters should be left to their default values (0.).

| VIMOS_mos_acq_Mask      |                          | VIMOS_mos_obs_Offset                    |                                          |
|-------------------------|--------------------------|-----------------------------------------|------------------------------------------|
| Exposure time (seconds) | 120                      | Exposure time (seconds)                 | 1800                                     |
| Filter                  | R                        | Filter                                  | OS-blue                                  |
| Catalogue Star from     | CATALOGUE                | Grism                                   | LR red                                   |
| Catalogue Star RA       | 0                        | Return to Origin P (T/F)                | <input checked="" type="checkbox"/> TRUE |
| Catalogue Star DEC      | 0                        | Number of Exposures P                   | 2                                        |
| Additional Velocity RA  | 0                        | Number of Telescope Offsets P           | 2                                        |
| Additional Velocity DEC | 0                        | List of offsets (arcsec) along the slit | 0 1.5                                    |
| ADP file 1              | fermax_user1_vmic_M1.adp | List of offsets (arcsec) perpendicular  | 0                                        |
| ADP file 2              | fermax_user2_vmic_M1.adp |                                         |                                          |
| ADP file 3              | fermax_user3_vmic_M1.adp |                                         |                                          |
| ADP file 4              | fermax_user4_vmic_M1.adp |                                         |                                          |

| Target          |               | Constraint Set    |         | Time Intervals |  | User Comments |  | Calibration Requirements |  |
|-----------------|---------------|-------------------|---------|----------------|--|---------------|--|--------------------------|--|
| Name            | Fulcrum       | Class             | Unknown |                |  |               |  |                          |  |
| Right Ascension | 05:30:55.541  | proper motion RA  | 0       |                |  |               |  |                          |  |
| Declination     | -38:48:12.056 | proper motion DEC | 0       |                |  |               |  |                          |  |
| Equinox         | J2000         | CI RA             | 0       |                |  |               |  |                          |  |
| Epoch           | J2000.0       | CI DEC            | 0       |                |  |               |  |                          |  |
| User Comments   |               |                   |         |                |  |               |  |                          |  |

Figure 4.2: P2PP: Example of a MOS OB

#### 4.5.2 VIMOS\_mos\_acq\_Mask

This template is used for the acquisition of a field to be observed in MOS mode. Note that the acquisition image covers only the central 2148x2444 pixels of each of the 4 CCDs whereas the spectroscopic images use the full CCD size (2148x4096 pixels).

Preparing a MOS acquisition template requires that masks have been prepared with VMMPs. The output of VMMPs are ‘Aperture Definition Files in Pixels’ (ADP files). Information concerning the target (coordinates and equinox) are taken from the ADP files and are automatically loaded in the P2PP target package when including the first ADP file in the template.

Figure 4.2 shows the P2PP window with an example of a MOS Acquisition.

The parameters specific to this template are:

- **Exposure time in (seconds).** This is the integration time for the acquisition image that will be taken through the mask to check the position of the reference stars and of the science target. A typical value is 60s, depending on the brightness of the reference stars defined with VMMPs.
- **Filter.** This filter is the one to be used for the acquisition image. **It is mandatory to choose the same filter as the one used for the pre-imaging.** Waiver is requested differently but see **Sect. 3.4.1** of this manual. The available filters are: U, B, V, R, I, z, OS-blue, OS-red, GG435, GG475 and no filter.
- **Optional additional velocity (moving targets) in alpha and delta (Additional Velocity RA and Additional Velocity DEC).** This velocity should be expressed in arcsec/sec. These parameters are mostly intended for imaging and possibly IFU programmes of moving targets and a priori unlikely to be of much use in MOS mode. If no additional velocity is needed, these parameters should be left to their default values.
- **ADP files 1 to ADP file 4.** These 4 parameters are the names of the ADP files produced by VMMPs. The files 1 to 4 correspond to the 4 quadrants of the instrument and should be set in the right order.



### 4.5.3 VIMOS\_mos\_acq\_Standard

The parameters of the template are:

- **Exposure time (seconds).** This is the integration time for the acquisition image that will be taken through the mask to check its positioning.
- **Filter.** This filter is the one to be used for the acquisition image. **It is recommended to choose**

the right combination filter/exptime to avoid the saturation of the star in the acquisition image

The available filters are: **U, B, V, R, I, z, OS-blue, OS-red, GG435, GG475** and no filter.

- Optional additional velocity (moving targets) in alpha and delta (**Additional Velocity RA** and **Additional Velocity DEC**). This velocity should be expressed in arcsec/sec. These parameters are mostly intended for imaging and possibly IFU programs of moving targets and a priori unlikely to be of much use in MOS mode. If no additional velocity is needed, these parameters should be left to their default values.

#### 4.5.4 VIMOS\_ifu\_acq\_Preset

This template is for the acquisition of a field to be observed with the Integral Field Unit. Unlike the MOS acquisition template, the IFU acquisition template requires to define the filter / grism combination that will be used in the following observation template. This is to optimize the setup of the instrument at acquisition time.

The coordinates of the field at the center of the IFU need to be entered in the Target Package at the bottom left of the main P2PP window.

The parameters for this template are:

- **Filter.** This should be the same filter as the one used in the first observation template following this acquisition template. **In order to produce a valid OB, the filter in the first observation template following this acquisition template should be the same as the one selected here.** The available filters are: **U, B, V, R, I, z, OS-blue, OS-red, GG435, GG475** and no filter.

**Only standard filter / grism combinations are allowed, as indicated in Table 2.3. Other combinations should be requested with a waiver.**

- **Grism.** This is the grism that is used in the first observation template following this acquisition template. **In order to produce a valid OB, the grism in the first observation template following this acquisition template should be the same as the one selected here.**

The available grisms are: **LR\_blue, LR\_red, MR, HR\_blue, HR\_orange, HR\_red.**

- **IFU Magnification.** This parameter defines the sampling, and correspondingly the field of view of the IFU. The 2 available samplings are 0.67 arcsec/fiber and 0.33 arcsec/fiber. **To produce a valid OB, the magnification should be the same in every template of the OB.** Observations requiring different samplings on the same field should resort to several OBs.
- **IFU Shutter.** This parameter is only meaningful when using the **LR** grisms. In **LR** spectroscopy, there are 4 rows of fibers along the dispersion direction, whereas in **MR** and **HR** modes only one 'slit' of fibers is used, hence the variation of field of view between **LR** and **MR** and **HR** modes. In case the object of interest is of limited size in **LR** mode, and if only the central  $1/4 \times 1/4$  of the field of view is needed, it is possible with the IFU Shutter option to mask the outer part of the IFU (always masked in **MR** and **HR** modes). This allows to avoid second order overlap due to sky and / or background objects between the 4 adjacent slits along the dispersion direction. **True** means that the shutter is in place and masks the external part of the IFU, **False** (the default value) means that the mask is not in place and the full IFU field of view is available. **To produce a valid OB, the mask shutter option should be the same in every template of the OB.** This parameter is ignored when **MR** or **HR** grisms are used.
- **Rotator on Sky.** This parameter defines the orientation of the instrument on sky. The IFU has a square geometry with the sides oriented N-S and E-W. With a rotator angle of zero the orientation of the IFU on the sky is with the North up and the East on the right. A rotation of +90 degrees will put East up and North left (see Fig. 2.2).
- Optional additional offsets in alpha and delta (**Alpha Offset** and **Delta offset**) with respect to the coordinates of the object. The coordinates at the center of the IFU are the coordinates specified in the Target Package plus the offsets here specified. This can be useful e.g. to define small offsets between the initial pointings on the same field over several OBs.
- **Get Guide Star From.** Two options are offered for this parameter:

- **CATALOGUE**. This is the default option for which the guide star will be searched through a catalog and will be selected at the telescope by the Telescope Operator. If this option is selected the subsequent fields (**Guide Star RA / DEC**) will be ignored. Given that the quasi-total Nasmyth field of view is accessible to choose a guide star in IFU mode, the users can safely use this option.
- **SETUPFILE**. If this option is selected, the guide star that will be used is the one for which the coordinates are defined in the (**Guide Star RA / DEC**) fields. As indicated above, it is unlikely that this option is required in IFU mode. Note that Guidecam tool is not intended to be used to select guide stars for the IFU mode, but only for the IMG and MOS mode of VIMOS.
- **Guide Star RA** and **Guide Star DEC**. Coordinates of the guide star when the **SETUPFILE** option of the **Get Guide Star From** field is selected.
- Optional additional velocity (moving targets) in alpha and delta (**Additional Velocity RA** and **Additional Velocity DEC**) This velocity should be expressed in arcsec/sec. If no additional velocity is needed, these parameters should be left to their default values.

## 4.6 Science templates

### 4.6.1 VIMOS\_img\_obs\_Offset

This template is used for Imaging mode observation both for **Science** and **Pre-Imaging**. The parameters of this template are the following:

- **Exposure time in (seconds)**. This is the integration time in seconds per exposure.
- **Filter**. The available filters are: **U, B, V, R, I, z**. **In order to produce a valid OB the filter in the first observation template following the acquisition template should be the same as the one selected in the acquisition template. For preimaging the use of R filter is mandatory.**
- **Return to Origin ? (T/F)**. This parameter defines if after a dithered set of observations the telescope will return at the position at beginning of template or not. As an example, this parameter should be set to 'T' when the aim of the OB is to execute multi-color imaging observations with the same pattern of offsets on the same field using different templates. This parameter is in general to be left to its default value in service mode. **If another value is desired in service mode a waiver should be requested.**
- **Number of Telescope Offsets?**. This parameter defines the number of telescope offsets during the template. The number of exposures taken at each telescope position is determined by the next parameter (**Number of exposures per telescope offset ?**). The pattern of offsets is defined in the parameters (**List of offsets (arcsec) in X or RA**) and (**List of offsets (arcsec) in Y or DEC**). Offsets are defined relative to the previous telescope position.
- **Number of exposures per telescope offset ?**. This parameter defines how many exposures are taken at each offset position.
- **Offset coordinates?**. This parameter defines the coordinate system used to define telescope offsets. If the **SKY** option is selected, telescope offsets are in RA and DEC, if **DETECTOR** is selected, telescope offsets are in X and Y on the detector. In both cases the offsets are given in **arcsec**.
- **List of offsets (arcsec) in X or RA** and **List of offsets (arcsec) in Y or DEC**. List of offsets between exposures. Offsets are defined relative to the previous telescope position. The offsets are defined in RA/DEC or X/Y on detector depending on which option is selected for the (**Offset coordinates**) field. If the first image is to be taken at the preset coordinates, the first set of offsets should be 0,0. If there are less offset values than the defined number of offset positions, the template will return at the beginning of the list of offsets and apply them again.

For example, if 5 telescopes offsets are defined (number of telescope offsets?) together with the following series of offsets:

|          |   |     |    |
|----------|---|-----|----|
| RA or X  | 0 | -10 | -5 |
| DEC or Y | 2 | 7   |    |

the offsets applied to the telescope will be:

|          |    |    |    |    |    |
|----------|----|----|----|----|----|
| RA or X  | 0  | 10 | 5  | 0  | 10 |
| DEC or Y | -2 | -7 | -2 | -7 | -2 |

(i.e., if SKY option for coordinates is used: after the preset the first offset will be applied to place the target at 2 arcsec North (0 offset in RA) and then the first exposure will be taken; after that the target will be moved further 10 arcsec West and 7 arcsec North to the position of the first exposure, etc. until the number of exposures is reached. If DETECTOR option is used the first exposure will be taken at pixels  $x = 0$   $y = 10$  (offset is 2", which is approximately 10 pixels), then the target will be moved further 10" West and 7" North to the position of the first exposure, etc., until the number of exposures is reached.)

- **Observation Category.** **PRE-IMAGE** is for a pre-imaging program, **SCIENCE** is for a regular science imaging program.

**It is critical that users define PRE-IMAGE for pre-imaging programs.** Failure to do so will result in slow processing and delivery (possibly till the end of the period) of the data.

#### 4.6.2 VIMOS\_mos\_obs\_Offset

The parameters of this template are:

- **Exposure time in (seconds).** This is the integration time in seconds per exposure.
- **Filter.** The available filters are: **U, B, V, R, I, z, OS-blue, OS-red, GG435, GG475** and no filter. The standard grism / filter configurations are given in Table 2.3.  
**Only standard filter / grism combinations are allowed, as indicated in Table 2.3. Other combinations should be requested with a waiver.**
- **Grism.** The available grisms are: **LR.blue, LR.red, MR, HR.blue, HR.orange, HR.red.**
- **Return to Origin ? (T/F).** This parameter defines if after a dithered set of observations the telescope will return at the position at beginning of template or not. **This parameter is in general to be left in the default value in service mode. If another value is desired in service mode a waiver should be requested.**
- **Number of Telescope Offsets?.** This parameter defines the number of telescope offsets during the template. The number of exposures taken at each telescope position is determined by the next parameter (**Number of exposures per telescope offset ?**). The pattern of offset is defined in the **List of offsets (arcsec) along the slit** and **List of offsets (arcsec) perpendicular to the slit** parameters.
- **Number of exposures per telescope offset ?.** This parameter defines how many exposures are taken at each offset position.
- **List of offsets (arcsec) along the slit and List of offsets (arcsec) perpendicular to the slit.** Offsets are defined relative to the previous telescope position. Offsets are defined in detector coordinates, along or perpendicular to the slits. Usually, only offsets along the slits need to be defined. If the first image is to be taken at the preset coordinates, the first set of offset should be 0,0. If there are less offset values than the required number of offset position, the template will return at the beginning of the list of offsets and apply them again.

For example, if 3 telescope offsets are defined (**Number of Telescope Offsets?**) together with the following series of offsets:



|                                                    |    |   |   |
|----------------------------------------------------|----|---|---|
| List of offsets (arcsec) along the slit            | -1 | 2 | 0 |
| List of offsets (arcsec) perpendicular to the slit | 0  |   |   |

the offsets applied to the telescope will be:

|                                            |   |    |   |
|--------------------------------------------|---|----|---|
| Offsets (arcsec) along the slit            | 1 | -2 | 0 |
| Offsets (arcsec) perpendicular to the slit | 0 | 0  | 0 |

(i.e., the first exposure is taken at  $-1''$  from the pointing position along the slit, then the second exposure is taken  $2''$  up along the slit with respect to the second exposure, and the third exposure is taken at the same position as the second.)

#### 4.6.3 VIMOS ifu\_obs\_Offset

This template is to be used for observation with the Integral Field Unit. The parameters of this template are:

- **Exposure time in (seconds).** This is the integration time in seconds per exposure.
- **Filter.** The available filters are: **U, B, V, R, I, z, OS-blue, OS-red, GG435, GG475** and no filter. **Only standard filter / grism combinations are allowed, as indicated in Table 2.3. Other combinations should be requested with a waiver.**
- **Grism.** The available grisms are: **LR\_blue, LR\_red, MR, HR\_blue, HR\_orange, HR\_red.**

**In order to produce a valid OB the filter and the grism in the first observation template following the acquisition template should be the same as the one selected in the acquisition template.**

- **IFU Magnification.** This parameter defines the sampling, and correspondingly the field of view of the IFU. The 2 available samplings are 0.67 arcsec/fiber and 0.33 arcsec/fiber.

**To produce a valid OB, the magnification should be the same in every template of the OB.**

- **IFU Shutter.** This parameter is only meaningful when using the **LR** grisms. In **LR** spectroscopy, there are 4 rows of fibers along the dispersion direction, whereas in **MR** and **HR** modes only one 'slit' of fibers is used, hence the variation of field of view between **LR** and **MR** and **HR** modes. In case the object of interest is of limited size in **LR** mode, and if only the central  $1/4 \times 1/4$  of the field of view is needed, it is possible with the **IFU Shutter** option to mask the outer part of the IFU (always masked in **MR** and **HR** modes). This allows to avoid second order overlap due to sky and / or background objects between the 4 adjacent slits along the dispersion direction. **True** means that the shutter is in place and masks the external part of the IFU, **False** (the default value) means that the mask is not in place and the full IFU field of view is available. **To produce a valid OB, the maks shutter option should be the same in every template of the OB.** This parameter is ignored when **MR** or **HR** grisms are used.
- **Return to Origin ? (T/F).** This parameter defines if after a dithered set of observations the telescope will return at the reference position or not. This parameter is in general to be left in the default value (T) in service mode. **If another value is desired in service mode a waiver should be requested.**
- **Number of Telescope Offsets?.** This parameter defines the number of telescope offsets during the template. The number of exposures taken at each telescope position is determined by the next parameter (Number of exposures per telescope offset ?). The pattern of offset is defined in the List of offsets (arcsec) CORRECT LABEL and List of offsets (arcsec) CORRECT LABEL parameters.
- **Number of exposures per telescope offset ?.** This parameter defines how many exposures are taken for each offset position.

- **List of offsets (arcsec)** List of offsets between exposures. Offsets are defined relative to the previous telescope position. The offsets are defined in RA/DEC.

If the first image is to be taken at the preset coordinates, the first set of offset should be 0,0. If there are less offset values than the defined number of offset positions, the template will return at the beginning of the list of offsets and apply them again.

For example, if 5 telescope offsets are defined (**Number of Telescope Offsets?**) together with the following series of offsets (in RA/DEC):

|     |   |     |    |
|-----|---|-----|----|
| RA  | 0 | -10 | -5 |
| DEC | 0 | 5   | 10 |

the offsets applied to the telescope will be:

|     |   |    |     |   |    |
|-----|---|----|-----|---|----|
| RA  | 0 | 10 | 5   | 0 | 10 |
| DEC | 0 | -5 | -10 | 0 | -5 |

## 4.7 Calibration Templates

### 4.7.1 Nighttime Imaging Calibrations

#### 4.7.1.1 VIMOS\_img\_cal\_Photom

This template is for acquisition of photometric data (zero points) in Imaging mode. The acquisition template to use with this template is VIMOS\_img\_acq\_Preset. The coordinates to define in the acquisition template are those of the standard star / field. Offsets from these coordinates to position the star / field in one or the four quadrants of the instrument will be applied automatically by the template.

The parameters of this template are:

- **Exposure time (seconds)**. This is the integration time in seconds per exposure.
- **Filter**. The available filters are: **U, B, V, R, I, z**. **To produce a valid calibration OB the filter in the first observation template following the acquisition template should be the same as the one in the acquisition template.**
- **Number of exposures per telescope offset ?**. This parameter defines how many exposures are taken at each offset position (1 or 4).
- **1 or 4 quadrants**. This parameter allows to define if the standard star image will be taken in 1 or 4 quadrants. If the **1 quadrant** option is selected, the standard star field will be positioned in only one quadrant by applying an appropriate offset from the coordinates of the star / field defined in the target package of the acquisition template.

If the **4 quadrant** option is selected, the standard star field will be positioned sequentially in the 4 quadrants by applying the appropriate offsets from the coordinates of the star / field defined in the target package of the acquisition template.

As many exposures as defined in (**Number of exposures per telescope offset ?**) will be taken at each telescope position (1 or 4).

At the end of the template the telescope will offset again to return to its position at the beginning of the template.

#### 4.7.1.2 VIMOS\_img\_cal\_Astrom

This template is for acquisition of astrometric data in Imaging mode. The acquisition template to use with this template is VIMOS\_img\_acq\_Preset. The coordinates to define in the acquisition template are those of the astrometric field. Offsets from these coordinates to position the field in one or the four quadrants of the instrument will be applied automatically by the template.

The parameters of the template are:

- **Exposure time (seconds).** This is the integration time in seconds per exposure.
- **Filter.** The available filters are: **U, B, V, R, I, z**. **To produce a valid calibration OB the filter in the first observation template following the acquisition template should be the same as the one in the acquisition template.**
- **Number of exposures per telescope offset ?.** This parameter defines how many exposures are taken at each offset position (1 or 4).
- **1 or 4 quadrants.** This parameter allows to define if the astrometric field image will be taken in 1 or 4 quadrants. If the **1 quadrant** option is selected, the astrometric field will be positioned in only one quadrant by applying an appropriate offset from the coordinates of the field defined in the target package of the acquisition template.

If the **4 quadrant** option is selected, the astrometric field will be positioned sequentially in the 4 quadrants by applying the appropriate offsets from the coordinates of the field defined in the target package of the acquisition template.

As many exposures as defined in (Number of exposures telescopes offset ?) will be taken at each telescope position (1 or 4).

At the end of the template the telescope will offset again to return to its position at the beginning of the template.

### 4.7.2 Nighttime MOS Calibrations

#### 4.7.2.1 VIMOS\_mos\_cal\_NightCalib

This template is to be used when screen flats or wavelength calibrations are needed during the night right after a science observation. **This template doesn't contain any setup definition and the data will be taken with exactly the same instrument setup as used in the previous science template.**

This template is to be used in a science OB **after** a science templates. It can not be used alone in an OB after an acquisition template. The exposure time parameters are automatically defined by the current instrument setup. The parameters of the template are:

- **Night Flat Field ? (T/F).** If set to 'T', 3 screen flats will be taken.
- **Night Arc ? (T/F).** If 'T' an arc spectrum will be taken.

#### 4.7.2.2 VIMOS\_mos\_cal\_Standard

This template is used for the observation of a standard star in mode MOS. This template can only be used following the VIMOS\_mos\_acq\_Standard acquisition template and **not** using the VIMOS\_mos\_acq\_Mask. This template can only use the specifically designed standard mask from the repository.

The parameters of this template are:

- **Exposure time (seconds).** This is the integration time in seconds per exposure.
- **Filter.** The available filters are: **U, B, V, R, I, z, OS-blue, OS-red, GG435, GG 475.**
- **Grism.** The available grisms are: **LR\_blue, LR\_red, MR, HR\_blue, HR\_orange, HR\_red.**

- **Number of exposures per telescope offset ?**. This parameter defines how many exposures are taken for each offset position.
- **1 or 4 quadrants?**. This determines whether the standar star will be positioned in one or four of the quadrants of the instrument. As many exposures as defined in (**Number of exposures per telescope offset ?**) will be taken at each telescope position (1 or 4).
- **Slit Width**. Define which of the slits available in the masks will be used. The width of the available slits are: **0.8, 1.0, 1.2, 1.5, 1.8, 2, 2.5, 5.0**. Once the user select the slit width the rest is totally automatic; the star will be offsetted from the acquisition box to the selected slit in the first quadrant and then moved to the same slit in the other quadrants.

### 4.7.3 Nighttime IFU Calibrations

#### 4.7.3.1 VIMOS\_ifu\_cal\_Specphot

This template is used for the observation of a spectrophotometric standard star in IFU mode. It is to be used in a OB with the VIMOS\_ifu\_acq\_Preset acquisition template.

The parameters of the template are:

- **Exposure time (seconds)**. This is the integration time in seconds per exposure.
- **Filter**. The available filters are: **U, B, V, R, I, z, OS-blue, OS-red, GG435, GG 475**.
- **Grism**. The available grisms are: **LR\_blue, LR\_red, MR, HR\_blue, HR\_orange, HR\_red**.  
**In order to produce a valid OB the filter and the grism in the first observation template following the acquisition template should be the same as the one selected in the acquisition template.**
- **IFU Magnification**. This parameter defines the sampling, and correspondingly the field of view of the IFU. The 2 available samplings are 0.67 arcsec/fiber and 0.33 arcsec/fiber.  
**To produce a valid OB, the magnification should be the same in every template of the OB.**
- **IFU Shutter**. This parameter is only meaningful when using the **LR** grisms. In **LR** spectroscopy, there are 4 rows of fibers along the dispersion direction, whereas in **MR** and **HR** modes only one 'slit' of fibers is used, hence the variation of field of view between **LR** and **MR** and **HR** modes. In case the object of interest is of limited size in **LR** mode, and if only the central  $1/4 \times 1/4$  of the field of view is needed, it is possible with the **IFU Shutter** option to mask the outer part of the IFU (always masked in **MR** and **HR** modes). This allows to avoid second order overlap due to sky and / or background objects between the 4 adjacent slits along the dispersion direction. **True** means that the shutter is in place and masks the external part of the IFU, **False** (the default value) means that the mask is not in place and the full IFU field of view is available. **To produce a valid OB, the maks shutter option should be the same in every template of the OB.** This parameter is ignored when **MR** or **HR** grisms are used.
- **Number of exposures per telescope offset ?**. This parameter defines how many exposures are taken for each offset position.
- **1 or 4 quadrants?**. This determines whether the standar star / field will be positioned in one or four of the quadrants of the IFU. As many exposures as defined in (**Number of exposures per telescope offset ?**) will be taken at each telescope position (1 or 4).

#### 4.7.3.2 VIMOS\_ifu\_cal\_NightCalib.

This template is to be used when screen flats or wavelength calibrations are needed during the night right after a science observation. **This template doesn't contain any setup definition and the data will be taken with exactly the same instrument setup as used in the previous observation template.**

This template is to be used in a science OB **after** observation templates. It can not be used alone in an OB after an acquisition template. The exposure time parameters are automatically defined by the instrument setup. The parameters of the template are:

- Night Flat Field ? (T/F). If set to 'T', 3 screen flats will be taken.
- Night Arc ? (T/F). If 'T' an arc spectrum will be taken.



# Appendix A

## VIMOS Filters

Table A.1 lists the characteristics of the VIMOS filters. The transmission curves are given in the following figures.

Table A.1: Characteristics of the VIMOS filters

| Filter  | $\lambda_0$ (nm) | FWHM (nm)               |
|---------|------------------|-------------------------|
| U'      | 370              | 50                      |
| B_BESS  | 430              | 97.0                    |
| V_BESS  | 546              | 89.0                    |
| R_BESS  | 648.5            | 130.0                   |
| I       | 830              | 80.0                    |
| z_GUNN  | 950              | 160                     |
| GG475   | edge filter      | high pass ( $> 475$ nm) |
| GG435   | edge filter      | high pass ( $> 435$ nm) |
| OS-red  | box filter       | band pass 550-950       |
| OS-blue | box filter       | band pass 370-670       |

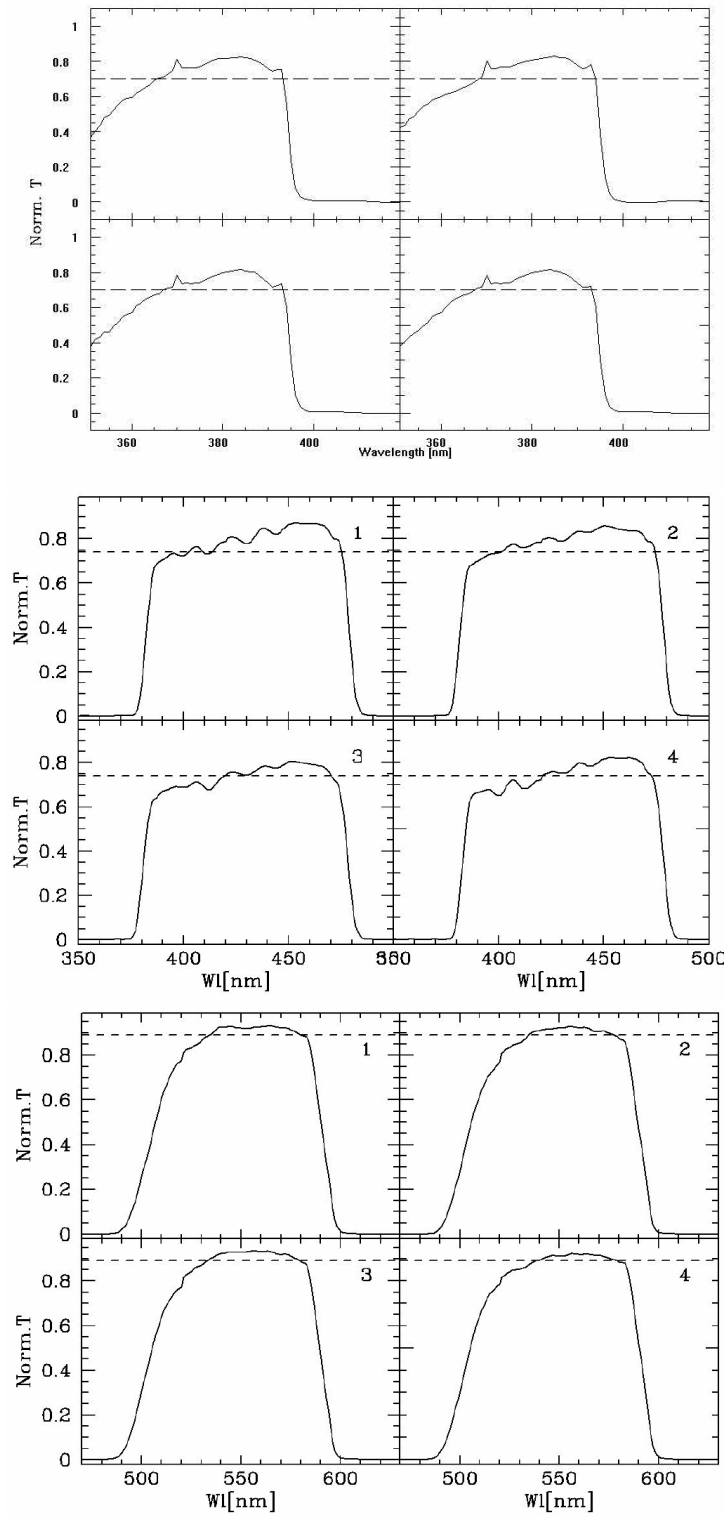


Figure A.1: U', B, V filter transmission curves



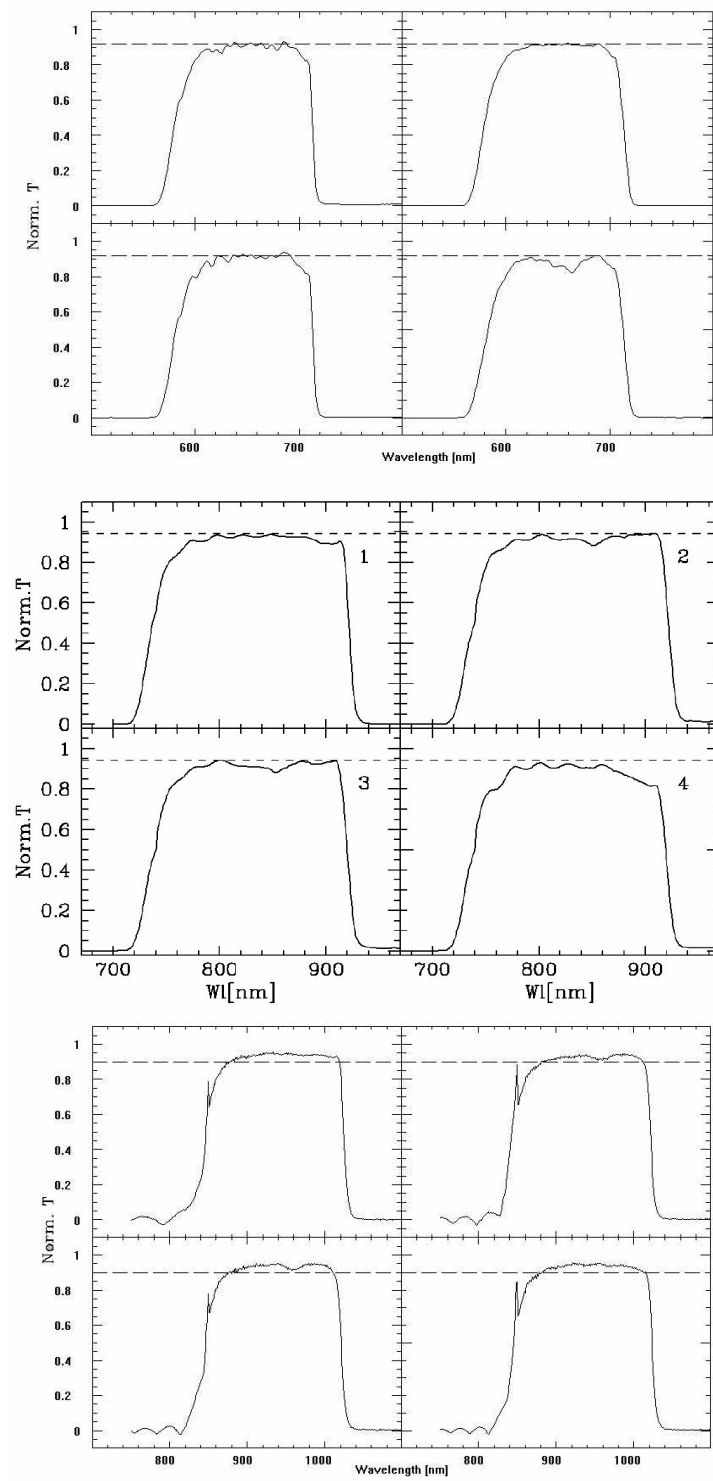


Figure A.2: R,I,z filter transmission curves

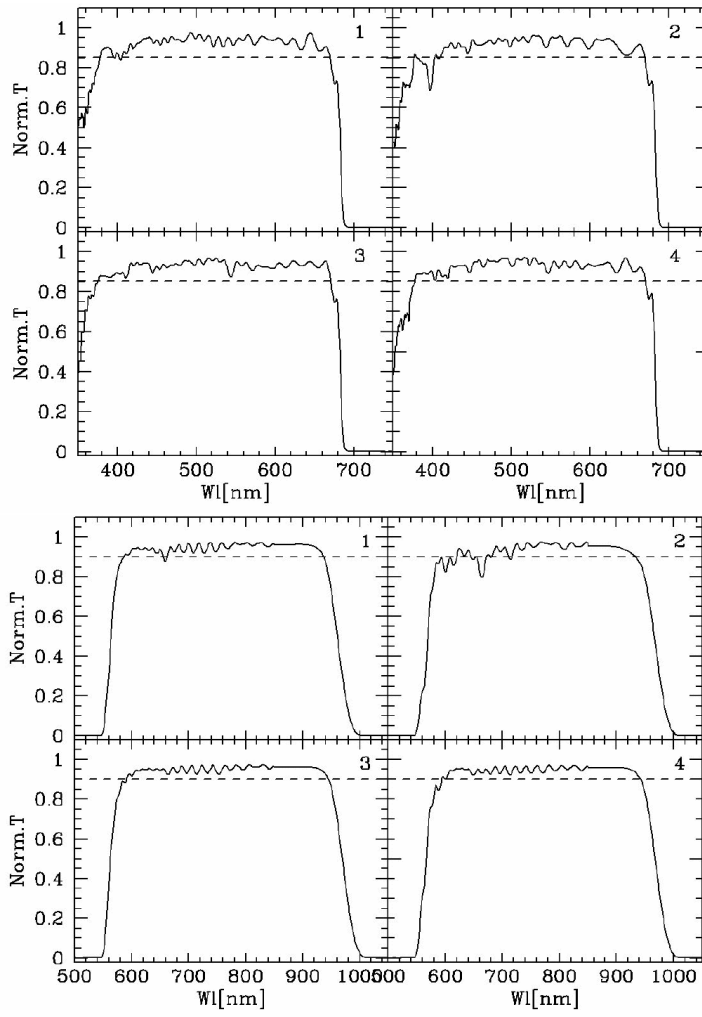


Figure A.3: OS-blue, OS-red transmission curves

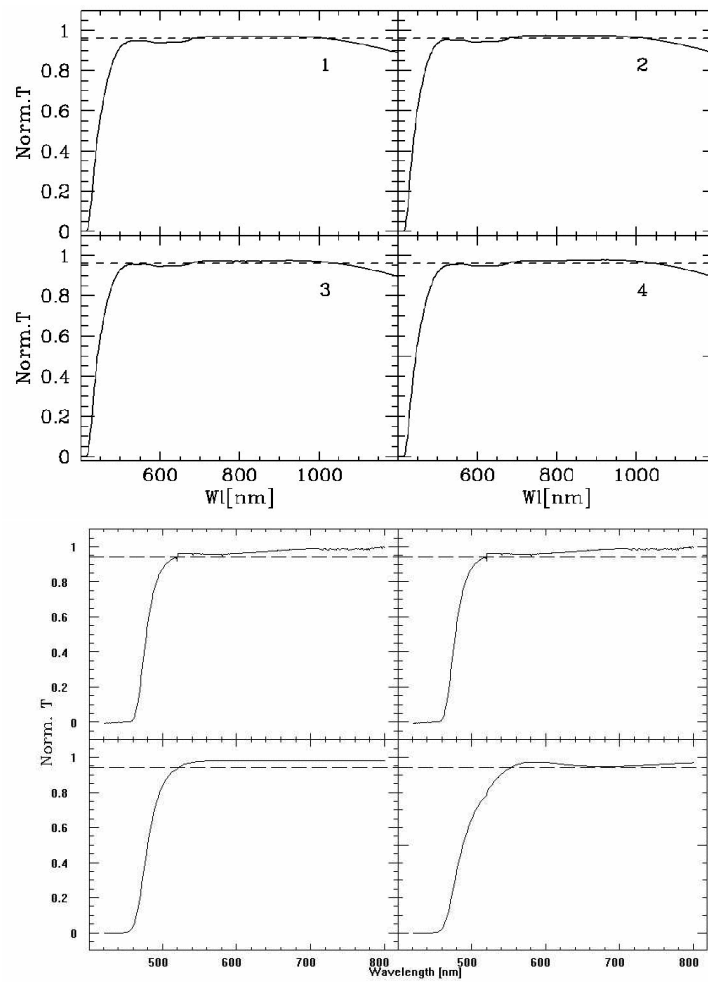


Figure A.4: GG435, GG475 transmission curves

## Appendix B

# VIMOS Grisms

This appendix contains the efficiency curves of all standard grisms available for VIMOS and the approximate wavelength range for a slit which is located in the field centre.

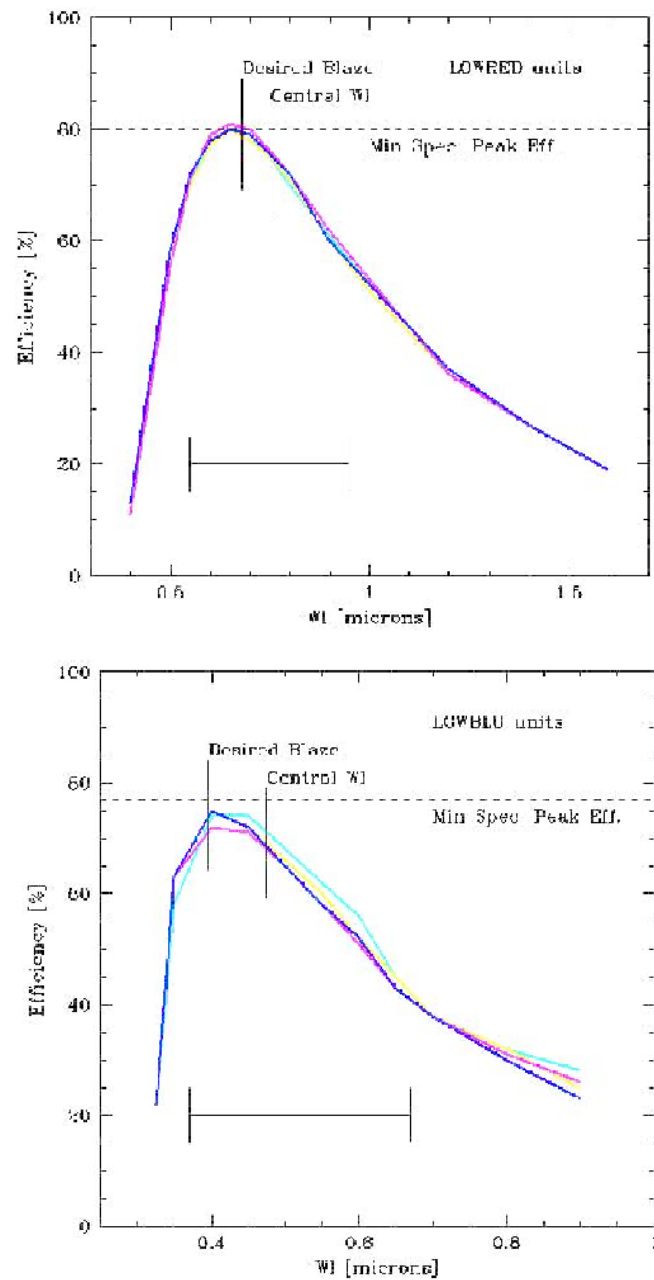


Figure B.1: Efficiency curve of the low resolution grisms LR-red, LR-blue. The vertical lines mark the desired blaze and the central wavelengths. The spectral range refers to the slit in the center of the field. The cutoff wavelength is in most cases given by the order separation filters, the red CCD limit or the limit of the VIMOS optics in the blue.

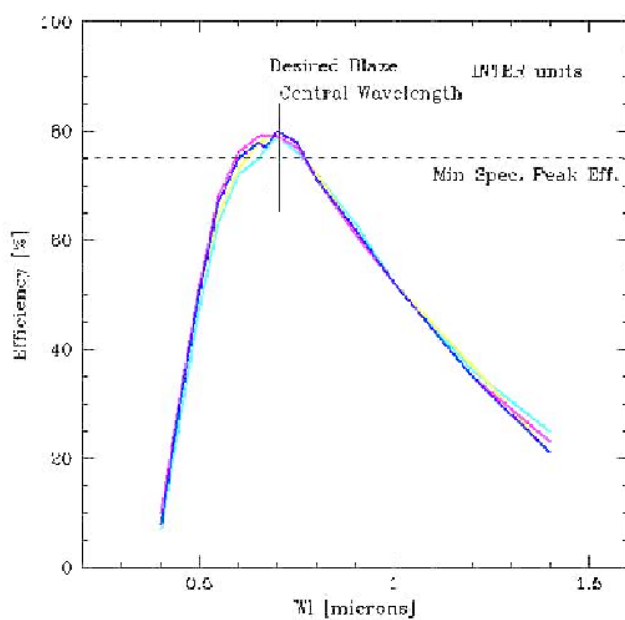


Figure B.2: Efficiency curve of the medium resolution grisms MR. The vertical lines mark the desired blaze and the central wavelengths. The spectral range refers to the slit in the center of the field. The cutoff wavelength is in most cases given by the order separation filters, the red CCD limit or the limit of the VIMOS optics in the blue

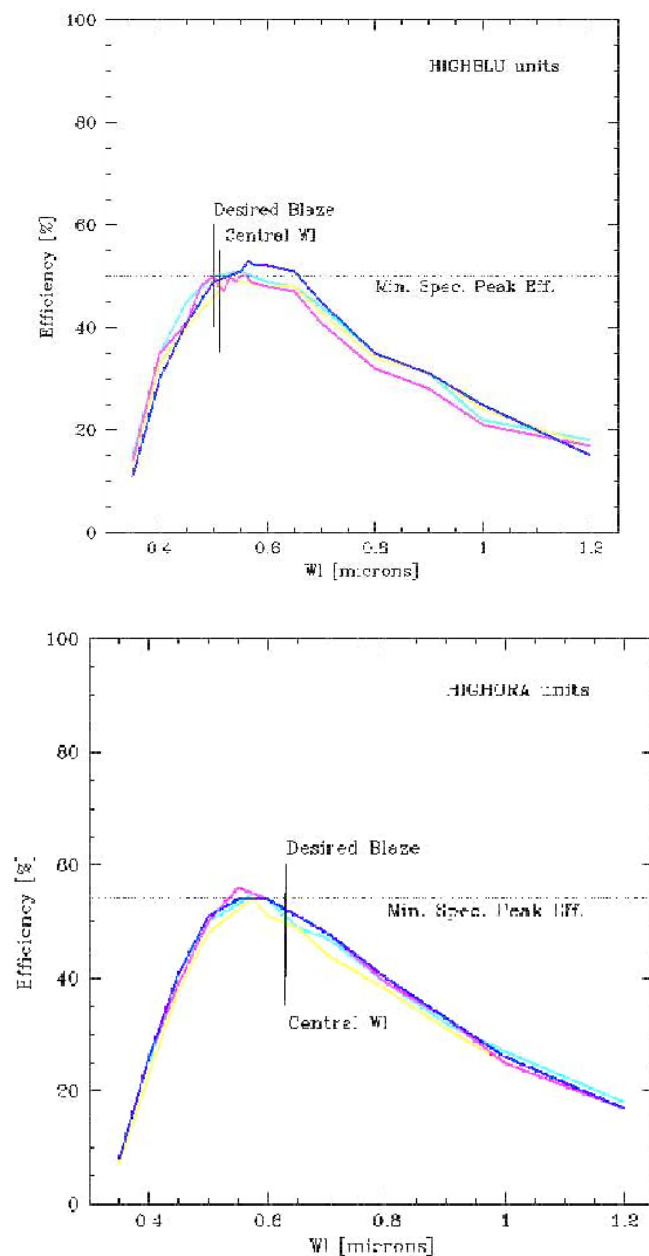


Figure B.3: Efficiency curve of the high resolution grisms HR-blue, HR-orange. The vertical lines mark the desired blaze and the central wavelengths. The spectral range refers to the slit in the center of the field. The cutoff wavelength is in most cases given by the order separation filters, the red CCD limit or the limit of the VIMOS optics in the blue

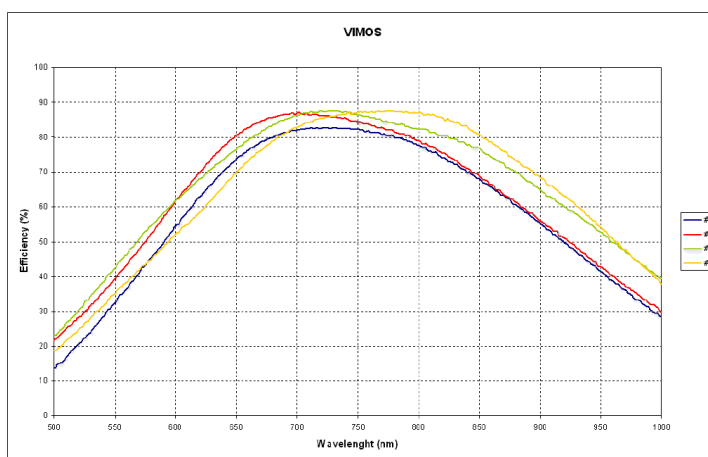


Figure B.4: Efficiency curve of the high resolution grisms HR-red VPHG. The spectral range and efficiency refers to the slit in the center of the field. The cutoff wavelength is in most cases given by the order separation filters, the red CCD limit or the limit of the VIMOS optics in the blue



## Appendix C

# Template Signature Files for Observation and Acquisition.

### C.1 Acquisition Template

The following tables list the parameters of the VIMOS acquisition templates.

| VIMOS_img_acq_Preset |                     |                         |
|----------------------|---------------------|-------------------------|
| Keyword              | Range               | Label in p2pp           |
| INS.FILT.NAME        | U B V R I z         | Filter                  |
| TEL.ROT.OFFANGLE     | -180 .. 180         | Rotator on sky          |
| TEL.TARG.OFFSETALPHA | -3600000 .. 3600000 | Alpha Offset            |
| TEL.TARG.OFFSETDELTA | -3600000 .. 3600000 | Delta Offset            |
| TEL.AG.GUIDESTAR     | SETUPFILE CATALOGUE | Get Guide Star From     |
| TEL.GS1.ALPHA        | 0..240000           | Guide Star RA           |
| TEL.GS1.DELTA        | -900000 .. 900000   | Guide Star Dec          |
| TEL.GS1.PPOS         | POS NEG             | Guide Probe Position    |
| TEL.TARG.ALPHA       |                     | RA of the Field Center  |
| TEL.TARG.DELTA       |                     | DEC of the Field Center |
| TEL.TARG.EQUINOX     | -2000...3000        | Equinox                 |
| TEL.TARG.EPOCH       | -2000...3000        | Epoch                   |
| TEL.TARG.PMA         | -10..10             | proper Motion Alpha     |
| TEL.TARG.PMD         | -10..10             | proper Motion Delta     |
| TEL.TARG.ADDVELALPHA | -15..15             | Additional Velocity RA  |
| TEL.TARG.ADDVELDELTA | -15..15             | Additional Velocity DEC |

### C.2 Observation Template.

### C.3 Template Signature Files for Calibrations

| VIMOS_mos_acq_Mask   |                     |                               |
|----------------------|---------------------|-------------------------------|
| Keyword              | Range               | Label in p2pp                 |
| DET.WIN1.UIT1        | 1..1000000          | Exposure time (seconds)       |
| INS.FILT.NAME        | U B V R I z         | Filter                        |
| TEL.AG.GUIDESTAR     | SETUPFILE CATALOGUE | Get Guide Star From           |
| TEL.GS1.ALPHA        | 0..240000           | Guide Star RA                 |
| TEL.GS1.DELTA        | -900000 .. 900000   | Guide Star Dec                |
| INS.ADP1             | file name           | Name of 1st quadrant ADP File |
| INS.ADP2             | file name           | Name of 2nd quadrant ADP File |
| INS.ADP3             | file name           | Name of 3rd quadrant ADP File |
| INS.ADP4             | file name           | Name of 4th quadrant ADP File |
| TEL.TARG.EQUINOX     | -2000...3000        | Equinox                       |
| TEL.TARG.EPOCH       | -2000...3000        | Epoch                         |
| TEL.TARG.PMA         | -10..10             | proper Motion Alpha           |
| TEL.TARG.PMD         | -10..10             | proper Motion Delta           |
| TEL.TARG.ADDVELALPHA | -15..15             | Additional Velocity RA        |
| TEL.TARG.ADDVELDELTA | -15..15             | Additional Velocity DEC       |

| VIMOS_ifu_acq_Preset |                                                 |                                     |
|----------------------|-------------------------------------------------|-------------------------------------|
| Keyword              | Range                                           | Label in p2pp                       |
| DET.WIN1.UIT1        | 1..1000000                                      | Exposure time (seconds)             |
| INS.FILT.NAME        | U B V R I z                                     | Filter                              |
| INS.GRIS.NAME        | OS-red OS-blue GG435 GG475                      | Grism                               |
| TEL.TARG.ALPHA       | LR_blue, LR_red, MR, HR_blue, HR_orange, HR_red | RA of the Field Center              |
| TEL.TARG.DELTA       |                                                 | DEC of the Field Center             |
| TEL.TARG.EQUINOX     | -2000.3000                                      | Equinox                             |
| TEL.TARG.EPOCH       | -2000...3000                                    | Epoch                               |
| TEL.TARG.PMA         | -10..10                                         | proper Motion Alpha                 |
| TEL.TARG.PMD         | -10..10                                         | proper Motion Delta                 |
| TEL.ROT.OFFANGLE     | -180 .. 180                                     | Rotator on sky                      |
| TEL.TARG.OFFSETALPHA | -3600000..3600000                               | Alpha Offset                        |
| TEL.TARG.OFFSETDELTA | -3600000..3600000                               | Delta Offset                        |
| TEL.AG.GUIDESTAR     | SETUPFILE CATALOGUE                             | Get Guide Star From                 |
| TEL.GS1.ALPHA        | 0..240000                                       | Guide Star RA                       |
| TEL.GS1.DELTA        | -900000 .. 900000                               | Guide Star Dec                      |
| INS.IFUE.NAME        | 0.33, 0.67                                      | IFU magnification in (arcsec/fiber) |
| INS.IFUS.MODE        | F T                                             | IFU Shutter ?                       |
| TEL.TARG.ADDVELALPHA | -15..15                                         | Additional Velocity RA              |
| TEL.TARG.ADDVELDELTA | -15..15                                         | Additional Velocity DEC             |

| VIMOS_img_obs_Offset |                   |                                         |
|----------------------|-------------------|-----------------------------------------|
| Keyword              | Range             | Label in p2pp                           |
| DET.WIN1.UIT1        | 1..100000         | Exposure time (seconds)                 |
| INS.FILT.NAME        | U B V R I z       | Filter                                  |
| SEQ.NEXPO            | 1..1000000        | Number of exposure per telescope offset |
| SEQ.NOFF             | 1..1000000        | Number of telescope offsets             |
| SEQ.RETURN           | F T               | Return to Origin ? (T/F)                |
| SEQ.OFFSET.COORDS    | SKY DETECTOR      | Offset Coordinates ?                    |
| SEQ.OFFSET1.LIST     |                   | List of offsets (arcsec) in X or RA     |
| SEQ.OFFSET2.LIST     |                   | List of offsets (arcsec) in Y or DEC    |
| SEQ.CATG             | PRE-IMAGE SCIENCE | Observation Category                    |

| VIMOS_mos_obs_Offset |                                                    |                                                    |
|----------------------|----------------------------------------------------|----------------------------------------------------|
| Keyword              | Range                                              | Label in p2pp                                      |
| DET.WIN1.UIT1        | 1..100000                                          | Exposure time (seconds)                            |
| INS.FILT.NAME        | U B V R I z<br>OS-red OS-blue GG435 GG475          | Filter                                             |
| INS.GRIS.NAME        | LR_blue, LR_red, MR, HR_blue,<br>HR_orange, HR_red | Grism                                              |
| SEQ.NEXPO            | 1..1000000                                         | Number of exposure per telescope offset            |
| SEQ.NOFF             | 1..1000000                                         | Number of telescope offsets                        |
| SEQ.RETURN           | F T                                                | Return to Origin ? (T/F)                           |
| SEQ.OFFSET1.LIST     |                                                    | List of offsets along the slit (arcsec)            |
| SEQ.OFFSET2.LIST     |                                                    | List of offsets perpendicular to the slit (arcsec) |

| VIMOS_ifu_obs_Offset |                                                    |                                         |
|----------------------|----------------------------------------------------|-----------------------------------------|
| Keyword              | Range                                              | Label in p2pp                           |
| DET.WIN1.UIT1        | 1..100000                                          | Exposure time (seconds)                 |
| INS.FILT.NAME        | U B V R I z<br>OS-red OS-blue GG435 GG475          | Filter                                  |
| INS.GRIS.NAME        | LR_blue, LR_red, MR, HR_blue,<br>HR_orange, HR_red | Grism                                   |
| INS.IFUE.NAME        | 0.33, 0.67                                         | IFU magnification in (arcsec/fiber)     |
| INS.IFUS.MODE        | F T                                                | IFU Shutter ?                           |
| SEQ.NEXPO            | 1..1000000                                         | Number of exposure per telescope offset |
| SEQ.NOFF             | 1..1000000                                         | Number of telescope offsets             |
| SEQ.RETURN           | F T                                                | Return to Origin ? (T/F)                |
| SEQ.OFFSET.COORDS    | SKY LENSLET                                        | Offset Coordinates ?                    |
| SEQ.OFFSET1.LIST     |                                                    | List of offsets (arcsec) in X or RA     |
| SEQ.OFFSET2.LIST     |                                                    | List of offsets (arcsec) in Y or DEC    |

| VIMOS_img_cal_Photom |             |                                           |
|----------------------|-------------|-------------------------------------------|
| Keyword              | Range       | Label in p2pp                             |
| DET.WIN1.UIT1        | 1..100000   | Exposure time (seconds)                   |
| INS.FILT.NAME        | U B V R I z | Filter                                    |
| SEQ.NEXPO            | 1..1000000  | Number of exposure per telescope offset ? |
| SEQ.CAL.FLAG         | 1 4         | 1 or 4 quadrants?                         |

| VIMOS_img_cal_Astrom |             |                                           |
|----------------------|-------------|-------------------------------------------|
| Keyword              | Range       | Label in p2pp                             |
| DET.WIN1.UIT1        | 1..100000   | Exposure time (seconds)                   |
| INS.FILT.NAME        | U B V R I z | Filter                                    |
| SEQ.NEXPO            | 1..1000000  | Number of exposure per telescope offset ? |
| SEQ.CAL.FLAG         | 1 4         | 1 or 4 quadrants?                         |

| VIMOS_mos_cal_NightCalib |            |                     |
|--------------------------|------------|---------------------|
| Keyword                  | Range      | Label in p2pp       |
| SEQ.FLATFIELD            | F T        | Night Flat Field?   |
| SEQ.ARC                  | F T        | Night Arc?          |
| SEQ.FLATFIELD.NEXPO      | 1..1000000 | No. of screenflats? |

| VIMOS_mos_acq_Standard |              |                         |
|------------------------|--------------|-------------------------|
| Keyword                | Range        | Label in p2pp           |
| DET.WIN1.UIT1          | 1..1000000   | Exposure time (seconds) |
| INS.FILT.NAME          | U B V R I z  | Filter                  |
| TEL.TARG.EQUINOX       | -2000...3000 | Equinox                 |
| TEL.TARG.EPOCH         | -2000...3000 | Epoch                   |
| TEL.TARG.PMA           | -10..10      | proper Motion Alpha     |
| TEL.TARG.PMD           | -10..10      | proper Motion Delta     |
| TEL.TARG.ADDVELALPHA   | -15..15      | Additional Velocity RA  |
| TEL.TARG.ADDVELDELTA   | -15..15      | Additional Velocity DEC |

| VIMOS_mos_cal_Standard |                                                    |                                           |
|------------------------|----------------------------------------------------|-------------------------------------------|
| Keyword                | Range                                              | Label in p2pp                             |
| DET.WIN1.UIT1          | 1..100000                                          | Exposure time (seconds)                   |
| INS.FILT.NAME          | U B V R I z<br>OS-red OS-blue GG435 GG475          | Filter                                    |
| INS.GRIS.NAME          | LR_blue, LR_red, MR, HR_blue,<br>HR_orange, HR_red | Grism                                     |
| INS.SLIT.NAID          | 0.8,1.0,1.2,1.5, 1.8, 2, 2.5, 5.0                  | 1.0                                       |
| SEQ.NEXPO              | 1..1000000                                         | Number of exposure per telescope offset ? |
| SEQ.CAL.FLAG           | 1 4                                                | 1 or 4 quadrants ?                        |

| VIMOS_ifu_cal_Specphot |                                                    |                                           |
|------------------------|----------------------------------------------------|-------------------------------------------|
| Keyword                | Range                                              | Label in p2pp                             |
| DET.WIN1.UIT1          | 1..100000                                          | Exposure time (seconds)                   |
| INS.FILT.NAME          | U B V R I z<br>OS-red OS-blue GG435 GG475          | Filter                                    |
| INS.GRIS.NAME          | LR_blue, LR_red, MR, HR_blue,<br>HR_orange, HR_red | Grism                                     |
| INS.IFUE.NAME          | 0.33, 0.67                                         | IFU magnification in (arcsec/fiber)       |
| INS.IFUS.MODE          | F T                                                | IFU Shutter ?                             |
| SEQ.NEXPO              | 1..1000000                                         | Number of exposure per telescope offset ? |
| SEQ.CAL.FLAG           | 1 4                                                | 1 or 4 quadrants ?                        |

| VIMOS_ifu_cal_NightCalib |            |                     |
|--------------------------|------------|---------------------|
| Keyword                  | Range      | Label in p2pp       |
| SEQ.FLATFIELD            | F T        | Night Flat Field ?  |
| SEQ.ARC                  | F T        | Night Arc ?         |
| SEQ.FLATFIELD.NEXPO      | 1..1000000 | No. of screenflats? |

## Appendix D

# Abbreviations and Acronyms

The following abbreviations and acronyms are used in this manual:

|               |                                               |
|---------------|-----------------------------------------------|
| ACQ           | Acquisition                                   |
| ADP           | Aperture Definition in Pixels                 |
| ADU           | Analogue-to-Digital Unite                     |
| CCD           | Charge Coupled Device                         |
| DEC           | Declination                                   |
| ESO           | European Southern Observatory                 |
| ETC           | Exposure Time Calculator                      |
| FEU           | Filter Exchange Unit                          |
| FIERA         | Fast Imager Electronic Readout Assembly       |
| FWHM          | Full Width Half Maximum                       |
| GEU           | Grism Exchange Unit                           |
| HR            | High Resolution                               |
| IFU           | Integral Field Unit                           |
| IMG           | Imaging                                       |
| MEU           | Mask Exchange Unit                            |
| MOS           | Multi Object Spectroscopy                     |
| OB            | Observation Block                             |
| OS            | Order Sorting                                 |
| PSF           | Point Spread Function                         |
| P2PP          | Phase 2 Proposal Preparation                  |
| RA            | Right Ascension                               |
| RMS           | Root Mean Square                              |
| RON           | Read Out Noise                                |
| TSF           | Template Signature File                       |
| VIMOS         | Visible Multi-Object Spectrograph             |
| VMMPS         | Vimos Mask-Manufacturing Preparation Software |
| VLT           | Very Large Telescope                          |
| Å             | Ångstrom                                      |
| $e^-$         | Electron                                      |
| cm            | Centimeter                                    |
| h             | Hour                                          |
| kpx           | KiloPixel                                     |
| min           | Minute                                        |
| mm            | Millimeter                                    |
| nm            | Nanometer                                     |
| px            | Pixel                                         |
| s             | Second                                        |
| $\mu\text{m}$ | Micrometer                                    |

